FLUIDIC PRODUCTS AND INFORMATION









For your information, our U.S. Federal Taxpayer Identification number is 01-0736657; European Taxpayer Identification number is DE 146581180.

ORDER PROCESSING/ SHIPPING

Orders placed on our website process and ship as specified during the checkout process. Orders that are placed via fax/email will be subject to a \$100 minimum order quantity and typically process within one business day and ship within two business days. Customers who have established contract pricing terms must place catalog orders directly with their Customer Service representative.

All shipments are FCA unless otherwise noted. For orders not placed via our website, shipping charges will be prepaid and added to the invoice. Or, if requested, we will charge your carrier account number.

TERMS

Credit cards are accepted for orders placed within the United States. Purchase orders may also be used by qualified customers. No international credit card payment is possible. No international shipments are possible. For terms and conditions of sale please visit www.idex-hs.com/legal-notices. Terms and prices are subject to change without notice.

RETURN OF MATERIALS

Returns must be authorized in advance. Please contact us within 30 days of purchase for your Returns Material Authorization (RMA) Number. A restocking fee may apply. To request a Return Materials Authorization for IDEX Health & Science products visit www.idex-hs.com/return.

Product purchased through a Distributor must be returned directly to that Distributor. Please contact them for Return Material Authorization procedures.

DIRECT FROM OUR GLOBAL DISTRIBUTION NETWORK

IDEX Health & Science has established relationships with a global network of Distributors that offer outstanding product support and assistance. In order to meet our customer's needs, and the demands of being a global supplier, we have come to rely on our authorized Distributors to provide value-added service and support of our products. While we are happy to take your order directly, we encourage you to contact your local IDEX Health & Science Distributor.

HOW TO FIND A LOCAL DISTRIBUTOR FOR FLUIDICS PRODUCTS

For a complete distributor listing please visit www.idex-hs.com/distributors or email:

US:

IHSOrders@idexcorp.com

Europe:

CustomerService.hsEurope@idexcorp.com

Asia:

CustomerService.hsAsia@idexcorp.com

Our customer service team will help you contact a distributor in your area.


HOW TO FIND A LOCAL DISTRIBUTOR FOR SEMROCK PRODUCTS

For a complete distributor listing please visit www.semrock.com/international-distributors.aspx

The products in this catalog are intended for use with analytical, biotechnology and diagnostic equipment and accessories. THE MANUFACTURERS OF THE PRODUCTS IN THIS CATALOG SHALL HAVE NO LIABILITY WHATSOEVER DUE TO ANY MISUSE OF THEIR PRODUCTS. Of course, the safe use of our products depends on our customers, since it is you who select and control the protective gear and safety procedures used, as well as the pressures, temperatures, solvents, samples, ventilation, and other variables. Product and material performance ratings are provided as guides only. Individual field tests should be performed by customers to determine safe operating parameters given your particular procedures and use. IDEX Health & Science facilities are certified 150 9001.

Present the state of the state

ORDERING INFORMATION www.biotech.se www.idex-hs.com

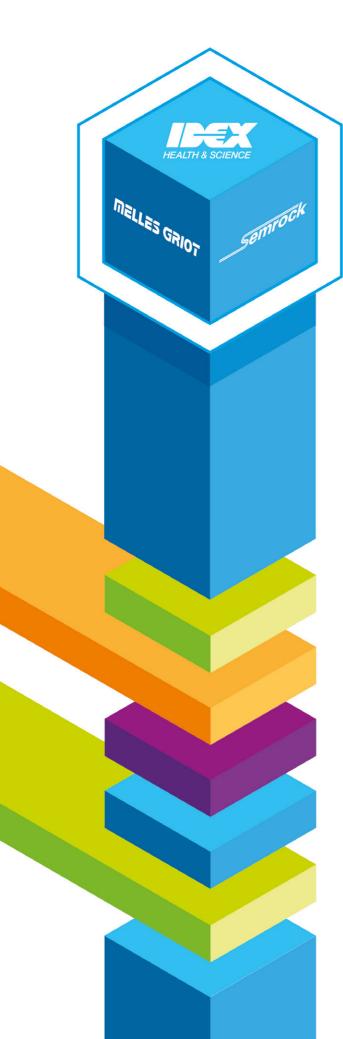

WE HAVE A NEW VISION FOR THE WORLD OF FLUIDICS & OPTICS


As your trusted partner in life science instrument development, we continue to invest in leading technologies and capabilities to solve your most demanding challenges.

For this reason, we have strategically integrated innovative optical companies, like Semrock and Melles Griot, to bring you a newly enhanced portfolio of state-of-the-art components and capabilities that are unrivaled in breadth, performance, quality and design.

Now, as the global authority in fluidics and optics, we deliver a complete spectrum of highly engineered solutions, and more powerful platforms than ever before.



Our vision of the complete path goes far beyond just meeting your needs — it anticipates them, with intelligent solutions for life.

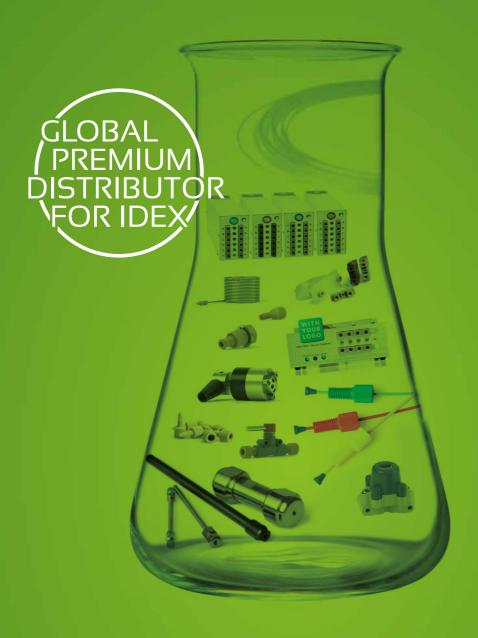

GLOBAL PREMIUM DISTRIBUTOR

Biotech AB in Sweden/Europe (Headquarter), Biotech USA LLC and Bionik Inc Japan are Global Premium Distributors of IDEX Health & Science.

FAST DELIVERY

With several weekly shipments from IDEX manufacturing sites in US, together with our extensive Swedish warehouse, we support the European & Global market with short delivery times as well as fast answers on technical questions and a high level of flexibility.

LOCAL DISTRIBUTORS

In almost every European country we work in close co-operation with local distributors that supports you with excellent daily service, long experience and highest level of technical assistance.

FULL IDEX ASSORTMENT

With the full IDEX assortment of components for fluidic systems along with our long experience and know-how in the field we can assure you to always do our very best in supplying you with excellent solutions and products. Customized OEM solutions, especially in the field of degassing, is something we offer to instrument manufacturers globally.



BIOTECH + IDEX =EVERYTHING FOR THE LAB

KNOW-HOW AND A WIDE ASSORTMENT

As Global Premium Distributor for IDEX we have special know-how about their assortment and can provide you with the right mix of products for a successful lab.

We have cutting edge knowledge to tailor the system for laboratories' specific needs for fluidic handling. You will find a wide range of products to suit every possible flow rate and type of fluid. In this field, we offer a number of well-known brands, such as Systec, Upchurch, Rheodyne and Isolation Technologies - all to give you the best solutions.

GREAT FLUIDIC SOLUTIONS

www.biotechfluidics.com



12 FLUIDIC CONNECTIONS

13	TUBING
16	PEEK Tubing
19	Stainless Steel Tubin
22	PEEKsil™ Tubing
24	PFA Tubing
26	FEP Tubing
27	ETFE Tubing
28	Tools

30	FITTINGS
32	One-Piece Fingertight Fittings
33	Stainless Steel Fittings
34	MicroTight® Fittings
35	Two-Piece Fingertight Fittings
36	Two-Piece SealTight [™] Fingertight Fittings
37	NanoTight [™] Fittings & Sleeves
38	LiteTouch® Fittings
39	$SuperFlangeless^{\scriptscriptstyleTM}Fittings$
42	VacuTight [™] Fittings
43	Metric Flangeless Fittings
44	Flanged Fittings
45	Flangeless Fittings
48	$MarvelX^{\scriptscriptstyleTM}$

50	Fittings Tools
52	Tubing Sleeves
55	Plugs & Caps
56	Large Bore Fittings
57	VHP Micro Fittings
58	VHP PK Fittings
59	Stainless Steel VHP Fittings
60	Reusable VHP Fittings
61	Assorted Fittings Kits

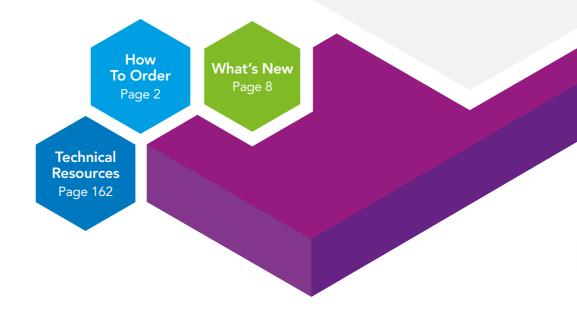


TABLE OF CONTENTS www.biotech.se www.idex-hs.com

62 CONNECTORS

- 63 Threaded Adapters
- 67 High Pressure Mixing Tees
- 69 High Pressure Unions
- 70 High Pressure PEEK
- Tees & Crosses
- 71 High Pressure
- MicroTight® Unions

 72 High Pressure Stainless Steel
 - Tees & Crosses
- 73 VHP Stainless Steel
- ZDV Unions
- 74 VHP Tees & Crosses
- 75 VHP MicroTight® Unions
- 77 MicroTight® Adapters
- 79 Accessories
- 80 NanoPort Assemblies
- 81 Low Pressure Manifolds
- 82 Low Pressure Bulkhead Unions
- 83 Low Pressure Y Connectors
- 84 Low Pressure Unions
- 86 Low Pressure Tees & Crosses
- 87 Quick Connect Luer Adapters
- 89 LuerTight™ Fittings
- 90 Barbed Connectors

94 FILTERS & FRITS

- 95 Stainless Steel Frits
- 97 PEEK Frits
- 99 Frit-in-a-Ferrule™
- 100 Bottom-of-the-Bottle™ Filters
- 102 Inlet Solvent Filters
- 103 Inline Filters
- 105 Precolumn Filters
- 107 Semi-Prep Filters
- 108 Inline MicroFilters
- 109 Mini MicroFilters
- 110 Precolumn MicroFilters
- 111 Bottle Caps
- 112 Bottle Cap Plugs

114 VALVES

- 119 Stand Alone Valve Products
- 120 Rapid Replacement Pods™
- 121 Manual Valves
- 124 Spare Parts
- 126 Valve accessories
- 127 Stainless Steel Sample Loops
- 128 PEEK Sample Loops
- 131 Suction Needle Adapter
- 132 Injection Port Adapters
- 133 Wrenches, Brackets, & Replacement Fittings

134 FLOW REGULATING VALVES

- 135 Inline Cartridge Check Valves
- 136 Standard 1/4-28 Inline Check Valves
- 137 Non-Metallic Check Valves
- 139 Quick-Stop Luer Inline Check Valve
- 140 Micro-Splitter Valves
- 141 Micro-Metering Valves
- 142 Shut-Off Valves

143 BACK PRESSURE REGULATORS

- 144 Ultra-Low Volume Back Pressure Regulators (BPRs)
- 145 Back Pressure Regulator Assemblies
- 146 Back Pressure
- Regulator Holders
- 147 Back Pressure Regulator Cartridges
- 148 Pressure Relief Valves

150 DEGASSERS

- 151 Debubblers
- 154 Full Stand Alone Degassing Systems

156 COLUMN HARDWARE

- 157 Guard Hardware
- 158 Guard Columns



GO WITH A HIGH FLOW INTRODUCING THE ALL NEW DEGASI PREP+

The World's first High Flow Rate, in-line Degasser, for Organic Solvents

Page 155





HIGH

ORGANIC SOLVENTS

CATALOG HINTS



Biocompatible Products

Look for components with the blue "Bio" icon, which designates products that use materials anticipated to maintain the integrity and improve the analysis of biological samples in your intended applications.



Navigation Breadcrumbs

Quickly pinpoint any product or product family within this catalog to the same navigation path each product is categorized on our website. Look to the bottom of every product page for easy-to-follow breadcrumbs to locate specific categories and sub-categories in sequential order, then visit: www.idex-hs.com.



Welcome to Fluidics



We are the premier provider of intelligently engineered fluidic components, assemblies, and integrated solutions for a wide range of life science applications that require precise control and measurement.

Perform precise sample preparation and analysis with our premium fluidic solutions.








Fluidic Connections

IDEX Health & Science has developed a comprehensive line of standard and custom tubing, connectors, fittings, and flow control devices that meet the increasingly demanding requirements of today's high performance analytical fluidic systems. We feature specialty, high-performance polymers and distinct materials designed to work with your system needs. We offer unique products such as biocompatible PEEK-lined stainless steel tubing as well as an assortment of high pressure and fluoropolymer tubing. All of our fittings, filters and frits and connectors come in a variety of materials and styles. We can provide micro and nano-scale dimensions and well as custom forming, assembly and kitting. We also offer our RI detector that provide high resolution and low dispersion detection for HPLC applications.

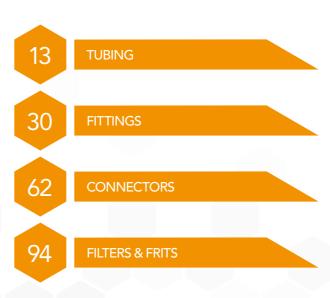





TUBING

Our high quality, versatile tubing is offered in a variety of materials and styles to meet your system requirements. Our high pressure tubing includes biocompatible PEEK selections and well as seamless, pre-cut stainless steel. Our flouropolymer tubing is constructed with genuine Teflon™ FEP and PFA resin, and our unique High Purity PFA. Many of our tubing options are color coded for easy detection and some are translucent making it easy to view the fluid pathway. Our tubing is manufactured to precise tight tolerances to ensure dependable product consistency.

- 16 HIGH PRESSURE TUBING
- **24** FLUOROPOLYMER TUBING
- 28 TOOLS

TUBING OVERVIEW & FUNCTIONS

HIGH PRESSURE TUBING











TUBING:	PEEK	Capillary PEEK	Fused Silica	Stainless Steel	PEEKsil™
Page	16	16	16	19	22
Description	Biocompatible, chemically inert to most commonly used solvents, PEEK tubing is flexible, offers a very smooth internal surface, and can be easily cut to desired lengths. • Great alternative for stainless steel tubing in high pressure applications • Many sizes available in color scheme to help identify ID	All the benefits of larger sized PEEK tubing, while serving as an excellent alternative to more traditional fused silica and stainless steel capillary tubing. Capillary PEEK tubing is available in a wide range of micro and nanoscale inner diameters. • Available in common capillary tubing sizes with tight tolerances on OD and ID • Tubing sleeves available for capillary tubing connections	Because of the tight tolerances of fused silica's inner diameters, this tubing is used for micro-scale analyses such as micro and nano-HPLC and capillary electrophoresis. Most commonly used OD and ID sizes available High quality, polyimide-clad fused silica Offered in convenient, two meter lengths	Seamless, pre-cut 316 stainless steel tubing meets the exacting requirements of today's analyses. Thorough preparation guarantees that the tubing is truly ready-to-use, with flat-bur-free ends and a clean finish. • Wide selection of outside and inside diameters and lengths • Pre-cut to ensure burr-free, flat connections • Many sizes feature a color-coded band for easy ID identification	PEEKsil is mechanically strong and has ideal characteristics for sealing with metal or polymer fittings. • Comprised of high quality fused silica sheathed by PEEK tubing • Excellent chemical compatibility • Very tight manufacturing tolerances • Good replacement for stainless steel, PEEK, or standard fused silica

HIGH PRESSURE T	UBING SPECIFICATION	NS			
OD (outside diameter)	1/32" (785 µm), 1/16" (1.55 mm), 1/8" (3.2 mm)	0.0145" (360 μm), 1/32" (785 μm), 0.020" (0.5 mm)	0.0145" (360 µm)	1/32" (785 μm), 1/16" (1.55 mm), 1/8" (3.2 mm)	0.0145" (360 μm), 1/32" (785 μm), 1/16" (1.55 mm)

ID (inside diameter)	0.001" (25 μm)– 0.080" (2.0 mm)	0.001" (25 μm)– 0.020" (0.50 mm)	0.0008" (20 μm)–0.006" (150 μm)	0.004" (100 μm)– 0.080" (2.0 mm)	0.001" (25 μm)– 0.012" (300 μm)
Operating Temp	-51 to 100 °C	-51 to 100 °C	-51 to 100 °C	-51 to 289 °C	-51 to 100 °C
Pressure Rating	500–10,000 psi (34–690 bar)	2,000–5,000 psi (138–345 bar)	N/A*	N/A*	10,000 psi (690 bar)
Typical Tolerances	±0.001" (25 μm) for 1/16" OD tubing; ±0.003" (75 μm) for 1/8" OD tubing	±0.0005" (12.5 μm)	±0.0004" (10 μm)	±0.001" (25 µm) for 1/16" OD tubing; ±0.003" (75 µm) for 1/8" OD tubing	±0.0004" (10 μm)
Refractive Index (Clarity)	Opaque	Opaque	1.78	Opaque	Opaque
pH Range	0–14	0–14	0–10	1–14	0–10
Sterilization Techniques	Gamma irradiation; ethylene oxide; thermal	Gamma irradiation; ethylene oxide; thermal	Ethylene oxide; thermal	Gamma irradiation; ethylene oxide; thermal	Ethylene oxide; thermal
Autoclavable?	Υ	Υ	Υ	Υ	Υ
*The manufacturer of this tubing or material does not publish this execution.					

TUBING OD SIZES

Please use this table as a reference tool to help quickly locate within this chapter the appropriate OD size tubing for your application.

1 mm OD

Page(s)

2 mm OD

Page(s)

3 mm OD



Page(s) 26

FLUOROPOLYMER TUBING



TUBING:	PFA	High Purity PFA	360 µm High Purity PFA	FEP	ETFE
Page	24	24	24	26	27
Description	Offers excellent chemical compatibility, plus due to its inner surface smoothness, PFA tubing tends to be more translucent than PTFE tubing. • Offers higher purity and enhanced translucence when compared with other fluoropolymer tubes • Great for more critical, low pressure applications	This polymer tubing is manufactured from a premium grade of PFA — one of the most contaminant-free polymers on the market. • Offers chemical stability, mechanical strength, and purity for applications such as medical, diagnostic, pharmaceutical, biotechnology, and semiconductor • Excellent replacement for PTFE where gas permeability and surface texture are issues • Clarity of tubing makes PFA an excellent choice for monitoring fluid movement	This tubing offers excellent chemical compatibility, transparency, very low contaminant levels and is available in the most commonly-used outside diameter for capillary tubing applications. • Replacement for capillary tubing in low pressure applications where excellent chemical compatibility is required. • Tubing sleeves available for capillary tubing connections	FEP tubing is a great alternative to traditional PTFE tubing, desirable for use because it is chemically inert to most solvents, easy to cut, and translucent for easy monitoring of solutions passing through. • Great for general, low pressure applications • Many sizes available in multiple colors for easy identification • Tight manufacturing tolerances to ensure product consistency	ETFE is chemically inert and more suitable for higher pressure applications (when using aqueous mobile phases) than PTFE, FEP, and PFA. Additionally, because ETFE is more rigid than PTFE, FEP, and PFA, this tubing better resists inner diameter collapse. • Excellent solvent resistance • More durable and less gas permeable than PTFE, FEP, and PFA • Operating temperatures up to 80 °C

FLUOROPOLYMER TUBING SPECIFICATIONS					
OD (outside diameter)	1/16" (1.55 mm), 1/8" (3.2 mm)	1/16" (1.55 mm), 1/8" (3.2 mm), 3/16" (4.8 mm), 1/4" (6.35 mm)	0.0145" (360 μm)	1/16" (1.55 mm), 0.080" (2.0 mm), 0.118" (3.0 mm), 1/8" (3.2 mm), 0.157" (4.0 mm), 3/16" (4.8 mm), 1/4" (6.35 mm), 5/16" (7.94 mm)	1/16" (1.6 mm), 1/8" (3.2 mm), 1/4" (6.35 mm)
ID (inside diameter)	0.020" (0.50 mm)– 0.062" (1.55 mm)	0.020" (0.50 mm)- 0.188" (4.80 mm)	0.002" (50 μm)– 0.006" (150 μm)	0.003" (0.075 mm) – 0.250" (6.35 mm)	0.010" (0.25 mm)– 0.188" (4.80 mm)
Operating Temp	-51 to 80 °C	-51 to 80 °C	-51 to 80 °C	-51 to 50 °C	-51 to 80 °C
Pressure Rating	500–2,000 psi (34–138 bar)	250–2,000 psi (17–138 bar)	1,750–3,500 psi (121–241 bar)	2,500–4,000 psi (172 - 276 bar)	250–4,000 psi (17–276 bar)
Typical Tolerances	±0.001" (25 µm) for 1/16" OD tubing; ±0.003" (75 µm) for 1/8" OD tubing	±0.001" (25 μm) or 1/16" OD tubing	±0.0005" (12.5 μm)	±0.001" (25 µm) for 1/16" OD tubing; ±0.003" (75 µm) for 1/8" OD tubing	±0.001" (25 µm) for 1/16" OD tubing; ±0.003" (75 µm) for 1/8" OD tubing
Refractive Index (Clarity)	1.34	1.34	1.34	1.338	1.4
pH Range	0–14	0–14	0–14	0–14	0–14
Sterilization Techniques	Ethylene oxide; thermal	Gamma irradiation; ethylene oxide; thermal	Gamma irradiation; ethylene oxide; thermal	Ethylene oxide; thermal	Ethylene oxide
Autoclavable?	Υ	Υ	Υ	Υ	Υ





PEEK Tubing

- ▶ 1/16" or 1/8" outside diameter available
- > Biocompatible, inert, and easily cut
- > Great for high pressure applications
- Maximum continuous use temperature: 100 °C

Our PEEK (polyetheretherketone) polymer tubing is biocompatible, chemically inert to most solvents, and can be used to replace stainless steel tubing in most liquid analytical systems. Unlike stainless steel tubing, PEEK tubing is flexible and can be easily cut to desired lengths.

PEEK tubing has a very smooth internal surface, which causes less turbulence than similarly sized metal tubing, contributing to improved resolution of sample bands. Of all our polymer tubing materials, PEEK is the least permeable to gas (see material properties on our website: www.idex-hs.com).

In addition, much of our 1/16" OD tubing is color-coded so different IDs are easily identified. Our proprietary extrusion process ensures color permanence in our tubing.

Our 5' length tubing is rough cut to approximately 5'1". A trim cut should be made before use, especially for smaller ID tubing. PEEK tubing can be cut easily with a razor blade. However for an improved cut, try our Tubing Cutters on page 28.



Capillary PEEK Tubing

- ▶ 360 µm or 1/32" outside diameter available
- IDs as small as 25 μm (0.001")

Capillary PEEK tubing offers all the benefits of larger sized PEEK tubing, while serving as an excellent alternative to more traditional fused silica and stainless steel capillary tubing (see Application Note, right). The capillary tubing can be coupled to many of the products in the Connectors chapter (starting on page 62) and to some of the valves in the Valves chapter (starting on page 114).



Fused Silica Tubing

- ightharpoonup Five inner diameters with most common capillary outside diameter, 360 μm
- > Cut in convenient lengths, up to 2 m

These products are manufactured from synthetic fused silica with a polyimide coating.



Because the thru-hole of our $25 \, \mu m$ ID PEEK tubing is very small, it is possible for some fittings to cause the ID to become occluded. Please use caution, especially with wrench-tightened fittings. For more information, please contact IDEX Health & Science or your local Distributor directly.



What Size PEEK Tubing Should I Use?

- ▶ It is usually safe to use 1/16" OD x 0.010" ID tubing throughout an analytical HPLC system. With a 0.010" ID, the pressure drop across most tubing lengths is negligible, and the ID is small enough to minimize band broadening.
- High pressure semi-prep LC systems will most likely use 1/8" OD tubing.
- Use our 1/32" OD tubing for the high pressure flow path of some microbore HPLC systems.
- ➤ Choose 360 µm OD tubing for most capillary systems.
- PEEK tubing is also available by the inch. Contact your local Distributor or IDEX Health & Science directly for pricing information.

APPLICATION NOTE

- ▶ An independent study conducted by a major pharmaceutical company indicated LC-MS chromatographic performance could be improved in some cases by switching the post-column transfer line from fused silica to PEEK polymer tubing. The switch dramatically reduced peak tailing and eliminated the degradation of peak symmetry as injection volume was reduced. For more information, please contact us or order the "Improved LC-MS Results Study" from the "Literature Request" section of our website at www.idex-hs.com.
- To straighten PEEK polymer tubing, first choose a piece of stainless steel tubing with an inner diameter slightly larger than the OD of your tubing and with an appropriate length for the PEEK tubing you wish to straighten. For instance, for 1/16" OD PEEK tubing with a length of 10", choose our U-825 tubing (stainless steel, 1/8" OD x 0.080" ID x 25 cm long, page 19. Slip your PEEK tubing into the stainless steel tubing. Place this "sleeved" tubing into an oven and bake at 425 °F (218 °C) for 30 minutes or 350 °F (177 °C) for 60 minutes. Allow the sleeved tubing to return to room temperature naturally (i.e., do not quench it with water). Once cooled, remove the PEEK tubing from the stainless steel sleeve and inspect it for straightness. If needed, repeat the process until the desired straightness is achieved.

SPECIFICATIONS & DETAILS

Tubing OD	Tubing ID	OD Tolerance	ID Tolerance
PEEK TUBING S	PECIFICATIONS		
1/16"	25 μm	±0.001" (25 μm)	±0.0005" (12.5 μm)
1/8"	All	±0.003" (75 μm)	±0.003" (75 μm)
CAPILLARY PEE	K TUBING SPECIFICATIONS		
360 µm	All	±0.0005" (12.5 μm)	±0.0005" (12.5 μm)
1/32"	All	±0.0005" (12.5 μm)	±0.0005" (12.5 μm)
FUSED SILICA T	UBING, 360 μm OD		
360 µm	20 μm (0.0008")	±0.0004" (10 μm)	±0.00008" (2 μm)
360 µm	50 μm (0.002") and 75 μm (0.003")	±0.0004" (10 μm)	±0.00012" (3 μm)
360 µm	$100~\mu m$ (0.004") and 150 μm (0.006")	±0.0004" (10 μm)	±0.00016" (4 µm)

PEEK Tubing (Cont.)

PEEK Tubing

N NI	ID.	C-I	Mary Dayson	0+
Part No. PEEK TUBING	ID : 1/14" OD	Color	Max. Pressure	Qt
PEER TUBING 1560	0.0025" (65 μm) ID x 5' (1.5 m)	Natural	7,000 psi (483 bar)	
1560L		Natural		ea.
1560XL	0.0025" (65 μm) ID x 50' (15 m)		7,000 psi (483 bar)	ea.
	0.0025" (65 μm) ID x 100' (30 m)	Natural	7,000 psi (483 bar)	ea.
560M	0.0025" (65 μm) ID x 1,000' (304 m)	Natural	7,000 psi (483 bar)	ea.
561	0.004" (0.10 mm) ID x 5' (1.5 m)	Black	7,000 psi (483 bar)	ea.
561L	0.004" (0.10 mm) ID x 50' (15 m)	Black	7,000 psi (483 bar)	ea.
561XL	0.004" (0.10 mm) ID x 100' (30 m)	Black	7,000 psi (483 bar)	ea.
561M	0.004" (0.10 mm) ID x 1,000' (304 m)	Black	7,000 psi (483 bar)	ea.
535	0.005" (0.125 mm) ID x 5' (1.5 m)	Red	7,000 psi (483 bar)	ea.
535L	0.005" (0.125 mm) ID x 50' (15 m)	Red	7,000 psi (483 bar)	ea.
535XL	0.005" (0.125 mm) ID x 100' (30 m)	Red	7,000 psi (483 bar)	ea.
535M	0.005" (0.125 mm) ID x 1,000' (304 m)	Red	7,000 psi (483 bar)	ea.
536	0.007" (0.175 mm) ID x 5' (1.5 m)	Yellow	7,000 psi (483 bar)	
				ea.
536L	0.007" (0.175 mm) ID x 50' (15 m)	Yellow	7,000 psi (483 bar)	ea.
536XL	0.007" (0.175 mm) ID x 100' (30 m)	Yellow	7,000 psi (483 bar)	ea.
536M	0.007" (0.175 mm) ID x 1,000' (304 m)	Yellow	7,000 psi (483 bar)	ea.
531	0.010" (0.25 mm) ID x 5' (1.5 m)	Natural	7,000 psi (483 bar)	ea.
531L	0.010" (0.25 mm) ID x 50' (15 m)	Natural	7,000 psi (483 bar)	ea.
531XL	0.010" (0.25 mm) ID x 100' (30 m)	Natural	7,000 psi (483 bar)	ea.
531M	0.010" (0.25 mm) ID x ID x 1,000' (304 m)	Natural	7,000 psi (483 bar)	ea.
531B	0.010" (0.25 mm) ID x 5' (1.5 m)	Blue	7,000 psi (483 bar)	ea.
531BL				
	0.010" (0.25 mm) ID x 50' (15 m)	Blue	7,000 psi (483 bar)	ea.
531BXL	0.010" (0.25 mm) ID x 100' (30 m)	Blue	7,000 psi (483 bar)	ea.
531BM	0.010" (0.25 mm) ID x 1,000' (304 m)	Blue	7,000 psi (483 bar)	ea
532	0.020" (0.50 mm) ID x 5' (1.5 m)	Orange	6,000 psi (414 bar)	ea.
532L	0.020" (0.50 mm) ID x 50' (15 m)	Orange	6,000 psi (414 bar)	ea.
532XL	0.020" (0.50 mm) ID x 100' (30 m)	Orange	6,000 psi (414 bar)	ea
532M	0.020" (0.50 mm) ID x 1,000' (304 m)	Orange	6,000 psi (414 bar)	ea
533	0.030" (0.75 mm) ID x 5' (1.5 m)	Green	4,000 psi (276 bar)	ea.
533L	0.030" (0.75 mm) ID x 50' (15 m)	Green	4,000 psi (276 bar)	ea
533XL	0.030" (0.75 mm) ID x 100' (30 m)	Green	4,000 psi (276 bar)	ea
533M	0.030" (0.75 mm) ID x 1,000' (304 m)	Green	4,000 psi (276 bar)	ea
538	0.040" (1.00 mm) ID x 5' (1.5 m)	Natural	3,000 psi (207 bar)	ea
538L	0.040" (1.00 mm) ID x 50' (15 m)	Natural	3,000 psi (207 bar)	ea.
538XL	0.040" (1.00 mm) ID x 100' (30 m)	Natural	3,000 psi (207 bar)	ea
538M	0.040" (1.00 mm) ID x 1,000' (304 m)	Natural	3,000 psi (207 bar)	ea.
PEEK TUBING				
534	•	Natural	4 000 mai (274 han)	
	0.062" (1.55 mm) ID x 5' (1.5 m)		4,000 psi (276 bar)	ea.
544	0.080" (2.00 mm) ID x 5' (1.5 m)	Natural	3,000 psi (207 bar)	ea.
CAPILLARY PE	EEK TUBING, 360 μm OD			
574	25 μm (0.001") ID x 5' (1.5 m)	Natural	5,000 psi (345 bar)	ea.
570	50 μm (0.002") ID x 5' (1.5 m)	Natural	2,000 psi (138 bar)	ea.
571	100 μm (0.004") ID x 5' (1.5 m)	Red	2,000 psi (138 bar)	ea.
572	150 µm (0.006") ID x 5' (1.5 m)	Yellow	2,000 psi (138 bar)	ea.
	EEK TUBING, 1/32" OD		_,	
	•		E 000 : (24E L)	
576	0.005" (0.125 mm) ID x 5' (1.5 m)	Red	5,000 psi (345 bar)	ea.
576L	0.005" (0.125 mm) ID x 50' (15 m)	Red	5,000 psi (345 bar)	ea
576XL	0.005" (0.125 mm) ID x 100' (30 m)	Red	5,000 psi (345 bar)	ea
576M	0.005" (0.125 mm) ID x 1,000' (304 m)	Red	5,000 psi (345 bar)	ea
577	0.007" (0.175 mm) ID x 5' (1.5 m)	Yellow	5,000 psi (345 bar)	ea
577L	0.007" (0.175 mm) ID x 50' (15 m)	Yellow	5,000 psi (345 bar)	ea
577L	0.007" (0.175 mm) ID x 100' (30 m)	Yellow	5,000 psi (345 bar)	ea
577M	0.007" (0.175 mm) ID x 1,000' (304 m)	Yellow	5,000 psi (345 bar)	ea
581	0.010" (0.25 mm) ID x 5' (1.5 m)	Blue	5,000 psi (345 bar)	ea
581L	0.010" (0.25 mm) ID x 50' (15 m)	Blue	5,000 psi (345 bar)	ea
581XL	0.010" (0.25 mm) ID x 100' (30 m)	Blue	5,000 psi (345 bar)	ea
581M	0.010" (0.25 mm) ID x 1,000' (304 m)	Blue	5,000 psi (345 bar)	ea
568	0.015" (0.40 mm) ID x 5' (1.5 m)	Natural	4,000 psi (276 bar)	ea
568L	0.015" (0.40 mm) ID x 50' (15 m)	Natural	4,000 psi (276 bar)	ea
568XL		Natural	4,000 psi (276 bar)	
	0.015" (0.40 mm) ID x 100' (30 m)			ea
568M	0.015" (0.40 mm) ID x 1,000' (304 m)	Natural	4,000 psi (276 bar)	ea
569	0.020" (0.50 mm) ID x 5' (1.5 m)	Orange	3,000 psi (207 bar)	ea
569L	0.020" (0.50 mm) ID x 50' (15 m)	Orange	3,000 psi (207 bar)	ea
569XL	0.020" (0.50 mm) ID x 100' (30 m)	Orange	3,000 psi (207 bar)	ea
569M	0.020" (0.50 mm) ID x 1,000' (304 m)	Orange	3,000 psi (207 bar)	ea
	TUBING, 360 µm OD			
	•	Natural	10,000 pci /400 bcs/	
S-120	20 μm (0.0008") ID x 6.4' (2 m)	Natural	10,000 psi (690 bar)	ea
S-150	50 μm (0.002") ID x 6.4' (2 m)	Natural	10,000 psi (690 bar)	ea.
S-175	75 μm (0.003") ID x 6.4' (2 m)	Natural	10,000 psi (690 bar)	ea.
S-110	100 μm (0.004") ID x 6.4' (2 m)	Natural	10,000 psi (690 bar)	ea.
			10,000 psi (690 bar)	



Stainless Steel Tubing

- > Precut 316 stainless steel
- Available ODs include 1/32", 1/16", and 1/8"
- Color-coded banding for easy identification of the inner diameter

IDEX Health & Science seamless, precut stainless steel tubing is designed to meet the exacting requirements of today's analyses. We machine cut and polish each end, deburr the inside and outside edges, and passivate the tubing (please see the passivation information on this page). Finally, we flush reagent-grade isopropanol through each piece.

Our thorough preparation and cleaning procedure guarantees tubing that is truly ready-to-use, with flat, burr-free ends and a clean finish. This care is important in achieving zero-dead-volume connections and good chromatographic results.

We offer a variety of precut lengths as well as longer lengths (5' and 25') of some sizes. Cutting the tubing disturbs and roughens the tubing's end surface, so we recommend using our precut tubing whenever possible. If you need to cut tubing to custom lengths, we suggest you then passivate the tubing.



- Our 1/32" OD tubing is designed for enhanced flexibility in high pressure applications.
- Standard 1/16" and 1/8" OD stainless steel tubing is suited for most analytical applications.



NOTE

PEEK polymer tubing can be used to replace stainless steel tubing in most liquid analytical systems. Unlike stainless steel tubing, PEEK tubing is biocompatible, flexible, and can easily be cut to desired lengths. See page 16. All Stainless Steel tubing longer than 1 m is coiled.

The Beauty of Precut Tubing







Precut tubing

Tubing cut by a commercially available tubing cutter

File cut tubing

Stainless Steel Tubing Passivation

Stainless steel is naturally self-passivating, forming an oxidized layer on newly created surfaces. IDEX Health & Science takes extra steps to ensure the chemical resistance of our stainless steel tubing by manually passivating before and after the tubing is cut into specified lengths (except in a few cases where size is prohibitive). In the precut stage, the internal wall is acid passivated and flushed. After the tubing is cut, deburred and polished, it is completely submerged in an acid passivation bath and again flushed clean. The table below summarizes the manual passivation steps performed for each size of our stainless steel tubing:

Tubing OD	Precut Passivation	Postcut Passivation
1/32"	All	All
1/16"	All	All, ex. 25' lengths
1/8"	None	All, ex. 3 & 5 m lengths

Stainless Steel Tubing (Cont.)




SPECIFICATIONS & DETAILS

- Maximum Recommended Operating Temperature: 750 °F (399 °C).
- Rockwell Hardness (B): Maximum of 95.
- Meets ASTM A269 and A213.

Tubing OD	OD Tolerance	Tubing ID	ID Tolerance
1/32"	+0.002"/-0.000" (+50 μm/-0 μm)	All, except 0.004" (0.10 mm)	+0.000"/-0.002" (+0 µm/-50 µm)
1/32"	+0.002"/-0.000" (+50 μm/-0 μm)	0.004" (0.10 mm)	+0.002"/-0.000" (+50 μm/-0 μm)
1/16"	+0.002"/-0.000" (+50 μm/-0 μm)	All	±0.001" (25 μm)
1/8"	±0.003" (75 μm)	All	±0.003" (75 μm)



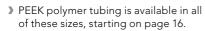



NOTE

Understanding the Maximum Pressure Value of Stainless Steel Tubing

Stainless steel is unique as a material. The Maximum Pressure value listed for each part number is the safe, continuous working pressure limit that IDEX Health & Science has assigned for the tubing. It reflects a safety margin before the tubing begins to "yield" — which is well below the tubing's "burst" pressure. For more information, contact IDEX Health & Science or your authorized Distributor.

Stainless Steel Tubing

Part No.	ID	Length	Color	Maximum Pressure	Qty.
STAINLESS STEEL,	. 1/32" OD				
U-1114	0.004" (0.10 mm)	2" (5 cm)	Red	19,300 psi (1,331 bar)	ea.
U-1115	0.004" (0.10 mm)	4" (10 cm)	Red	19,300 psi (1,331 bar)	ea.
U-1116	0.004" (0.10 mm)	8" (20 cm)	Red	19,300 psi (1,331 bar)	ea.
U-1117	0.004" (0.10 mm)	12" (30 cm)	Red	19,300 psi (1,331 bar)	ea.
U-1120	0.006" (0.15 mm)	2" (5 cm)	Yellow	19,300 psi (1,331 bar)	ea.
U-1122	0.006" (0.15 mm)	8" (20 cm)	Yellow	19,300 psi (1,331 bar)	ea.
U-1125	0.008" (0.20 mm)	2" (5 cm)	Clear	17,800 psi (1,227 bar)	ea.
U-1126	0.008" (0.20 mm)	4" (10 cm)	Clear	17,800 psi (1,227 bar)	ea.
U-1128	0.008" (0.20 mm)	12" (30 cm)	Clear	17,800 psi (1,227 bar)	ea.
U-1130	0.010" (0.25 mm)	2" (5 cm)	Blue	16,200 psi (1,117 bar)	ea.
U-1131	0.010" (0.25 mm)	4" (10 cm)	Blue	16,200 psi (1,117 bar)	ea.
U-1132	0.010" (0.25 mm)	8" (20 cm)	Blue	16,200 psi (1,117 bar)	ea.
U-1133	0.010" (0.25 mm)	12" (30 cm)	Blue	16,200 psi (1,117 bar)	ea.
U-1140	0.015" (0.40 mm)	2" (5 cm)	Green	12,300 psi (848 bar)	ea.
U-1141	0.015" (0.40 mm)	4" (10 cm)	Green	12,300 psi (848 bar)	ea.
U-1142	0.015" (0.40 mm)	8" (20 cm)	Green	12,300 psi (848 bar)	ea.
U-1143	0.015" (0.40 mm)	12" (30 cm)	Green	12,300 psi (848 bar)	ea.
U-1145	0.018" (0.45 mm)	2" (5 cm)	Black	10,000 psi (689 bar)	ea.
U-1146	0.018" (0.45 mm)	4" (10 cm)	Black	10,000 psi (689 bar)	ea.
U-1148	0.018" (0.45 mm)	12" (30 cm)	Black	10,000 psi (689 bar)	ea.

Stainless Steel Tubing (Cont.) Part No. ID

Part No.	ID	Length	Color	Maximum Pressure	Qty.
STAINLESS STEEL,	1/16" OD				
U-152	0.005" (0.125 mm)	2" (5 cm)	Red	21,600 psi (1,489 bar)	ea.
U-153	0.005" (0.125 mm)	4" (10 cm)	Red	21,600 psi (1,489 bar)	ea.
U-154	0.005" (0.125 mm)	8" (20 cm)	Red	21,600 psi (1,489 bar)	ea.
U-155	0.005" (0.125 mm)	12" (30 cm)	Red	21,600 psi (1,489 bar)	ea.
U-156	0.005" (0.125 mm)	1.6' (0.5 m)	Red	21,600 psi (1,489 bar)	ea.
U-157	0.005" (0.125 mm)	3.2' (1 m)	Red	21,600 psi (1,489 bar)	ea.
U-158	0.005" (0.125 mm)	5' (1.5 m)	Red	21,600 psi (1,489 bar)	ea.
U-160	0.005" (0.125 mm)	25' (7.6 m)	Red	21,600 psi (1,489 bar)	ea.
U-126	0.007" (0.175 mm)	2" (5 cm)	Black	20,900 psi (1,441 bar)	ea.
U-127	0.007" (0.175 mm)	4" (10 cm)	Black	20,900 psi (1,441 bar)	ea.
U-128	0.007" (0.175 mm)	8" (20 cm)	Black	20,900 psi (1,441 bar)	ea.
U-129	0.007" (0.175 mm)	12" (30 cm)	Black	20,900 psi (1,441 bar)	ea.
U-130	0.007" (0.175 mm)	1.6' (0.5 m)	Black	20,900 psi (1,441 bar)	ea.
U-131	0.007" (0.175 mm)	3.2' (1 m)	Black	20,900 psi (1,441 bar)	ea.
U-108	0.007" (0.175 mm)	5' (1.5 m)	Black	20,900 psi (1,441 bar)	ea.
U-161	0.007" (0.175 mm)	25' (7.6 m)	Black	20,900 psi (1,441 bar)	ea.
U-111	0.010" (0.25 mm)	2" (5 cm)	Blue	19,700 psi (1,358 bar)	ea.
U-112	0.010" (0.25 mm)	4" (10 cm)	Blue	19,700 psi (1,358 bar)	ea.
U-113	0.010" (0.25 mm)	8" (20 cm)	Blue	19,700 psi (1,358 bar)	ea.
U-114	0.010" (0.25 mm)	12" (30 cm)	Blue	19,700 psi (1,358 bar)	ea.
U-132	0.010" (0.25 mm)	1.6' (0.5 m)	Blue	19,700 psi (1,358 bar)	ea.
U-133	0.010" (0.25 mm)	3.2' (1 m)	Blue	19,700 psi (1,358 bar)	ea.
U-106	0.010" (0.25 mm)	5′ (1.5 m)	Blue	19,700 psi (1,358 bar)	ea.
U-162	0.010" (0.25 mm)	25' (7.6 m)	Blue	19,700 psi (1,358 bar)	ea.
U-101	0.010 (0.23 mm)	2" (5 cm)	Yellow	15,800 psi (1,089 bar)	
U-102	0.020" (0.5 mm)	4" (10 cm)	Yellow	15,800 psi (1,089 bar)	ea.
U-103	0.020" (0.5 mm)	8" (20 cm)	Yellow	15,800 psi (1,089 bar)	
U-104		, ,	Yellow		ea.
U-134	0.020" (0.5 mm)	12" (30 cm)	Yellow	15,800 psi (1,089 bar)	ea.
U-135	0.020" (0.5 mm)	1.6' (0.5 m) 3.2' (1 m)		15,800 psi (1,089 bar)	ea.
	0.020" (0.5 mm)		Yellow	15,800 psi (1,089 bar)	ea.
U-105 U-163	0.020" (0.5 mm)	5′ (1.5 m)	Yellow	15,800 psi (1,089 bar)	ea.
	0.020" (0.5 mm)	25' (7.6 m)	Yellow	15,800 psi (1,089 bar)	ea.
U-115	0.030" (0.75 mm)	2" (5 cm)	White	12,000 psi (827 bar)	ea.
U-116	0.030" (0.75 mm)	4" (10 cm)	White	12,000 psi (827 bar)	ea.
U-117	0.030" (0.75 mm)	8" (20 cm)	White	12,000 psi (827 bar)	ea.
U-118	0.030" (0.75 mm)	12" (30 cm)	White	12,000 psi (827 bar)	ea.
U-136	0.030" (0.75 mm)	1.6' (0.5 m)	White	12,000 psi (827 bar)	ea.
U-137	0.030" (0.75 mm)	3.2′ (1 m)	White	12,000 psi (827 bar)	ea.
U-107	0.030" (0.75 mm)	5′ (1.5 m)	White	12,000 psi (827 bar)	ea.
U-164	0.030" (0.75 mm)	25' (7.6 m)	White	12,000 psi (827 bar)	ea.
U-138	0.040" (1.0 mm)	2" (5 cm)	N/A	8,100 psi (558 bar)	ea.
U-139	0.040" (1.0 mm)	4" (10 cm)	N/A	8,100 psi (558 bar)	ea.
U-140	0.040" (1.0 mm)	8" (20 cm)	N/A	8,100 psi (558 bar)	ea.
U-141	0.040" (1.0 mm)	12" (30 cm)	N/A	8,100 psi (558 bar)	ea.
U-142	0.040" (1.0 mm)	1.6' (0.5 m)	N/A	8,100 psi (558 bar)	ea.
U-143	0.040" (1.0 mm)	3.2' (1 m)	N/A	8,100 psi (558 bar)	ea.
U-144	0.040" (1.0 mm)	5′ (1.5 m)	N/A	8,100 psi (558 bar)	ea.
U-165	0.040" (1.0 mm)	25' (7.6 m)	N/A	8,100 psi (558 bar)	ea.
U-145	0.046" (1.15 mm)	2" (5 cm)	N/A	5,800 psi (400 bar)	ea.
U-146	0.046" (1.15 mm)	4" (10 cm)	N/A	5,800 psi (400 bar)	ea.
U-147	0.046" (1.15 mm)	8" (20 cm)	N/A	5,800 psi (400 bar)	ea.
U-148	0.046" (1.15 mm)	12" (30 cm)	N/A	5,800 psi (400 bar)	ea.
U-149	0.046" (1.15 mm)	1.6' (0.5 m)	N/A	5,800 psi (400 bar)	ea.
U-150	0.046" (1.15 mm)	3.2' (1 m)	N/A	5,800 psi (400 bar)	ea.
U-151	0.046" (1.15 mm)	5′ (1.5 m)	N/A	5,800 psi (400 bar)	ea.
STAINLESS STEEL,	1/8" OD				
U-825	0.080" (2.0 mm)	10" (25 cm)	N/A	7,600 psi (524 bar)	ea.
	0.080" (2.0 mm)	3.2' (1 m)	N/A	7,600 psi (524 bar)	ea.
U-800	0.000 (2.0 11111)	. ,			
U-800	0.080" (2.0 mm)	9.8' (3 m)	N/A	7,600 psi (524 bar)	ea.




PEEKsil Tubing

- PEEK covered fused silica
- ▶ 1/32" and 1/16" outside diameters with a wide variety of inside diameters
- > Precut to numerous standard lengths

PEEKsil's sheathing is mechanically strong and has ideal characteristics for sealing with many styles of fittings. The fused silica core provides a consistent and rigid fluid pathway with very tight tolerances and industry-accepted chemical properties. Together, this makes PEEKsil tubing ideal for numerous applications. In fact, PEEKsil can be used as a direct replacement for conventional stainless steel or PEEK tubing in many analytical systems.

Like traditional fused silica tubing, PEEKsil has excellent chemical compatibility and extremely low adsorption characteristics, especially when compared with stainless steel.

Please Note: Do not cut this tubing. It should be used at its precut lengths because of permanent damage caused by conventional cutters.

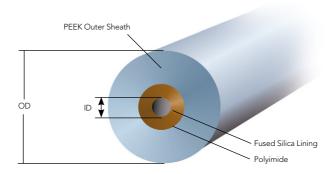



SPECIFICATIONS & DETAILS

Tubing OD	OD Tolerance	Tubing ID	ID Tolerance
		25 μm	±0.00004" (1 μm)
1/32"	±0.0008" (20 μm)	50–100 μm	±0.00012" (3 μm)
1/16"	±0.0012" (30 μm)	0.15-0.30 mm	±0.0002" (5 μm)



SPECIFICATIONS & DETAILS

Because PEEKsil tubing has fused silica tubing at its core, the pressure rating for this tubing is determined by the inner diameter of the tubing. The following chart highlights the Maximum Pressure values for this tubing, as determined by SGE International Pty., Ltd., the manufacturer of this tubing:

Tubing ID	Maximum Pressure
25 μm	25,000 psi (1,723 bar)
50 μm	20,000 psi (1,379 bar)
75–100 μm	15,000 psi (1,034 bar)
150–175 μm	8,500 psi (586 bar)
200–300 μm	6,000 psi (414 bar)

The pressure ratings provided are indicative of the performance capabilities of the tubing. The actual pressure limits achievable will depend upon the fittings used, the quality of the receiving port, and other factors. Contact IDEX Health & Science or your authorized Distributor for more information

PEEKsil[™] Tubing

" (10 cm)	ID 22" OD	Length	Color	Q
PEEKSIL TUBING, 1/3 255	32" OD 0.001" (25 μm)	2" (5 cm)	Orongo	2-
2510			Orange	
	0.001" (25 μm)	4" (10 cm)	Orange	2-
515	0.001" (25 μm)	6" (15 cm)	Orange	2-
520	0.001" (25 μm)	8" (20 cm)	Orange	2-
550	0.001" (25 μm)	1.6' (50 cm)	Orange	2-
05	0.002" (50 μm)	2" (5 cm)	Natural	2-
010	0.002" (50 μm)	4" (10 cm)	Natural	2-
015	0.002" (50 μm)	6" (15 cm)	Natural	2-
020	0.002" (50 μm)	8" (20 cm)	Natural	2-
EEKSIL TUBING, 1/3	32" OD			
050	0.002" (50 μm)	1.6' (50 cm)	Natural	2-
55	0.003" (75 μm)	2" (5 cm)	Black	2-
510	0.003" (75 μm)	4" (10 cm)	Black	2-
515	0.003" (75 μm)	6" (15 cm)	Black	2-
520	0.003" (75 µm)	8" (20 cm)	Black	2-
550	0.003" (75 μm)	1.6' (50 cm)	Black	2-
005			Red	
	0.004" (100 µm)	2" (5 cm)		2-
0010	0.004" (100 μm)	4" (10 cm)	Red	2-
0015	0.004" (100 μm)	6" (15 cm)	Red	2-
0020	0.004" (100 μm)	8" (20 cm)	Red	2-
0050	0.004" (100 μm)	1.6' (50 cm)	Red	2-
505	0.006" (150 μm)	2" (5 cm)	Purple	2-
5010	0.006" (150 μm)	4" (10 cm)	Purple	2-
5015	0.006" (150 μm)	6" (15 cm)	Purple	2-
5020	0.006" (150 μm)	8" (20 cm)	Purple	2-
5050	0.006" (150 µm)	1.6' (50 cm)	Purple	2-
EEKSIL TUBING, 1/	·		·	
255	0.001" (25 μm)	2" (5 cm)	Orange	5-
510	0.001" (25 μm)	4" (10 cm)		5-
			Orange	
2515	0.001" (25 μm)	6" (15 cm)	Orange	5-
2520	0.001" (25 µm)	8" (20 cm)	Orange	5-
2550	0.001" (25 μm)	1.6' (50 cm)	Orange	2-
05	0.002" (50 μm)	2" (5 cm)	Natural	5-
5010	0.002" (50 μm)	4" (10 cm)	Natural	5-
015	0.002" (50 μm)	6" (15 cm)	Natural	5-
6020	0.002" (50 μm)	8" (20 cm)	Natural	5-
5050	0.002" (50 μm)	1.6' (50 cm)	Natural	2-
' 55	0.003" (75 µm)	2" (5 cm)	Black	5-
7510	0.003" (75 μm)	4" (10 cm)	Black	5-
515	0.003" (75 µm)	6" (15 cm)	Black	5-
7520	0.003 (75 μm)	8" (20 cm)	Black	5-
550	0.003" (75 μm)	1.6' (50 cm)	Black	2-
1005	0.004" (100 μm)	2" (5 cm)	Red	5-
0010	0.004" (100 μm)	4" (10 cm)	Red	5-
0015	0.004" (100 μm)	6" (15 cm)	Red	5-
0020	0.004" (100 μm)	8" (20 cm)	Red	5-
0050	0.004" (100 μm)	1.6' (50 cm)	Red	2-
505	0.006" (150 μm)	2" (5 cm)	Purple	5-
5010	0.006" (150 μm)	4" (10 cm)	Purple	5-
5015	0.006" (150 µm)	6" (15 cm)	Purple	5-
5020	0.006" (150 μm)	8" (20 cm)	Purple	5-
5050	0.006" (150 µm)	1.6' (50 cm)	Purple	2-
7515	0.007" (175 µm)	6" (15 cm)	Yellow	5-
7513 7520	0.007 (175 µm)	8" (20 cm)	Yellow	5-
7550 7550				
	0.007" (175 μm)	1.6' (50 cm)	Yellow	2-
005	0.008" (200 µm)	2" (5 cm)	Blue	5-
0015	0.008" (200 μm)	6" (15 cm)	Blue	5-
0020	0.008" (200 μm)	8" (20 cm)	Blue	5-
0050	0.008" (200 μm)	1.6' (50 cm)	Blue	2-
8005	0.012" (300 μm)	2" (5 cm)	Gray	5-
80010	0.012" (300 µm)	4" (10 cm)	Gray	5-
80015	0.012" (300 μm)	6" (15 cm)	Gray	5-
30020	0.012" (300 µm)	8" (20 cm)	Gray	5-
	0.0 12 (300 pm)	0 (20 011)	U.U,	J-

PFA Tubing



- ▶ 1/16" and 1/8" ODs available
- Excellent solvent resistance and low gas permeability
- ➤ Constructed with genuine Teflon[™] PFA resin

PFA (perfluoroalkoxyalkane) tubing offers excellent solvent resistance (virtually identical to FEP and PTFE) while adding several advantages. These include smoother surface texture, higher continuous service temperature and superior polymer purity. The recommended maximum operating temperature for our PFA tubing is 80 °C.



-) 360 μ m, 1/16", 1/8", 3/16", and 1/4" outside diameters available
- > PFA HP and PFA HP Plus Grades available
- > Virtually contaminant free
- Constructed with genuine Teflon™ High Purity PFA resin

PFA High Purity (HP) tubing offers all of the benefits of standard PFA tubing, with the additional benefit of being manufactured from a premium grade of PFA that is one of the most contaminant-free polymers available. In PFA HP, we offer tubing with the following outer diameters: 1/16", 1/8", 3/16", and 1/4".

PFA High Purity (HP) Plus tubing carries all of the benefits of PFA HP tubing, with the additional benefits of increased ability to withstand repeated flexing; improved resistance to stress cracking when exposed to aggressive fluorosurfactants; and smoother, clearer walls. In PFA HP Plus, we offer tubing with the following outer diameters: 360 μ m, 1/16", and 1/8".

(Please Note: Due to the physical nature of the 360 µm OD tubing, we recommend using our A-350 Polymer Tubing Cutter from page 28 when cutting this tubing to length. Additionally, extra care should be taken to ensure fittings are not overtightened and to ensure the tubing is not stretched once secured in place, to ensure the dimensional stability of the tubing.)



SPECIFICATIONS & DETAILS

Tubing OD	OD Tolerances	Tubing ID	ID Tolerance		
PFA TUBING SPEC	CIFICATIONS				
1/16"	±0.001" (25 μm)	All	±0.001" (25 μm)		
1/8"	±0.003" (75 μm)	All	±0.003" (75 μm)		
HIGH PURITY PFA TUBING SPECIFICATIONS					
1/16"	±0.001" (25 μm)	All	±0.001" (25 μm)		
1/8"	±0.003" (75 μm)	All	±0.003" (75 μm)		
3/16"	±0.003" (75 μm)	All	±0.003" (75 μm)		
1/4"	±0.004" (100 μm)	All	±0.004" (100 μm)		
360 μM OD PFA H	IP TUBING SPECIFICATIONS				
360 µm	±0.0005" (12.5 μm)	All	±0.0005" (12.5 μm)		

PFA Tubing

Part No.	ID	Length	Color	Max. Pressure	Qty.
PFA TUBING, 1/1	6" OD				
1500	0.020"(0.50 mm)	5' (1.5 m)	Natural	2,000 psi (138 bar)	ea.
1512L	0.020"(0.50 mm)	50' (15 m)	Natural	2,000 psi (138 bar)	ea.
1512M	0.020"(0.50 mm)	1,000' (304 m)	Natural	2,000 psi (138 bar)	ea.
1502	0.030" (0.75 mm)	5' (1.5 m)	Natural	1,000 psi (69 bar)	ea.
1514L	0.030" (0.75 mm)	50' (15 m)	Natural	1,000 psi (69 bar)	ea.
1514M	0.030" (0.75 mm)	1,000' (304 m)	Natural	1,000 psi (69 bar)	ea.
1503	0.040" (1.0 mm)	5' (1.5 m)	Natural	500 psi (34 bar)	ea.
1507L	0.040" (1.0 mm)	50' (15 m)	Natural	500 psi (34 bar)	ea.
1507M	0.040" (1.0 mm)	1,000' (304 m)	Natural	500 psi (34 bar)	ea.
PFA TUBING, 1/8	" OD				
1509-5	0.062" (1.55 mm)	5' (1.5 m)	Natural	500 psi (34 bar)	ea.
1509L	0.062" (1.55 mm)	50' (15 m)	Natural	500 psi (34 bar)	ea.
PFA HP TUBING,	1/16" OD				
1622-5	0.020"(0.50 mm)	5' (1.5 m)	Natural	2,000 psi (138 bar)	ea.
1622L	0.020"(0.50 mm)	50' (15 m)	Natural	2,000 psi (138 bar)	ea.
1622M	0.020"(0.50 mm)	1,000' (304 m)	Natural	2,000 psi (138 bar)	ea.
1632-5	0.030" (0.75 mm)	5' (1.5 m)	Natural	1,000 psi (69 bar)	ea.
1632L	0.030" (0.75 mm)	50' (15 m)	Natural	1,000 psi (69 bar)	ea.
1632M	0.030" (0.75 mm)	1,000' (304 m)	Natural	1,000 psi (69 bar)	ea.
PFA HP TUBING,	1/8" OD				
1641-5	0.062" (1.55 mm)	5' (1.5 m)	Natural	500 psi (34 bar)	ea.
1641L	0.062" (1.55 mm)	50' (15 m)	Natural	500 psi (34 bar)	ea.
PFA HP PLUS TUE	BING, 1/16" OD				
1902-5	0.010 (0.25 mm)	5' (1.5 m)	Natural	3,000 psi (207 bar)	ea.
1902L	0.010 (0.25 mm)	50' (15 m)	Natural	3,000 psi (207 bar)	ea.
1902M	0.010 (0.25 mm)	1,000' (304 m)	Natural	3,000 psi (207 bar)	ea.
1907-5	0.020"(0.50 mm)	5' (1.5 m)	Natural	2,000 psi (138 bar)	ea.
1907L	0.020"(0.50 mm)	50' (15 m)	Natural	2,000 psi (138 bar)	ea.
1907M	0.020"(0.50 mm)	1,000' (304 m)	Natural	2,000 psi (138 bar)	ea.
1912-5	0.030" (0.75 mm)	5' (1.5 m)	Natural	1,000 psi (69 bar)	ea.
1912L	0.030" (0.75 mm)	50' (15 m)	Natural	1,000 psi (69 bar)	ea.
1912M	0.030" (0.75 mm)	1,000' (304 m)	Natural	1,000 psi (69 bar)	ea.
PFA HP PLUS TUE	BING, 1/8" OD				
1921-5	0.062" (1.55 mm)	5' (1.5 m)	Natural	500 psi (34 bar)	ea.
1921L	0.062" (1.55 mm)	50' (15 m)	Natural	500 psi (34 bar)	ea.



FEP Tubing

- Great for moderate-to-low pressure applicationse
- **▶** 1/16", 1/8", 3/16", 1/4", or 5/16" outside diameter available
- ▶ 1 mm, 2 mm, or 3 mm outside diameter available
- Maximum continuous use temperature: 50 °C
- Constructed with genuine Teflon™ FEP resin

With virtually identical chemical resistance to PFA at a lower price, FEP tubing is great for general, low pressure applications. Compared to PTFE, FEP (fluorinated ethylene propylene) tubing is held to tighter tolerances and has lower gas permeability (see material properties on our website: www.idex-hs.com).

Much of our FEP Tubing — even the color-tinted options — is translucent, making it possible to watch fluid flow. Using different colored tubing can help identify transfer lines in multisolvent systems. Color coding also allows easy identification of the tubing thru-hole size. Black FEP tubing is available for light-sensitive applications (such as enzymatic and chemi-luminescent reactions) and entering/exiting flow cells.

SPECIFICATIONS & DETAILS

Tubing Size	OD Tolerances	ID Tolerances
1/16" OD	±0.001" (25 μm)	±0.001" (25 μm)
1/8" OD	±0.003" (75 μm)	±0.003" (75 μm)
3/16" OD	±0.004" (0.10 mm)	±0.004" (0.10 mm)
5/16" OD	±0.004" (0.10 mm)	±0.004" (0.10 mm)
1 mm OD	±0.001" (25 μm)	±0.001" (25 μm)
2 mm OD	±0.003" (75 μm)	±0.003" (75 μm)
3 mm OD	±0.003" (75 µm)	±0.003" (75 µm)

Part No.	ID	Length	Color	Max. Pressure	Qty.
FEP TUBING, 1/	/16" OD				
1527-5	0.010" (0.25 mm)	5' (1.5 m)	Natural	3,000 psi (207 bar)	ea.
1527L	0.010" (0.25 mm)	50' (15 m)	Natural	3,000 psi (207 bar)	ea.
1527XL	0.010" (0.25 mm)	100' (30 m)	Natural	3,000 psi (207 bar)	ea.
1527M	0.010" (0.25 mm)	1,000' (304 m)	Natural	3,000 psi (207 bar)	ea.
1548-5	0.020" (0.50 mm)	5' (1.5 m)	Natural	2,000 psi (138 bar)	ea.
1548L	0.020" (0.50 mm)	50' (15 m)	Natural	2,000 psi (138 bar)	ea.
1548XL	0.020" (0.50 mm)	100' (30 m)	Natural	2,000 psi (138 bar)	ea.
1548M	0.020" (0.50 mm)	1,000' (304 m)	Natural	2,000 psi (138 bar)	ea.
1520-5	0.030" (0.75 mm)	5' (1.5 m)	Natural	1,000 psi (69 bar)	ea.
1520L	0.030" (0.75 mm)	50' (15 m)	Natural	1,000 psi (69 bar)	ea.
1520XL	0.030" (0.75 mm)	100' (30 m)	Natural	1,000 psi (69 bar)	ea.
1520M	0.030" (0.75 mm)	1,000' (304 m)	Natural	1,000 psi (69 bar)	ea.
FEP TUBING, 1/	/8" OD				
1521-5	0.062" (1.55 mm)	5' (1.5 m)	Natural	500 psi (34 bar)	ea.
1521L	0.062" (1.55 mm)	50' (15 m)	Natural	500 psi (34 bar)	ea.
1521XL	0.062" (1.55 mm)	100' (30 m)	Natural	500 psi (34 bar)	ea.
FEP TUBING, 3/	/16" OD				
1524L	0.125" (3.2 mm)	50' (15 m)	Natural	500 psi (34 bar)	ea.
1524XL	0.125" (3.2 mm)	100' (30 m)	Natural	500 psi (34 bar)	ea.
FEP TUBING, 1/	⁄4" OD				
1650L	0.188" (4.8 mm)	50' (15 m)	Natural	250 psi (17 bar)	ea.
1650XL	0.188" (4.8 mm)	100' (30 m)	Natural	250 psi (17 bar)	ea.
FEP TUBING, 1.	0 MM OD				
1671L	0.020" (0.50 mm)	50' (15 m)	Natural	500 psi (34 bar)	ea.
1671XL	0.020" (0.50 mm)	100' (30 m)	Natural	500 psi (34 bar)	ea.
FEP TUBING, 2.	.0 MM OD				
1673L	0.40" (1.0 mm)	50' (15 m)	Natural	500 psi (34 bar)	ea.
1673XL	0.40" (1.0 mm)	100' (30 m)	Natural	500 psi (34 bar)	ea.
FEP TUBING, 3.					
1677L	0.080" (2.0 mm)	50' (15 m)	Natural	500 psi (34 bar)	ea.
1677XL	0.080" (2.0 mm)	100′ (30 m)	Natural	500 psi (34 bar)	ea.



ETFE Tubing

- > Excellent chemical resistance
- Constructed with genuine Tefzel[™] resin
- > Holds pressure up to 4,000 psi (276 bar)
- ▶ 1/16" or 1/8" outside diameter available
- Maximum continuous operating temperature: 80 °C

ETFE (ethylene-tetrafluoroethylene) tubing is an excellent fluoropolymer product that offers several benefits over tubing manufactured from PTFE, FEP, or PFA. These benefits include enhanced pressure holding capabilities, increased mechanical stability and lower gas permeability.



APPLICATION NOTE

ETFE tubing is an ideal choice for the fluid pathway between the vacuum degasser and the system's pump. Its low gas permeability will help ensure the mobile phase solvents do not regas while in transit.



RELATED PRODUCTS

Other tubing materials and dimensions may be available. Please contact IDEX Health & Science or your local representative directly.



SPECIFICATIONS & DETAILS

Tubing OD	Tubing ID	OD/ID Tolerances
1/16" OD	0.010" (0.25 mm), 0.020" (0.50 mm), 0.030" (0.75 mm)	±0.001" (25 μm)
1/16" OD	0.040" (1.0 mm)	±0.002" (50 μm)
1/8" OD	All	±0.003" (75 μm)

Part No.	ID	Length	Color	Max. Pressure	Qty.
ETFE TUBING, 1	I/16" OD				
1529	0.010 (0.25 mm)	5' (1.5 m)	Natural	4,000 psi (276 bar)	ea.
1529L	0.010 (0.25 mm)	50' (15 m)	Natural	4,000 psi (276 bar)	ea.
1529XL	0.010 (0.25 mm)	100' (30 m)	Natural	4,000 psi (276 bar)	ea.
1529M	0.010 (0.25 mm)	1,000' (304 m)	Natural	4,000 psi (276 bar)	ea.
1516	0.020"(0.50 mm)	5' (1.5 m)	Natural	3,000 psi (207 bar)	ea.
1516L	0.020"(0.50 mm)	50' (15 m)	Natural	3,000 psi (207 bar)	ea.
1516XL	0.020"(0.50 mm)	100' (30 m)	Natural	3,000 psi (207 bar)	ea.
1516M	0.020"(0.50 mm)	1,000' (304 m)	Natural	3,000 psi (207 bar)	ea.
1528	0.030" (0.75 mm)	5' (1.5 m)	Natural	2,000 psi (138 bar)	ea.
1528L	0.030" (0.75 mm)	50' (15 m)	Natural	2,000 psi (138 bar)	ea.
1528XL	0.030" (0.75 mm)	100' (30 m)	Natural	2,000 psi (138 bar)	ea.
1528M	0.030" (0.75 mm)	1,000' (304 m)	Natural	2,000 psi (138 bar)	ea.
1517	0.040" (1.0 mm)	5' (1.5 m)	Natural	500 psi (34 bar)	ea.
1517L	0.040" (1.0 mm)	50' (15 m)	Natural	500 psi (34 bar)	ea.
1517XL	0.040" (1.0 mm)	100' (30 m)	Natural	500 psi (34 bar)	ea.
1517M	0.040" (1.0 mm)	1,000' (304 m)	Natural	500 psi (34 bar)	ea.
ETFE TUBING, 1	I/8" OD				
1530	0.062" (1.55 mm)	5' (1.5 m)	Natural	1,000 psi (69 bar)	ea.
1530L	0.062" (1.55 mm)	50' (15 m)	Natural	1,000 psi (69 bar)	ea.
1530XL	0.062" (1.55 mm)	100' (30 m)	Natural	1,000 psi (69 bar)	ea.
1648	0.093" (2.4 mm)	5' (1.5 m)	Natural	500 psi (34 bar)	ea.
1648L	0.093" (2.4 mm)	50' (15 m)	Natural	500 psi (34 bar)	ea.
1648XL	0.093" (2.4 mm)	100' (30 m)	Natural	500 psi (34 bar)	ea.



loos



Fused Silica Tubing Cutters

We offer a precision cutter for fused silica tubing — SGT's Shortix™ Cutter (FS-315). This cutter ensures clean, trouble-free cutting of fused silica tubing, providing better cuts than any other product on the market. It also includes a built-in magnifying glass to examine the cut tubing ends. Order the FS-315-02 Maintenance Kit, as needed, to replace a worn or damaged cutting wheel.

When using traditional fused silica tubing cutters, only a small part of the tubing wall is scratched, then the tubing is snapped or pulled in two, often resulting in a jagged, uneven cut. With a Shortix Cutter, a clean cut is made every time, regardless of skill or experience, as the cut is made by rotating a diamond blade around the entire circumference of the tubing.

Please Note: The FS-315 Fused Silica Tubing Cutters are designed to cut only tubing with ODs of 350 μ m-780 μ m and IDs of 100 μ m-350 μ m.



Polymer Tubing Cutters

For 1/16", 1/8", 3/16", 1/4", and 5/16" OD tubing

A flat, 90°, burr-free cut is difficult to obtain with most commercial polymer tubing cutters. Our experts have designed several tubing cutters specifically to cut polymer tubing. This line of tubing cutters includes a standard cutter for 1/16" and 1/8" OD tubing (A-327), and another for large bore tubing (A-329). Each has guide holes to ensure precise cutting. These cutters are durable, reliable, and easy to operate. Five replacement blades are included with each tool.



Capillary Polymer Tubing Cutters

Our A-350 Cutter is designed to cut capillary-sized polymer tubing. The cutter makes clean, perpendicular cuts without collapsing thin capillary walls. A set of ten tubing sleeves, required for cutting, are included with each cutter, along with five replacement blades. The included tubing sleeves are for cutting 360 µm OD polymer capillary tubing. Alternative sleeves are available for cutting 1/32" OD tubing. All tubing sleeves are 2" long and constructed with genuine Teflon™ FEP resin.



NOTE

- ▶ The A-350 Capillary Polymer Tubing Cutter can be used to cut tubing $\overline{\text{OD}}$ sizes other than 360 µm and 1/32". Simply use the proper NanoTight™ Tubing Sleeve found on page 52. Please note, however, that these sleeves are shorter than those listed on this page, and therefore will last through fewer cuts.
- Our tubing cutters are material specific: the A-327, A-329, and A-350 should only be used to cut polymer tubing, where as the FS-315 should only be used to cut fused silica tubing.

Part No.	Description	Qty.
FUSED	SILICA TUBING CUTTERS	
FS-315	Shortix Fused Silica Tubing Cutter	ea.
CAPILL	ARY POLYMER TUBING CUTTER	
A-350	Capillary Polymer Tubing Cutter* for 360 µm-1/32" OD tubing Includes (1) F-262x 10-pack of sleeves and (1) M-438-03 wrench	ea.
F-262x	Replacement Sleeves for A-350, 0.0155" ID, Green, for cutting 360 μm OD tubing	10-pk
A-327	Standard Polymer Tubing Cutter* for 1/16" and 1/8" OD tubing	ea.
A-329	Large Bore Polymer Tubing Cutter* for 3/16" – 5/16" OD tubing	ea.
A-328	Replacement Blades for A-350, A-327, and A-329	5-pk
* Include	s (1) A-328 5-pack of replacement blades.	

-	





FITTINGS

We offer a wide and diverse selection of fittings to meet your system requirements. A "fitting" refers to a complete product ready to assemble and connect tubing into a part. This could be a one-piece connector or a nut and ferrule packaged together. A "Nut" indicates the male or female threaded product sold separately, and a "Ferrule" is sold separately when indicated in the description. For your convenience we ship most Fittings and Ferrules in 10-packs. We are dedicated to providing the most reliable, proven products on the market. We have implemented more stringent testing protocols and a generous safety margin to our ratings to ensure your safety.

- 32 CONED FITTINGS
- 39 FLAT-BOTTOM FITTINGS
- 48 MARVELXACT™

- 50 FITTINGS TOOLS
- 56 LARGE BORE FITTINGS
- 57 VHP FITTINGS
- 61 FITTINGS KITS



All testing is performed with water at room temperature unless otherwise specified. Please contact IDEX Health & Science directly for further details. Results may vary depending on the material of the receiving port and tubing, actual tubing diameters (with stated tolerances), temperature and solvents used. If a pressure range is listed for a product's specification, the pressure rating depends on the tubing material used. The lower end of the range will represent testing performed on softer tubing such as FEP, and the higher end of the range will represent testing performed on harder tubing such as Stainless Steel. For more detail, please see the product specification sheets on our website, www.idex-hs.com, or contact us directly.



What Threads Do I Have?

Hold your fitting over the thread silhouettes below to identify the threads.

U.S. CUSTOMARY THREADS 6-40 6-32 10-32 1/4-28
10-32
1/4-28
1/4-28

5/16-24
••••••
4/2.22
1/2-20
METRIC THREADS
M4 x 0.7
M6 x 1





One-Piece Fingertight Fittings

- > The original One-Piece Fingertight Fitting
- All polymer construction
- Versions available for 1/16", 1/32" or 1/8" OD tubing
- Convenient and easy to use
- Reusable




NOTE

- Some of the fittings on this page are available in additional colors. Please contact your distributor or us for more information.
- ➤ Fingertight is generally equal to 3–4 in-lbs (0.34–0.45 N·m).

Our One-Piece Fingertight Fittings provide convenience and ease of use because the ferrule will not stick in a receiving port and the fitting is more easily found if dropped. The fittings for 1/16" OD tubing and 10-32 coned ports are available in a variety of colors, materials and lengths to suit virtually every application.

Beyond the standard 10-32 fittings, also featured in this product family are specialty fittings for specific applications. Our M-645 Fitting is a direct replacement for the 6-40 threaded VICI® (Valco) fitting. The P-100 can be used in 1/4-28 coned ports for 1/8" OD tubing including some of the inlet filters starting on page 102.

RheFlex® One-Piece Fittings are included in many of the manual valves, starting on page 121. The One-Piece RheFlex M4 Fittings, for use with our MX Nano-Scale Modules, are listed on page 61.

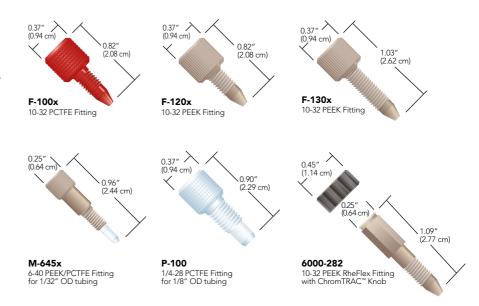

Part No.	Description	Port	Pressure Rating	Head Style	Material	Qty.	
ONE-PIECE FINGERTIGHT FITTINGS							
6000-282	Fingertight Fitting for 1/16" OD Tubing	10-32 Coned	5,000 psi (345 bar)	ChromTRAC knob	PEEK, Natural	10-pk	
F-100x	Fingertight Fitting for 1/16" OD Tubing	10-32 Coned	4,000 psi (276 bar)	Diamond Knurl	PCTFE, Red	10-pk	
F-120x	Fingertight Fitting for 1/16" OD Tubing	10-32 Coned	5,000 psi (345 bar)	Standard Knurl	PEEK, Natural	10-pk	
F-130x	Fingertight Fitting for 1/16" OD Tubing, Long	10-32 Coned	5,000 psi (345 bar)	Standard Knurl	PEEK, Natural	10-pk	
M-645x	Fingertight Fitting for 1/32" OD Tubing	6-40 Coned	1,750-3,250 psi (121-224 bar)	Headless Knurl	PEEK, Natural/PCTFE, Natural	10-pk	
P-100	Fingertight Fitting for 1/8" OD Tubing	1/4-28 Coned	1,000 psi (69 bar)	Diamond Knurl	PCTFE, Natural	ea.	



Stainless Steel Fittings

NOTE

- Do not use metal fittings in plastic ports, as this can damage the port.
- The recommended torque to tighten these fittings is 20 in-lbs (2.25 N·m).

These 316 Stainless Steel Fittings are rated to 20,000 psi (1,380 bar) when wrench tightened. Choose IDEX Health & Science standard fittings, or select from the Common Valve Fittings or other manufacturer-compatible offerings.

Standard Stainless Steel Fittings



Common Valve Fittings



Beckman® Compatible Fittings



VICI® (Valco) Compatible Fittings



Part No.	Description	Port	Pressure Rating	Head Style	Material	Qty.
STANDAR	D STAINLESS STEEL FITTINGS					
C-235x	Nut for 1/8" OD Tubing	1/4-28 Coned	20,000 psi (1,380 bar)	5/16" Hex	SST	10-pk
C-236x	Ferrule for 1/8" OD Tubing	1/4-28 Coned	20,000 psi (1,380 bar)	_	SST	10-pk
U-400x	Nut for 1/16" OD Tubing	10-32 Coned	20,000 psi (1,380 bar)	1/4" Hex	SST	10-pk
U-401x	Ferrule for 1/16" OD Tubing	10-32 or M6 Coned	20,000 psi (1,380 bar)	_	SST	10-pk
U-450x	Nut for 1/16" OD Tubing	M6 Coned	20,000 psi (1,380 bar)	5/16" Hex	SST	10-pk
COMMON	VALVE FITTINGS					
6000-082	Fitting for 1/8" OD Tubing	5/16-24 Coned	20,000 psi (1,380 bar)	5/16" Hex	SST	ea.
6000-083	Ferrule for 1/8" OD Tubing	5/16-24 Coned	20,000 psi (1,380 bar)	_	SST	5-pk
6000-209	Fitting for 1/16" OD Tubing	10-32 Coned	20,000 psi (1,380 bar)	1/4" Hex	SST	10-pk
6000-210	Ferrule for 1/16" OD Tubing	10-32 Coned	20,000 psi (1,380 bar)	_	SST	10-pk
6000-211	Fitting for 1/16" OD Tubing, Long	10-32 Coned	20,000 psi (1,380 bar)	1/4" Hex	SST	10-pk
6000-262	Fitting for 1/16" OD Tubing, Extra Long	10-32 Coned	20,000 psi (1,380 bar)	1/4" Hex	SST	10-pk
MANUFAC	CTURER COMPATIBLE FITTINGS					
U-320x	Nuts for 1/16" OD Tubing, Valco/VICI Compatible	10-32 Coned	20,000 psi (1,380 bar)	1/4" Hex	SST	10-pk
U-321x	Ferrule for 1/16" OD Tubing, Valco/VICI Compatible	10-32 Coned	20,000 psi (1,380 bar)	_	SST	10-pk
U-410x	Nuts for 1/16" OD Tubing, Waters Compatible	10-32 Coned	20,000 psi (1,380 bar)	5/16" Hex	SST	10-pk



MicroTight Fittings

- Comprehensive Fitting System for Connecting Capillary Tubing
- Made from PEEK Polymer

MicroTight® One-Piece Fittings are designed for use with the NanoPort™ and MicroTight Unions, Adapters and Inline MicroFilters. Specifically made for 360 µm OD tubing, 1/32" OD tubing, or our MicroTight Tubing Sleeves (see page 52), these fittings make superior fingertight connections with capillary tubing. MicroTight Fittings withstand temperatures up to 125 °C.

The MicroTight family also includes a female nut matched with one of five dedicated ferrules for connecting specific tubing ODs.

Use the P-277 Extender Tool to tighten standard micro knurl 6-32 fittings in hardto-reach places. Tighten micro headless 6-32 fittings with our N-290 Tool. See page 50 for more information.




MicroTight fittings and MicroFerrules

While the MicroTight Female Nuts may be used with any of the separate MicroFerrules, the MicroFerrules themselves are port-specific and are thus not interchangeable. Additionally, the one-piece MicroTight fittings are also portspecific and should not be exchanged.



Standard Head Fitting for use with 360 µm OD tubing

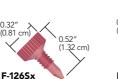

Standard Head Fitting for use with 1/32" OD tubing



F-124Hx

Headless Fitting

P-555x Standard Head Plug



Standard Head Fitting for use with MicroTight Sleeves







Headless Fitting for use with MicroTight Sleeves




MicroFerrule for 1/16" OD tubing



F-172 MicroFerrule for 0.025" OD tubing



P-116 MicroFerrule Plug

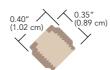

P-416 5/16-24 internal threads



P-416BLK Female Nut 5/16-24 internal threads



P-416G Female Nut 5/16-24 internal threads



F-152

MicroFerrule for 360 μm OD tubing

RELATED PRODUCTS

- Connectors for Capillary Tubing can be found on page 75.
- > Very High Pressure fittings for capillary tubing can be found on page 58.
- Capillary tubing is featured on page 16.

Part No.	Description	Port	Pressure Rating	Head Style	Material	Qty.
MICROTIG	HT FITTINGS					
F-124Hx	MicroTight Fitting for 360 μm OD Tubing	6-32 Coned	5,000 psi (345 bar)	Micro Headless Knurl	PEEK Blue	10-pk
F-124Sx	MicroTight Fitting for 360 μm OD Tubing	6-32 Coned	5,000 psi (345 bar)	Standard Micro Knurl	PEEK Blue	10-pk
F-125Hx	MicroTight Fitting for MicroTight Tubing Sleeves	6-32 Coned	4,000 psi (276 bar)	Micro Headless Knurl	PEEK Natural	10-pk
F-125x	MicroTight Fitting for MicroTight Tubing Sleeves	6-32 Coned	4,000 psi (276 bar)	Standard Micro Knurl	PEEK Natural	10-pk
F-126Sx	MicroTight Fitting for 1/32" OD Tubing	6-32 Coned	5,000 psi (345 bar)	Standard Micro Knurl	PEEK Red	10-pk
P-555x	MicroTight Plug	6-32 Coned	5,000 psi (345 bar)	Standard Micro Knurl	PEEK Natural	10-pk
MICROFE	RRULES AND FEMALE NUTS					
F-112	MicroFerrule for 1/32" OD Tubing	5/16-24 Coned	5,000 psi (345 bar)	_	PEEK Natural	ea.
F-132	MicroFerrule for 1/16" OD Tubing	5/16-24 Coned	5,000 psi (345 bar)	_	PEEK Natural	ea.
F-152	MicroFerrule for 360 µm OD Tubing	5/16-24 Coned	5,000 psi (345 bar)	_	PEEK Natural	ea.
F-152BLK	MicroFerrule for 360 µm OD Tubing	5/16-24 Coned	5,000 psi (345 bar)	_	PEEK Black	ea.
F-172	MicroFerrule for MicroTight Tubing Sleeves	5/16-24 Coned	4,000 psi (276 bar)	_	PEEK Black	ea.
P-116	MicroFerrule Plug	5/16-24 Coned	5,000 psi (345 bar)	_	PEEK Black	ea.
P-416	MicroTight Female Nut	5/16-24 Coned	4,000-5,000 psi (276-345 bar)	Female Knurl	PEEK Natural	ea.
P-416BLK	MicroTight Female Nut	5/16-24 Coned	4,000-5,000 psi (276-345 bar)	Female Knurl	PEEK Black	ea.
P-416G	MicroTight Female Nut	5/16-24 Coned	4,000-5,000 psi (276-345 bar)	Female Knurl	PEEK Green	ea.



Two-Piece Fingertight Fittings

- Designed to connect tubing to 10-32 coned ports
- > Ferrules available for directly connecting 1/16", 1/32", 360 µm, or 190 µm OD tubing
- > Economical, replace only the ferrule



APPLICATION NOTE

Some Fingertight Nuts feature wings in addition to a knurled head, which provide more leverage when tightening the fitting into a receiving port. Choose our single or double-winged design.

Please Note: customers can use the standard knurl head fittings with our tightening tools found on page 50.



Two-Piece Fingertight Fittings feature a separate ferrule. Use a standard knurled head fitting for traditional fingertight applications, or use a fitting with wings built into the head for extra tightening leverage. A stainless steel hex headed fitting can be used for applications where a wrench may be needed for added tightening torque.

The M-215 Conductive Perfluoroelastomer Ferrule is designed for mass spectrometer electrospray applications. Unlike most graphite ferrules, the elastomeric properties of this ferrule let you use it through many tightening/retightening cycles. It also eliminates any possibility of graphite contamination in your system. Like graphite ferrules, you can apply voltage through a metallic port block or metallic nut, allowing voltage to translate to the flow path through the ferrule.







PCTFE Ferrule for 360 µm OD tubing




F-113 PEEK Ferrule for 1/32" OD tubing




F-142N ETFE Ferrule for 1/16" OD tubing



Conductive Perfluoroelastometer Ferrule for 360 µm OD tubing

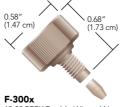

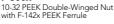





Part No.	Description	Port	Pressure Rating	Head Style	Material	Qty.		
TWO-PIECE FINGERTIGHT FITTINGS (INCLUDES F-142 FERRULES)								
F-300x	Fingertight Fitting for 1/16" OD Tubing	10-32 Coned	6,000 psi (414 bar)	Double Wing	PEEK Natural	10-pk		
F-330x	Fingertight Fitting for 1/16" OD Tubing, Long	10-32 Coned	6,000 psi (414 bar)	Standard Knurl	PEEK Natural	10-pk		
F-331x	Fingertight Fitting for 1/16" OD Tubing	10-32 Coned	6,000 psi (414 bar)	Standard Knurl	PEEK Natural	10-pk		
REPLACEMENT FERRULES								
F-113	Ferrule for 1/32" OD Tubing	10-32 Coned	6,000 psi (414 bar)	_	PEEK Natural	ea.		
F-142	Ferrule for 1/16" OD Tubing	10-32 Coned	6,000 psi (414 bar)	_	PEEK Natural	ea.		
F-142N	Ferrule for 1/16" OD Tubing	10-32 Coned	4,000 psi (276 bar)	_	ETFE Natural	ea.		
F-148	Ferrule for 190 µm OD tubing	10-32 Coned	6,000 psi (414 bar)	_	PCTFE Natural	ea.		
F-151	Ferrule for 360 µm OD Tubing	10-32 Coned	6,000 psi (414 bar)	_	PCTFE Natural	ea.		
M-215	Conductive Ferrule for 360 µm OD tubing	10-32 Coned	1,500 psi (103 bar)	_	Conductive Perfluoroelastomer	ea.		



Two-Piece SealTight Fingertight Fittings

- Several nut lengths and head styles to fit into a variety of applications
- Designed to connect 1/16" OD tubing to 10-32 coned ports
- ▶ Hold up to 9,000 psi (620 bar)





Overtightening these fittings on fluoropolymer (e.g., FEP, PFA, and ETFE) tubing can cause the ID of your tubing to collapse.



- ➤ Find tightening tools for these fittings on page 50.
- Try the F-350x FlushNut[™] for the ultimate streamline design.

The dual compression created by the specially designed nut and ferrule enables our SealTight™ Fittings system to outperform standard finger tightened fittings. The forward cone of the SealTight Ferrule provides gripping power and a leak-free seal via conventional compression by the receiving port. The slotted end creates the second compression zone in conjunction with a SealTight Nut. All SealTight Nuts are for use with 1/16″ OD tubing and are designed to be used with the F-192x Ferrule. A wide variety of fitting head styles are available for various space constraints.



F-193x 10-32 Short PEEK Hex Head Nut, with F-192x Ferrule







F-192x SealTight Ferrule, for 1/16" OD tubing



Part No.	Description	Port	Pressure Rating	Head Style	Material	Qty.		
SEALTIGHT TWO-PIECE FITTINGS (INCLUDES F-192 FERRULES)								
F-193x	SealTight Fitting for 1/16" OD Tubing, Short	10-32 Coned	7,000–9,000 psi (483–620 bar)	1/4" Hex	PEEK Black	10-pk		
F-195x	SealTight Fitting for 1/16" OD Tubing, Short	10-32 Coned	7,000–9,000 psi (483–620 bar)	Standard Knurl	PEEK Black	10-pk		
F-196x	SealTight Fitting for 1/16" OD Tubing, Long	10-32 Coned	7,000-9,000 psi (483-620 bar)	Standard Knurl	PEEK Black	10-pk		
F-287x	SealTight Fitting for 1/16" OD Tubing, Long	10-32 Coned	7,000-9,000 psi (483-620 bar)	Knurl-1/4" Hex	PEEK Black	10-pk		
F-350x	SealTight Fitting for 1/16" OD Tubing, FlushNut	10-32 Coned	7,000-9,000 psi (483-620 bar)	FlushNut	SST	10-pk		
REPLACEN	REPLACEMENT FERRULES							
F-192x	SealTight Ferrule for 1/16" OD Tubing	10-32 or M6 Coned	7,000-9,000 (483-620 bar)	_	PEEK/Black	10-pk		



NanoTight[™] Fittings & Sleeves

- For connecting 1/16" OD or capillary tubing using tubing sleeves to standard 10-32 coned ports
- > Multiple nut styles available
- Nuts manufactured from PEEK polymer, ferrules manufactured from ETFE



- Find tightening tools for these head styles on page 50.
- NanoTight Tubing sleeves start on page 52.

NanoTight™ Fittings and Sleeves are designed to connect 70 µm−1 mm OD capillary tubing to any standard 10-32 coned port normally intended for 1/16″ OD tubing using the NanoTight Tubing Sleeves on page 52. The fittings can also be used to connect any 1/16″ OD tubing. The ETFE ferrule material is softer than PEEK, making it a good candidate for connecting thin walled semi-rigid tubing such as FEP and ETFE into 10-32 ports with minimal constricting to the inner diameter.

Select from our expansive line of PEEK NanoTight Fittings, featuring several head style and length options. Each 10-pack of nuts includes ten ETFE F-142Nx ferrules.



F-333Nx Short Headless Nut with F-142Nx Ferrule


Part No.	Description	Port	Pressure Rating	Head Style	Material (Nut/Ferrule)	Qty.
NANOTI	GHT FITTINGS (INCLUDES F-142N FERRULES)					
F-333Nx	NanoTight Fitting for 1/16" OD Tubing and NanoTight Sleeves, Short	10-32 Coned	4,000 psi (276 bar)	Headless Knurl	PEEK Natural/ETFE Natural	10-pk
REPLACE	EMENT FERRULES					
F-142Nx	NanoTight Ferrule for 1/16" OD Tubing and NanoTight Sleeves	10-32 Coned	4,000 psi (276 bar)	_	ETFE Natural	10-pk



LiteTouch® Fittings

- Helps prevent twisting of polymer tubing
- > High pressure with fingertight convenience
- Options available for 1/32", 1/16", or 1/8" OD tubing



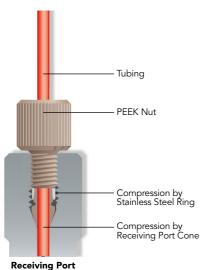



The stainless steel nuts on page 59 can also be used with the LiteTouch ferrules on this page. The LiteTouch® Fittings System grips tubing at two compression points (see diagram), holding to high pressures with Fingertight convenience. It also prevents polymer tubing from twisting, a potential problem when using standard Fingertight fittings. LiteTouch Fittings are available for use with 1/32", 1/16", or 1/8" OD tubing sizes, and for 10-32 or 1/4-28 coned ports.

For those space-limited applications where nut heads interfere with each other, try the FlushNut™ Fittings. (FlushNut Fittings require a tightening tool. Please see page 50 for more information about these products.)

To avoid collapsing the ID of your tubing, the LiteTouch system can be used on hard tubing only, such as stainless steel and PEEK polymer tubing. The LiteTouch Ferrule System is not recommended for repeated use in plastic ports.



F-354x 10-32 Stainless Steel FlushNut for 1/32" and 1/16" OD tubing




LT-110x 10-32 PEEK Nut for 1/32" and 1/16" OD tubing



C-235x 1/4-28 Stainless Steel Nut for 1/8" OD tubing



F-364x 1/4-28 Stainless Steel FlushNut for 1/8" OD tubing

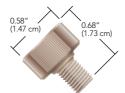

LT-210x 1/4-28 PEEK Double-Winged Nut for 1/8" OD tubing



LT-132x PEEK Ferrule with Stainless Steel Lock Ring for 1/32" OD tubing



LT-100x PEEK Ferrule with Stainless Steel Lock Ring for 1/16" OD tubing



LT-200x PEEK Ferrule with Stainless Steel Lock Ring for 1/8" OD tubing

Part No.	Description	Port	Pressure Rating	Head Style	Material	Qty.
LITETOU	CH NUTS					
C-235x	LiteTouch Nut for 1/8" OD Tubing	1/4-28 Coned	4,500 psi (310 bar)	5/16" Hex	SST	10-pk
F-354x	LiteTouch Nut for 1/16" or 1/32" OD Tubing, FlushNut	10-32 Coned	5,000 psi (345 bar)	FlushNut	SST	10-pk
F-364x	LiteTouch Nut for 1/8" OD Tubing, FlushNut	1/4-28 Coned	4,500 psi (310 bar)	FlushNut	SST	10-pk
LT-110x	LiteTouch Nut for 1/16" or 1/32" OD Tubing	10-32 Coned	5,000 psi (345 bar)	Standard Knurl	PEEK Natural	10-pk
LT-210x	LiteTouch Nut for 1/8" OD Tubing	1/4-28 Coned	4,500 psi (310 bar)	Double Wing	PEEK Natural	10-pk
LITETOU	CH FERRULES					
LT-100x	LiteTouch Ferrule for 1/16" OD Tubing	10-32 Coned	5,000 psi (345 bar)	_	PEEK Natural/SST	10-pk
LT-132x	LiteTouch Ferrule for 1/32" OD Tubing	10-32 Coned	5,000 psi (345 bar)	_	PEEK Natural/SST	10-pk
LT-200x	LiteTouch Ferrule for 1/8" OD Tubing	1/4-28 Coned	4,500 psi (310 bar)	_	PEEK Natural/SST	10-pk

^{*} When used with a stainless steel 10-32 nut from page 59.



Highest pressure holding flat-bottom fitting system we offer

- Eliminates loosening of fittings due to tubing twist
- Excellent for Tubing Assemblies
- > Holds tight even through vibration




ASSEMBLY HINT

Make sure the locking ring is oriented correctly! The flattened end of the ring should face towards the nut with the narrow end of the ferrule towards the ring.



Enlarged to show detail

SUPER FLANGELESS FITTINGS SYSTEM

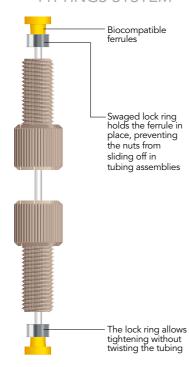

Super Flangeless Fittings

Our Super Flangeless™ Fittings provide the highest pressure holding capability in a flat-bottom fitting system that we offer. Our unique design eliminates loosening of fittings due to tubing twist and holds tight even through vibration. Our high pressure fittings are excellent for tubing assemblies and those times when connections need to be broken frequently.

6-40 & 6-32 options (for 1-16" OD tubing)



10-32 options (for 1-16" OD tubing)



Ferrules for M6x1, 1/4-28, 5/16-24



P-140x

for 3/16" OD tubing



P-360x

for 1/8" OD tubing

Super Flangeless[™] Fittings (Cont.)

M6x1, 1/4-28, 5/16-24 Options for 1/32"-3/16" OD Tubing



(1/4-28, ≤ 1/16" OD Tubing)

(1/4-28, > 1/16" OD Tubing)

F-356x

F-364x



LT-105x (1/4-28, ≤ 1/16" OD Tubing) **C-235x** (1/4-28, > 1/16" OD Tubing)



P-213x (M6X1, ≤ 1/16" OD Tubing) P-337x (M6X1, > 1/16" OD Tubing) P-232x

P-232x (1/4-28, ≤ 1/16" OD Tubing) **P-336x** (1/4-28, > 1/16" OD Tubing)



P-287x (1/4-28, ≤ 1/16" OD Tubing) **P-387x** (1/4-28, > 1/16" OD Tubing)



P-219x (M6X1, ≤ 1/16" OD Tubing) P-319x (M6X1, > 1/16" OD Tubing) LT-115x (1/4-28, ≤ 1/16" OD Tubing)




P-217x (M6X1, ≤ 1/16" OD Tubing) **P-317x**

P-317x (M6X1, > 1/16" OD Tubing) P-246x, P-255x, P-281x (1/4-28, < 1/16" OD Tubing) P-331x, P-381x (1/4-28, > 1/16" OD Tubing)



LT-210x P-4 (1/4-28, > 1/16" OD Tubing) (1/4

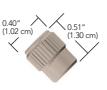

P-420x (1/4-28, ≤ 1/16" OD Tubing) **F-156x** (1/4-28, > 1/16" OD Tubing)



(5/16-24, 3/16" OD Tubing)

P-137x

One-Piece Super Flangeless™ Fittings

- ➤ All-PEEK construction
- **>** For 1/16" OD and 1/8" OD tubing
- ➤ M6x1 and 1/4-28 options
- ➤ Finger tight (2–3 in-lbs / 0.23–0.34 N·m)
- > Extremely easy to use
- > Reusable one piece design that requires no swaging



P-329x M6X1, for 1/8" OD Tubing P-249x 1/4-28, for 1/16" OD Tubing P-349x 1/4-28, for 1/8" OD Tubing



P-229x M6X1, for 1/16" OD Tubing

Super Flangeless Tubing OD / Thread Comparison

, ,	1/32"	1/16"	1.8 MM	2.5 MM	1/8″	3/16"
6-40		✓				
6-32		✓				
10-32		✓				
M6x1	✓	✓	✓	✓	✓	
1/4-28	✓	✓	✓	✓	✓	
5/16-24						~

Qty.

Head Style

Material

Super Flangeless $^{\scriptscriptstyle\mathsf{TM}}$ and One-Piece Super Flangeless Fittings

Description

CLIDED EL	ANGELESS FERRULES FOR 1/32", 1/16", 1/8", 3/16",	1 9MM 2 0MM 2 5MM		Tiona Ctyle		,·
M-250x	Super Flangeless Ferrule for 1/16" OD Tubing	10-32 Flat-Bottom			PEEK Natural/SST	10
M-650x	Super Flangeless Ferrule for 1/16 OD Tubing Super Flangeless Ferrule for 1/16" OD Tubing	6-32 or 6-40 Flat Bottom	1,000–5,000 psi (69–345 bar)		PEEK Natural/SST	10-pk 10-pk
P-248x		10-32 Flat-Bottom	750–3,750 psi (52–259 bar)	_	ETFE Green/SST	10-pk
P-250x	Super Flangeless Ferrule for 1/32" OD Tubing Super Flangeless Ferrule for 1/16" OD Tubing		2,500 psi (172 bar)	_		
	, ,	1/4-28 or M6 Flat Bottom	2,500 psi (172 bar)	_	PEEK Natural/SST	10-pk
P-259x	Super Flangeless Ferrule for 1/16" OD Tubing	1/4-28 or M6 Flat Bottom	1,350 psi (93 bar)	_	ETFE Yellow/SST	10-pk
P-260x	Super Flangeless Ferrule for 1/16" OD Tubing	1/4-28 or M6 Flat Bottom	1,850 psi (128 bar)	_	PEEK Natural/SST	10-pk
P-350x	Super Flangeless Ferrule for 1/8" OD Tubing	1/4-28 Flat Bottom	2,500 psi (172 bar)	_	PEEK Natural/SST	10-pk
P-352x	Super Flangeless Ferrule for 1/16" OD Tubing	1/4-28 or M6 Flat Bottom	2,500 psi (172 bar)	_	PEEK Black/SST	10-pk
P-355x	Super Flangeless Ferrule for 1.8 mm OD Tubing	1/4-28 or M6 Flat Bottom	2,500 psi (172 bar)	_	PCTFE Green/SST	10-pk
P-357-2x	Super Flangeless Ferrule for 2.0 mm OD Tubing	M6 Flat Bottom 1/4-28 Flat Bottom	5,000 psi (345 bar)	_	PEEK Natural/SST	10-pk
P-359x	Super Flangeless Ferrule for 1/8" OD Tubing		1,000 psi (69 bar) 1,500 psi (102 bar)	_	ETFE Yellow/SST	10-pk
P-360x	Super Flangeless Ferrule for 1/8" OD Tubing	1/4-28 Flat Bottom		_	PEEK Natural/SST	10-pk
P-366x P-140x	Super Flangeless Ferrule for 2.5" OD Tubing	1/4-28 Flat Bottom	1,000 psi (69 bar)	_	PEEK Natural/SST	10-pk
	Super Flangeless Ferrule for 3/16" OD Tubing	5/16-24 Flat Bottom	500 psi (34 bar)	_	ETFE Green/SST	10-pk
	6-32 FITTINGS FOR 1/16" OD TUBING					
M-660x	Super Flangeless Nut for 1/16" OD Tubing	6-32 Flat Bottom	750–3,750 psi (52–259 bar)	Micro Headless	PEEK Natural	10-pk
	Super Flangeless Nut for 1/16" OD Tubing	6-40 Flat Bottom	750-3,750 psi (52-259 bar)	Micro Headless	PEEK Green	10-pk
	TINGS FOR 1/16" OD TUBING					
M-652x	Super Flangeless Nut for 1/16" OD Tubing	10-32 Flat Bottom	1,000–5,000 psi (69–345 bar)	1/4" Hex	PEEK Green	10-pk
M-653x	Super Flangeless Nut for 1/16" OD Tubing	10-32 Flat Bottom	1,000–5,000 psi (69–345 bar)	Headless Knurl	PEEK Green	10-pk
M-655x	Super Flangeless Nut for 1/16" OD Tubing, Long	10-32 Flat Bottom	1,000–5,000 psi (69–345 bar)	1/4" Hex	PEEK Green	10-pk
M6X1 FIT	TINGS FOR 1/16" AND 1/32" OD TUBING					
P-213x	Super Flangeless Nut for 1/16" or 1/32" OD Tubing, Short	M6 Flat Bottom	*	Headless Knurl	PEEK Black	10-pk
P-217x	Super Flangeless Nut for 1/16" or 1/32" OD Tubing	M6 Flat Bottom	*	Standard Knurl	PPS Black	10-pk
P-219x	Super Flangeless Nut for 1/16" or 1/32" OD Tubing, Short	M6 Flat Bottom	*	Standard Knurl	PEEK Black	10-pk
M6X1 FIT	TINGS FOR 1.8 MM, 20. MM, 2.5 MM, 1/8" OD TUBI	NG				
P-317x	Super Flangeless For >1/16"-≤ 1/8" OD Tubing	M6 Flat Bottom	*	Standard Knurl	PPS Black	10-pk
P-319x	Super Flangeless Nut for 1/8" OD Tubing, Short	M6 Flat Bottom	*	Standard Knurl	PEEK Black	10-pk
P-337x	Super Flangeless For >1/16"-≤ 1/8" OD Tubing, Short	M6 Flat Bottom	*	Headless Knurl	PEEK Black	10-pk
P-357x	Super Flangeless Fitting for 2.0 mm OD Tubing	M6 Flat Bottom	*	Standard Knurl	PEEK Black, Natural/SST	10-pk
1/4-28 FIT	TTINGS FOR 1/16" AND 1/32" OD TUBING					
F-356x	Super Flangeless Nut for 1/16" or 1/32" OD Tubing, FlushNut	1/4-28 Flat Bottom	*	FlushNut	SST	10-pk
LT-105x	Super Flangeless Nut for 1/16" or 1/32" OD Tubing, Short	1/4-28 Flat Bottom	*	1/4" Hex	SST	10-pk
LT-115x	Super Flangeless Nut for 1/16" or 1/32" OD Tubing, Short	1/4-28 Flat Bottom	*	Standard Knurl	PEEK Natural	10-pk
P-232x	Super Flangeless Nut for 1/16" or 1/32" OD Tubing, Short	1/4-28 Flat Bottom	*	Headless Knurl	PEEK Natural	10-pk
P-246x	Super Flangeless Nut for 1/16" or 1/32" OD Tubing	1/4-28 Flat Bottom	*	Standard Knurl	PFA Natural	10-pk
P-255x	Super Flangeless Nut for 1/16" or 1/32" OD Tubing	1/4-28 Flat Bottom	*	Standard Knurl	PEEK Natural	10-pk
P-281x	Super Flangeless Nut for 1/16" or 1/32" OD Tubing	1/4-28 Flat Bottom	*	Standard Knurl	PPS Natural	10-pk
P-287x	Super Flangeless Nut for 1/16" or 1/32" OD Tubing	1/4-28 Flat Bottom	*	Headless Knurl	PPS Natural	10-pk
P-420x	Super Flangeless Nut for 1/16" or 1/32" OD Tubing, Female	1/4-28 Flat Bottom	*	Female Knurl	PEEK Natural	10-pk
	TTINGS FOR 1.8MM, 2.5 MM, 1/8" OD TUBING					
C-235x	Super Flangeless Nut for 1/8" OD Tubing	1/4-28 Flat Bottom	*	1/4" Hex	SST	10-pk
F-156x	Super Flangeless Nut for 1/8" OD Tubing, Female	1/4-28 Flat Bottom	*	Female Knurl	PEEK Black	10-pk
F-364x	Super Flangeless Nut for 1/8" OD Tubing, FlushNut™	1/4-28 Flat Bottom	*	FlushNut	SST	10-pk
LT-210x	Super Flangeless Nut for 1/8" OD Tubing	1/4-28 Flat Bottom	*	Double Wings	PEEK Natural	10-pk
P-331x	Super Flangeless Nut for 1/8" OD Tubing	1/4-28 Flat Bottom	*	Standard Knurl	PEEK Natural	10-pk
P-336x	Super Flangeless Nut for 1/8" OD Tubing, Short	1/4-28 Flat Bottom	*	Headless Knurl	PEEK Natural	10-pk
P-381x	Super Flangeless Nut for 1/8" OD Tubing	1/4-28 Flat Bottom	*	Standard Knurl	PPS Natural	10-pk
P-387x	Super Flangeless Nut for 1/8" OD Tubing	1/4-28 Flat Bottom	*	Standard Knurl	PPS Natural	10-pk
	ITTINGS FOR 1/16", 1/8", 3/16" OD TUBING	., . Lo . ide Dottoiii		Tanaa a Kilali		.5 pk
P-137x	Super Flangeless Fitting for 3/16" OD Tubing	5/16-24 Flat Bottom	*	Standard Knurl	PEEK Black	10-pk
P-141x	Super Flangeless Fitting for 1/16" OD Tubing	5/16-24 Flat Bottom	*	Standard Knurl	PEEK Natural	10-pk
	CE SUPER FLANGELESS FITTINGS FOR 1/16" AND 1/			Standard Kriuff	I ELIX INGCUI di	10-рк
			1 000 mai (40 k)	Ctondon-I V	DEEK	10 1
P-229x	One Piece Super Flangeless Fitting for 1/16" OD Tubing	M6 Flat Bottom	1,000 psi (69 bar)	Standard Knurl	PEEK	10-pk
P-249x	One Piece Super Flangeless Fitting for 1/16" OD Tubing	1/4-28 Flat Bottom	1,000 psi (69 bar)	Standard Knurl	PEEK	10-pk
P-329x	One Piece Super Flangeless Fitting for 1/16" OD Tubing One Piece Super Flangeless Fitting for 1/16" OD Tubing	M6 Flat Bottom	1,000 psi (69 bar)	Standard Knurl	PEEK	10-pk
P-349x		1/4-28 Flat Bottom	1,000 psi (69 bar)	Standard Knurl	PEEK	10-pk



VacuTight Fittings


- > For 1/16" or 1/8" OD tubing connections into 10-32, 1/4-28, or M6 flat-bottom ports
- > Vacuum Rated to 25 in-Hg (84 kPa)
- > Improve transfer volume consistency

VacuTight Fittings are designed to provide airtight, dependable connections under vacuum and low pressure conditions. Many of the VacuTight Nuts have streamlined profiles for use in systems requiring a large number of connections in a small area. Furthermore, the VacuTight Ferrule's small size ensures sufficient nut/thread engagement, even in shallow ports. These features make VacuTight Fittings ideal for "combichem" high throughput screening, clinical diagnostic, and other automated liquid handling applications.

The configuration of the VacuTight flat-bottom ferrules prevents overcompression and tubing ID reduction that can occur with many coned fittings. The result is more consistent aspirating and dispensing volumes across all system connections.

The VacuTight fittings can also work well in some positive pressure applications. The pressure range for each fitting is listed below and depends upon the tubing used for the connection. Please contact your distributor or IDEX Health & Science for more information. Additionally, please note that some of the VacuTight fittings have changed in color from red to black; however, this color change does not affect product performance.

All VacuTight Nuts must be used exclusively with VacuTight Ferrules.

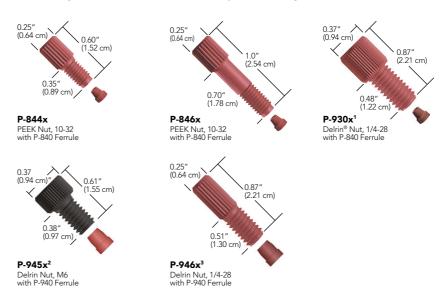

- ¹ The dimensions shown apply to P-930x, P-931x, P-938x, P-942x, and P-948x.
- ² The dimensions shown apply to P-945x ³ The dimensions shown apply to P-946x

Description	Port	Pressure Rating	Head Style	Material	Qty.
T FITTINGS (INCLUDES P-840 OR P-940 FERI	RULES)				
VacuTight Fitting for 1/16" OD Tubing, Short	10-32 Flat-Bottom	400-800 psi (27-55 bar)	1/4" Hex	PEEK Red	10-pk
VacuTight Fitting for 1/16" OD Tubing, Short	10-32 Flat-Bottom	400–800 psi (27–55 bar)	Headless Knurl	PEEK Red	10-pk
VacuTight Fitting for 1/16" OD Tubing, Long	10-32 Flat-Bottom	400–800 psi (27–55 bar)	Headless Knurl	PEEK Red	10-pk
VacuTight Fitting for 1/16" OD Tubing	1/4-28 Flat-Bottom	400–800 psi (27–55 bar)	Standard Knurl	Delrin Red	10-pk
VacuTight Fitting for 1/16" OD Tubing	M6 Flat-Bottom	400–800 psi (27–55 bar)	Standard Knurl	Delrin Red	10-pk
VacuTight Fitting for 1/16" OD Tubing	1/4-28 Flat-Bottom	400–800 psi (27–55 bar)	Standard Knurl	PEEK Natural	10-pk
VacuTight Fitting for 1/8" OD Tubing	1/4-28 Flat-Bottom	500-1,000 psi (34-69 bar)	Standard Knurl	Delrin Red	10-pk
VacuTight Fitting for 1/8" OD Tubing, Short	M6 Flat-Bottom	500-1,000 psi (34-69 bar)	Standard Knurl	Delrin Black	10-pk
VacuTight Fitting for 1/8" OD Tubing	1/4-28 Flat-Bottom	500-1,000 psi (34-69 bar)	Headless Knurl	Delrin Red	10-pk
VacuTight Fitting for 1/8" OD Tubing	1/4-28 Flat-Bottom	500-1,000 psi (34-69 bar)	Standard Knurl	PEEK Natural	10-pk
ENT FERRULES					
VacuTight Ferrule for 1/16" OD Tubing	M6 or 1/4-28 Flat-Bottom	400-800 psi (27-55 bar)	_	ETFE Red	ea.
VacuTight Ferrule for 1/8" OD Tubing	M6 or 1/4-28 Flat-Bottom	500-1,000 psi (34-69 bar)	_	ETFE Red	10-pk
	FITTINGS (INCLUDES P-840 OR P-940 FER! VacuTight Fitting for 1/16" OD Tubing, Short VacuTight Fitting for 1/16" OD Tubing, Short VacuTight Fitting for 1/16" OD Tubing, Long VacuTight Fitting for 1/16" OD Tubing VacuTight Fitting for 1/16" OD Tubing VacuTight Fitting for 1/16" OD Tubing VacuTight Fitting for 1/8" OD Tubing VacuTight Fitting for 1/8" OD Tubing, Short VacuTight Fitting for 1/8" OD Tubing ENT FERRULES VacuTight Ferrule for 1/16" OD Tubing	VacuTight Fitting for 1/16" OD Tubing, Short VacuTight Fitting for 1/16" OD Tubing, Short VacuTight Fitting for 1/16" OD Tubing, Short VacuTight Fitting for 1/16" OD Tubing, Long VacuTight Fitting for 1/16" OD Tubing VacuTight Fitting for 1/8" OD Tubing VacuTight Fitting for 1/8" OD Tubing VacuTight Fitting for 1/8" OD Tubing, Short VacuTight Fitting for 1/8" OD Tubing N6 or 1/4-28 Flat-Bottom M6 or 1/4-28 Flat-Bottom	VacuTight Fitting for 1/16" OD Tubing	VacuTight Fitting for 1/16" OD Tubing, Short 10-32 Flat-Bottom 400-800 psi (27-55 bar) 1/4" Hex	VacuTight Fitting for 1/16" OD Tubing, Short 10-32 Flat-Bottom 400-800 psi (27-55 bar) 1/4" Hex PEK Red VacuTight Fitting for 1/16" OD Tubing, Short 10-32 Flat-Bottom 400-800 psi (27-55 bar) Headless Knurl PEK Red VacuTight Fitting for 1/16" OD Tubing, Long 10-32 Flat-Bottom 400-800 psi (27-55 bar) Headless Knurl PEK Red VacuTight Fitting for 1/16" OD Tubing 1/4-28 Flat-Bottom 400-800 psi (27-55 bar) Headless Knurl PEK Red VacuTight Fitting for 1/16" OD Tubing 1/4-28 Flat-Bottom 400-800 psi (27-55 bar) Standard Knurl Delrin Red VacuTight Fitting for 1/16" OD Tubing 1/4-28 Flat-Bottom 400-800 psi (27-55 bar) Standard Knurl PEK Natural VacuTight Fitting for 1/8" OD Tubing 1/4-28 Flat-Bottom 400-800 psi (34-69 bar) Standard Knurl Delrin Red VacuTight Fitting for 1/8" OD Tubing, Short M6 Flat-Bottom 500-1,000 psi (34-69 bar) Standard Knurl Delrin Black VacuTight Fitting for 1/8" OD Tubing 1/4-28 Flat-Bottom 500-1,000 psi (34-69 bar) Standard Knurl Delrin Red VacuTight Fitting for 1/8" OD Tubing 1/4-28 Flat-Bottom 500-1,000 psi (34-69 bar) Standard Knurl Delrin Red VacuTight Fitting for 1/8" OD Tubing 1/4-28 Flat-Bottom 500-1,000 psi (34-69 bar) Standard Knurl PEK Natural VacuTight Fitting for 1/8" OD Tubing 1/4-28 Flat-Bottom 500-1,000 psi (34-69 bar) Standard Knurl PEK Natural VacuTight Fitting for 1/8" OD Tubing 1/4-28 Flat-Bottom 500-1,000 psi (34-69 bar) Standard Knurl PEK Natural VacuTight Fitting for 1/8" OD Tubing 1/4-28 Flat-Bottom 500-1,000 psi (34-69 bar) Standard Knurl PEK Natural PEK PEK Natural PEK Natural PEK Natural PEK Natural PEK PEK N



Metric Flangeless Fittings

Metric Flangeless Ferrules are designed to connect 1.8, 2.0, 2.5, 3.0, or 4.0 mm OD tubing to flat-bottom ports when paired with the appropriate M6, 1/4-28, or 5/16-24 Flangeless Nuts. We also offer M6-threaded nuts to connect 1/16" or 1/8" OD tubing, plus a tubing sleeve to facilitate 1.0 mm OD tubing connections. Please refer to the "Metric Connections" chart on this page for information regarding which nuts and ferrules to use with your tubing.

- For 1/16", 1.8 mm, 2.0 mm, 2.5 mm, 3.0 mm, 4.0 mm, or 1/8" OD tubing
- Convenience of flangeless fittings for metric tubing sizes and M6 flat-bottom ports












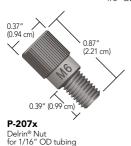

METRIC CONNECTIONS

Use this chart to determine the low pressure fittings needed to connect metric and English-sized tubing into the indicated ports.

Tubing Size	Port	Ferrules	Nuts
1.0 mm	M6	P-200x (w/F-252 sleeve, not included)	P-207x, P-207Sx, P-247x
	1/4-28	P-200x (w/F-252 sleeve, not included)	Any 1/4-28 nut for 1/16" OD tubing from page 47
1.8 mm	M6	P-342x	P-307x, P-307Sx, P-347x
	1/4-28	P-342x	Any 1/4-28 nut for 1/8" OD tubing from page 47
2.0 mm	M6	P-363Rx	P-307x, P-307Sx, P-347x
	1/4-28	P-363Rx	Any 1/4-28 nut for 1/8" OD tubing from page 47
2.5 mm	M6	P-353x	P-307x, P-307Sx, P-347x
	1/4-28	P-353x	Any 1/4-28 nut for 1/8" OD tubing from page 47
3.0 mm	M6	P-343x	P-307x, P-307Sx, P-347x
	1/4-28	P-343x	Any 1/4-28 nut for 1/8" OD tubing from page 47
4.0 mm	5/16-24	P-344x	XP-132x from page 56
1/16"	M6 M6	P-200x P-840	P-207x, P-207Sx, P-247x, P-931, page 42
1/8"	M6 M6	P-300x P-940x	P-307x, P-307Sx, P-347x, P-945x, page 42

	Page		Page
MORE METRIC-SIZED PRODUCTS			
High Pressure Polymer Fittings	58	Low Pressure Unions	84
High Pressure Stainless Steel Fittings	59	Bulkhead Unions	82
Luer Adapters	87	PEEK (1.8 mm OD and Capillary) and Fused Silica Tubing	16
Metric Threaded Adapters	65	PEEKsil™ Tubing	22
External NPT Adapters	66	FEP Tubing (1.0–4.0 mm OD) and PFA Capillary Tubing	26
VacuTight™ Fittings	42	Frit-In-A-Ferrule™	99
Plugs and Cans	55		

In addition, many of our 1/4-28 threaded Filters, Valves and Flow Control Accessories can be converted to accept 1.8, 2.0, 2.5 and 3.0 mm tubing, using the ferrules listed for 1/4-28 ports in the "Metric Connections" table, this page.


Part No.	Description	Port	Pressure Rating	Head Style	Material	Qty.
METRIC F	FLANGELESS NUTS					
P-207x	Flangeless Nut for 1/16" OD Tubing	M6 Flat-Bottom	2,000 psi (138 bar)	Standard Knurl	Delrin Black	10-pk
P-207Sx	Flangeless Nut for 1/16" OD Tubing, Short	M6 Flat-Bottom	2,000 psi (138 bar)	Standard Knurl	Delrin Black	10-pk
P-247x	Flangeless Nut for 1/16" OD Tubing, Short	M6 Flat-Bottom	2,000 psi (138 bar)	Standard Knurl	PEEK Black	10-pk
P-307x	Flangeless Nut for 1.8 mm, 2.0 mm, 3.0 mm, 1/8" OD Tubing	M6 Flat-Bottom	500 psi (34 bar)	Standard Knurl	Delrin Black	10-pk
P-307Sx	Flangeless Nut for 1.8 mm, 2.0 mm, 3.0 mm, 1/8" OD Tubing	M6 Flat-Bottom	500 psi (34 bar)	Standard Knurl	Delrin Black	10-pk
P-347x	Flangeless Nut for 1.8 mm, 2.0 mm, 3.0 mm, 1/8" OD Tubing	M6 Flat-Bottom	500 psi (34 bar)	Standard Knurl	PEEK Black	10-pk
FLANGE	LESS FERRULES					
F-252x	1/16" OD Tubing Sleeve for 1.0 mm ID Tubing	M6 or 1/4-28 Flat-Bottom	500 psi (34 bar)	_	FEP Purple	10-pk
P-200x	Flangeless Ferrule for 1/16" OD Tubing	M6 or 1/4-28 Flat-Bottom	2,000 psi (138 bar)	_	ETFE Blue	10-pk
P-300x	Flangeless Ferrule for 1/8" OD Tubing	M6 or 1/4-28 Flat-Bottom	500 psi (34 bar)	_	ETFE Yellow	10-pk
P-342x	Flangeless Ferrule for 1.8 mm OD Tubing	M6 or 1/4-28 Flat-Bottom	500 psi (34 bar)	_	ETFE Green	10-pk
P-343x	Flangeless Ferrule for 3.0 mm OD Tubing	M6 or 1/4-28 Flat-Bottom	500 psi (34 bar)	_	ETFE Orange	10-pk
P-344x	Flangeless Ferrule for 4.0 mm OD Tubing	5/16-24	250 psi (17 bar)	_	ETFE Natural	10-pk
P-353x	Flangeless Ferrule for 2.5 mm OD Tubing	M6 or 1/4-28 Flat-Bottom	500 psi (34 bar)	_	ETFE Natural	10-pk
P-363Rx	Flangeless Ferrule for 2.0 mm OD Tubing	M6 or 1/4-28 Flat-Bottom	500 psi (34 bar)	_	ETFE Red	10-pk



Flanged Fittings

- Fittings for 1/16" or 1/8" OD tubing, supplied with nut and 316 stainless steel washer
- Multiple head styles and materials available; contact IDEX Health & Science for more information
- For 1/4-28 and M6 flat-bottom ports
- Some color options available; call for more information


Flanged Fittings are compatible with most standard 1/4-28 or M6 Flat-Bottom flanged fittings. The Delrin® (acetal resin) nut resists cross threading or loosening during use.

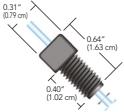

P-401x1

Flanged Fitting for 1/16" OD tubing*

- ¹ The dimensions shown apply to all square-head Flanged Fittings
- * Flanged tubing not included

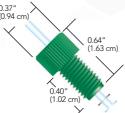

P-482G3

Flanged Fitting for 1/16" OD tubing*

- ³ The dimensions shown apply to all knurled-head Flanged Fittings
- * Flanged tubing not included

RELATED PRODUCTS

For an alternative to flanging tubing, we highly recommend the Flangeless Fittings found on page 45, the Super Flangeless™ Fittings found on page 39, or the VacuTight™ Fittings on page 42.

Part No.	Description	Port Geometry	Head Style	Material (Nut/Washer)	Qty.
FLANGED FITT	INGS (INCLUDES STAINLESS STEEL WASHE	ERS)			
P-401x	Flanged Fitting for 1/16" OD Tubing	1/4-28 Flat-Bottom	5/16" Square	Delrin Black/SST	10-pk
P-482BLK	Flanged Fitting for 1/16" OD Tubing	1/4-28 Flat-Bottom	Standard Knurl	Delrin Black/SST	ea.
P-501x	Flanged Fitting for 1/8" OD Tubing	1/4-28 Flat-Bottom	5/16" Square	Delrin Black/SST	10-pk
P-982BLKx	Flanged Fitting for 1/16" OD Tubing	M6 Flat-Bottom	Standard Knurl	Delrin Black/SST	10-pk
P-1082BLKx	Flanged Fitting for 1/8" OD Tubing	M6 Flat-Bottom	Standard Knurl	Delrin Black/SST	10-pk
REPLACEMENT	WASHERS				
P-407x	Washer for 1/16" OD Tubing	1/4-28 Flat-Bottom	_	SST	10-pk
P-507x	Washer for 1/8" OD Tubing	1/4-28 Flat-Bottom	_	SST	10-pk
P-1087x	Washer for 1/8" OD Tubing	M6 Flat-Bottom	_	SST	10-pk



THE CONVENIENCE OF FLANGELESS FITTINGS



Flangeless Fittings

Flangeless Fittings eliminate the need to flange tubing. This removable and reusable system provides several benefits:

Convenience: Flangeless Fittings are easy to use. Just slip the nut and ferrule over the tubing and finger tighten the assembly into your receiving port. In tests, it is shown that the ideal amount of torque to achieve expected part performance should be approximately 3–4 in-lbs (0.34–0.45 N·m). Check out the line of special tightening tools designed to adapt to many standard torque wrenches, on page 50 and the adjustable torque driver, VHP-4000 on page 51.

Minimal Down-Time: Component replacement is quick, taking only a few seconds — unlike the significant time required to flange tubing.

Cost-Effectiveness: Repairing a flanged tubing assembly requires a costly flanging tool or the purchase of a complete replacement assembly, including a new length of tubing and a set of fittings. The Flangeless Fittings system typically requires only one new ferrule at minimal cost when repairing a connection.

The 1/4-28 and M6 Flangeless Fittings for 1/16", 1/8", and metric sized OD tubing are summarized on the following page and listed on page 47.

1/4-28 Flangeless Fittings – Nuts



Flangeless Fittings (Cont.)



NOTE

The XP-340x ferrule is designed for use with shallow receiving ports, such as those used on some low pressure valves.



XP-340x

ETFE Small Valve Ferrule

The XF-368x FlushNut is an excellent choice for applications where port-to-port spacing is limited; see page 31 for more information on this innovative product line. As an alternative, consider one of the "headless" fittings shown on this page.

Ferrules



Standard 1/16" P-200x XP-200Nx



Small Valve 1/8" XP-340x



Small Valve 1/16" XP-240x



Standard 4.0 mm



Standard 1/8" and Metric Ferrules XP-300x



RELATED PRODUCTS

- ➤ For the Large Bore Flangeless Fittings, please refer to page 56.
- Nuts for M6 threaded ports are on page 43; nuts for 5/16-24 threaded ports are on page 56.

Flangeless Fittings for 1/16" and 1/8"OD Tubing

- Wide variety of materials and geometries to fit most applications
- Excellent replacement for flanged fittings
- Convenient and easy to use
- > Fittings and ferrules packaged together for easy ordering convenience

Flangeless Fittings for 1/16" OD Tubing, and for 1/8" OD Tubing are available in a variety of materials. The replacement ferrules are manufactured from inert ETFE, and are sold in a colored version or ETFE's natural color as the N option. The smaller ferrules, XP-240x and XP-340x are designed for shallow ports.

The fittings shown on this page are packaged together with the ferrules in convenient 10-packs. The 1/16" version include the XP-200x ferrule, and the 1/8" version include the XP-300x ferrule. Please visit our website, www.idex-hs.com, for single packaging options.

For higher pressure and temperature applications consider our Super Flangeless™ found on page 39.

Lock Nut

The XP-312x Lock Nut is for use with any 1/4-28 male Flangeless Fitting. Use this product in applications where vibrations can loosen fittings.

To Use: Thread the lock nut onto the male fitting. When the male fitting is firmly seated into the receiving port, tighten the lock nut down against the receiving port to securely hold the male fitting in place.



XP-312x Lock Nut White Delrin

Flangeless Fittings

Part No.	Description	Port	Pressure Rating	Head Style	Material	Qty.
FLANGE	LESS FITTINGS FOR 1/16" OD TUBING					
XF-358x	Flangeless Fitting for 1/16" OD Tubing, FlushNut	1/4-28 Flat-Bottom	2,000 psi (138 bar)	FlushNut	SST	10-pk
XLT-111x	Flangeless Fitting for 1/16" OD Tubing	10-32 Flat-Bottom	2,500 psi (172 bar)	Standard Knurl	PEEK Natural	10-pk
XP-201x	Flangeless Fitting for 1/16" OD Tubing	1/4-28 Flat-Bottom	2,000 psi (138 bar)	Standard Knurl	Delrin Black	10-pk
XP-202x	Flangeless Fitting for 1/16" OD Tubing	1/4-28 Flat-Bottom	2,000 psi (138 bar)	Standard Knurl	Delrin Red	10-pk
XP-218x	Flangeless Fitting for 1/16" OD Tubing	1/4-28 Flat-Bottom	2,000 psi (138 bar)	Standard Knurl	ETFE Natural	10-pk
XP-230x	Flangeless Fitting for 1/16" OD Tubing	1/4-28 Flat-Bottom	2,000 psi (138 bar)	Standard Knurl	PEEK Natural	10-pk
XP-235x	Flangeless Fitting for 1/16" OD Tubing, Short	1/4-28 Flat-Bottom	2,000 psi (138 bar)	Standard Knurl	PEEK Natural	10-pk
XP-238x	Flangeless Fitting for 1/16" OD Tubing	1/4-28 Flat-Bottom	2,000 psi (138 bar)	Standard Knurl	Delrin Purple	10-pk
XP-245x	Flangeless Fitting for 1/16" OD Tubing	1/4-28 Flat-Bottom	2,000 psi (138 bar)	Standard Knurl	PFA Natural	10-pk
XP-286x	Flangeless Fitting for 1/16" OD Tubing	1/4-28 Flat-Bottom	2,000 psi (138 bar)	Headless Knurl	PPS Natural	10-pk
REPLACE	MENT FERRULES FOR 1/16" OD TUBING					
XP-200x	Flangeless Ferrule for 1/16" OD Tubing	1/4-28 Flat-Bottom	2,000 psi (138 bar)	_	ETFE Blue	10-pk
XP-200Nx	Flangeless Ferrule for 1/16" OD Tubing	1/4-28 Flat-Bottom	2,000 psi (138 bar)	_	ETFE Natural	10-pk
XP-240x	Flangeless Ferrule for 1/16" OD Tubing, Small Valve	1/4-28 or 10-32 Flat-Bottom	2,500 psi (172 bar)	_	ETFE Natural	10-pk
FLANGE	LESS FITTINGS FOR 1/8" OD TUBING (INCLUDE	S P-300 FERRULES)				
XF-368x	Flangeless Fitting for 1/8" OD Tubing, FlushNut	1/4-28 Flat-Bottom	500 psi (34 bar)	FlushNut	SST	10-pk
XP-301x	Flangeless Fitting for 1/8" OD Tubing	1/4-28 Flat-Bottom	500 psi (34 bar)	Standard Knurl	Delrin Black	10-pk
XP-302x	Flangeless Fitting for 1/8" OD Tubing	1/4-28 Flat-Bottom	500 psi (34 bar)	Standard Knurl	Delrin Red	10-pk
XP-305x	Flangeless Fitting for 1/8" OD Tubing	1/4-28 Flat-Bottom	500 psi (34 bar)	Standard Knurl	Delrin Green	10-pk
XP-308x	Flangeless Fitting for 1/8" OD Tubing, Short	1/4-28 Flat-Bottom	500 psi (34 bar)	Standard Knurl	Delrin Black	10-pk
XP-315x	Flangeless Fitting for 1/8" OD Tubing	1/4-28 Flat-Bottom	500 psi (34 bar)	Standard Knurl	ETFE Natural	10-pk
XP-330x	Flangeless Fitting for 1/8" OD Tubing	1/4-28 Flat-Bottom	500 psi (34 bar)	Standard Knurl	PEEK Natural	10-pk
XP-335x	Flangeless Fitting for 1/8" OD Tubing, Short	1/4-28 Flat-Bottom	500 psi (34 bar)	Standard Knurl	PEEK Natural	10-pk
REPLACE	MENT FERRULES FOR 1/8" OD TUBING					
XP-300x	Flangeless Ferrule for 1/8" OD Tubing	1/4-28 Flat-Bottom	500 psi (34 bar)	_	ETFE Yellow	10-pk
XP-300Nx	Flangeless Ferrule for 1/8" OD Tubing	1/4-28 Flat-Bottom	500 psi (34 bar)	_	ETFE Natural	10-pk
XP-340x	Flangeless Ferrule for 1/8" OD Tubing, Small Valve	1/4-28 Flat-Bottom	500 psi (34 bar)	_	ETFE Natural	10-pk
XP-312x	Lock Nut for Flangeless Nuts	1/4-28 Flat-Bottom	_	_	Delrin White	10-pk
	-					




- Finger-tight to 19,000 psi
- Reusable more than 100 times
- Zero Dead Volume
- Biocompatible
- Limit Wear
- Increase Product Life
- Many IDs and Lengths Available



Learn more about MarvelXact at www.biotechfluidics.com/ products/marvelxact/

NOTE

MarvelXACT tubing includes a sleeve that assists in product identification, with ID, length and part number information.

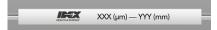


Minimum recommended bend-radius with MarvelXACT tubing is 1/4" (~6.35 mm).

MarvelXACTTM

Our New MarvelXACT™ Fitting System with a "click" feedback

Ensure a perfect connection every time with our new MarvelXACT™ fitting system for trouble-free liquid chromatography. Our MMarvelXACT™ connection systems have been expertly designed to eliminate the risk of under- or over-tightening with our patented torque limiting mechanism. This unique feature emits a haptic "click" feedback when it reaches the optimum torque, assuring a perfect installation every time. MarvelXACT™ incorporates our advanced MarvelX™ Sealing Technology to deliver precise face sealing (sealing at the port bottom), which eliminates additional internal volume, and minimizes carryover risk, peak tailing, and peak broadening.

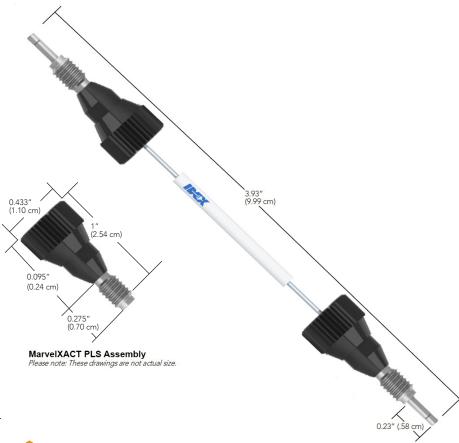

SPECIFICATIONS & DETAILS

Pressure Capability	19,000 psi (~1,310 bar, 131 MPa) for routine use
Installation Method	Finger-tighten until the first "click" feedback is received
Tubing Type	1/32" OD flexible 316 Stainless Steel with 1/16" OD rigid tube ends
Fitting Type	10-32 threaded, PEEK fittings with 316 Stainless Steel threads
Wetted Materials PEEK- Lined versions	PEEK
Stainless Steel versions	PEEK and 316 Stainless Steel
Maximum Use Temperature	120 °C

NOTE: The above performance specifications apply to use with appropriately-designed receiving ports under optimal conditions, using water at up to 120 °C for the testing process. If different conditions are used, the expected pressure threshold will received.


first "click" feedback is




MARVELXACT™ UHPLC FITTINGS VS CONVENTIONAL CONED FITTINGS

Conventional coned fittings require a ferrule in conjunction with a fitting for proper sealing. They depend on complex techniques, including tools, to improve sealing performance, which significantly increases probability of extra internal volume and poor chromatography results. The excessive force needed for tightening increases wear of expensive components and the likelihood of replacement, adding to overall costs.

MarvelXACTTM fittings do not depend on ferrules. They seal at the bottom of the port, without complex techniques, which significantly reduces required torque and enables many more connects and disconnects. MarvelXACT significantly reduces wear on your hardware, increasing product life. An enhanced proprietary tip design also ensures zero dead volume (ZDV) and better chromatography results.

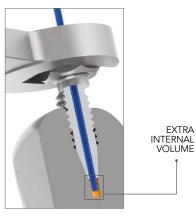

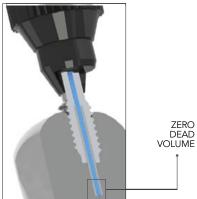


TECHNICAL SPECIFICATIONS


Peek-lined Stainless Steel Assemblies* 25 μm ID UPFP-7025150 UPFP-7025250 UPFP-7025350 UPFP-7025500 UPFP-7025 50 μm ID UPFP-7050150 UPFP-7050250 UPFP-7050350 UPFP-7050500 UPFP-7050	m								
	Peek-lined Stainless Steel Assemblies*								
50 μm ID UPFP-7050150 UPFP-7050250 UPFP-7050350 UPFP-7050500 UPFP-7050	5600								
	0600								
75 μm ID UPFP-7075150 UPFP-7075250 UPFP-7075350 UPFP-7075500 UPFP-7075	5600								
100 μm ID UPFP-7100150 UPFP-7100250 UPFP-7100350 UPFP-7100500 UPFP-7100	0600								
Stainless Steel Assemblies*									
100 μm ID UPFS-7100150 UPFS-7100250 UPFS-7100350 UPFS-7100500 UPFS-7100	0600								
125 μm ID UPFS-7125150 UPFS-7125250 UPFS-7125350 UPFS-7125500 UPFS-7125	5600								
254 μm ID UPFS-7254150 UPFS-7254250 UPFS-7254350 UPFS-7254500 UPFS-7254	4600								


 $[*]Product\ availability\ and\ lead\ times\ may\ vary\ depending\ on\ the\ configuration.$



VHP-1000 VHP Wrench 1/4 in Hex 10 in-lbs (1.13 N·m)

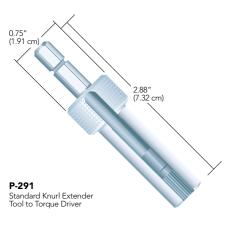

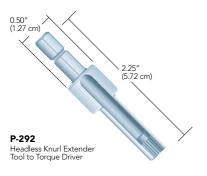

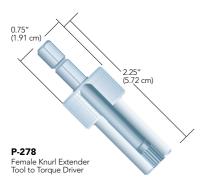

Fittings Tools

Tightening Tools for VHP & Other Fittings

- Configured for the optimum torque to provide assurance of a strong connection
- Prolongs the lifetime of reusable fittings by not overtightening
- Available for multiple fitting head styles

This new line of tightening tools is designed for the VHP fittings and can also be used with any fitting in this chapter described to have a corresponding head style to the tool listed below. There are three styles of tightening tools available for various applications. The Torque Tools (VHP-1000 and VHP-2000) are breakaway torque wrenches designed to deliver a precise amount of torque to the fitting system. These torque wrenches come calibrated according to ISO 6789:2003 (± 6% of setting) and have been tested extensively with the reusable VHP fittings on page 60. Choose the appropriate torque delivered and the proper head style to work with the VHP fittings, increasing the ease of use with these fittings.

The VHP-4000 Torque Driver couples with the specially designed Extender Tools listed below and provides an externally adjustable torque setting. This tool along with the appropriate Extender Tools will tighten any IDEX Health & Science knurled polymer fitting in your system. Reference the head style found in the tables at the bottom of each page for information on the proper Extender Tool to select.

Because of the small hex-head on the M4 fittings (VHP-900 and VHP-920), a custom wrench, the VHP-9000, is available in the table at the bottom of the page 51.

Extender Tools

These tools can be used to tighten most of our knurled nuts in hard to reach places. See the application note on this page for knurl size and corresponding extender tool.

For precise tightening, the extender tools listed with 1/4" hex drives are designed to adapt to any torque wrench with a female 1/4" socket, such as the VHP-4000 Torque Driver on page 51. The tools featured on this page also include the FlushNut™ wrenches, used to tighten the FlushNuts found throughout this chapter and described in detail on page 51.



The drawings represent actual size of the various knurled head designs of the IDEX Health & Science nuts featured in this chapter. Select the appropriate extender tool for the knurl pattern of the nut you've selected.







Removal Tool

Use the LT-300 Removal Tool to detach LiteTouch® and Super Flangeless™ Ferrules from tubing. Simply slide the appropriate tool blade slot between the lock ring and the ferrule body. With a slight twist, the ring will pop off, releasing the ferrule from the tubing. Please Note: This Removal Tool will not work with the LT-135x Ferrule System.

Wrenches

For your convenience, we offer wrenches in three standard sizes. You will need two A-304 wrenches to tighten most nuts into unions found on page 36 (for union 1593, you need one A-304 and one A-320 wrench).

The IDEX Wrench is slotted to fit over 1/16" and 1/8" OD tubing, and has 1/4" and 5/16" internal hex ends, to engage with the heads of the hex-head fittings most commonly used with IDEX Health & Science valves and the stainless steel fittings listed on page 59.

Fittings Tools

Part No.	Description	Use With Head Style	Torque Delivered	Qty.
VHP TIGHTENI	ING TOOLS			
F-347	Extender Tool to Torque Driver	FlushNut (1/4-28)	_	ea.
N-291	Extender Tool to Torque Driver	Micro Headless	_	ea.
P-268	Extender Tool to Torque Driver	1/4" Hex	_	ea.
P-278	Extender Tool to Torque Driver	Female Nut Knurl	_	ea.
P-279	Extender Tool to Torque Driver	Micro Nut Knurl	_	ea.
P-291	Extender Tool to Torque Driver	Standard Nut Knurl	_	ea.
P-292	Extender Tool to Torque Driver	Headless Nut Knurl	_	ea.
P-1000	Standard Knurl Torque Tool	Standard Knurl	4 in-lbs (0.45 N·m)	ea.
VHP-1000	VHP Torque Tool	1/4" Hex	10 in-lbs (1.13 N·m)	ea.
VHP-2000	VHP Torque Tool	1/4" Hex	14 in-lbs (1.58 N·m)	ea.
VHP-4000	VHP Torque Driver	Extender Tool 1/4" Drive	Adjustable between 2–12 in-lbs (0.23–1.35 N·m)	ea.
VHP-9000	4 mm Wrench	4 mm Hex	_	ea.
EXTENDER TO	OLS			
Part No.	Description	Material		Qty.
P-291	Extender Tool for Standard Head Nuts, with 1/4" Hex Drive	Aluminum		ea.
P-298	Extender Tool for Standard Head Nuts	Delrin®		ea.
P-299	Extender Tool for Standard Head Nuts	Aluminum		ea.
P-399	Extender Tool for Standard Head Nuts, Short	Aluminum		ea.
P-297	Extender Tool for Headless Nuts	Aluminum		ea.
P-292	Extender Tool for Headless Nuts, with 1/4" Hex Drive	Aluminum		ea.
P-277	Extender Tool for Standard Micro Nuts	Aluminum		ea.
N-290	Extender Tool for Micro Headless Nuts	Aluminum		ea.
P-278	Extender Tool for Female Nuts, with 1/4" Hex Drive	Aluminum		ea.
MISCELLANEC	OUS TOOLS			
A-304	Wrench, 1/4" x 5/16"	Steel		ea.
A-305	Wrench, 1/2" x 9/16"	Steel		ea.
A-320	Wrench, 3/8" x 7/16"	Steel		ea.
5810	IDEX Wrench, 1/4" x 5/16"	Steel		ea.
F-345	FlushNut Wrench for 10-32 Threaded Fittings	Steel/Plastic Handle		ea.
F-346	FlushNut Wrench for 1/4-28 Threaded Fittings	Steel/Plastic Handle		ea.
LT-300	Removal Tool for LiteTouch and Super Flangeless Ferrules	Steel/Plastic Handle		ea.
M-150	Swaging Tool for TinyTight Fittings, for 6–40 MINSTAC Port	SST		ea.





Tubing Sleeves



MicroTight® Tubing Sleeves

- Manufactured from PEEK polymer
- Pressure rated to 4,000 psi (276 bar)
- Color-coded for easy inner diameter identification

IDEX Health & Science MicroTight Tubing Sleeves feature an outer diameter of 0.025" and offer a wide assortment of inner diameters to help facilitate capillary tubing connections with our MicroTight accessories. Because the sleeves are manufactured from PEEK polymer, they carry an upper temperature threshold of 125 °C.

To use these sleeves properly, choose a sleeve with an inner diameter 0.001"-0.002" (25–50 µm) larger than the outer diameter of your capillary tubing. Then, slip the sleeve over your flow path tubing, such that your tubing extends all the way through the sleeve, but not beyond the end of the sleeve. Choose the correct fitting that corresponds with your receiving port, slide it over the sleeved flow path tubing and connect as normal.



NanoTight[™] Tubing Sleeves

- Manufactured from FEP fluoropolymer
- Pressure rated to 4,000 psi (276 bar)
- ▶ Outer diameter of 1/16" the most popular size used on most instrumentation

NanoTight Tubing Sleeves are manufactured using FEP fluoropolymer and precisely cut to a 1.6" length. A wide assortment of sleeves is available, ensuring the availability of a NanoTight sleeve for most applications. Many of the sleeves feature a light color tint that can help more easily identify the inner diameter for future orders. Because FEP is the base polymer for these sleeves, there is a maximum recommended continuous operating temperature of 50 °C.

Our NanoTight sleeves were designed primarily for use with the NanoTight fittings, found on page 37 and also work well with the Super Flangeless™ fittings for 1/16" OD tubing on page 39. For tubing sleeves that can be used effectively with stainless steel fittings and at higher temperatures, consider using the PEEK Tubing Sleeves, found below.



1/16" OD PEEK Tubing Sleeves

- For connecting capillary tubing to standard 10-32 ports
- > Require the use of wrench tightened stainless steel nuts
- Pressure rated to 6,000 psi (414 bar)

Like the NanoTight™ FEP Sleeves on the previous page, these PEEK Tubing Sleeves are designed to be used with 1/16" OD, 10-32 threaded fittings to adapt capillary tubing to standard coned ports. Made of PEEK polymer, these 1.3" long sleeves can be used up to

These sleeves require a wrench tightened nut to achieve proper sealing. We recommend our SealTight™ fittings on page 36. Many researchers also use a stainless steel nut and ferrule with these sleeves, such as our U-400 and U-401 (page 33).





1/32" OD PEEK Tubing Sleeves

These 1.6" long 1/32" OD PEEK Tubing Sleeves can be used with any fitting designed for 1/32" OD tubing when smaller tubing must be connected. Select the appropriate sleeve from the product listing for your capillary tubing OD size. The 1/32" OD PEEK Tubing Sleeves have

a maximum recommended temperature of 125 °C and have a pressure rating of 5,000 psi (345 bar).

1/32" OD FEP Tubing Sleeves

These 1.6" long sleeves facilitate connecting capillary tubing into ports designed for 1/32" OD tubing. Please refer to the product listing below to select the appropriate sleeve for your capillary OD size. These sleeves can be used at up to 50 °C and have a pressure rating of 1,750 psi (121 bar).

Clockwise, starting at top:

- ▶ 1/16" OD PEEK Tubing Sleeves, use with U-400 and U-401 (not shown)
- ▶ 1/32" OD FEP Tubing Sleeves, shown with F-126Sx Fitting
- Fittings and tubing only shown to highlight how sleeves are designed to be used; they are not included with the sleeves





Why use Sleeves?

Because most capillary tubing connections are made into coned receiving ports, where the port is not designed to be used with capillary tubing directly, special care must be used to ensure a good connection. While custom ferrules can help make these connections, they only offer a fixed-length nose — and because most tubing pockets will vary slightly in length, this can lead to leaking or dead volume.

To help save overall expense while maintaining a concentric connection with minimal dead volume, IDEX Health & Science recommends the use of sleeves. Because sleeves are not permanently attached to a ferrule, they can easily adapt to varying tubing pocket depths. Additionally, because they are manufactured using extruded polymer tubing, you are assured of the concentricity of the resultant connection.

Tubing Sleeves (Cont.)



Use 1/32" OD PEEK or FEP Sleeves to connect capillary tubing with the following:

- The F-113 Ferrule and Two-Piece Fingertight Fittings for 10-32 ports (page 35).
- The F-112 and P-416BLK MicroTight® Fittings (page 34) 1/32" OD PEEK Tubing Sleeves only.
- ▶ The 1/32" OD MicroTight Fittings on page 34.
- The RheFlex M4 Fitting (page 61) for MX Module applications; the M-645 Valco®-Compatible Fitting (page 33) for Valco Nanovolume® valve applications.

Tubing Sleeves


Part No.	ID	For Tubing OD Size	Color	Qty.
MICROTIC	HT PEEK TUBING SLEEVES AND KITS, 0.025" OD			
F-180	125 µm (0.005")	70–110 μm	Red	ea.
F-181	180 µm (0.007")	125–165 μm	Yellow	ea.
F-182	230 µm (0.009")	175–215 μm	Natural	ea.
F-183	280 μm (0.011")	225–265 μm	Blue	ea.
F-184	330 µm (0.013")	275–315 μm	Orange	ea.
F-185	395 μm (0.0155")	340–380 μm	Green	ea.
F-186	455 µm (0.018")	400–440 μm	Black	ea.
F-187	535 µm (0.021")	480–520 μm	Natural	ea.
F-188	152 µm (0.006")	95–135 μm	Purple	ea.
1328	MicroTight Tubing Sleeve Kit, contains (6) each of the sleeve sizes listed above	_	_	ea.
1356	MicroTight Connector Kit, contains: a 10-pack of each MicroTight Tubing Sleeve (F-180–F-187); (2) P-770 MicroTight Adapters; and (2) MicroTight P-720 Unions	_	_	ea.
NANOTIG	HT FEP TUBING SLEEVES, 1/16" OD			
F-237	125 µm (0.005")	70–110 μm	Red	ea.
F-238	180 µm (0.007")	125–165 µm	Yellow	ea.
F-239	215 µm (0.0085")	160-200 µm	Natural	ea.
F-240	280 μm (0.011")	225–265 µm	Blue	ea.
F-241	330 µm (0.013")	275–315 μm	Orange	ea.
F-242	395 μm (0.0155")	340–380 μm	Green	ea.
F-243	455 µm (0.018")	400–440 μm	Black	ea.
F-244	535 µm (0.021")	480–520 μm	Natural	ea.
F-245	610 µm (0.024")	555–595 μm	Red	ea.
F-246	685 µm (0.027")	630–670 µm	Yellow	ea.
F-247	840 µm (0.033")	785–825 μm	Green	ea.
F-252	1.07 mm (0.042")	1 mm	Purple	ea.
PEEK TUB	ING SLEEVES FOR 1/16" OD FITTINGS			
F-225	125 µm (0.005")	70–110 μm	Red	ea.
F-226	180 µm (0.007")	125–165 µm	Yellow	ea.
F-227	230 µm (0.009")	175–215 µm	Yellow	ea.
F-228	250 µm (0.011")	225–265 µm	Blue	ea.
F-229	330 µm (0.013")	275–315 μm	Natural	ea.
F-230	405 μm (0.016")	350–390 µm	Orange	ea.
F-231	560 µm (0.022")	505–545 μm	Natural	ea.
F-232	785 µm (0.031")	730–770 µm	Natural	ea.
F-233	865 µm (0.034")	785–825 µm	Blue	ea.
F-234	685 µm (0.027")	630–670 µm	Yellow	ea.
PEEK TUB	ING SLEEVES FOR 1/32" OD FITTINGS	·		
F-381	180 µm (0.007")	125–165 µm	Yellow	ea.
F-382	205 µm (0.008")	150–190 µm	Natural	ea.
F-384	255 µm (0.010")	200–240 µm	Blue	ea.
F-385	380 µm (0.015")	325–365 µm	Natural	ea.
F-386	510 µm (0.020")	455–495 µm	Orange	ea.
F-387	250 µm (0.011")	225–265 µm	Red	ea.
F-388	330 µm (0.013")	275–315 μm	Black	ea.
FEP TUBIN	NG SLEEVES FOR 1/32" OD FITTINGS			
F-374	280 µm (0.011")	225–265 µm	Blue	ea.
F-375	330 µm (0.013")	275–315 μm	Orange	ea.
F-376	395 µm (0.0155")	340–380 µm	Green	ea.
1-570	373 μπ (0.0133)	340–300 μπ	Green	C



Plugs & Caps

Seal 6-32, 6-40, 10-32, 1/4-28, M6, or 5/16-24 threaded ports or fittings Use our plugs to close off unused ports in valves and multi-port connectors. Our color-coded 10-32 threaded plugs are perfect for identifying stored columns that have different packing materials, or in which different mobile phases have been utilized. Cap off tubing with one of the PEEK or ETFE caps presented on this page and the appropriate fittings from this chapter.


To help determine which plug or cap is best suited for your application, please visit www.idex-hs.com for detailed chemical compatibility data.

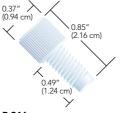

P-316 PFA Plug for 1/4-28 flat-bottom ports



U-467RDelrin® Column Plug for 10-32 coned ports

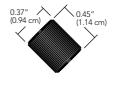

P-755 ETFE Cap for 1/4-28 flat-bottom fittings



Part No.	Description	Head Style	Material	Qty.
PLUGS				
P-120	Plug for 1/4-28 Coned Ports for 1/8" OD Tubing	Standard Knurl	PCTFE Natural	ea.
P-123	Plug for 1/4-28 Flat-Bottom Ports	5/16" Hex	ETFE Natural	ea.
P-309	Plug for 1/4-28 Flat-Bottom Ports	Standard Knurl	Delrin Black	ea.
P-311	Plug for 1/4-28 Flat-Bottom Ports	Standard Knurl	ETFE Natural	ea.
P-314	Plug for M6 Flat-Bottom Ports	Standard Knurl	ETFE Black	ea.
P-316	Plug for 1/4-28 Flat-Bottom Ports	Standard Knurl	PFA Natural	ea.
P-321	Plug for 1/4-28 Flat-Bottom Ports, FlushNut™	FlushNut	PEEK Natural	ea.
P-520	Plug for 10-32 Coned Ports	5/16" Hex	SST	ea.
P-550	Plug for 10-32 Coned Ports, Extra Long	Standard Knurl	PEEK Natural	ea.
P-551	Plug for 10-32 Coned Ports	Standard Knurl	PEEK Natural	ea.
P-555	Plug for 6-32 Coned Ports	Standard Micro Knurl	PEEK Natural	ea.
P-556	Plug for 5/16-24 Flat-Bottom Ports	Standard Knurl	PEEK Natural	ea.
P-558	Plug for 6-40 Flat-Botton Ports	Micro Headless Knurl	PEEK Green	ea.
P-849	Plug for 10-32 Flat-Bottom Ports	Standard Knurl	Delrin Black	ea.
U-467R	Plug for 10-32 Coned Ports	Standard Knurl	Delrin Red	ea.
VHP-600	VHP Plug for 10-32 Coned Ports	3/8" Hex	PK-SST	ea.
CAPS				
P-754	Cap for 10-32 Coned Ports	Standard Knurl	ETFE Yellow	ea.
P-755	Cap for 1/4-28 Flat-Bottom Ports	Standard Knurl	ETFE Black	ea.
P-756	Cap for M6 Flat-Bottom Ports	Standard Knurl	ETFE Blue	ea.

Large Bore Fittings

0.88

- **>** 5/16-24 or 1/2-20 threads
- For use with 1/16", 1/8", 3/16", 1/4", 5/16", 3.0 mm, or 4.0 mm OD tubing

NOTE

Each of the Large Bore Fittings shown on this page comes in a convenient 10-pack and is packaged with the most popularly chosen Ferrule option. The Fittings can be ordered separately by removing the preceding letter "X" from the part number. Additionally, to connect metric-sized tubing with outer diameters less than 4.0 mm to 5/16-24 threaded ports, reference the chart on page 43 to choose the correct nut/ ferrule combination.

RELATED PRODUCTS

	Page
MORE LARGE BORE PRODUCTS	
5/16-24 Coned Fittings	41
Barbed Adapters	90
Threaded Adapters	63
Plugs	55, 112
Y Connector	90
ETFE Tubing	27
FEP, PFA Tubing	26, 24
Tubing Cutter	28
Pressure Relief Valve	148
Solvent Inlet Filters	102
Bottle Caps	55
Semi-Prep Inline Filters	107



XP-130x

PEEK Nut, for 1/8" OD tubing shown with P-300 Flangeless Ferrule (page 43)

5/16-24

PEEK Nut, for 1/16" OD tubing shown with P-200 Flangeless Ferrule (included and found on page 43)

XP-136x

0.50



XP-137x

XP-131x

PEEK Nut, for 3/16" OD tubing shown with P-140 Super Flangeless Ferrule (included and found on page 56)



XP-132x

PEEK Nut, for 3/16" and 4.0 mm OD tubing shown with P-133 Flangeless Ferrule (included and found on this page)



XP-141x

PEEK Nut, for 1/16" OD tubing shown with P-259 Super Flangeless Ferrule (included and found on page 41)



XU-620x

PEEK Nut, for 1/4" OD tubing shown with inverted U-650 Flangeless Ferrule for coned ports (included and found on this page)

PEEK Nut, for 1/4" OD tubing shown with U-650 Flangeless Ferrule (included and found on this page)

XU-662x PEEK Nut, for 5/16" OD tubing shown with U-660 Flangeless Ferrule (included and found on this page)


Part No.	Description	Port	Pressure Rating	Head Style	Material (Nut/Washer)	Qty.
LARGE BO	DRE FITTINGS					
XP-130x	Flangeless Fitting for 1/8" OD tubing	5/16-24 Flat-Bottom	500 psi (34 bar)	Standard Knurl	PEEK Natural/ETFE Yellow	10-pk
XP-131x	Super Flangeless Fitting for 1/8" OD tubing	5/16-24 Flat-Bottom	1,000 psi (69 bar)	Standard Knurl	PEEK Natural/ETFE Yellow/SST	10-pk
XP-132x	Flangeless Fitting for 3/16" OD tubing	5/16-24 Flat-Bottom	500 psi (34 bar)	Standard Knurl	PEEK Natural/ETFE Blue	10-pk
XP-136x	Flangeless Fitting for 1/16" OD tubing	5/16-24 Flat-Bottom	2,000 psi (138 bar)	Standard Knurl	PEEK Natural/ETFE Blue	10-pk
XP-137x	Super Flangeless Fitting for 3/16" OD tubing	5/16-24 Flat-Bottom	500 psi (34 bar)	Standard Knurl	PEEK Black/ETFE Green/SST	10-pk
XP-141x	Super Flangeless Fitting for 1/16" OD tubing	5/16-24 Flat-Bottom	1,350 psi (93 bar)	Standard Knurl	PEEK Natural/ETFE Yellow/SST	10-pk
XU-620x	Flangeless Fitting for 1/4" OD tubing	1/2-20 Coned	250 psi (17 bar)	Large Knurl	PEEK Red/ETFE Natural	10-pk
XU-655x	Flangeless Fitting for 1/4" OD tubing	1/2-20 Flat-Bottom	250 psi (17 bar)	Large Knurl	PEEK Black/ETFE Natural	10-pk
XU-662x	Flangeless Fitting for 5/16" OD tubing	1/2-20 Flat-Bottom	250 psi (17 bar)	Large Knurl	PEEK Black/ETFE Natural	10-pk
REPLACE	MENT FERRULES					
P-133x	Flangeless Ferrule for 3/16" OD tubing	5/16-24 Flat-Bottom	500 psi (34 bar)	_	ETFE Blue	10-pk
P-133Nx	Flangeless Ferrule for 3/16" OD tubing	5/16-24 Flat-Bottom	500 psi (34 bar)	_	ETFE Natural	10-pk
P-140x	Super Flangeless Ferrule for 3/16" OD tubing	5/16-24 Flat-Bottom	500 psi (34 bar)	_	ETFE Green	10-pk
U-650x	Flangeless Ferrule for 1/4" OD tubing	1/2-20 Flat-Bottom	250 psi (17 bar)	_	ETFE Natural	10-pk
U-660x	Flangeless Ferrule for 5/16" OD tubing	1/2-20 Flat-Bottom	250 psi (17 bar)	_	ETFE Natural	10-pk



VHP Micro Fittings

Micro Fittings are specifically designed for use with microferrules. They are manufactured from a proprietary PEEK blend (PK) which allow them to be used at higher temperatures (up to 200° C) and higher pressures ideal for UHPLC applications.

VHP MicroFerrules and Fittings are made from a proprietary high performance PEEK polymer blend, a material which is unique in its ability to enable the use of capillary tubing in UHPLC environments. The new high pressure MicroFerrules are available for use with 1/32" or 360 µm OD tubing, and they are incorporated into several of our VHP products for capillary tubing.

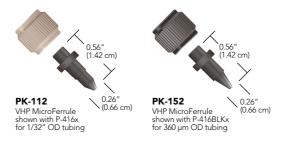



Caution: While the proprietary blend of the PK fittings will allow a fitting to attain a higher pressure and minimal cold flow properties relative to pure PEEK, some fittings molded of PK are known to be conductive. Use caution when employing PK fittings in high voltage applications.



MicroTight fittings and MicroFerrules

While the MicroTight Female Nuts may be used with any of the separate MicroFerrules, the MicroFerrules themselves are port-specific and are thus not interchangeable. Additionally, the one-piece MicroTight fittings are also port-specific and should not be exchanged.

Part No.	Description	Port	Pressure Rating	Required Torque	Head Style	Material	Qty.
PK MICRO	FERRULES AND FEMALE NUTS						
P-416	Female Nut for Microferrule	5/16-24 Coned	15,000 psi (1,035 bar)	4.0 in-lbs (0.45 N·m)	Female Knurl	PEEK, Natural	ea.
P-416BLK	Female Nut for Microferrule	5/16-24 Coned	15,000 psi (1,035 bar)	4.0 in-lbs (0.45 N·m)	Female Knurl	PEEK, Black	ea.
PK-112	VHP MicroFerrule for 1/32" OD Tubing	5/16-24 Coned	15,000 psi (1,035 bar)	_	_	PK	ea.
PK-152	VHP MicroFerrule for 360 µm OD Tubing	5/16-24 Coned	15,000 psi (1,035 bar)	_	_	PK	ea.





VHP PK Fittings



> Find unions, tees and crosses for VHP applications on page 75, and 74.

Ultra High Performance fittings are manufactured from a proprietary PEEK blend (PK) which allow them to be used at higher temperatures (up to 200 °C) and higher pressures.

The VHP PK One-Piece fittings are available for 10-32 coned, 6-32 coned, or M4 coned ports, and Two-Piece fittings are available to connect either 1/16" or 1/32" OD tubing into 10-32 coned ports in multiple styles.













PK-126Hx 6-32 PK MicroTight® Fitting for 1/32" OD tubing



PK-126x 6-32 PK MicroTight® Fitting for 1/32" OD tubing



UH-904x M4, 1/32" Fitting for IDEX Health & Science

Part No.	Description	Port	Pressure Rating	Required Torque	Head Style	Material	Qty.
PK VHP ON	E-PIECE FITTINGS						
PK-120BLKx	PK One-Piece Fitting for 1/16" OD Tubing	10-32 Coned	12,000 psi (827 bar)	8.0 in-lbs (0.90 N·m)	Standard Knurl	PK	10-pk
PK-126Hx	PK One-Piece Headless Fitting for 1/32" OD Tubing	6-32 Coned	15,000 psi (1,035 bar)	3.0 in-lbs (0.34 N·m)	Headless Micro Knurl	PK	10-pk
PK-126x	PK One-Piece Fitting for 1/32" OD Tubing	6-32 Coned	15,000 psi (1,035 bar)	3.0 in-lbs (0.34 N·m)	Standard Micro Knurl	PK	10-pk
UH-904x	PK One-Piece Fitting for 1/32" OD Tubing	M4 Coned	15,000 psi (1,035 bar)	4.0 in-lbs (0.45 N·m)	Headless Knurl	PK	10-pk
PK VHP FIT	TINGS (LITETOUCH® STYLE, NUTS AND FER	RULES SOLD SE	PARATELY)				
PK-100x	PK Ferrule for 1/16" OD Tubing	10-32 Coned	16,500 psi (1,140 bar)	_	_	PK	10-pk
PK-110x	PK Nut for 1/16" OD Tubing	10-32 Coned	16,500 psi (1,140 bar)	8.0 in-lbs (0.90 N·m)	Standard Knurl	PK	10-pk
PK-132x	PK Ferrule for 1/32" OD Tubing	10-32 Coned	16,500 psi (1,140 bar)	_	_	PK	10-pk



Stainless Steel VHP Fittings

- > Pressure rated to 30,000 psi (2,070 bar)
- > Double compression ferrule design
- Available with 10-32 threads for 1/16" OD tubing and M4 threads for 1/32" OD tubing



APPLICATION NOTE

In order to seal up to the stated pressure rating, the VHP-200-01 ferrule requires 20 in-lbs ($2.25~N\cdot m$) of torque. Similar ferrules on the market require tightening torque of at least 30 in-lbs ($3.3~N\cdot m$), which can result in a restricted tubing passage, as shown in the picture below. This restriction can increase turbulence and add a 'throttling' effect to the fluid pathway, resulting in mixing and other potential chromatographic problems.

IDEX Health & Science VHP-200



Conventional Two Piece Ferrule Design



Uniform Tubing Passage

Constricted Tubing Passage

The all Stainless-Steel VHP Fittings include a unique ferrule system with two compression points to provide twice the grip of a standard ferrule. This design also allows the bite on the tubing to be less concentrated and does not restrict the inner diameter, as discussed in the Application Note. The ferrules for 1/16" OD tubing and 10-32 coned ports are two pieces, while the grooved ferrule for 1/32" OD tubing and M4 coned ports is a one-piece design for easier handling, but it will act as two pieces with double compression on the tubing as it is tightened down.



VHP-200x VHP 10-32 Fitting for 1/16" OD tubing



VHP-700x VHP 6-40 Fitting for 1/32" OD tubing

Part No.	Description	Port	Pressure Rating	Required Torque	Head Style	Material	Qty.	
STAINLESS S	STAINLESS STEEL VHP FITTINGS (INCLUDES NUT AND FERRULE)							
VHP-200x	VHP Fitting for 1/16" OD	10-32 Coned	30,000 psi (2,070 bar)	20 in-lbs (2.25 N·m)	1/4" Hex	SST	10-pk	
VHP-700x	VHP Fitting for 1/32" OD	6-40 Coned	30,000 psi (2,070 bar)	20 in-lbs (2.25 N·m)	4 mm Hex	SST	10-pk	
STAINLESS STEEL VHP FERRULES								
VHP-200-01x	VHP Ferrule for 1/16" OD	10-32 Coned	30,000 psi (2,070 bar)	20 in-lbs (2.25 N·m)	_	SST	10-pk	





Reusable VHP Fittings

- > Pressure rated up to 25,000 psi (1,720 bar)
- > Patented innovative design
- Capable of up to ten repeat assembly cycles with no impact on pressure holding ability or carry-over
- ▶ Available in 10-32 threads for 1/16" OD tubing and M4 threads for 1/32" OD tubing
- > Materials of construction: stainless steel and proprietary PEEK polymer blend (PK)
- Quick component replacement, minimal downtime



> Find tightening tools on page 50 designed to deliver the torque necessary for these fittings.

IDEX Health & Science introduces an innovative line of Very High Pressure (VHP) fittings, designed to withstand extreme pressures. This patented line of ground-breaking fitting systems is perfect for use within the increasingly demanding requirements of today's high performance analytical systems.

The Reusable VHP fittings can be reused when following the tightening torque specification listed below. With a polymer front ferrule, there is no damage to the tubing or receiving port, also increasing the life of these components.





for 1/16" OD tubing



10-32 VHP Fitting, Long for 1/16" OD tubing







Part No.	Description	Port	Pressure Rating	Required Torque	Head Style	Material	Qty.
REUSABLE V		1011	r ressure reading	required forque	ricua style	Macciai	a.y.
VHP-320x	VHP Fitting for 1/16" OD	10-32 Coned	25,000 psi (1,720 bar)	10 in-lbs (1.10 N·m)	1/4" Hex	SST/PK	10-pk
VHP-325x	VHP Fitting for 1/16" OD, Long	10-32 Coned	25,000 psi (1,720 bar)	10 in-lbs (1.10 N·m)	1/4" Hex	SST/PK	10-pk
VHP-920x	VHP Fitting for 1/32" OD	M4 Coned	25,000 psi (1,720 bar)	8 in-lbs (0.90 N·m)	4 mm Hex	SST/PK	10-pk
VHP-3200x	VHP Fitting for 1/16" OD	10-32 Coned	11,000 psi (760 bar)	3.5 in-lbs (0.40 N·m)	1/2" Knurl	SST/PK	10-pk
VHP-1001	VHP Fingertight 1/4" Hex Tool	_	_	_	_	PPS	ea.



Assorted Fittings Kits



0.45" (1.14 cm) 0.25" (0.64 cm) 0.95" (2.41 cm) 6000-254 10-32 PEEK Nut with 6000-251 PEEK Ferrule (RD ChromTRAC suffix)



- For PEEK tubing sleeves that can be used with these M4 RheFlex fittings, see page 52.
- For reusable fittings that both work in UHPLC applications and can help ensure the tubing is fully inserted into the receiving port, see page 60.

RheFlex® M4 Fittings

- Incorporates M4 coned threads for 1/32" OD tubing
- Pressure rated to 5,000 psi (345 bar)

Our RheFlex M4 Fitting is designed to connect 1/32" OD tubing in MX Series II™ valves (see Actuated Valves, starting on page 119). This PEEK fitting has a one piece design, which eliminates the need for a separate nut and ferrule. The M4 Fitting design provides dependable zero dead volume connections for micro and nano applications. Due to the unique RheFlex gripping design, the M4 Fitting will hold to 5,000 psi (345 bar) on PEEK or with a PEEK tubing sleeve on fused silica tubing. A PEEK M4 Plug is also available.

Use ChromTRAC™ knobs with the RheFlex M4 Fitting for fingertight convenience and to color-code connections.

Two-Piece RheFlex Fingertight Fittings

The RheFlex Precision Two-Piece PEEK Fittings sets provide inert, biocompatible connections for instrumentation. These fittings have a reliable, time-tested design. Each 1/16" fittings set contains a 10-32 threaded nut and a specially-designed PEEK ferrule. Three lengths of the 1/16" nut are available: Standard, Short, and Extra Long. RheFlex Fingertight Fittings are rated for use up to 7,000 psi (483 bar). Also offered in this product line is the 6000-078 fitting, designed to connect 1/8" OD tubing into our manual preparative-scale injection valves. (See page 115 for more information on these valves.) View the online product bulletin at: www.idex-hs.com.

ChromTRAC™

▶ Brightly colored knobs to easily track inlets and outlets of valves, columns, and detectors

All ChromTRAC-compatible RheFlex fittings offer the ChromTRAC knob option. Specify the ChromTRAC two letter suffix for the color choice when ordering. Please see the ChromTRAC Suffix Codes table below. For example, to order red ChromTRAC knobs with the RheFlex One-Piece Fitting on this page, specify 6000-282RD. No suffix indicates black knobs. View the online product bulletin for RheFlex fittings at: www.idex-hs.com.

ChromTRAC Sufficx Codes

CODE	COLOR	CODE	COLOR
BL	Blue	WH	White
GN	Green	YL	Yellow
GY	Gray	MC	Multi-color (two each of blue, green, gray, red, and yellow)
RD	Red		

Add these letter suffixes to the end of the seven-digit part numbers of the 10-32 and M4 threaded RheFlex Fittings listed below.

Part No.	Description	Port	Pressure Rating	Head Style	Material	Qty.			
RHEFLEX ONE-PIECE	FITTINGS								
6000-360	RheFlex Fitting for 1/32" OD Tubing	M4 Coned	5,000 (345 bar)	1/4" Hex	PEEK, Natural	10-pk			
RHEFLEX TWO-PIECE	FITTINGS (INCLUDES FERRULES)								
6000-078	RheFlex Fitting for 1/8" OD Tubing	5/16-24 Coned	5,000 psi (345 bar)	5/16" Hex	PEEK, Natural	ea.			
6000-254	RheFlex Fitting for 1/16" OD Tubing	10-32 Coned	7,000 psi (483 bar)	ChromTRAC knob	PEEK, Natural	10-pk			
6000-255	RheFlex Fitting for 1/16" OD Tubing, Short	10-32 Coned	7,000 psi (483 bar)	1/4" Hex	PEEK, Natural	10-pk			
REPLACEMENT FERR	REPLACEMENT FERRULES								
6000-079	RheFlex Ferrule for 1/8" OD Tubing	5/16-24 Coned	7,000 psi (483 bar)	ChromTRAC knob	PEEK, Natural	5-pk			
6000-251	RheFlex Ferrule for 1/16" OD Tubing	10-32 Coned	7,000 psi (483 bar)	ChromTRAC knob	PEEK, Natural	10-pk			





CONNECTORS


Connectors are designed to securely join tubing together or to facilitate the joining of tubing to other fluid pathway components. We offer multiport connectors with different thread and port configurations to meet your system requirements and connection needs. Some of our connectors feature a True ZDV (Zero Dead Volume) internal configuration that helps minimize the formation of dead volume in your fluidic pathway. Our versatile adapters help bring two connectors with different configurations together. Connectors are manufactured from 316 stainless steel or from inert polymers to ensure chemical compatibility with the fluid passing through. Peristaltic tube connectors are ideal for making connections with soft-walled, peristaltic tubing. Our extensive line of connectors includes tees, crosses, Luer Adapters, barbed and threaded adapters, and a variety of other options.

63	THREADED ADAPTERS	79	ACCESSORIES
67	HIGH PRESSURE MULTIPORT	80	NANOPORT ASSEMBLIES
72	CONNECTORS ULTRA HIGH PRESSURE MULTIPORT	81	LOW PRESSURE MULTIPORT CONNECTORS
	CONNECTORS	90	LUER ADAPTERS
77	MICROTIGHT® ADAPTERS	91	PERISTALTIC TUBE CONNECTORS



Threaded Adapters

- Threaded adapters in a variety of configurations
- English, Metric, and NPT threaded adapters offered
- Bring together connectors with different threads
- Manufactured from inert polymers PEEK, PCTFE, ETFE, and PTFE

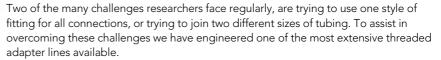


Threaded Adapters come in a wide variety of configurations to meet your system requirements. They are designed to effectively bring together connectors with different threads. We offer them in English, Metric, and NPT versions. Manufactured from inert polymers and stainless steel they deliver excellent chemical resistance.



English Threaded Adapters

Our versatile English Threaded Adapters are used specifically to securely attach connectors with different threads. We designed these adapters to work with English to English threaded geometries. Manufactured from Stainless Steel, PEEK, or Tefzel™ (ETFE), they deliver excellent solvent resistance.

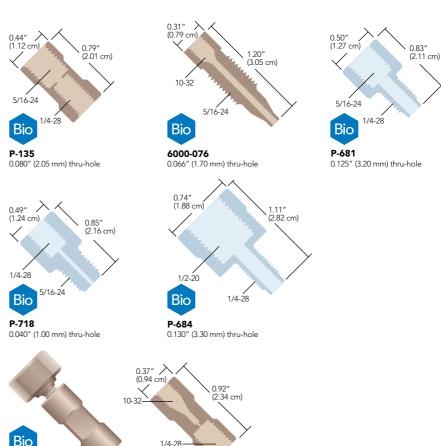

P-627

0.020" (0.50 mm) thru-hole Includes (1) F-300 for 1/16" OD tubing


Threaded Adapters (Cont.)



NOTE

When using an adapter with male (external) threads, we recommend you first attach the adapter body into the receiving port, and then connect your tubing and fitting into the head of the adapter body.

- Use the 6000-076 Adapter to connect 1/16" OD tubing to the Preparative-Scale Injector Valve (page 66).
- You may not need an adapter to connect 1/16" OD tubing into your flat-bottom port. A less expensive alternative is to use a Flangeless Nut and Ferrule starting on page 45 or a Super Flangeless™ Nut and Ferrule starting on page 39.

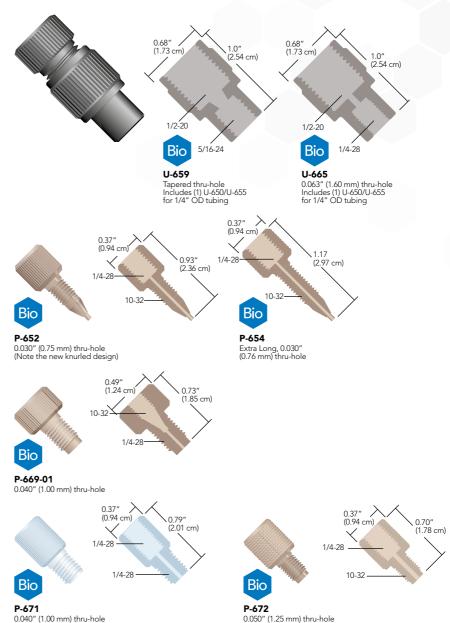




Here are application ideas using two of our popular adapters:

- Many injection valves used in HPLC systems have 10-32 coned ports designed to accept 1/16" OD tubing. However, this may be a problem if large injection volumes are required (in excess of 10 mL). The most popular loops for large volume samples are made from 1/8" OD tubing, making it impossible to connect these larger volume loops to your injection valve. The solution: use our P-654 Adapter and the appropriate fittings for your sample loop. This set-up allows connection of 1/8" OD sample loop leads to your injection valve.
- Another potential application is connecting tubing to low-pressure solenoid valves with 1/4-28 flat-ottom ports. Most low-pressure valves of this type have very shallow threaded ports, which typically preclude the use of our Flangeless Fittings. However, by first threading our P-671 Adapter into the valve port(s), you can effectively use standard 1/4-28 fittings to connect your tubing into the backside of the adapter body. This also saves "wear and tear" on the threads in the valve ports.



Metric Threaded Adapters

Our versatile Metric Threaded Adapters are used specifically to effectively attach connectors with different threads. We designed these adapters to work with English to Metric threaded geometries. Manufactured from Stainless Steel, PEEK or Tefzel™ (ETFE) they deliver excellent solvent resistance.

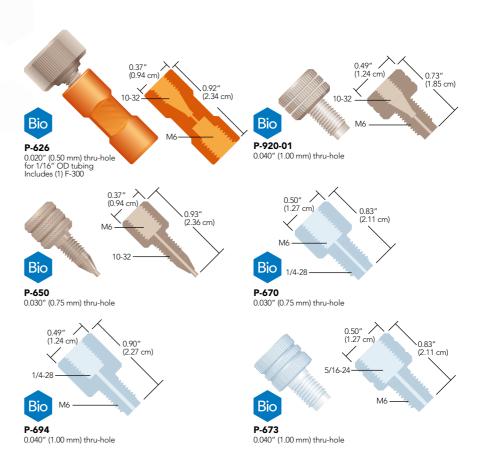



- For an alternative to the Female M6 Adapters presented above, try a P-602 or P-622 Low Pressure Metric Union from page 84, along with the appropriate Metric Flangeless Fittings on page 43.
- To direct connect your tubing into a flat-bottom port, find the appropriate Flangeless or Super Flangeless™ Fittings on page 45 and page 39 respectively.
- Need metric fittings for your connections? See page 43.

Threaded Adapters (Cont.)



National Pipe Thread Adapters

These adapters make connections to female 1/8" and 1/4" National Pipe Thread (NPT) ports.

Manufactured from PEEK polymer, our NPT Adapters are durable and chemically resistant. We provide versions with either 1/4-28 or 5/16-24 flat-bottom threads, suitable for most low pressure applications.

Please Note: Wrap the threads on the NPT side of these adapters with thread seal tape (plumber's tape) to ensure a leak-free seal.



Replacement fittings for these adapters are located on the pages indicated below:

	Page(s)
1/4-28 for 1/8" OD tubing	41
5/16-24 for 1/8" OD tubing	41, 56
5/16-24 for 3/16" OD tubing	41

Other tubing/fitting combinations are available. For more information, please contact your local Distributor or IDEX Health & Science directly.



1/8" NPT to 1/4-28 Flat-Bottom Female Adapter for 1/8" OD tubing Includes (1) XP-308 Fitting



1/8" NPT to 5/16-24 Flat-Bottom Female Adapter for 3/16" OD tubing Includes (1) XP-132 Fitting



Our U-500 and U-510 NPT Adapters are great for attaching 1/8" OD fluoropolymer sparging lines to sparging gas tank regulating valves. Simply thread the appropriatelysized NPT Adapter into the valve's receiving port and then attach your sparging tubing to the adapter body using the fittings provided.

Threaded Adapters

Part No.	Description			Includes	Thru-hole	Swept Volume	Pressure Rating	Qty.
ENGLISH	THREADED ADAPTERS							
6000-076	PEEK Adapter, 5/16-24 C, M to 10-32 C,	F		N/A	0.066" (1.70 mm)	49.8 µL	3,000 psi (207 bar)	ea.
P-135	PEEK Adapter, 5/16-24 FB, F to 1/4-28 F		N/A	0.080" (2.05 mm)	4.1 μL	1,000 psi (69 bar)	ea.	
P-627	PEEK Adapter, 10-32 C, F to 1/4-28 FB, F			(1) F-300	0.020" (0.50 mm)	0.30 µL	1,000 psi (69 bar)	ea.
P-681	PCTFE Adapter, 5/16-24 FB, F to 1/4-28 FB, M			N/A	0.125" (3.20 mm)	96.6 µL	1,000 psi (69 bar)	ea.
P-684	PCTFE Adapter, 1/2-20 FB, F to 1/4-28 FB, M		N/A	0.130" (3.30 mm)	121.7 μL	250 psi (17 bar)	ea.	
P-718	PCTFE Adapter, 5/16-24 FB, M to 1/4-28 FB, F		N/A	0.040" (1.00 mm)	10.3 μL	1,000 psi (69 bar)	ea.	
U-659	PEEK Adapter, 5/16-24 FB, F to 1/2-20 FB, F		(1) XU-655	Tapered**	42.0 µL	250 psi (17 bar)	ea.	
U-665	PEEK Adapter, 1/2-20 FB, F to 1/4-28 FB, F			(1) XU-655	0.063" (1.60 mm)	6.6 µL	250 psi (17 bar)	ea.
P-652	PEEK Adapter, 1/4-28 FB, F to 10-32 C, M		N/A	0.030" (0.75 mm)	6.7 µL	1,000 psi (69 bar)	ea.	
P-654	PEEK Adapter, 1/4-28 FB, F to 10-32 C, M, Extra Long			N/A	0.030" (0.75 mm)	9.5 µL	1,000 psi (69 bar)	ea.
P-669-01	PEEK Adapter, 10-32 C, F to 1/4-28 FB, M		N/A	0.040" (1.00 mm)	6.6 µL	1,000 psi (69 bar)	ea.	
P-671	PTFE Adapter, 1/4-28 FB, F to 1/4-28 FB, M		N/A	0.040" (1.00 mm)	8.0 µL	1,000 psi (69 bar)	ea.	
P-672	PEEK Adapter, 1/4-28 FB, F to 10-32 FB, M			N/A	0.050" (1.25 mm)	11.4 µL	1,000 psi (69 bar)	ea.
METRIC	M6 THREADED ADAPTERS							
P-626	PEEK Adapter, 10-32 C, F to M6 FB, F			(1) F-300	0.020" (0.50 mm)	0.3 μL	1,000 psi (69 bar)	ea.
P-650	PEEK Adapter, M6 FB, F to 10-32 C, M S	tandard		N/A	0.030" (0.75 mm)	6.7 µL	1,000 psi (69 bar)	ea.
P-670	PCTFE Adapter, M6 FB, F to 1/4-28 FB, I	M		N/A	0.030" (0.75 mm)	2.6 µL	1,000 psi (69 bar)	ea.
P-673	PCTFE Adapter, 5/16-24 FB, F to M6 FB,	М		N/A	0.040" (1.00 mm)	9.9 µL	1,000 psi (69 bar)	ea.
P-694	PCTFE Adapter, 1/4-28 FB, F to M6 FB, I	M		N/A	0.040" (1.00 mm)	11.3 µL	1,000 psi (69 bar)	ea.
P-920-01	PEEK Adapter, 10-32 C, F to M6 FB, M			N/A	0.040" (1.00 mm)	8.0 µL	1,000 psi (69 bar)	ea.
1/8" MA	LE NATIONAL PIPE THREAD ADAPTERS							
Part No.	Description	Color	Tubing OD	Includes	Thru-hole	Swept Volume	Pressure Rating	Qty.
U-510	PEEK 1/8" NPT, M to 1/4-28 FB, F Adapter	Red	1/8"	(1) XP-308	0.062" (1.60 mm)	17.3 μL	500 psi (34 bar)	ea.
U-514	PEEK 1/8" NPT, M to 5/16-24 FB, F Adapter	Natural	3/16"	(1) XP-132	0.125" (3.2 mm)	70.4 μL	500 psi (34 bar)	ea.
1/4" MA	LE NATIONAL PIPE THREAD ADAPTERS					·		
U-500	PEEK 1/4" NPT, M to 1/4-28 FB, F Adapter	Red	1/8"	(1) XP-308	0.062" (1.60 mm)	17.3 μL	500 psi (34 bar)	ea.
U-504	PEEK 1/4" NPT, M to 5/16-24 FB, F Adapter	Natural	3/16"	(1) XP-132	0.125" (3.2 mm)	70.4 uL	500 psi (34 bar)	ea.

F = Female (internal) threads; M = Male (external) threads; XL = extra long; C = Coned; FB = Flat-Bottom * The pressure rating of this adapter exceeds the pressure holding ability of the fittings and tubing used with it. ** Thru-hole tapers from 0.188" (4.80 mm) to 0.125" (3.20 mm).



High Pressure Mixing Tees

Mixing Tees utilize a specifically engineered internal geometry to efficiently mix two fluid streams into one combined stream. Mixing Tees are ideal for microbore or analytical gradient HPLC. These mixing tees are specifically designed for high pressure applications.

Static Mixing Tees

- PEEK body with two-piece fingertight fittings
- Low swept volume

Static Mixing Tees are ideal for microbore or analytical gradient HPLC. They have a low swept volume of 2.2 μL (includes frit volume) and are designed for flow rates of 0.5 to 3 mL/min and a maximum pressure of 5,000 psi (345 bar). The back pressure caused by the tee is typically only 10 to 20 psi (0.7 to 1.4 bar) at these flow rates. The thru-holes are 0.020" (0.50 mm) and the center port features a 10 μm UHMWPE or stainless steel frit that aids mixing.



NOTE


- Turbulent mixing of solvents often increases outgassing. To maintain a bubble-free fluid pathway, we recommend solvent degassing when using this product.
- The frit incorporated into our U-466 and U-466S Static Mixing Tees is not replaceable. If it becomes clogged, the Mixing Tee must be replaced.

U-466 and U-466S Static Mixing Tees with F-300 Fingertight Fittings for 1/16" OD tubing



- Constructed of inert PEEK and PCTFE
- **)** Low swept volume of 0.95 μL
- Designed for flow rates of 20–250 μL/min

Our Micro Static Mixing Tee utilizes a specifically engineered internal geometry to efficiently mix two fluid streams into one combined stream. The center port also features a 0.5 µm porosity PEEK polymer frit to aid in mixing. This frit adds a maximum of 20 psi (1.4 bar) back pressure to most systems (within the stated flow rate range). The Mixing Tee handles a maximum pressure of 5,000 psi (345 bar) when directly connecting 1/16" OD tubing, or up to 4,000 psi (276 bar) with capillary tubing when using our NanoTight[™] Fittings and Tubing Sleeves (page 37).

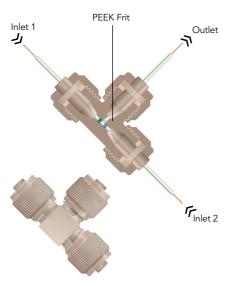

M-540

Micro Static Mixing Tee 0.010" thru-hole with fittings included (tubing and tubing sleeves not included)



- See our Vacuum Degassing Systems on page 154.
- Our standard Static Mixing Tees are designed for flow rates from 0.5 mL/min to 3 mL/min.

High Pressure Mixing Tees (Cont.)



APPLICATION NOTE

Several researchers use our PEEK MicroTee to introduce ionizing voltage to their fluid stream just prior to a Mass Spectrometer¹. MicroTees are well suited for this application due to advantageous internal geometry and PEEK polymer's electrical resistance. The materials required for this setup are as follows: one gold or platinum conducting wire, one P-775 or P-875 MicroTee (this page), one MicroTight Tubing Sleeve (page 52) for the conducting wire (as needed to accommodate wire diameter), and at least two more MicroTight Tubing Sleeves (page 52) to connect your capillary tubing.

To set up a similar connection, first thread your wire through the appropriate tubing sleeve, if necessary, with the wire extending beyond both ends of the sleeve. Slip the female nut included with the MicroTee over the wire or sleeved wire, followed by the ferrule — ensuring the wire (and its sleeve) extends well past the end of the ferrule tip. Align the tip of the wire with the thru-hole of the MicroTee and gently insert the wire until it bottoms out. Now finger tighten the female nut into place. Attach your flow path tubing to the MicroTee's two other available ports, following the instructions provided with the MicroTee.

Begin fluid flow through the tee and apply voltage to the conducting wire lead. This setup typically provides effective electrospray ionization in applications having a flow rate of 100 µL/min or greater.

¹One such paper describing pioneering electrospray work: Protein Identification at the Low Femtomole Level from Silver-Stained Gels Using a New Fritless Electrospray Interface for Liquid Chromatography-Microspray and Nanospray Mass Spectrometry. Christine L. Gatlin, Gerd R. Kleemann, Lara G. Hays, Andrew J. Link, John R. Yates III (1998) Analytical Biochemistry 263, 93-101.

MicroTee & Cross for Capillary Tubing


- Direct connect 1/16", 1/32", 360 μm OD tubing, plus other capillary tubing
- Low swept volume

Use our MicroTees and MicroCrosses to join capillary tubing. All of these products are made entirely of PEEK and have 0.006" (0.150 mm) thru-holes, with resulting swept volumes ranging from 29 to 81 nL.



NOTE

Use only the ferrules supplied with each connector — they are not interchangeable. Replacement ferrules and female nuts are available on page 35. For MicroUnions, MicroTees, and MicroCrosses for UHPLC applications, see page 72.

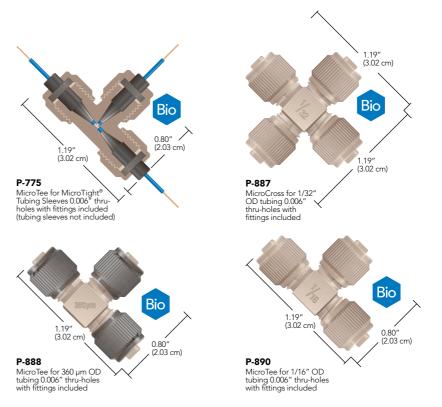

High Pressure Mixing Tees

Part No.	Description	Threads	Includes	Thru-hole	Swept Volume	Pressure Rating	Qty.
STATIC	MIXING TEE						
U-466	PEEK Static Mixing Tee for 1/16" OD Tubing, 10 μm UHMWPE Frit	10-32 Coned	(3) F-300	0.020" (0.50 mm)	2.2 μL	5,000 psi (345 bar)	ea.
U-466S	PEEK Static Mixing Tee for 1/16" OD Tubing, 10 µm SST Frit	10-32 Coned	(3) F-300	0.020" (0.50 mm)	2.2 μL	5,000 psi (345 bar)	ea.
MICRO:	STATIC MIXING TEE						
M-540	PEEK Micro Static Mixing Tee, for 1/16" OD Tubing	5/16-24 Coned	(3) F-132/P-416	0.010" (0.250 mm)	0.95 µL	5,000 psi (345 bar)	ea,
MICROT	EE, MICROCROSS AND MICROELBOW						
P-775	PEEK MicroTee for MicroTight Sleeves	5/16-24 Coned	(3) F-172, (3) P-416	0.006" (0.150 mm)	29 nL	4,000 psi (276 bar)	ea.
P-777	PEEK MicroCross for MicroTight Sleeves	5/16-24 Coned	(4) F-172, (4) P-416	0.006" (0.150 mm)	38 nL	4,000 psi (276 bar)	ea.
P-875	PEEK MicroTee with Mounting Hole, for MicroTight Sleeves	5/16-24 Coned	(3) F-172, (3) P-416	0.006" (0.150 mm)	29 nL	4,000 psi (276 bar)	ea.
P-885	PEEK MicroTee for 1/32" OD Tubing	5/16-24 Coned	(3) F-112, (3) P-416	0.006" (0.150 mm)	29 nL	5,000 psi (345 bar)	ea.
P-887	PEEK MicroCross for 1/32" OD Tubing	5/16-24 Coned	(4) F-112, (4) P-416	0.006" (0.150 mm)	38 nL	5,000 psi (345 bar)	ea.
P-888	PEEK MicroTee for 360 µm OD Tubing	5/16-24 Coned	(3) F-152, (3) P-416BLK	0.006" (0.150 mm)	29 nL	5,000 psi (345 bar)	ea.
P-889	PEEK MicroCross for 360 µm OD Tubing	5/16-24 Coned	(4) F-152, (4) P-416BLK	0.006" (0.150 mm)	38 nL	5,000 psi (345 bar)	ea.
P-890	PEEK MicroTee for 1/16" OD Tubing	5/16-24 Coned	(3) F-132, (3) P-416	0.006" (0.150 mm)	58 nL	5,000 psi (345 bar)	ea.
P-891	PEEK MicroCross for 1/16" OD Tubing	5/16-24 Coned	(4) F-132, (4) P-416	0.006" (0.150 mm)	81 nL	5,000 psi (345 bar)	ea.



High Pressure Unions

PEEK ZDV Unions

Our PEEK zero-dead-volume (ZDV) Unions come complete with two F-300 Fingertight Fittings for 1/16" OD tubing and are pressure rated to 5,000 psi (344 bar).







Bio-Inert UHPLC Unions

- Unique, Patent-Pending Process allows a fully-PEEK fluid contact area combined with the strength of stainless steel
- > Pressure rated to 17,400 psi (1,200 bar)
- Two inner diameters available: 0.008" and 0.016"


These unions are specifically engineered for Bio-Inert UHPLC applications. Combining the physical strength of 316 stainless steel with the inertness and biocompatibility of an all-PEEK fluid pathway, these unions will work well in applications where pressures reach up to 17,400 psi (1,200 bar) — without allowing metal contact by the fluid.

Neither union comes with fittings, but can be paired successfully with any 10-32 coned fitting that uses a polymer nose or ferrule.

Note: All-stainless steel fittings should NOT be used with these unions, as they will damage the internal conical seat.

NanoTight[™] Union

NanoTight Unions improve capillary tubing connections in several ways. The internal design of the union greatly reduces the incidence of tubing misalignment. When using 1/16" OD tubing sleeves (found on page 52) to connect capillary tubing, the webbed thru-hole minimizes breaking of fused silica while adding only miniscule swept volume. The results are fewer blockages, fewer flow rate reductions and fewer back pressure problems.

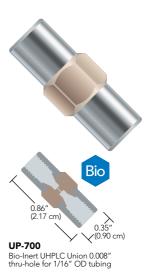

10-32 0.83"
(2.11 cm)

P-779
NanoTight Union
0.005" thru-hole with
F-331N Fittings


Part No.	Description	Threads	Includes	Thru-hole	Swept Volume	Pressure Rating	Qty.		
BIO-INER	BIO-INERT UHPLC UNIONS								
UP-700	Bio-Inert UHPLC Union for 1/16" OD Tubing, Natural (Tan)	10-32 Coned	N/A	0.008" (0.20 mm)	0.05 μL	17,400 psi (1,200 bar)	ea.		
PEEK ZDV	UNIONS								
P-704	PEEK Union for 1/16" OD Tubing	10-32 Coned	(2) F-300	0.020" (0.50 mm)	0.28 μL	5,000 psi (344 bar)	ea.		
P-742	PEEK Union for 1/16" OD Tubing	10-32 Coned	(2) F-300	0.010" (0.25 mm)	0.07 μL	5,000 psi (344 bar)	ea.		
P-760	PEEK Union for 1/16" OD Tubing	10-32 Coned	(2) F-300	0.050" (1.25 mm)	1.2 µL	5,000 psi (344 bar)	ea.		
NANOTIG	NANOTIGHT UNION								
P-779	PEEK NanoTight Union for 1/16" OD Tubing and Tubing Sleeves	10-32 Coned	(2) F-331N	0.005" (125 µm)	8 nL	5,000 psi (344 bar)	ea.		



High Pressure PEEK Tees & Crosses

- > Highest pressure holding flat-bottom fitting system we offer
- > Eliminates loosening of fittings due to tubing twist
- > Excellent for Tubing Assemblies
- > Holds tight even through vibration


Our PEEK Tees and Crosses include high pressure F-300 PEEK Fingertight Fittings allowing maximum operating pressures to 3,500 psi (241 bar) when used with 1/16" OD PEEK or stainless steel tubing.

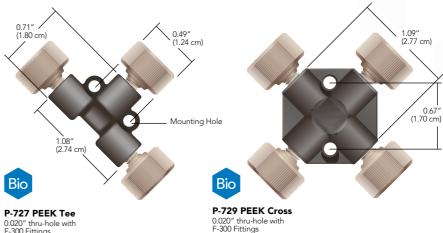

0.020" thru-hole with F-300 Fittings

PEEK 7-Port Manifold

Combine several streams into one or split one fluid stream into several. This PEEK 7-Port Manifold comes complete with F-331 Fingertight Fittings for 1/16" OD tubing and offers a pressure rating of 5,000 psi (345 bar). Seal unused ports with any of our polymer 10-32 coned plugs on page 55.

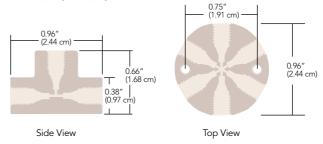



Part No.	Description	Threads	Includes	Thru-hole	Swept Volume	Pressure Rating	Qty.
PEEK TEE	ES AND CROSSES						
P-727	PEEK Tee for 1/16" OD Tubing	10-32 Coned	(3) F-300	0.020" (0.50 mm)	0.57 μL	3,500 psi (241 bar)	ea.
P-728	PEEK Tee for 1/16" OD Tubing	10-32 Coned	(3) F-300	0.050" (1.25 mm)	3.0 µL	3,500 psi (241 bar)	ea.
P-729	PEEK Cross for 1/16" OD Tubing	10-32 Coned	(4) F-300	0.020" (0.50 mm)	0.72 μL	3,500 psi (241 bar)	ea.
PEEK MA	NIFOLD						
P-170	PEEK 7-Port Manifold for 1/16" OD Tubing	10-32 Coned	(7) F-331	0.020" (0.50 mm)	2.2 μL	5,000 psi (345 bar)	ea.



High Pressure MicroTight® Unions

Conductive MicroTight Union

The Conductive MicroTight Union manufactured by IDEX Health & Science provides an excellent opportunity to introduce voltage into an electrospray or capillary electrophoresis system. With an extremely low internal volume of 16 nL, this union can be placed inline with 360 µm OD capillary tubing. Mount and apply voltage to these unions using our Insulating Mounting Bracket below.



Conductive MicroTight Union for 360 µm OD tubing with fittings and Capsule Union included



APPLICATION NOTE

For an example of using a Conductive MicroTight Union in a pressure driven ion preconcentration application see: "Self-Sealed Vertical Polymeric Nanoporous Junctions for High Throughput Nanofluidic Applications."


Sun Jae Kim and Jong Yoon Han. Analytical Chem. 2008, 80: 3507-3511.



Easily integrate the Conductive MicroTight Union into your system with our Insulating Mounting Bracket, shown on page 79.

MicroTight® Connectors for Capillary Tubing

Connect two pieces of capillary tubing with our PEEK MicroTight Connectors. The True ZDV Unions allow two pieces of tubing to connect directly to each other — using the included gauge plug to ensure proper alignment. The standard union and elbow both feature a 0.006" (0.150 mm) thru-hole, adding only a small amount of additional flow-path volume to help ensure proper chromatographic results.

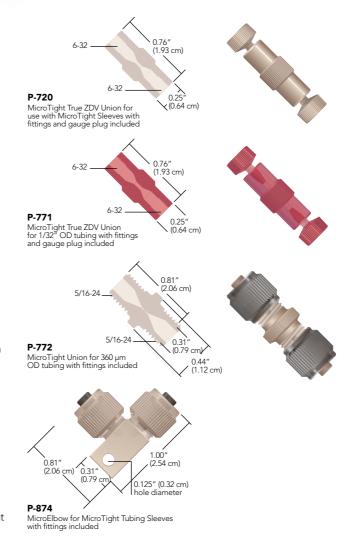


Part No.	Description	Threads	Includes	Thru-hole	Swept Volume	Pressure Rating	Qty.			
MICRO	MICROTIGHT UNIONS									
P-720	PEEK True ZDV Union for MicroTight Sleeves	6-32 Coned	(2) F-125, (1) P-553	N/A	N/A	4,000 psi (276 bar)	ea.			
P-771	PEEK True ZDV Union for 1/32" OD Tubing	6-32 Coned	(2) F-126S, (1) P-553	N/A	N/A	5,000 psi (345 bar)	ea.			
P-772	PEEK Union for 360 µm OD Tubing	5/16-24 Coned	(2) F-152, (2) P-416BLK	0.006" (0.150 mm)	5 nL	5,000 psi (345 bar)	ea.			
P-874	PEEK MicroElbow for MicroTight Sleeves	5/16-24 Coned	(2) F-172, (2) P-416	0.006" (0.150 mm)	20 nL	4,000 psi (276 bar)	ea.			
REPLAC	EMENT GAUGE PLUGS (TO ACHIEVE TRUE	ZDV CONNEC	TIONS WITH OUR P-720 AND	P-771 UNIONS)						
P-553	Gauge Plug, Delrin®	6-32 Coned	N/A	N/A	N/A	N/A	ea.			
CONDU	CONDUCTIVE MICROTIGHT UNIONS									
M-572	Conductive Union for 360 µm OD Tubing, PEEK/SST	5/16-24 Coned	(2) F-152, (2) P-416BLK, (1) M-128NF	0.011" (0.279 mm)	16 nL	5,000 psi (345 bar)	ea.			




High Pressure Stainless Steel Tees & Crosses

These 316 stainless steel connectors come complete with 10-32 stainless steel fittings for use with 1/16" OD tubing and are rated to 20,000 psi (1,380 bar). They are compatible with any 10-32 coned threaded fittings.

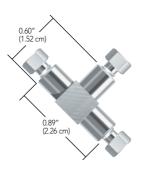

U-428 Stainless Steel Tee 0.020" thru-hole with U-400 and U-401 Fittings

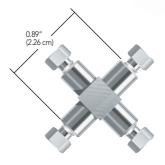

U-430 Stainless Steel Cross 0.020" thru-hole with U-400 and U-401 Fittings

Part No.	Description	Threads	Includes	Thru-hole	Swept Volume	Pressure Rating	Qty.
VHP TEE	FOR 1/16" OD TUBING						
U-428	Stainless Steel Tee for 1/16" OD Tubing	10-32 Coned	(3) U-400, (3) U-401	0.020" (0.50 mm)	0.57 μL	20,000 psi (1,380 bar)	ea.
U-429	Stainless Steel Tee for 1/16" OD Tubing	10-32 Coned	(3) U-400, (3) U-401	0.040" (1.00 mm)	2.1 µL	20,000 psi (1,380 bar)	ea.
U-430	Stainless Steel Cross for 1/16" OD Tubing	10-32 Coned	(4) U-400, (4) U-401	0.020" (0.50 mm)	0.72 μL	20,000 psi (1,380 bar)	ea.
U-431	Stainless Steel Cross for 1/16" OD tubing	10-32 Coned	(4) U-400, (4) U-401	0.040" (1.00 mm)	2.5 µL	20,000 psi (1,380 bar)	ea.



VHP Stainless Steel ZDV Unions


- Supplied with fittings for 1/16" OD or 1/8" OD tubing
- Manufactured from 316 stainless steel
- All union assemblies rated to 20,000 psi (1,380 bar) or higher



NOTE

It is possible to order the products on this page without the fittings. Simply use a -01 at the end of the product number to order the union body without fittings.

Our high pressure, zero-dead-volume (ZDV) unions are precision machined from 316 stainless steel, carefully passivated, then thoroughly rinsed. Each comes complete with stainless steel nuts and ferrules.

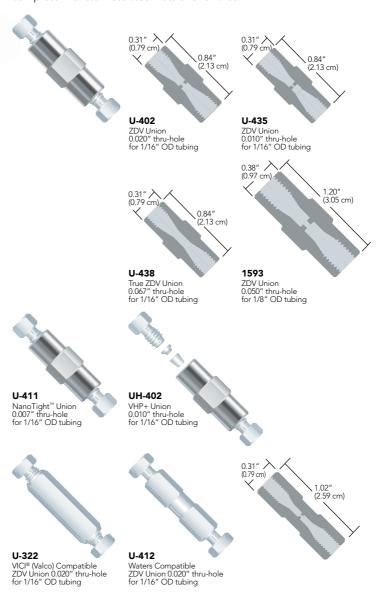

Part No.	Description	Threads	Includes	Thru-hole	Swept Volume	Pressure Rating	Qty.
VHP STA	AINLESS STEEL ZDV UNIONS						
1593	Stainless Steel Union for 1/8" OD Tubing	1/4-28 Coned	(2) C-235/C-236	0.050" (1.25 mm)	1.48 µL	20,000 psi (1,380 bar)	ea.
U-402	Stainless Steel Union for 1/16" OD Tubing	10-32 Coned	(2) U-400/U-401	0.020" (0.50 mm)	0.13 μL	20,000 psi (1,380 bar)	ea.
U-411	Stainless Steel Union for 1/16" OD Tubing	10-32 Coned	(2) U-400/U-401	0.007" (178 µm)	13 nL	20,000 psi (1,380 bar)	ea.
U-435	Stainless Steel Union for 1/16" OD Tubing	10-32 Coned	(2) U-400/U-401	0.010" (0.25 mm)	20 nL	20,000 psi (1,380 bar)	ea.
U-438	Stainless Steel Union for 1/16" OD Tubing	10-32 Coned	(2) U-400/U-401, (1) P-554 Gauge Plug	0.067" (1.70 mm)	Near 0 µL	20,000 psi (1,380 bar)	ea.
UH-402	VHP+ Stainless Steel Union for 1/16" OD Tubing	10-32 Coned	(2) VHP-200	0.010" (0.25 mm)	20 nL	30,000 psi (2,070 bar)	ea.
VICI (VA	LCO) COMPATIBLE ZDV UNION						
U-322	Stainless Steel Union for 1/16" OD Tubing	10-32 Coned	(2) U-320/U-321	0.020" (0.50 mm)	0.15 μL	20,000 psi (1,380 bar)	ea.
WATERS	© COMPATIBLE ZDV UNION						
U-412	Stainless Steel Union for 1/16" OD Tubing	10-32 Coned	(2) U-410/U-401	0.020" (0.50 mm)	0.10 μL	20,000 psi (1,380 bar)	ea.



VHP Tees & Crosses



1.19" (3.02 cm) 0.80" (2.03 cm) UH-750 VHP MicroTee for 360 µm OD tubing with fittings included (3.02 cm) (3.02 cm)UH-752



VHP MicroCross

for 360 µm OD tubing with fittings included

Our VHP Stainless Steel Tees and Crosses are precision machined from durable stainless steel. It is mechanically designed for bringing together three or four pieces of tubing. Our VHP Tees & Crosses have an extremely high pressure rating of 30,000 psi (2070 bar).

VHP Tees & Crosses for Capillary Tubing

- Direct-connect either 360 μm or 1/32" OD tubing no sleeves required!
- Available in both tee and cross configurations
- Pressure rated to 15,000 psi (1,034 bar)

To help facilitate multi-port connections in UHPLC applications, our experts have developed a line of MicroTees and MicroCrosses, manufactured from stainless steel and featuring small thru-holes and very low internal volume. Additionally, the stainless steel construction allows these products to be used in applications where electrical conductivity is desired.

Included with the MicroTees and MicroCrosses are the VHP MicroFerrules found on page 59. The P-278 Extender Tool on page 33 can be used to tighten the female nuts that are included with these connectors.



Why 1/32" OD Tubing and 360 µm OD Tubing?

IDEX Health & Science has focused strongly on the development of a variety of connectors and accessories for 1/32" OD tubing and 360 µm OD tubing. We have focused on these specific sizes due to their overwhelming popularity in analytical instruments, especially where micro and nano-scale analyses are being performed. By creating products designed for these popular sizes, the overall connection is easier to make and generally holds to increased pressures over connections where tubing sleeves are involved.

VHP Tee for 1/16" OD Tubing

IDEX Health & Science offers this Very High Pressure (VHP) Tee Connector, designed to bring three pieces of tubing together. The all-316 stainless steel connector is designed for 1/16" OD tubing and is pressure rated to 30,000 psi (2,070 bar).

Part No.	Description	Threads	Includes	Thru-hole	Swept Volume	Pressure Rating	Qty.
VHP TEE F	OR 1/16" OD TUBING						
UH-427	VHP Tee for 1/16" OD Tubing, SST	10-32 Coned	(3) VHP-200	0.020" (0.50 mm)	0.57 μL	30,000 psi (2,070 bar)	ea.
VHP TEES	& CROSSES FOR CAPILLARY TUBING						
UH-700	VHP MicroTee for 1/32" OD Tubing, PEEK/SST	5/16-24 Coned	(3) PK-112, (3) P-416	0.010" (0.25 mm)	84 nL	15,000 psi (1,034 bar)	ea.
UH-750	VHP MicroTee for 360 µm OD Tubing, PEEK/SST	5/16-24 Coned	(3) PK-152, (3) P-416BLK	0.010" (0.25 mm)	84 nL	15,000 psi (1,034 bar)	ea.
UH-752	VHP MicroCross for 360 µm OD Tubing, PEEK/SST	5/16-24 Coned	(4) PK-152, (4) P-416BLK	0.010" (0.25 mm)	101 nL	15,000 psi (1,034 bar)	ea.



VHP MicroTight® Unions

- Featuring stainless steel bodies and PK/PEEK fittings
- > Pressure rated up to 15,000 psi (1,034 bar)
- Options to direct-connect both 1/32" OD tubing and 360 µm OD tubing

VHP Unions for Capillary Tubing

IDEX Health & Science has expanded its line of specialized fittings and connectors for UHPLC applications to include several innovative unions and adapters.

Two of these products — the UH-432 and UH-436 — follow the design of our popular Mini MicroFilters (see page 109) and allow a convenient union between either 1/32" OD tubing or 360 μ m OD tubing. Each features a stainless steel union body and a unique stainless steel union capsule, enabling both excellent chemical compatibility as well as conductivity, making these a great choice for electrical interfacing in certain LC-MS applications. Each is also coupled with direct-connect ferrules made from our proprietary PEEK polymer blend (PK), allowing tubing connections up to 15,000 psi (1,034 bar). (*Please Note: While these connectors can be used at elevated pressures, they are not recommended for applications above 100 °C.*)

The UH-632 is a more traditionally designed connector, incorporating internally threaded ports. The union (UH-632) features a true ZDV (zero dead volume) connection between both tubes. This unique product is coupled with our one-piece Ultra-High Performance Fingertight fittings manufactured from our proprietary PEEK polymer blend, allowing them to be used in high temperature applications (up to 200 °C) at pressures up to 6,000 psi (414 bar) — or use these connectors at room temperature up to 15,000 psi (1,034 bar)!


The 1959-01 is a new VHP union designed to accept the popular M4x0.7 threaded fittings for 1/32" OD tubing. These unions will work nicely with both the VHP-900 fittings (found on page 59) as well as the reusable VHP-920 (found on page 60).

UHP MicroTight Union for 1/32" OD tubing with fittings included



UHP MicroTight® Union for 1/32" OD tubing

VHP MicroTight® Unions (Cont.)

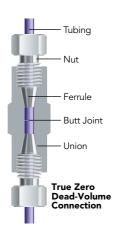



What is a True ZDV Union?

True zero dead volume (ZDV) unions are designed so that the two joined pieces of tubing butt perfectly together as shown in the image to the right. These products have no swept volume contained within the union body. The fluid moves directly from one tube into another in this type of connector.

When using true ZDV unions, it is important to take care to ensure connecting tubing has burr-free 90 degree ends. Find tubing cutters on page 28 to assist with cleanly cutting polymer and fused silica tubing. Gauge plugs are supplied with True ZDV Unions to assist with assembly. With the gauge plug inserted into one side of the union, a hard stop is created for the tubing to bottom out against as it is connected to the opposite port. The gauge plug is removed and then the second piece of tubing is connected, using the first piece of tubing to bottom out against resulting in the two tubes joined together in the center of the union.



- > Find replacement VHP fittings on page 59.
- > Find Fused Silica tubing on page 16.
- Find 1/32" OD Stainless Steel tubing on page 19.
- To achieve 15,000 psi (1,034 bar) with the female threaded fittings used with some of these products, use the P-278 extender tool found on page 50.

VHP MicroTight Unions

Part No.	Description	Threads	Includes	Thru-hole	Volume	Pressure Rating	Qty.
VHP UN	IIONS FOR CAPILLARY TUBING						
UH-432	VHP Union for 1/32" OD Tubing, PEEK/SST	5/16-24 Coned	(2) PK-112, (2) P-416	0.006" (0.150 mm)	5 nL	15,000 psi (1,034 bar)	ea.
UH-436	VHP Union for 360 µm OD Tubing, PEEK/SST	5/16-24 Coned	(2) PK-152, (2) P-416BLK	0.006" (0.150 mm)	5 nL	15,000 psi (1,034 bar)	ea.
UH-632	VHP True ZDV Union for 1/32" OD Tubing, PEEK/SST	6-32 Coned	(2) PK-126, (1) P-553 Gauge Plug	N/A	N/A	15,000 psi (1,034 bar)	ea.
1959-01	VHP Union for 1/32" OD Tubing, SST	M4x0.7	N/A (Fittings must be ordered separately)	0.007" (178 µm)	16 nL	30,000 psi (2,070 bar)	ea.



Convenient adapters for common 1/16" OD to capillary tubing

- ▶ Direct connect to 1/32" OD or 360 µm OD tubing options available
- > VHP adapters pressure rated to 12,000 psi (828 bar)



NOTE


While many 10-32 coned fittings are interchangeable, coned fittings using different threads are generally not interchangeable. As such, IDEX Health & Science recommends that only the style of coned fittings that accompanies these connectors be used for replacements.

MicroTight® Adapters

Create a true zero dead volume (ZDV) connection between 1/16" OD tubing and capillary tubing with our MicroTight Adapters.

For Very High Pressure applications the UH-630 will connect 1/16" OD to 1/32" OD tubing in an inline true ZDV connection with the ability to withstand 12,000 psi (828 bar)! The materials of construction also allow this product to be used up to 200 °C, which reduces the pressure rating to 8,000 psi (552 bar). For more information on the fittings used with the VHP adapter, please see page 59.



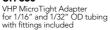



VHP MicroTight Adapting Cross 10-32 Coned for 1/16" OD tubing and 5/16-24 Coned for 360 µm OD tubing



1958-01 VHP MicroTight Adapter 10-32 Coned for 1/16" OD tubing and M4x0.7 for 1/32" OD tubing





VHP MicroTight Adapter for 1/16" and 360 µm OD tubing with fittings included



UH-753 VHP MicroTight Adapting Tee 360 µm (2 ports) to 10-32 C for 1/16" OD tubing (1 port)



UH-631-01 VHP MicroTight Adapter 10-32 Coned for 1/16" OD tubing and 6-40 Coned for 1/32" OD tubing fittings not included



MicroTight ZDV Adapter for 1/16" OD to MicroTight Sleeves with fittings included



MicroTight ZDV Adapter for 1/16" to 360 µm OD tubing with fittings included

MicroTight® Adapters (Cont.)



- > Replacement 6-32 fittings are on page 34.
- Replacement F-120 style nuts are on page 32 (when ordering, replace the "x" with an "R" or "B" to order either red or blue fittings).
- **)** Use this list to find micro flow products outside this chapter.

	Page
$360~\mu m,510~\mu m$ (0.020"), and 1/32" OD PEEK Tubing	16
360 μm OD Fused Silica Tubing	16
1/16" and 1/32" OD PEEKsil™ Tubing	22
1/32" OD FEP Tubing	26
360 µm OD High Purity PFA Tubing	24
510 µm (0.020") and 1/32" OD Stainless Steel Tubing	19
Polymer Capillary and Fused Silica Tubing Cutters	28
MX Series II™ Injection and Switching Valves	119
Manual Injection Valves	121
Micro Injection Port Adapters	132
Micro-Splitter Valves	140
Micro-Metering Valves	141
Microbore Guard Column	158
Ultra-Low Volume Back-Pressure Regulators	144
Nonmetallic 10-32 Micro-Volume Inline Check Valve	137

MicroTight® Adapters

Part No.	Description	Threads	Includes	Color	Swept Volume	Pressure Rating	Qty.
MICROTI	GHT ADAPTERS						
P-770	PEEK Micro Adapter, True ZDV, for 1/16" OD Tubing to MicroTight Tubing Sleeve	10-32 C to 6-32 C	(1) F-120, (1) F-125, (1) P-554	Natural	N/A	4,000 psi (276 bar)	ea.
P-881	PEEK Micro Adapter, True ZDV, for 1/16" to 1/32" OD Tubing	10-32 C to 6-32 C	(1) F-120R, (1) F-126S, (1) P-554	Red	N/A	5,000 psi (345 bar)	ea.
P-882	PEEK Micro Adapter, True ZDV, for 1/16" to 360 µm OD Tubing	10-32 C to 6-32 C	(1) F-120B, (1) F-124S, (1) P-554	Blue	N/A	5,000 psi (345 bar)	ea.
UH-630	Stainless Steel VHP Micro Adapter, for 1/16" to 1/32" OD Tubing	10-32 C to 6-32 C	(1) PK-120BLK, (1) PK-126, (1) P-554	SST/Black	N/A	12,000 psi (827 bar)	ea.
UH-634	Stainless Steel VHP Micro Adapter, for 1/16" to 360 µm OD Tubing	10-32 C to 6-32 C	(1) PK-120BLK, (1) PK-124, (1) P-554	SST/Black	N/A	12,000 psi (827 bar)	ea.
UH-753	Stainless Steel VHP Micro Adapting Tee, for 1/16" to 360 µm OD Tubing	10-32 C to 5/16-24 C	(2) P-416BLK, (2) PK-152	SST/Black	152 nL	15,000 psi (1,035 bar)*	ea.
1958-01	Stainless Steel VHP Micro Adapter, for 1/16" to 1/32" OD Tubing	10-32 C to M4x0.7 C	N/A	SST	16 nL	30,000 psi (2,070 bar)*	ea.
UH-631-01	Stainless Steel VHP Micro Adapter, for 1/16" to 1/32" OD Tubing	10-32 C to 6-40 C	N/A	SST	13 nL	30,000 psi (2,070 bar)*	ea.
UH-906	Stainless Steel VHP Micro Adapting Cross, for 1/16" to 360 µm OD Tubing	10-32 C to 5/16-24 C	(2) PK-120BLK, (2) P-416BLK, (2) PK-152	SST/Black	0.11 μL	15,000 psi (1,035 bar)*	ea.
REPLACE	MENT GAUGE PLUGS (TO ACHIEVE	TRUE ZDV CONNE	CTIONS WITH THE ABOVE AD	APTERS)			
P-554	Delrin® Gauge Plug	10-32 C		White	N/A	N/A	ea.
C = Coned * Pressure i	l rating depends upon the fitting used.						



Accessories



Insulating Mounting Bracket, shown with lead wire and Conductive Micro Tight Union, not included.

Insulating Mounting Bracket

Use our Insulating Mounting Bracket to easily integrate the Conductive MicroTight Union (shown on page 71) or our Conductive Mini MicroFilters (on page 109) into your system or lab.

The product snaps into place. Voltage from your lead wire is conducted through the attaching stainless steel nut and screw (included), then onto the mounted product via the stainless steel clip.

The bracket's base includes two holes (#2 screw clearance) for easy mounting onto any lab surface. Dimensions are 1.25" L x 0.45" W x 0.63" H.

Part No.	Description	Threads	Includes	Thru-hole	Swept Volume	Pressure Rating	Qty.
INSULA	TING MOUNTING BRACKET						
M-447	Insulating Mounting Bracket	N/A	N/A	N/A	N/A	N/A	ea.



- > For lab-on-a-chip applications
- Doptions to connect 1/16" OD Tubing directly, or 360µm and 1/32" OD Tubing with tubing sleeves
- Wetted materials: PEEK and perfluoroelastomer

NanoPort Assemblies

NanoPort Assemblies provide consistent fluid connections for chip-based analyses. NanoPort connections will bond to a variety of substrate materials with the use of Loctite.™

All NanoPort components are made of inert, biocompatible PEEK polymer (nuts and ports), Perlast® perfluoroelastomer (gaskets), and ETFE (ferrules). Their unique design also prevents adhesive contamination of the fluid path. And NanoPort connections add no additional volume to the fluid path, virtually eliminating dead volume traditionally associated with chip-based fluid connections.





Our NanoPort Assembly will readily connect 1/16" OD tubing with the included fittings. To connect 1/32" OD or 360µm OD, tubing sleeves for each size are included in each assembly.

Adhesive is not included in the N-333 NanoPort Assembly. Please contact IDEX Health & Science for bonding information or use common bonding adhesives such as Loctite.

Part No.	Description	Threads	For Chip Hole	Tubing OD	Qty.			
NANOPORT AS	SSEMBLIES							
10-32 Coned NanoPort Assembly								
N-333	F-333N	F-142N	Up to 0.063" (1.6 mm)	1/16"	ea.			
NANOPORT RE	PLACEMENT PARTS							
F-333Nx	Headless Fittings	10-32 C	Up to 0.063" (1.6 mm)	1/16"	10-pk			
F-142Nx	Ferrules	10-32 C	Up to 0.063" (1.6 mm)	1/16"	10-pk			
Gaskets								
N-123-02	Gasket, For all assemblies	except 6-32 Coned Assemblies	N/A	N/A	ea.			




Low Pressure Manifolds


Choose a 5, 7, or 9 Port Manifold to combine several streams into one, or split one fluid stream into several. Each PEEK manifold comes complete with 1/4-28 Super Flangeless™ Fittings for either 1/16" or 1/8" OD tubing, with pressure ratings of 2,000 psi (138 bar) and 500 psi (34 bar), respectively.

A few useful applications include:

- Multiport mixing chamber
- Gas sparging splitting union
- Sample injection onto multi-well plates or a multiple direction flow path union

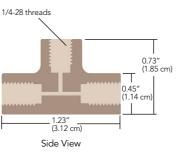



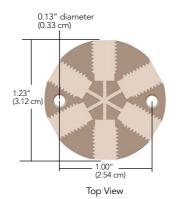


Part No.	Description	Threads	Includes	Thru-hole	Swept Volume	Pressure Rating	Qty.	
MANIFO	MANIFOLDS CONTROL OF THE PROPERTY OF THE PROPE							
Standard	d .							
P-150	PEEK 7-Port Manifold for 1/16" OD Tubing	1/4-28 FB	(7) P-255, (7) P-250	0.040" (1.00 mm)	42.0 µL	1,000 psi (69 bar)	ea.	
P-154	PEEK 5-Port Manifold for 1/16" OD Tubing	1/4-28 FB	(5) P-255, (5) P-250	0.040" (1.00 mm)	22.3 µL	1,000 psi (69 bar)	ea.	
P-155	PEEK 5-Port Manifold for 1/8" OD Tubing	1/4-28 FB	(5) P-331, (5) P-359	0.062" (1.60 mm)	53.8 µL	500 psi (34 bar)	ea.	
P-190	PEEK 9-Port Manifold for 1/8" OD Tubing	1/4-28 FB	(9) P-331, (9) P-359	0.062" (1.60 mm)	160 μL	500 psi (34 bar)	ea.	
P-191	PEEK 9-Port Manifold for 1/16" OD Tubing	1/4-28 FB	(9) P-255, (9) P-250	0.040" (1.00 mm)	139 µL	1,000 psi (69 bar)	ea.	

FB = Flat Bottom



Low Pressure Bulkhead Unions

- Designed for plumbing tubing through equipment housing
- For use with standard 10-32 coned or 1/4-28 flat-bottom threaded fittings

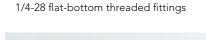





Thread PEEK Bulkhead Unions directly through your equipment housing to connect internal tubing to the outside. Each union has unique 3/8-24 external threads and comes complete with a stainless steel nut and lock washer to hold it in place. Requires a 3/8" hole to mount. The recommended torque limit for these unions is 15 in.— lbs (1.7 N·m).



Bulkhead Union includes stainless steel nut/lock washer

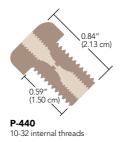

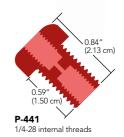

P-430 PEEK Elbow comes with Flangeless Fittings



Elbow Connectors

Use these Elbow Connectors to easily navigate tight corners. One Elbow is designed for use with 1/16" OD tubing and has a 0.020" (0.50 mm) thru-hole. Use 1/8" OD tubing with the other Elbow, which has a 0.062" (1.6 mm) thru-hole. Both come complete with 1/4-28 PEEK nuts and ETFE ferrules, and are pressure rated to 1,000 psi (69 bar).

Large Bore Union

> 5/16-24 flat-bottom threads

Use any of the 5/16-24 fittings on page 55 and the appropriate ferrule to create a true zero dead volume (ZDV) connection with the P-134 Union.



- > Stainless Steel Bulkhead Unions are also available. Please contact us for more information.
- To use Elbows in higher pressure applications, simply replace the provided fittings with Super Flangeless™ Nuts and Ferrules, found on page 39.


Part No.	Description	Threads	Color	Includes	Thru-hole	Swept Volume	Qty.			
BULKHEA	BULKHEAD UNIONS									
P-440	PEEK Bulkhead Union	10-32 Coned	Natural	(1) SST Nut/Washer	0.020" (0.50 mm)	1.9 μL	ea.			
P-441	PEEK Bulkhead Union	1/4-28 Flat-Bottom	Red	(1) SST Nut/Washer	0.040" (1.00 mm)	2.9 μL	ea.			
P-441N	PEEK Bulkhead Union	1/4-28 Flat-Bottom	Natural	(1) SST Nut/Washer	0.040" (1.00 mm)	2.9 µL	ea.			
ELBOW CO	ONNECTORS									
P-430	PEEK Elbow for 1/16" OD Tubing	1/4-28 Flat-Bottom	Natural	(2) XP-235	0.020" (0.50 mm)	1.4 µL	ea.			
P-432	PEEK Elbow for 1/8" OD Tubing	1/4-28 Flat-Bottom	Natural	(2) XP-335	0.062" (1.60 mm)	13.6 µL	ea.			
LARGE BC	LARGE BORE UNION									
P-134	PEEK True ZDV Union	5/16-24 Flat-Bottom	Natural	N/A	N/A	N/A	ea.			



Low Pressure Y Connectors

PEEK Y Connectors are designed to split a stream or join two streams together, just like a tee. However, the configuration of a tee can lead to turbulent flow and solvent outgassing, which increases baseline noise and reduces sensitivity. The geometry of a Y connector creates less turbulence and thus can improve analytical results.

All of these Y Connectors use 1/4-28 Flangeless fittings, except P-515 which uses 5/16-24 fittings (to accommodate larger tubing).

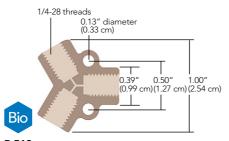

P-512PEEK Y comes with Flangeless Fittings

Part No.	Description	Threads	Includes	Thru-hole	Swept Volume	Pressure Rating	Qty.
Y CONN	IECTORS						
P-512	PEEK Y for 1/16" OD Tubing	1/4-28 FB	(3) XP-235	0.020" (0.50 mm)	1.7 µL	1,000 psi (69 bar)	ea.
P-513	PEEK Y for 1/8" OD Tubing	1/4-28 FB	(3) XP-335	0.040" (1.00 mm)	6.0 µL	500 psi (34 bar)	ea.
P-514	PEEK Y for 1/8" OD Tubing	1/4-28 FB	(3) XP-335	0.060" (1.50 mm)	13.6 µL	500 psi (34 bar)	ea.
P-515	PEEK Y for 3/16" OD Tubing	5/16-24 FB	(3) XP-132	0.125" (3.20 mm)	47.7 μL	500 psi (34 bar)	ea.

FB = Flat-Bottom




- Manufactured from PEEK, ETFE, Delrin®, polypropylene, or PCTFE
- Available with 1/4-28, M6, or 10-32 flat-bottom threads

Low Pressure Unions


Our Low Pressure Unions are available in a variety of polymers, providing several low-cost and chemically-resistant options. The union assemblies below include fittings as shown in the table. The unions in the right column do not include fittings, allowing for customizing the fitting selection. In some cases, a union can be configured to connect two different tubing sizes—for example, if 1/4-28 Flangeless fittings for 1/16" and 1/8" OD tubing were selected from page 45 they can be used with the P-603 union to connect the two different tubing sizes.

Low Pressure PEEK Union Assemblies

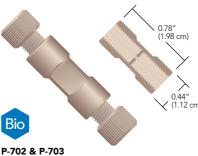


P-702 & P-703
PEEK Unions with Flangeless Fittings (1/4-28 internal threads)


Low Pressure ETFE Union Assemblies

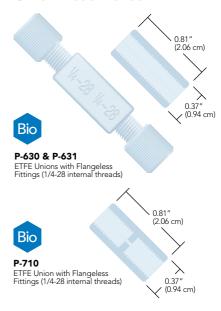

Low Pressure Standard Unions

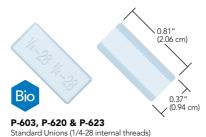

Low Pressure Metric Unions

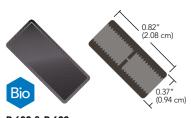

P-602 & P-622 Metric Unions (M6 internal threads)

Low Pressure Male Union



VacuTight[™] Union





- To use connectors in higher pressure applications, simply replace the provided fittings with Super Flangeless™ Nuts and Ferrules, found on page 39.
- Use any of the 10-32 flat-bottom fittings on 39 and 42 to make an inline connection with our VacuTight Union. This product is designed for use with 1/16" OD tubing.

Low Pressure Unions

Part No.	Description	Color	Threads	Includes	Thru-hole	Swept Volume	Pressure Rating	Qty.
PEEK UN	IION ASSEMBLIES							
P-702	PEEK Union for 1/16" OD Tubing	Natural	1/4-28 FB	(2) XP-235	0.020" (0.50 mm)	0.41 μL	1,000 psi (69 bar)	ea.
P-703	PEEK Union for 1/8" OD Tubing	Natural	1/4-28 FB	(2) XP-335	0.050" (1.25 mm)	2.57 μL	1,000 psi (69 bar)	ea.
ETFE UN	ION ASSEMBLIES							
P-630	ETFE True ZDV Union for 1/16" OD Tubing	Natural	1/4-28 FB	(2) P-200N/P-245	N/A	N/A	1,000 psi (69 bar)	ea.
P-631	ETFE True ZDV Union for 1/8" OD Tubing	Natural	1/4-28 FB	(2) P-300N/P-345	N/A	N/A	1,000 psi (69 bar)	ea.
P-710	ETFE Union for 1/16" OD Tubing	Natural	1/4-28 FB	(2) XP-245	0.030" (0.75 mm)	0.93 μL	1,000 psi (69 bar)	ea.
STANDA	RD UNIONS							
P-603	Delrin True ZDV Standard Union	Natural	1/4-28 FB	N/A	N/A	N/A	N/A*	ea.
P-620	Polypropylene True ZDV Standard Union	Natural	1/4-28 FB	N/A	N/A	N/A	N/A*	ea.
P-623	ETFE True ZDV Standard Union	Natural	1/4-28 FB	N/A	N/A	N/A	N/A*	ea.
METRIC	UNIONS							
P-602	Delrin Metric Union	Black	M6 FB	N/A	0.020" (0.50 mm)	0.41 μL	N/A*	ea.
P-622	ETFE Metric Union	Blue	M6 FB	N/A	0.020" (0.50 mm)	0.41 μL	N/A*	ea.
MALE UI	NION							
P-645	PCTFE Male Union	Natural	1/4-28 FB	N/A	0.062" (1.60 mm)	61.3 µL	500 psi (34 bar)	ea.
VACUTIO	SHT UNION							
P-845-01	PEEK Union for 1/16" OD Tubing	Red	10-32 FB	N/A	0.020" (0.50 mm)	0.20 μL	N/A*	ea.
* Pressure FB = Flat-E	Rating depends on Fittings selected. See pressu Bottom	re rating for fitting	gs on appropriate	e page.				



Low Pressure Tees & Crosses

0.49" (1.24 cm) Mounting Holes 0.13" diameter (0.33 cm) Bio P-632 & P-633 0.37" (0.94 cm) ETFE Tees come with Flangeless Fittings 1/4-28 threads 1.08 0.13" diameter P-722

Our Low Pressure Tees and Crosses are available in two inert polymers and can handle pressures to 500 psi (34 bar) or 1,000 psi (69 bar), depending upon the configuration of the products. Each is designed with handy mounting holes. All ETFE Tees and Crosses ship complete with 1/4-28 PFA Flangeless nuts and ETFE ferrules, while their PEEK polymer counterparts ship with 1/4-28 PEEK nuts and ETFE ferrules. Replacement fittings are located on page 47.





comes with Flangeless Fittings

PEEK Cross

To order just the body of one of our tees and crosses without fittings, simply add a '-01' to the part number — e.g., P-632-01.

RELATED PRODUCTS

- > Seal off unused ports with any of our 1/4–28 flat-bottom plugs found on page 55.
- To use the PEEK polymer versions of our Tees and Crosses in higher pressure and Ferrules, found on page 39.
- ▶ High Pressure Tees, Crosses, and a 7-Port Manifold (all with 10-32 threaded ports) are on page 74.

Part No.	Description	Threads	Includes	Thru-hole	Swept Volume	Pressure Rating	Qty.
LOW PRESS	URE TEES AND CROSSES						
P-632	ETFE Tee for 1/16" OD Tubing	1/4-28 Flat-Bottom	(3) P-245, (3) P-200N	0.020" (0.50 mm)	2.9 μL	1,000 psi (69 bar)	ea.
P-633	ETFE Tee for 1/8" OD Tubing	1/4-28 Flat-Bottom	(3) P-345, (3) P-300N	0.050" (1.25 mm)	17.5 μL	500 psi (34 bar)	ea.
P-634	ETFE Cross for 1/16" OD Tubing	1/4-28 Flat-Bottom	(4) P-245, (4) P-200N	0.020" (0.50 mm)	3.8 µL	1,000 psi (69 bar)	ea.
P-635	ETFE Cross for 1/8" OD Tubing	1/4-28 Flat-Bottom	(4) P-345, (4) P-300N	0.050" (1.25 mm)	22.8 μL	500 psi (34 bar)	ea.
P-712	PEEK Tee for 1/16" OD Tubing	1/4-28 Flat-Bottom	(3) XP-235	0.020" (0.50 mm)	2.9 μL	1,000 psi (69 bar)	ea.
P-713	PEEK Tee for 1/8" OD Tubing	1/4-28 Flat-Bottom	(3) XP-335	0.050" (1.25 mm)	17.5 μL	500 psi (34 bar)	ea.
P-714	PEEK Tee for 1/16" OD Tubing	1/4-28 Flat-Bottom	(3) XP-235	0.040" (1.00 mm)	11.4 μL	1,000 psi (69 bar)	ea.
P-722	PEEK Cross for 1/16" OD Tubing	1/4-28 Flat-Bottom	(4) XP-235	0.020" (0.50 mm)	3.8 µL	1,000 psi (69 bar)	ea.
P-723	PEEK Cross for 1/8" OD Tubing	1/4-28 Flat-Bottom	(4) XP-335	0.050" (1.25 mm)	22.8 µL	500 psi (34 bar)	ea.



Quick Connect Luer Adapters

Delrin®, polypropylene, ETFE, or PEEK Versions

Adapts luers to 1/4-28, 10-32, 5/16-24, or M6 threaded ports

These luer adapters were designed to work in a variety of applications. By connecting any male luer to any female luer, you can create your own quick connect union or adapter. Each Quick Connect Luer Adapter conforms to ISO requirements for medical luer taper configuration and performance (45 psi/3.1 bar).

Find fittings to connect tubing to the threaded ports of these adapters in the Fittings chapter, starting on page 30.

Please Note: Our Female Quick Connect Luer Adapters can be used with any of the Male Luers on this page, i.e., those with and without lock hubs.




- Our A-626 Bottle Cap Plug (page 55) can be used to plug any of the female luer adapters on this page.
- > To prevent a chemical spill when disconnecting your solvent reservoir tubing from the pump, try our Quick-Stop Luer Check Valve on page 139.
- > To economically prime an HPLC pump, simply remove the 10-32 fitting on the outlet check valve (standard on most models), insert a P-642 luer adapter, attach a syringe (such as our B-310) and draw the mobile phase through the pump head.




P-604, P-618, P-624 Female Luer to 1/4-28 Male (luer end of P-604 different than shown)



P-619, P-625 Male Luer to 1/4-28 Male

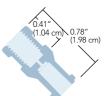

P-628 Female Luer to 1/4-28 Female

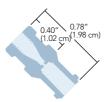

P-629 Female Luer to 10-32 Female



P-642 Female Luer to 10-32 Male




P-719 Female Luer to 10-32 Male





Male Luer to 10-32 Female

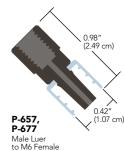

0.75" (1.91 cm)

Quick Connect Luer Adapters (Cont.)


Luer-To-MicroTight® Adapter

> Easily connect 360 μm OD tubing to a syringe



P-662 Luer-to-MicroTight Adapter, shown with a B-310 Syringe (see table below) and PEEK capillary tubing (page 16), not included.

The Luer-to-MicroTight Adapter is ideal for infusing sample into lab-on-a-chip devices. This product is made entirely of biocompatible PEEK polymer and introduces only 14 nL of additional volume to the flow path. Use it to directly connect a luer-tip syringe or other product that terminates with a standard male luer to 360 μm OD capillary tubing without tubing sleeves (see photo). MicroTight Fittings are included.

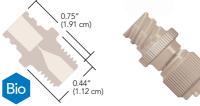

P-662 Luer-To-MicroTight Adapter for Luer to 360 μm OD tubing with fittings included

Quick Connect Luer Adapters

Part No.	Description	Body Material	Lock Hub Material	Thru-hole	Qty.
QUICK CONI	NECT LUER ADAPTERS				
P-604	F Luer to 1/4-28 FB, M	Nat. Delrin	N/A	0.05" (1.3 mm)	ea.
P-618	F Luer to 1/4-28 FB, M	Nat. PP	N/A	0.05" (1.3 mm)	ea.
P-619	M Luer to 1/4-28 FB, M	Nat. PP	None *	0.05" (1.3 mm)	ea.
P-624	F Luer to 1/4-28 FB, M	Nat. ETFE	N/A	0.05" (1.3 mm)	ea.
P-625	M Luer to 1/4-28 FB, M	Nat. ETFE	None *	0.04" (1.0 mm)	ea.
P-628	F Luer to 1/4-28 FB, F	Nat. ETFE	N/A	0.04" (1.0 mm)	ea.
P-629	F Luer to 10-32 C, F	Nat. ETFE	N/A	0.04" (1.0 mm)	ea.
P-642	F Luer to 10-32 C, M	Nat. ETFE	N/A	0.05" (1.3 mm)	ea.
P-655	M Luer to 1/4-28 FB, F	Red PEEK	Black PEEK	0.04" (1.3 mm)	ea.
P-656	M Luer to 10-32 C, F	Nat. PEEK	Black PEEK	0.05" (1.3 mm)	ea.
P-657	M Luer to M6 FB, F	Black PEEK	Black PEEK	0.05" (1.3 mm)	ea.
P-658	F Luer to 1/4-28 FB, F	Red PEEK	N/A	0.05" (1.3 mm)	ea.
P-659	F Luer to 10-32 C, F	Nat. PEEK	N/A	0.05" (1.3 mm)	ea.
P-660	F Luer to M6 FB, F	Black PEEK	N/A	0.05" (1.3 mm)	ea.
P-661	F Luer to 5/16-24 FB, M	Nat. ETFE	N/A	0.05" (1.3 mm)	ea.
P-675	M Luer to 1/4-28 FB, F	Red ETFE	Natural PP	0.05" (1.3 mm)	ea.
P-677	M Luer to M6 FB, F	Black ETFE	Natural PP	0.05" (1.3 mm)	ea.
P-678	F Luer to 1/4-28 FB, F	Red ETFE	N/A	0.05" (1.3 mm)	ea.
P-680	F Luer to M6 FB, F	Black ETFE	N/A	0.05" (1.3 mm)	ea.
P-683	M Luer to 1/4-28 FB, M	Nat. PEEK	Black PEEK	0.04" (1.0 mm)	ea.
P-686	F Luer to M6 FB, M	Black ETFE	N/A	0.05" (1.3 mm)	ea.
P-719	F Luer to 10-32 C, M	Nat. PEEK	N/A	0.05" (1.3 mm)	ea.
SYRINGE WIT	TH MALE LUER LOCK				
B-310	10 cc Disposable Luer-Lock Syringe. For	use with any Female Luer Adapter		0.05" (1.3 mm)	ea.
LUER-TO-MIC	ROTIGHT ADAPTER				
P-662	Luer-to-MicroTight Adapter	(1) F-152, (1) P-416	0.006" (0.150 mm)	45 psi (2.4 bar)	ea.
	rnal) threads; M = Male (external) threads; Nat. = ene; FB = Flat-Bottom; C = Coned luer.	= Natural; N/A = Not Applicable;			




LuerTight Fittings


- Luer fittings for fluoropolymer tubing
- Quick disconnect and barbless
- For 1/16" and 1/8" OD tubing

Our LuerTight fittings are specifically designed to connect fluoropolymer tubing without barbs or nuts! By integrating ferrules into the luer bodies, LuerTights will reliably hold your tubing in place while giving you the convenience of a luer connection. An inline set of LuerTight fittings provides a quick and easy disconnection option. LuerTight connections are also less bulky and more economical than nut-to-luer style fittings.

The bodies of these products are manufactured from polypropylene and the ferrules, where used, are made of ETFE.

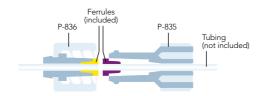

P-837 LuerTight Fittings System for 1/16" OD tubing

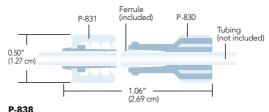

LuerTight Fittings System for 1/8" OD tubing



LuerTight fittings are designed to be used exclusively within the LuerTight family. Combining LuerTight fittings with non-LuerTight luer products may result in a poor connection.

Part No.	Description	Includes	Thru-hole	Pressure Rating	Qty.
LUERTIGH'	T FITTINGS SYSTEMS				
P-837	LuerTight System for 1/16" OD Tubing	(1) P-835, (1) P-836, (1) P-830T	N/A	100 psi (7 bar)	ea.
P-838	LuerTight System for 1/8" OD Tubing	(1) P-830, (1) P-831, (1) P-830T	N/A	100 psi (7 bar)	ea.
LUERTIGH"	T FITTING COMPONENTS				
P-830	Female Fitting for 1/8" OD Tubing	(1) Ferrule	N/A	N/A	ea.
P-830T	Set Plug to swage Ferrules into P-835 and P-830	N/A	N/A	N/A	ea.
P-831	Male Fitting for 1/8" OD Tubing	No Ferrule Required	N/A	N/A	ea.
P-835	Female Fitting for 1/16" OD Tubing	(1) Ferrule	N/A	N/A	ea.
P-836	Male Fitting for 1/16" OD Tubing	(1) Ferrule	N/A	N/A	ea.
Female = inte	ernal receiving luer pocket; Male = external luer nose (surr	ounded by internally-threaded locking ring)			



Barbed Connectors

Our Type 1 Barbed Unions have been engineered to effectively join two pieces of soft-walled tubing together. This type of connector is typically the connector of choice for joining two peristaltic tubes with similar inner diameters together. Our unions are manufactured from either polypropylene or nylon.



Barbed to Barbed Adapters

- Adapters on this page feature various luer to barb adaptations
- Adapters on the next page feature a variety of barb-to-barb connectors

Use these barbed adapters to connect peristaltic-type flexible tubing for general, low pressure applications. The polypropylene used to manufacture the majority of these products is a Class VI material. Due to the low melt point of polypropylene (PP), these adapters are not autoclavable, however, they can be sterilized via gamma radiation. There are also Barb to Female Luer-Lock connectors available from ETFE, which has superior solvent resistance and a higher temperature rating (80 °C).




Barbed "Y" Adapters

Our Barbed "Y" Adapters, manufactured from polypropylene, are engineered to effectively join three pieces of soft-walled tubing together in a Y configuration, offering less turbulence and gentler mixing of fluids than a traditional Tee Connector. This type of connector works well for joining three peristaltic tubes with similar inner diameters together.



- Three barb sizes, for 1/16", 1/8", and 3/16" ID flexible tubing
- Adapt to 1/4-28 flat-bottom, 5/16-24 flat-bottom, or 10-32 coned receiving ports

These adapters make it easy to connect flexible tubing to any standard 1/4-28 flat-bottom or 10-32 coned receiving port. Simply thread the adapter into a receiving port and slip tubing over the barbed stem to create a reliable low pressure connection.

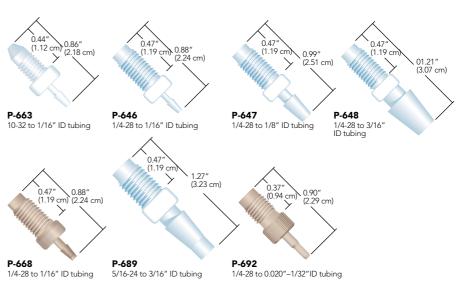





- To connect low pressure fluoropolymer tubing, try the LuerTight™ Adapters on page 89.
- To connect peristaltic tubing to low pressure fluoropolymer tubing, see page 92.
- ▶ For peristaltic tubing, see page 15.



Swivel Barb Adapters


- ▶ Barb connection spins freely from the nut to prevent twist during installation
- Manufactured from polypropylene

The Swivel Barb Adapters from IDEX Health & Science are made up of two captive pieces acting as a one-piece fitting for ease of use. Manufactured from polypropylene and available in three barb sizes, the Swivel Barb will facilitate connection between flexible tubing to a 1/4-28 flat-bottom port. The barbed insert spins freely from the threaded nut in order to prevent the tubing from twisting during installation.



Luer to Barbed Adapters

Our Luer to Barbed Adapters are an excellent choice when connecting between soft-walled tubing and luer-based products, such as a syringe or a low-pressure filter, for example. We offer several different configurations, allowing you to connect various sizes of soft-walled tubing to receiving ports that accept a male luer "slip" style connection; a male luer "lock" style connection; and a female-luer style of connector.

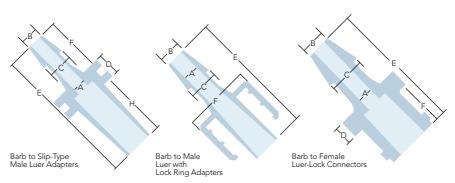

Barbed Connectors (Cont.)



Conical Adapters

- Direct connect 1/16" and 1/8" OD rigid and semi-rigid tubing to peristaltic tubing
- ▶ Accept 0.020"-1/8" (0.50-3.2 mm) ID peristaltic tubing
- Biocompatible flow path with excellent chemical compatibility

Conical Adapters provide a reliable connection between rigid/semi-rigid tubing and peristaltic-type flexible tubing, such as Tygon® and PharMed®. These adapters are composed of a PEEK polymer female nut, our Super Flangeless™ ferrule system and an ETFE or PEEK conical adapter body. The narrow coned end of the adapter body allows peristaltic tubing to slide on more easily than it does onto conventional barbed adapters. Peristaltic tubing is also easier to remove from our Conical Adapters, since no cutting or excessive pulling is required.



APPLICATION NOTE

To help secure peristaltic tubing more firmly to the Conical Adapters, simply attach a cable tie to the outside of the peristaltic tubing once it has been placed onto the Adapter body.








Peristaltic Tubing Adapters

These unique adapters connect peristaltic tubing to standard 1/16" or 1/8" OD tubing. A specially-designed nose allows the peristaltic tubing to simply press fit over the nose and then be held tightly in place by the retainer sleeve. Your 1/16" OD tubing may then be connected with the Flangeless Fittings supplied with the adapter. To connect your peristaltic tubing to tubing with a different OD, simply replace the supplied fittings with your choice of Flangeless Fittings from page 45.

One popular application for these adapters is to use them as "stops" for your peristaltic pump. By doing so, you can reduce the amount of peristaltic tubing required for your flow path, thus reducing cost.

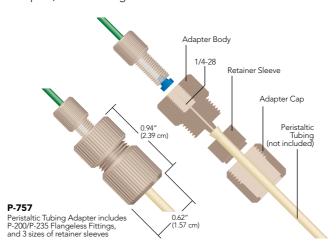

COLUMN HARDWARE

Barbed Connectors

Part No.	Tubing ID		Material			Qty.
BARBED T	TO BARBED ADAPTERS					
P-801	0.06" (1.5 mm)		Polypropylene			ea.
P-802	0.12" (3.0 mm)		Polypropylene			ea.
BARBED '	"Y" CONNECTORS					
P-860	0.06" (1.5 mm)		Polypropylene			ea.
P-861	0.10" (2.5 mm)		Polypropylene			ea.
P-862	0.12" (3.0 mm)		Polypropylene			ea.
P-863	0.18" (4.8 mm)		Polypropylene			ea.
P-864	0.25" (6.4 mm)		Polypropylene			ea.
THREAD '	TO BARBED ADAPTERS					
Part No.	Description		Material	Threads	Thru-hole	Qty.
P-663	Barb Adapter, 1/16" (1.55 mm) ID Tubing		ETFE	10-32 Coned	0.04" (1.0 mm)	ea.
P-646	Barb Adapter, 1/16" (1.55 mm) ID Tubing		ETFE	1/4-28 Flat-Bottom	0.04" (1.0 mm)	ea.
P-647	Barb Adapter, 1/8" (3.20 mm) ID Tubing		ETFE	1/4-28 Flat-Bottom	0.08" (2.0 mm)	ea.
P-648	Barb Adapter, 3/16" (4.75 mm) ID Tubing		ETFE	1/4-28 Flat-Bottom	0.10" (2.5 mm)	ea.
P-668	Barb Adapter, 1/16" (1.55 mm) ID Tubing		PEEK	1/4-28 Flat-Bottom	0.04" (1.0 mm)	ea.
P-689	Barb Adapter, 3/16" (4.75 mm) ID Tubing		ETFE	5/16-24 Flat-Bottom	0.10" (2.5 mm)	ea.
P-692	Barb Adapter, 0.020" to 1/32" (0.50 to 0.8	0 mm) ID Tubing	PEEK	1/4-28 Flat-Bottom	0.02" (0.5 mm)	ea.
SWIVEL B	ARB ADAPTERS					
D-646	Swivel Barb Adapter, 1/16" (1.55 mm) ID	Tubing	Polypropylene	1/4-28 Flat-Bottom	0.03" (0.75 mm)	ea.
D-647	Swivel Barb Adapter, 3/32" (2.40 mm) ID	Tubing	Polypropylene	1/4-28 Flat-Bottom	0.056" (1.5 mm)	ea.
D-648	Swivel Barb Adapter, 1/8" (3.20 mm) ID To	ubing	Polypropylene	1/4-28 Flat-Bottom	0.08" (2.0 mm)	ea.
BARB TO	SLIP-TYPE MALE LUER ADAPTERS					
Part No.	Description		Material			Qty.
P-854	Male Luers (Slip-type) for use with 1/16" I	D (1.55 mm) Tubing	PP			_
P-834	A=0.046" B=0.064" C=0.090" D=0.129" I	E=0.711" F=0.198" H=0.384"	rr			ea.
	ype male luer fittings are for use in systems					
BARB TO	MALE LUER WITH LOCK RING ADA	NPTERS				
P-850	Male Luers with Lock Ring for use with 1/ A=0.049" B=0.065" C=0.090" E=0.583" F		PP			ea.
P-851	Male Luers with Lock Ring for use with 3/ A=0.071" B=0.100" C=0.139" E=0.681" F	32" ID (2.40 mm) Tubing	PP			ea.
P-852	Male Luers with Lock Ring for use with 1/A=0.099" B=0.132" C=0.184" E=0.777" F	8" ID (3.20 mm) Tubing	PP			ea.
BARB TO	FEMALE LUER-LOCK CONNECTOR					
	Female Luer Connectors for use with 1/10					
P-857	A=0.030" B=0.063" C=0.106" D=0.100" I		PP			ea.
P-858	Female Luer Connectors for use with 3/32 A=0.056" B=0.102" C=0.145" D=0.100" I		PP			ea.
P-859	Female Luer Connectors for use with 1/8' A=0.080" B=0.135" C=0.187" D=0.100" I		PP			ea.
P-870	For use with 1/16" (1.55 mm) ID Tubing A=0.030" B=0.063" C=0.106" D=0.100" I	E=0.598" F=0.253"	ETFE			ea.
P-872	For use with 1/8" (3.20 mm) ID Tubing A=0.080" B=0.137" C=0.187" D=0.100" I	E=0.733" F=0.253"	ETFE			ea.
PERISTAL	TIC TUBING ADAPTERS					
Part No.	Description	Tubing OD	Peristaltic Tubing ID		Thru-Hole	Qty.
P-757	Standard Adapter	up to 0.180" (4.55 mm)	0.048"- 0.110" (1.20-	-2.80 mm)	0.030" (0.75 mm)	ea.
P-767	Large Bore Adapter	up to 0.250" (6.35 mm)	0.100"- 0.150" (2.55-		0.070" (1.78 mm)	ea.
CONICAL	. ADAPTER ASSEMBLIES					
Part No.	Description	Rigid or Semi-Rigid Tubing OD	Peristaltic Tubing ID		Thru-Hole	
P-794	Conical Adapter	1/16"	0.020"-0.030" (0.50 r	mm-0.75 mm)	0.020" (0.50 mm)	ea.
P-797	Conical Adapter	1/16"	1/16"-3/32" (1.55 mr		0.040" (1.0 mm)	ea.
P-798	Conical Adapter	1/8"	1/16"-3/32" (1.55 mr		0.040" (1.0 mm)	ea.
2-799	Conical Adapter	1/8"	3/32"-1/8" (2.40 mm		0.060" (1.5 mm)	ea.
	ADAPTER REPLACEMENT PARTS					
Part No.	Description	For Use With	Material			
F-156	Female Nut, 1/8", 1/4-28	P-798, P-799	Black PEEK			ea.
P-420	Female Nut, 1/16", 1/4-28	P-794, P-797	Natural PEEK			ea.
P-259	Super Flangeless Ferrule, 1/16"	P-794, P-797	Yellow ETFE/SST			ea.
,		P-798, P-799	Yellow ETFE/SST			ea.
2-359						ca.
P-359 P-691	Super Flangeless Ferrule, 1/8" Conical Adapter Body	P-799	Natural ETFE			ea.





FILTERS & FRITS

Our Filters offer an optimal way to filter your solvents, preventing pump cavitation and system damage. We offer different style filters for specific system specifications. Our filters protect your system from particulate matter from the solvent that may otherwise damage expensive hardware.

We offer a complete line of Frits manufactured from two different materials: PEEK and stainless steel. Both materials offer a variety of sizes of frit discs, as well as being available in numerous porosities. All our frits are designed with exceptional uniform porosity and a long filtration life.

95 FRITS

100 FILTERS

111 BOTTLE CAPS & PLUGS



Stainless Steel Frits

Our Analytical-scale 316 Stainless Steel Frits are available in 0.5 µm or 2 µm porosity—the most common HPLC filtration ratings. Each frit includes a PCTFE or PEEK polymer sealing ring.


Many of the frits shown have the common 0.250" (0.64 cm) and 0.254" (0.64 cm) ODs, which allow them to be used in many of the Precolumn and Inline Filters found starting on page 103. Choose the larger diameter faces and/or larger porosity frits for faster flow rates. Choose frits with a smaller diameter face and/or smaller porosity for applications sensitive to extra flow path volume.

APPLICATION NOTE


To Clean Or Not To Clean?

It is rarely worth the time and effort to clean frits, given the relatively low cost of replacements. Furthermore, cleaning may leave some debris embedded in the frit pores. If the washed frit is accidently returned to your instrument in a reverse orientation, any remaining debris could be flushed out and deposited further down the fluid path. If this frit is being used as a column head frit, the debris may be washed directly onto the column bed.

0.5 µm Stainless Steel Frits

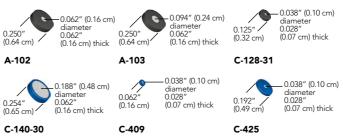

2 µm Stainless Steel Frits

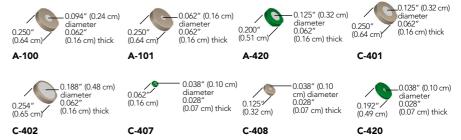

1 NOTE

Frits without the polymer rings cannot be used with our standard Precolumn and Inline Filter assemblies.

Semi-Prep Stainless Steel Frits

Many of these frits come complete with a PCTFE, ETFE, or PTFE sealing ring. Choose from 2 μm , 5 μm , 10 μm , and 20 μm filtration porosities and a range of diameters to match your intended flow rate and filtration requirements.

2 µm Semi-Prep Stainless Steel Frits




Stainless Steel Frits (Cont.)


5 µm Semi-Prep Stainless Steel Frits



10 µm Semi-Prep Stainless Steel Frits

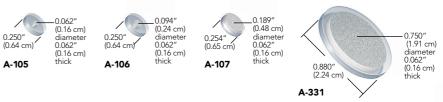

20 µm Semi-Prep Stainless Steel Frits

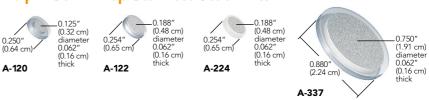

Stainless Steel Frits

Part No.	Porosity	Disc Diameter	Disc Thickness	Ring OD	Ring Material	Frit Volume	Qty.
STAINLESS ST	TEEL FRITS						
A-100	2 µm	0.094" (0.24 cm)	0.062" (0.16 cm)	0.250" (0.64 cm)	PEEK	1.7 µL	ea.
A-101	2 µm	0.062" (0.16 cm)	0.062" (0.16 cm)	0.250" (0.64 cm)	PEEK	0.7 µL	ea.
A-102	0.5 µm	0.062" (0.16 cm)	0.062" (0.16 cm)	0.250" (0.64 cm)	PEEK	0.6 μL	ea.
A-103	0.5 µm	0.094" (0.24 cm)	0.062" (0.16 cm)	0.250" (0.64 cm)	PEEK	1.4 µL	ea.
A-420	2 µm	0.125" (0.32 cm)	0.062" (0.16 cm)	0.200" (0.51 cm)	PCTFE	3.0 µL	ea.
C-128-31	0.5 µm	0.038" (0.10 cm)	0.028" (0.07 cm)	0.125" (0.32 cm)	PEEK	0.1 μL	ea.
C-140-30	0.5 µm	0.188" (0.48 cm)	0.062" (0.16 cm)	0.254" (0.65 cm)	PCTFE	6.5 µL	ea.
C-401	2 µm	0.125" (0.32 cm)	0.062" (0.16 cm)	0.250" (0.64 cm)	PEEK	3.0 µL	ea.
C-402	2 µm	0.188" (0.48 cm)	0.062" (0.16 cm)	0.254" (0.65 cm)	PEEK	7.8 µL	ea.
C-407	2 µm	0.038" (0.10 cm)	0.028" (0.07 cm)	0.062" (0.16 cm)	PCTFE	0.1 μL	ea.
C-408	2 µm	0.038" (0.10 cm)	0.028" (0.07 cm)	0.125" (0.32 cm)	PEEK	0.1 μL	ea.
C-409	0.5 µm	0.038" (0.10 cm)	0.028" (0.07 cm)	0.062" (0.16 cm)	PCTFE	0.1 μL	ea.
C-420	2 µm	0.038" (0.10 cm)	0.028" (0.07 cm)	0.192" (0.49 cm)	PCTFE	0.1 μL	ea.
C-425	0.5 µm	0.038" (0.10 cm)	0.028" (0.07 cm)	0.192" (0.49 cm)	PCTFE	0.1 μL	ea.
SEMI-PREP ST	TAINLESS STEEL FRI	TS					
A-105	10 µm	0.062" (0.16 cm)	0.062" (0.16 cm)	0.250" (0.64 cm)	PCTFE	0.9 μL	ea.
A-106	10 µm	0.094" (0.24 cm)	0.062" (0.16 cm)	0.250" (0.64 cm)	PCTFE	2.0 µL	ea.
A-107	10 µm	0.189" (0.48 cm)	0.062" (0.16 cm)	0.254" (0.65 cm)	PCTFE	9.1 µL	ea.
A-120	20 µm	0.125" (0.32 cm)	0.062" (0.16 cm)	0.250" (0.64 cm)	PCTFE	3.7 µL	ea.
A-122	20 µm	0.188" (0.48 cm)	0.062" (0.16 cm)	0.254" (0.65 cm)	PCTFE	9.7 µL	ea.
A-224	20 μm	0.188" (0.48 cm)	0.062" (0.16 cm)	0.254" (0.65 cm)	PTFE	9.7 µL	ea.
A-331	10 μm	0.750" (1.91 cm)	0.062" (0.16 cm)	0.880" (2.24 cm)	ETFE	141.9 μL	ea.
A-332	2 μm	0.750" (1.91 cm)	0.062" (0.16 cm)	0.880" (2.24 cm)	ETFE	141.9 µL	ea.
A-337	20 μm	0.750" (1.91 cm)	0.062" (0.16 cm)	0.880" (2.24 cm)	ETFE	152 μL	ea.
A-343	2 µm	0.625" (1.59 cm)	0.062" (0.16 cm)	0.750" (1.91 cm)	PCTFE	112.6 μL	ea.



PEEK Frits

- > Inert, biocompatible, and metal-free
- > Uniform porosity, longer filtration life
- > Sealing rings manufactured from PCTFE


Patented IDEX Health & Science PEEK Frits offer exceptionally uniform porosity. This property ensures longer filtration life and consistent frit-to-frit swept volumes. The PEEK polymer frit discs are biocompatible and inert to most solvents, making them well-suited for bioanalytical applications. PEEK's robust properties make these products suitable for low and high pressure applications.

Disc rings, included on all PEEK frits, are made of PCTFE and are slightly thicker than the frit disc, providing enhanced sealing and excellent chemical resistance. PCTFE surrounded PEEK frits can be used up to 80 °C.

0.5 µm PEEK Frits



_0.062" (0.16 cm) diameter 0.062" (0.16 cm) thick

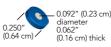

0.195" (0.5 cm) diameter 0.062" (0.64 cm) (0.16 cm) thick

A-701

A-703

A-707

2 µm PEEK Frits



_0.062" (0.16 cm) diameter 0.062" (0.16 cm) thick



A-702

0.200" / (0.51 cm)

A-708

_ 0.091" (0.23 cm) diameter 0.062" (0.16 cm) thick



0.125" (0.32 cm) diameter 0.062" (0.16 cm) thick

A-700



A-706


_0.188" (0.48 cm) diameter 0.062" (0.16 cm) thick

___ 0.062" (0.16 cm) diameter 0.062" (0.16 cm) thick A-704


0.125" (0.32 cm) diameter 0.062" (0.16 cm) thick

A-710

2 µm Semi-Prep PEEK Frits

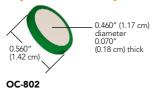

$5 \, \mu m$ and $10 \, \mu m$ PEEK Frits

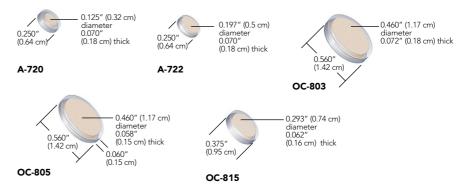

PEEK Frits (Cont.) NOTE The thickness dimension in the part

- The thickness dimension in the part drawings and the pricing tables represents the thickness of the frit disc not the frit ring. Frit rings are often slightly thicker to ensure a proper seal. When tightened into a filter holder the ring compresses to nearly match the thickness of the frit disc.
- The manufacturing process may cause some slight color variance in our PEEK frits. This does not affect their quality or performance. Frit dimensions are approximate. Actual batch-to-batch frit dimensions may vary slightly.



RELATED PRODUCTS


Any 0.247" to 0.254" diameter frit (including polymer ring) can be used with the Standard HPLC Inline Solvent Filters on page 102 and the Standard Precolumn Filters on page 105.



Frit Volume

The term "frit volume" refers to the volume of the various fluid pathways that comprise the matrix of a frit. A standard frit is a mass of small particles fused together through a controlled process of compression and heat. Because of their shape, there are gaps between the fused particles. Fluid makes its way through these gaps, creating a pathway from one side of the frit to the other (see the diagram, below, where the white circles represent frit particles, and the black area represents the void between the particles.)

Generally, when the frit particles increase in size, the frit's porosity increases as well. The larger the particles, the larger the gaps between particles. Cumulatively, these gaps comprise what is known as "frit volume." Using gravimetric determination, it has been experimentally shown that the total volume of any given frit may range from 18%–30%, depending upon the porosity of the frit.

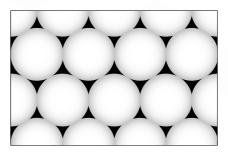

Frit volume is calculated by determining what the mass of the frit would be if it were a solid block of material of equal size. Then the solid mass of the frit is multiplied by the percentage assigned to the porosity to determine the theoretical frit volume.

20% for 0.5 μ m frits 24% for 2 μ m frits

26% for 5 µm frits 28% for 10 µm frits

30% for 20 µm frits

From a chromatographic perspective, it's important to know the volume of the frit used in your system. It is possible for a frit to negatively impact your chromatography if the total frit volume is too large and if it is placed in an area through which the sample will pass. To avoid frit-related problems like band broadening and loss of resolution, most inline filters placed after the sample introduction point (e.g., between the injection valve and the column) are smaller in size and porosity than inline filters that are placed in areas before the sample is introduced into the flow path (e.g., between the pump and the injection valve).

PEEK Frits

Part No.	Porosity	Disc Diameter	Disc Thickness	Ring OD	Ring Material	Frit Volume	Qty.
PEEK FRITS							
A-700	2 µm	0.062" (0.16 cm)	0.062" (0.16 cm)	0.250" (0.64 cm)	PCTFE	0.7 μL	ea.
A-701	0.5 µm	0.062" (0.16 cm)	0.062" (0.16 cm)	0.250" (0.64 cm)	PCTFE	0.6 μL	ea.
A-702	2 µm	0.091" (0.23 cm)	0.062" (0.16 cm)	0.250" (0.64 cm)	PCTFE	1.7 µL	ea.
A-703	0.5 µm	0.092" (0.23 cm)	0.062" (0.16 cm)	0.250" (0.64 cm)	PCTFE	1.4 µL	ea.
A-704	2 µm	0.125" (0.32 cm)	0.062" (0.16 cm)	0.250" (0.64 cm)	PCTFE	3.0 µL	ea.
A-706	2 µm	0.188" (0.48 cm)	0.062" (0.16 cm)	0.250" (0.64 cm)	PCTFE	7.1 µL	ea.
A-707	0.5 µm	0.195" (0.5 cm)	0.062" (0.16 cm)	0.250" (0.64 cm)	PCTFE	6.1 µL	ea.
A-708	2 µm	0.062" (0.16 cm)	0.062" (0.16 cm)	0.200" (0.51 cm)	PCTFE	0.7 μL	ea.
A-710	2 μm	0.125" (0.32 cm)	0.062" (0.16 cm)	0.200" (0.51 cm)	PCTFE	3.0 µL	ea.
SEMI-PREP PI	EEK FRITS						
A-720	10 µm	0.125" (0.32 cm)	0.070" (0.18 cm)	0.250" (0.64 cm)	PCTFE	4.2 µL	ea.
A-722	10 µm	0.197" (0.5 cm)	0.070" (0.18 cm)	0.250" (0.64 cm)	PCTFE	9.9 µL	ea.
OC-802	2 μm	0.460" (1.17 cm)	0.070" (0.18 cm)	0.560" (1.42 cm)	PCTFE	46.4 µL	ea.
OC-803	10 µm	0.460" (1.17 cm)	0.072" (0.18 cm)	0.560" (1.42 cm)	PCTFE	57.2 μL	ea.
OC-805	5 μm	0.460" (1.17 cm)	0.058" (0.15 cm)	0.560" (1.42 cm)	PCTFE	41.1 µL	ea.
OC-815	5 μm	0.293" (0.74 cm)	0.062" (0.16 cm)	0.375" (0.95 cm)	PCTFE	17.8 µL	ea.



Frit-in-a-Ferrule

- > Seals and filters simultaneously
- Less expensive and more convenient than traditional inline filter systems
- Available in both Flangeless and Super Flangeless™ versions

Now you can filter at any point in your system where 1/16" or 1/8" OD tubing is used in a flat-bottom 1/4-28, M6 or 5/16-24 connection.

Our Frit-In-A-Ferrule product line is designed to seal and filter simultaneously by incorporating a frit into the body of a flat-bottom ferrule. This simple design allows you to eliminate traditional inline filters and reduce the number of additional connections in your system.





P-372 Flangeless Frit-In-A-Ferrule for 1/8" OD tubing

P-276Super Flangeless Frit-In-A-Ferrule for 1/16" OD tubing



Part No.	Description	Porosity	Frit Material	Frit Diameter	Frit Thickness	Swept Volume	Maximum Pressure	Qty.
FRIT-IN-	A-FERRULE FOR 1/16" OD TUBING							
P-270	Super Flangeless, Natural PEEK, SST lock ring	2 µm	SST	0.062"	0.062"	0.74 μL	2,500 psi (172 bar)	ea.
P-272	Flangeless, Green PCTFE	2 µm	SST	0.062"	0.062"	0.74 μL	2,000 psi (138 bar)	ea.
P-273	Flangeless, Blue PCTFE	0.5 µm	SST	0.062"	0.062"	0.61 μL	2,000 psi (138 bar)	ea.
P-274	Super Flangeless, Natural PEEK, SST lock ring	2 µm	PEEK	0.046"	0.030"	0.20 μL	2,500 psi (172 bar)	ea.
P-275	Super Flangeless, Black PEEK, SST lock ring	$0.5\mu m$	PEEK	0.046"	0.030"	0.16 μL	2,500 psi (172 bar)	ea.
P-276	Super Flangeless, Stainless Steel, Natural ETFE, SST lock ring	10 µm	SST	0.062"	0.062"	0.90 μL	2,500 psi (172 bar)	ea.
FRIT-IN-	A-FERRULE FOR 1/8" OD TUBING							
P-372	Flangeless, Green PCTFE	2 µm	SST	0.094"	0.062"	1.69 µL	500 psi (34 bar)	ea.
P-373	Flangeless, Blue PCTFE	0.5 µm	SST	0.094"	0.062"	1.41 µL	500 psi (34 bar)	ea.
P-374	Super Flangeless**, Natural PEEK, SST lock ring	2 µm	PEEK	0.094"	0.042"	1.15 μL	2,500 psi (172 bar)	ea.
	olumes include/reflect theoretical frit volume values. " Super Flangeless versions cannot be used in M6 ports.							





Bottom-of-the-Bottle[™] Filters

Our uniquely designed Bottom-of-the-Bottle™ Filters effectively protect your system by filtering out particulate matter that my otherwise damage expensive hardware.

1.70" (4.32 cm) 0.50" (1.27 cm)

Stainless Steel Bottom-of-the-Bottle Solvent Filters

- Draws solvent from within 1/8" of the bottom of the bottle
- Replaceable stainless steel filter cups
- > Versions for 1/8" and 3/16" OD tubing
- Materials of construction: PEEK, ETFE, and 316 Stainless Steel

Patented Stainless Steel Bottom-of-the-Bottle Solvent Filter Assemblies feature a 2 µm or 10 µm replaceable stainless steel filter cup and a design that allows solvent to be drawn from within 1/8" of the bottom of your solvent bottle. The filter cups are inexpensive and easy to replace, making this an economical, trouble-free choice.

Port for solvent intake line Port for your sparging line (optional use)

Bottom-of-the-Bottle Inlet Solvent Filter

Maximum Flow Rates: 2 μm-up to 10 mL/min. 10 μm-up to 40 mL/min




All-PEEK Bottom-of-the-Bottle Solvent Filters

- Most recommended filtering unit
- > 100% PEEK polymer construction
- > Easy operation no fittings required

These biocompatible filters are made from 100% PEEK polymer, including the two built-in PEEK frits. The bottom frit (2 µm or 10 µm) will draw solvents from within 0.080" (2.0 mm) of the bottom of the solvent bottle. The 2 μ m frit on the side may be used for a 1/8" OD helium sparging line.

To use, simply press fit your appropriately sized fluoropolymer tubing firmly into the top holes. That's it!

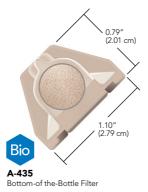

Maximum Flow Rate: up to 30 mL/min




Bottom-of the-Bottle Filter (for small neck bottles)

UHMWPE Bottom-of-the-Bottle[™] Solvent Filters

- > Replaceable filter cup
- Economical
- Materials of construction: UHMWPE, ETFE
- > Versions for 1/16" and 1/8" OD tubing

The design of the UHMWPE solvent filters allows tubing to pass through to the bottom of the filter cup, enabling the filter to draw solvent from within 0.10" (2.5 mm) of the bottom of your solvent bottle.

Please Note: UHMWPE is a hydrophobic material. To establish proper surface wetting, you may need to prime the filter with methanol or acetonitrile.

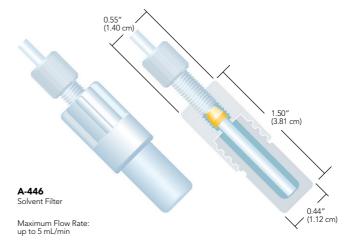

Bottom-of the-Bottle Filters

Part No.	Description	Porosity	For Tubing Size	Includes	Qty.
STAINLESS	STEEL BOTTOM-OF-THE-BOTTLE SOLVENT	FILTERS			
A-550	SST Filter Assembly, with A-520 filter cup	10 μm	1/8" OD	(1) XP-130	ea.
A-551	SST Filter Assembly, with A-522 filter cup	2 µm	1/8" OD	(1) XP-130	ea.
A-520x	SST Replacement Solvent Filter Cups, 10-pk	10 µm	_	_	ea.
A-522x	SST Replacement Solvent Filter Cups, 10-pk	2 μm	_	_	ea.
ALL-PEEK I	BIOCOMPATIBLE BOTTOM-OF-THE-BOTTLE S	SOLVENT FILTERS			
A-435	PEEK Filter	2 µm	1/8" OD	_	ea.
A-437	PEEK Filter, for small-neck (GL-38) bottles	2 µm	1/8" OD	_	ea.
A-438	PEEK Filter, for small-neck (GL-38) bottles	10 µm	1/8" OD	_	ea.
A-440	PEEK Filter	10 µm	1/8" OD	_	ea.
A-441	PEEK Filter	10 μm	3/16" OD	_	ea.
A-451	PEEK Filter	10 μm	1/16" OD	_	ea.
UHMWPE I	BIOCOMPATIBLE BOTTOM-OF-THE-BOTTLE S	SOLVENT FILTERS			
A-445	UHMWPE Filter Assembly	10 μm	1/16" OD	(1) XP-245	ea.
A-446	UHMWPE Filter Assembly	10 µm	1/8" OD	(1) XP-345	ea.
A-427	UHMWPE Replacement Solvent Filter Cups, 5-pk	10 µm	_	_	ea.



Inlet Solvent Filters

- Large surface areas prevent pump cavitation
- Disposable
- 2 μm, 10 μm, and 20 μm pore sizes available
- General use and prep filters for higher flow applications




APPLICATION NOTE

Why Use An Inlet Solvent Filter?

- To filter out particulate matter from the solvent that may otherwise damage expensive hardware. (Use a 10 µm or 20 µm version for this purpose. The A-309 and A-230A filters have an added "Bottom of the Bottle™" feature to help draw solvent to within 1/8" of the bottom of your solvent bottle.)
- To prevent particulates originating from the sparging system from entering the mobile phase reservoir and to help disperse the sparging gas efficiently. (Use a 2 µm filter for this purpose.)
- To hold your tubing in place at the bottom of the bottle. (Most stainless steel filter options work best for this purpose.)

Note: It is usually a good idea to change the inlet filter as part of your semi-annual or annual preventative maintenance program. It is good practice to filter your solvents to prevent pump damage. Our 316 stainless steel filters provide that protection.

Because filters should be changed periodically, we make it easy to replace them without tools. For those filters using a plastic fitting, the tubing can be reconnected by finger tightening the fitting into the new filter. The filters with stems allow easy insertion into the inlet tubing.

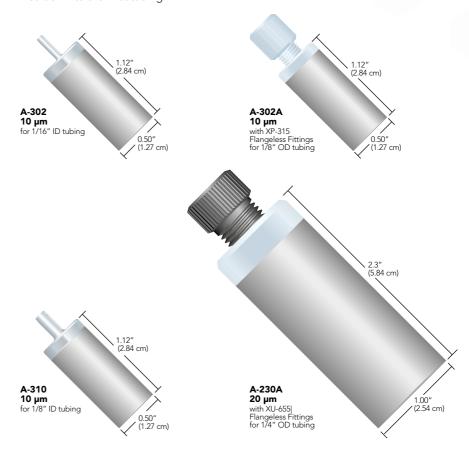

Part No.	Description	Porosity	Material	For Tubing Size	Includes	Max. Suggested Flow Rate*	Qty.
INLET SO	LVENT FILTERS FOR ANALYTICAL HPLO						
A-242	Inlet Solvent Filter with One-Piece Fitting	2 µm	PCTFE, SST	1/8" OD	(1) P-100	10 mL/min	ea.
A-243	A-242, 5-pack	2 µm	PCTFE, SST	1/8" OD	(5) P-100	10 mL/min	ea.
A-228	Inlet Solvent Filter with stem	2 µm	SST	1/8" ID	_	80 mL/min	ea.
A-302	Inlet Solvent Filter with stem	10 µm	SST	1/16" ID	_	40 mL/min	ea.
A-302A	Inlet Solvent Filter with Flangeless Fittings	10 µm	PCTFE, SST	1/8" OD	(1) XP-315	40 mL/min	ea.
A-309	Inlet Solvent Filter with stem	10 µm	SST	1/16" ID	_	40 mL/min	ea.
A-231A	Inlet Solvent Filter with Flangeless Fittings	20 µm	PCTFE, SST	3/16" OD	(1) XP-132	100 mL/min	ea.
A-310	Inlet Solvent Filter with stem	10 µm	SST	1/8" ID	_	40 mL/min	ea.
INLET SO	LVENT FILTERS FOR PREPARATIVE HPL	C SYSTEMS					
A-225	Inlet Solvent Filter with stem	20 µm	SST	1/16" ID	_	100 mL/min	ea.
A-225A	Inlet Solvent Filter with Flangeless Fittings	20 µm	PCTFE, SST	1/8" OD	(1) P-315, (1) P-300N	100 mL/min	ea.
A-227A	Inlet Solvent Filter with Flangeless Fittings	10 µm	PCTFE, SST	1/4" OD	(1) XU-655	100 mL/min	ea.
A-230A	Inlet Solvent Filter with Flangeless Fittings	20 µm	PCTFE, SST	1/4" OD	(1) XU-655	100 mL/min	ea.
A-311	Inlet Solvent Filter with stem	10 µm	SST	1/16" ID	_	100 mL/min	ea.
A-311A	Inlet Solvent Filter with Flangeless Fittings	10 µm	PCTFE, SST	1/8" OD	(1) XP-315	100 mL/min	ea.
+ 1 4	suggested flow rates are determined by perseit		_				

^{*} Maximum suggested flow rates are determined by porosity and surface area.



Inline Filters

- > Specially engineered for inline filtration
- Versions include Micro, Standard, and Semi-Preparative
- Bio-inert and stainless steel options offered
- Variety of porosities, application appropriate




Fittings

All Standard Inline Solvent Filters have 10-32 threads for 1/16" OD tubing, allowing the use of most standard chromatography high pressure fittings. Our Inline Filters are specially engineered for inline filtration. It is specifically designed to help prevent particulate contamination from clogging sensitive equipment. It is ideally suited for placement along the flow path line between the pump and injection valve/autosampler. We offer a variety of porosities for your application.

Standard Inline Solvent Filters

- For 1/16" OD tubing
- Versions for Standard HPLC (6,000 psi/414 bar) and UHPLC (25,000 psi/1,725 bar)
- Replacement frits availableVersions for Standard HPLC (6,000 psi/414 bar) and UHPLC (25,000 psi/1,725 bar)
- > Help prevent particulate contamination from clogging sensitive equipment
- Ideally suited for placement along the flow path line between the pump and injection valve/autosampler

Inline filter assemblies that begin with the letter "A" are engineered for standard HPLC applications (up to 6,000 psi/414 bar). Inline Filter Assemblies that begin with the "VHP" prefix are suitable for use in UHPLC systems, where pressures can reach 25,000 psi (1,725 bar).

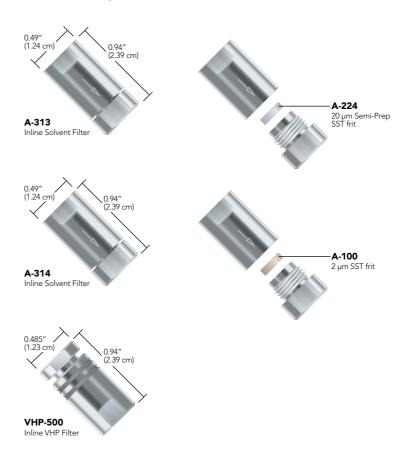


Inline Filters (Cont.)

Biocompatible Standard Inline Filters

- ▶ 0.5 µm and 2 µm versions available
- > Features 100% PEEK flow path

Our A-430 and A-431 Inline Filters consist of a stainless steel body and two PEEK end fittings. Maximum recommended flow rate is 25 mL/min for the A-430 Filter and 10 mL/min for the A-431 Filter. And, you get the added benefit of biocompatibility since all wetted surfaces are PEEK. When you need to replace the frit, simply dispose of the end fitting that contains the frit and replace it with a new one.

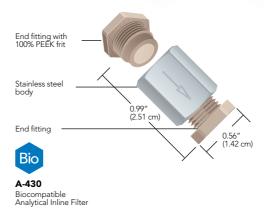

Inline Filters

Part No.	Description	Porosity	For Tubing Size	Threads	Includes	Swept Volume	Pressure Rating	Qty.
STANDA	RD INLINE SOLVENT FILTERS							
A-313	Solvent Filter Assembly	20 µm	1/16" OD	10-32 Coned	(1) A-224	12.3 µL	6,000 psi (414 bar)	ea.
A-314	Solvent Filter Assembly	2 µm	1/16" OD	10-32 Coned	(1) A-100	4 μL	6,000 psi (414 bar)	ea.
A-100	Replacement Frits, Stainless Steel, ea.	2 µm	N/A	_	_	1.4 µL	N/A	ea.
A-224	Replacement Frits, Stainless Steel, ea.	20 µm	N/A	_	_	9.7 µL	N/A	ea.
VHP-500	Inline VHP Filter	0.5 µm	1/16" OD	10-32 Coned	(5) VHP-501	1.2 µL	25,000 psi (1,725 bar)	ea.
VHP-505	Inline VHP Filter	0.2 µm	1/16" OD	10-32 Coned	(5) VHP-506	1.1 µL	25,000 psi (1,725 bar)	ea.
VHP-501	Replacement Inline VHP Frit	0.5 µm	N/A	N/A	N/A	0.60 μL	N/A	ea.
VHP-506	Replacement Inline VHP Frit	0.2 µm	N/A	N/A	N/A	0.54 μL	N/A	ea.
BIOCOM	IPATIBLE INLINE FILTERS							
A-430	Biocompatible Filter Assembly	2 µm		10-32 Coned	(1) A-429	7.1 µL	6,000 psi (414 bar)	ea.
A-431	Biocompatible Filter Assembly	0.5 µm		10-32 Coned	(1) A-428	5.9 µL	6,000 psi (414 bar)	ea.
A-428x	PEEK Filter End Fittings, Black PEEK body, 10-pk	0.5 µm		10-32 Coned	_	5.7 µL	N/A	10-pk
A-429x	PEEK Filter End Fittings, Natural PEEK body, 10-pk	2 μm		10-32 Coned	_	6.9 µL	N/A	10-pk
*Swent vo	lumes include/reflect theoretical frit volume values							



Precolumn Filters


Our economical Precolumn Filters offer secure protection for analytical columns in HPLC and UHPLC. We offer traditional versions that can successfully connect tubing on both sides and our direct-connect versions attach to the inlet port of most standard columns. All versions feature a 10-32 coned ports for 1/16" OD tubing.

Standard Precolumn Filters

- > Economical protection for larger columns and injections
- Traditional versions connect tubing on both sides
- Direct-connect versions attach to the inlet port of most standard columns
- All versions feature 10-32 coned ports for 1/16" OD tubing

These are designed to protect columns by filtering out particulate matter originating from the sample or from rotor seal wear.

- Assemblies that begin with the letter "A" are traditional versions for standard HPLC
- Assemblies that begin with "VHP" are direct-connect versions for UHPLC applications
- Versions that begin with "9085" are direct-connect for standard HPLC and must be used with polymer fittings

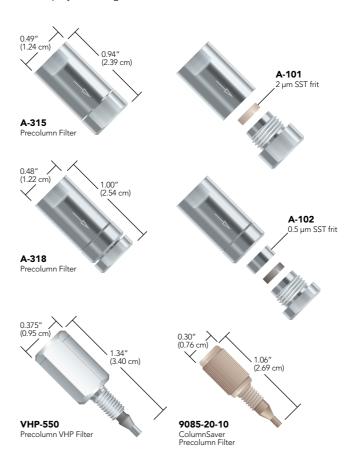


Precolumn Filters (Cont.)

Biocompatible Precolumn Filters

- > Pre-assembled with either 0.5 μm or 2 μm porosity frits
- Great column protection
- > Feature PEEK bodies and PCTFE-surrounded PEEK frits

Biocompatible Precolumn Filters have 0.020" (0.50 mm) diameter thru-holes and 8° distribution cones for minimal band spreading and mixing. The bodies of these filters are manufactured from biocompatible PEEK polymer and are pressure rated to 5,000 psi (345 bar). These filters are designed for use with 1/16" OD tubing, which can be connected to these filters using standard Fingertight fittings.

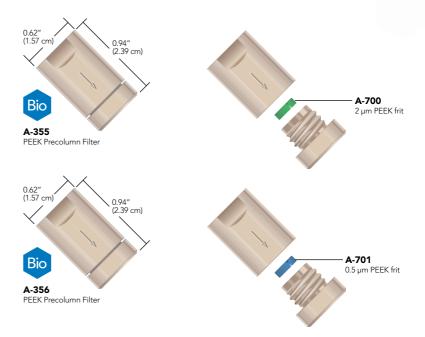

Precolumn Filters

Part No.	Description	Porosity	For Tubing Size	Threads	Includes	Swept Volume*	Pressure Rating	Qty.
STANDAR	RD PRECOLUMN FILTERS							
A-315	Solvent Filter Assembly	2 µm	1/16" OD	10-32 Coned	(1) A-101	1.4 µL	6,000 psi (414 bar)	ea.
A-318	Solvent Filter Assembly	0.5 µm	1/16" OD	10-32 Coned	(1) A-102	0.84 μL	6,000 psi (414 bar)	ea.
A-101	Replacement Frits, Stainless Steel, ea.	2 µm	N/A	_	_	0.74 µL	N/A	ea.
A-102	Replacement Frits, Stainless Steel, ea.	0.5 µm	N/A	_	_	0.61 μL	N/A	ea.
VHP-550	Precolumn VHP Filter	0.5 µm	1/16" OD	10-32 Coned	(5) VHP-551	1.9 µL	20,000 psi (1,380 bar)	ea.
VHP-555	Precolumn VHP Filter	0.2 µm	1/16" OD	10-32 Coned	(5) VHP-556	1.8 µL	20,000 psi (1,380 bar)	ea.
VHP-551	Replacement Precolumn VHP Frit Assembly	0.5 µm	N/A	N/A	N/A	1.9 µL	N/A	ea.
VHP-556	Replacement Precolumn VHP Frit Assembly	0.2 µm	N/A	N/A	N/A	1.8 µL	N/A	ea.
9085-05-10	ColumnSaver Precolumn Filter, with SST frit	0.5 µm	1/16" OD	10-32 Coned	N/A	3.1 µL	6,000 psi (414 bar)	10-pk
9085-20-10	ColumnSaver Precolumn Filter, with SST frit	2 µm	1/16" OD	10-32 Coned	N/A	3.1 µL	6,000 psi (414 bar)	10-pk
BIOCOME	PATIBLE PRECOLUMN FILTERS							
A-355	Solvent Filter Assembly, Biocompatible	2 µm		10-32 Coned	(1) A-700	1.4 μL	5,000 psi (345 bar)	ea.
A-356	Solvent Filter Assembly, Biocompatible	0.5 µm		10-32 Coned	(1) A-701	1.3 µL	5,000 psi (345 bar)	ea.
A-700	Replacement Frit, PEEK Polymer	2 µm		_	_	0.74 µL	N/A	ea.
A-701	Replacement Frit, PEEK Polymer	0.5 µm		_	_	0.61 µL	N/A	ea.

SST = Stainless Stee

^{*}Swept volumes include/reflect theoretical frit volume values.



Semi-Prep Filters


Biocompatible Semi-Prep Filters consist of a stainless steel body, two PEEK end fittings, and a separate PEEK frit. These filters are ideal for many higher flow analytical, semi-prep and preparative applications. Best of all, if the filter becomes clogged, simply unscrew the assembly, remove the frit and replace it. The frits are interchangeable.

Semi-Prep Inline Filters

- Designed for high-flow applications
- > Economical protection for larger columns and injections
- > SFC and HPLC compatible





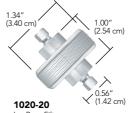

Iso-Prep Filter Shown with standard 10-32 stainless steel nuts and ferrules (not included)



Biocompatible Semi-Prep Inline Filters

- Versions for 1/16", 1/8", 3/16", 1/4", and 5/16" OD tubing
- ▶ 100% PEEK flow path

Biocompatible Semi-Prep Filters consist of a stainless steel body, two PEEK end fittings, and a separate PEEK frit. These filters are ideal for many higher flow analytical, semi-prep and preparative applications. Best of all, if the filter becomes clogged, simply unscrew the assembly, remove the frit and replace it. The frits are interchangeable.

Part No.	Description	Porosity	Threads	Includes	Swept Volume*	Pressure Rating	Qty.
SEMI-PR	REP INLINE FILTERS						
A-330	Semi-Prep Filter Assembly	10 µm	10-32 Coned	(1) A-331	223 μL	7,500 psi (517 bar)	ea.
A-360	Semi-Prep Filter Assembly	10 µm	5/16-24 Flat Bottom	(1) A-331	235 μL	3,500 psi (207 bar)	ea.
A-331	Stainless Steel Frits, Natural ETFE ring	10 µm	N/A	N/A	142 µL	N/A	ea.
A-332	Stainless Steel Frits, Natural ETFE ring	2 µm	N/A	N/A	122 µL	N/A	ea.
A-337	Stainless Steel Frits, Natural ETFE ring	20 µm	N/A	N/A	152 µL	N/A	ea.
ISO-PRE	P FILTERS						
1020-05	21.2 mm Filter Holder	0.5 µm	10-32 Coned	(1) 7031-05	203 uL	8,000 psi (552 bar)	ea.
1020-20	21.2 mm Filter Holder	2 µm	10-32 Coned	(1) 7031-20	196 uL	8,000 psi (552 bar)	ea.
7031-05	21.2 mm Replacement Filter	0.5 µm	N/A	N/A	122 uL	8,000 psi (552 bar)	ea.
7031-20	21.2 mm Replacement Filter	2 µm	N/A	N/A	115 uL	8,000 psi (552 bar)	ea.
BIOCON	MPATIBLE SEMI-PREP INLINE FILTERS						
A-410	Biocompatible Filter Assembly	2 µm	10-32 Coned	(1) OC-802	89 µL	6,000 psi (414 bar)	ea.
A-411	Biocompatible Filter Assembly	10 µm	10-32 Coned	(1) OC-803	103 μL	6,000 psi (414 bar)	ea.
A-510	Biocompatible Filter Assembly	5 µm	5/16-24 Flat Bottom	(1) OC-805	89 µL	500 psi (34 bar)	ea.
OC-802	PEEK Frit, Green PCTFE ring	2 µm	N/A	N/A	46 µL	N/A	ea.
OC-803	PEEK Frit, Natural PCTFE ring	10 µm	N/A	N/A	57 μL	N/A	ea.
OC-805	PEEK Frit, Natural PCTFE ring	5 µm	N/A	N/A	50 μL	N/A	ea.
*Swept vo	lumes include/reflect theoretical frit volume values.						

^{*}Swept volumes include/reflect theoretical frit volume values.




Inline MicroFilters

- ▶ 100% biocompatible PEEK polymer option available
- Miniscule 240 nL void volume
- Two versions: direct connect 1/32" OD tubing or use MicroTight® tubing sleeves for 70-520 µm OD capillary tubing

Our Inline MicroFilters protect your column from particles originating in the mobile phase or sample, or from pump seal and sample injection valve wear. These filters have a 0.006" (150 μ m) thru-hole. Choose the M-520 with a $0.5~\mu$ m 100% PEEK frit to connect to capillary tubing using the MicroTight tubing sleeves (page 52). You may also directly connect 1/32" OD tubing using the M-525 which contains a 0.5 µm PEEK frit.



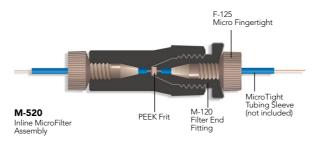



Part No.	Description	Porosity	For Tubing Size	Threads	Includes	Swept Volume	Pressure Rating	Qty.
INLINE I	MICROFILTERS							
M-520	Inline MicroFilter Assembly, PEEK Frit	0.5 µm	MicroTight Tubing Sleeve	MicroTight Tubing Sleeve	(5) M-120, (2) F-125	240 nL	4,000 psi (276 bar)	ea.
M-525	Inline MicroFilter Assembly, PEEK Frit	0.5 µm	1/32" OD	1/32" OD	(5) M-140, (2) F-126	240 nL	4,000 psi (276 bar)	ea.
REPLAC	EMENT INLINE MICROFILTER EN	ID-FITTIN	GS					
M-120x	End-Fittings, Black, with PEEK Frit	0.5 µm	MicroTight Tubing Sleeve	MicroTight Tubing Sleeve	N/A	216 nL	N/A	10-pk
M-140x	End-Fittings, Natural, with PEEK Frit	0.5 µm	1/32" OD	1/32" OD	N/A	216 nL	N/A	10-pk



Mini MicroFilters

- Total volume as low as 10 nL
- Conductive version for CEC and mass spectrometry applications



APPLICATION NOTE

The Mini MicroFilters can be used to pack capillary tubing. Simply place one of these filters on the effluent side of the capillary tubing, then slurry pack. Once packed, place a filter at the head of the tubing. This creates a reliable capillary column without fusing the silica to make frits or pressing filter paper inside the capillary tubing.

Increase the Life of Your Column

Why use a Precolumn Filter when there is a frit at the head of the column itself? Changing the column frit is extremely difficult to do without disturbing the column packing. A Precolumn Filter provides relatively inexpensive insurance against column damage, and changing its frit is easy. A Precolumn Filter placed between the sample injection valve and the HPLC column protects the column from particles originating in the sample and from pump and valve seal wear.

Our Inline Mini MicroFilter Assemblies filter effectively with internal volumes low enough to ensure reliable chromatographic results — even at nanoliter per minute flow rates! Internal volumes of these encapsulated filters are as low as 85 nL with the micro-screen and 10 nL to 22 nL with the frit disc option.





SPECIFICATIONS & DETAILS

Because of the size-specific nature of the ferrules included with each Mini MicroFilter assembly, please note that these ferrules are not interchangeable with other MicroFerrules for different tubing sizes.

Filter Capsule Color Identification











What's the Difference Between Precolumn & Inline Filters?

You may have noticed that the bodies of Precolumn and Inline Filters look similar, and as such, you may have wondered what the differences are. Because Precolumn Filters, by definition, are typically placed in a volume-sensitive area immediately preceding the column, these filters usually feature smaller thru-holes and smaller frit diameters. In contrast, Inline Filters are often placed where the internal volume is not as critical and where longer life and less fluid restriction is more important.

Part No.	Description	Porosity	Frit Type	For use with Tubing	Includes	Swept Volume	Pressure Rating	Qty.
MINI MI	CROFILTER ASSEMBLY							
M-547	Mini MicroFilter Assembly	1 µm	SST Frit	1/32" (790 µm) OD	(5) M-133, (2) F-112, (2) P-416	22 nL	4,000 psi (276 bar)	ea.
M-548	Mini MicroFilter Assembly	1 µm	Ti Frit	1/32" (790 µm) OD	(5) M-134, (2) F-112, (2) P-416	22 nL	4,000 psi (276 bar)	ea.
REPLAC	EMENT MINI MICROFILTE	R CAPSUL	ES					
Part No.	Description	Porosity	Frit Type	For Use With	Material	Swept Volume		Qty.
M-121	Filter Capsule	1 µm	SST Screen	M-530 and M-531	PEEK	85 nL		2-pk
M-125	NanoFilter Capsule	1 µm	SST Frit	M-537 and M-538	PEEK	10 nL		2-pk
M-126	NanoFilter Capsule	1 µm	Ti Frit	M-537 and M-538	PEEK	10 nL		2-pk
M-131	Filter Capsule	1 µm	SST Screen	M-543	PEEK	85 nL		2-pk
M-133	NanoFilter Capsule	1 µm	SST Frit	M-547 and M-548	PEEK	10 nL		2-pk
M-134	NanoFilter Capsule	1 µm	Ti Frit	M-547 and M-548	PEEK	10 nL		2-pk
M-128	Conductive NanoFilter Capsule	1 µm	SST Frit	M-534	SST/PEEK	10 nL		2-pk


SST = Stainless Steel; Ti = Titanium



Precolumn MicroFilters

- Direct connects to columns with 10-32 threads
- > Total void volume of 0.5 μL
- > Two versions: direct connect 1/16" OD tubing or use MicroTight® tubing sleeves for 70-520 µm OD capillary tubing

The Precolumn MicroFilters directly connect into your microbore or analytical column. Total theoretical void volume is only $0.5\,\mu L$ (includes frit volume) and the PEEK tubing used in the assembly of these units has a 0.005" (125 μ m) ID, virtually eliminating any mixing of the sample with the mobile phase.

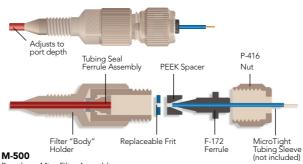


Precolumn MicroFilter Assembly (Includes indicated products)

Part No.	Description	Porosity	For Tubing Size	Threads	Includes	Swept Volume*	Pressure Rating	Qty.
PRECOL	LUMN MICROFILTER ASSEMBLIES							
M-500	Precolumn MicroFilter Assembly, SST Frit	0.5 µm	MicroTight Tubing Sleeve	10-32 Coned	(5) C-425, (1) F-172, (1) P-416	0.5 μL	4,000 psi (276 bar)	ea.
M-510	Precolumn MicroFilter Assembly, PEEK Frit	0.5 µm	MicroTight Tubing Sleeve	10-32 Coned	(5) A-735, (1) F-172, (1) P-416	0.5 μL	4,000 psi (276 bar)	ea.
M-550	Precolumn MicroFilter Assembly, SST Frit	0.5 µm	1/16" OD	10-32 Coned	(5) C-425, (1) F-132, (1) P-416	0.5 μL	4,000 psi (276 bar)	ea.
M-560	Precolumn MicroFilter Assembly, PEEK Frit	0.5 µm	1/16" OD	10-32 Coned	(5) A-735, (1) F-132, (1) P-416	0.5 μL	4,000 psi (276 bar)	ea.
REPLAC	EMENT PRECOLUMN MICROFILTER	FRITS (I	FRIT DIAMETER X FRIT	THICKNESS	X OVERALL DIAMETER)			
A-735	PEEK Frits, 0.045" x 0.031" x 0.192"	0.5 µm	N/A	N/A	N/A	216 nL	N/A	ea.
C-420	SST Frits, 0.038" x 0.028" x 0.192"	2 µm	N/A	N/A	N/A	101 nL	N/A	ea.
C-425	SST Frits, 0.038" x 0.028" x 0.192"	$0.5~\mu m$	N/A	N/A	N/A	101 nL	N/A	ea.



Bottle Caps

- Extremely simple no threaded ports or fittings
- Manufactured from ETFE and Polypropylene

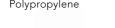



APPLICATION NOTE


- A self-regulating sparging system can help reduce helium consumption and improve pump performance. Set this up by pressing your tubing through the appropriate holes in your bottle cap and attaching each line to a filter. Sparge your mobile phase with an inert gas (preferably helium) for 15-20 minutes. Then reduce the outlet pressure of the sparging gas to a maximum of 5 psi (0.34 bar) and insert a plug (A-626 or A-628) into the remaining port of the cap. The sparging gas will shut off once the incoming pressure equals the pressure inside the reservoir. As the mobile phase is consumed and the internal pressure lowers, sparging gas will enter to keep the system pressurized and degassed. Please Note: If gas leaks while pressurizing the bottle, try removing the sealing ring from the bottle, as it sometimes interferes with the sealing of these bottle caps.
- Done concern with sparging systems is the possibility of solvent backing up the sparging inlet line. This can occur if the gas tank completely evacuates with the regulating valves open, creating a vacuum in the tubing. Solvent backup may damage sparging system components and cause cross-contamination of mobile phase reservoirs. To help prevent solvent backup, install the CV-3010 Inline Check Valve (page 135) along the tubing line that runs between the gas supply and the solvent bottle.
- For a more efficient degassing system, please see the HPLC Vacuum Degassing Systems on page 154.
- Please see the Quick-Stop Luer Check Valve on page 139 for another solvent inlet Application Note.

If you are looking for a bottle cap that is quick and easy to use, but still allows many connect ion options, we have just what you need! The Bottle Caps fit standard GL-45 (1 L) or smaller-neck GL-38 (4 L) glass bottles.

Each cap has three holes. With two of the holes you simply push your tubing straight through. The third hole, with a luer taper, can be used for a number of options. Any male luer (such as a luer-lock syringe) will fit snugly in this hole, or you can use the A-626 or A-627 Plug. Exceptions are the A-610 Bottle Caps. Please see the note below.





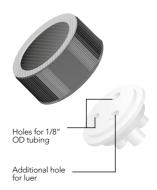

A-622 GL-38 Bottle Cap for 1/8" tubing



NOTE

The A-610 Bottle Cap has a slightly different configuration than other caps. One hole accepts 3/16" OD tubing, the typical size used with some Waters® systems. The remaining two holes accept 1/8" OD tubing. Unlike the other caps, the A-610 does not have a tapered luer hole. If desired, use our A-628 Plug or A-629 Filter Plug for one of the 1/8" holes.



RELATED PRODUCTS

To ensure a tight seal, use fluoropolymer tubing with these bottle caps (page 55).

Part No.	Description	Qty.
BOTTLE CAPS	FOR GL-45, 1 L BOTTLES	
A-610	for 3/16" OD tubing, Red	ea.
A-620	for 1/8" OD tubing, Red	ea.
A-630	for 1/16" OD tubing, Red	ea.
BOTTLE CAPS	FOR GL-38, 4 L BOTTLES	
A-622	for 1/8" OD tubing. Black	ea.



Bottle Cap Plugs

Use the A-626 Bottle Cap Plug to seal the third "tapered" luer hole found in most IDEX Health & Science Bottle Caps. Or, use the A-628 Plug to seal any unused 1/16" or 1/8" bottle cap holes.

Alternatively, try the A-627 or A-629 Filter Bottle Cap Plug to cap an unused hole in your bottle cap. The 20 µm stainless steel frit in these products prevents foreign matter from contaminating your solvent while leaving the bottle open to the atmosphere, thus allowing fluid to be pulled out without creating a vacuum (generally not used with sparging applications). All plug bodies are manufactured from ultra-high molecular weight polyethylene (UHMWPE).




20 µm stainless steel frit

A-626 Bottle Cap Plug

A-629 Filter Bottle Cap Plug

Part No.	Description	Qty.
BOTTLE CAP PI	LUGS	
A-626	Bottle Cap Plug for luer hole, UHMWPE	ea.
A-627	Filter Bottle Cap Plug for luer hole, UHMWPE with 20 µm stainless steel frit	ea.
A-628	Bottle Cap Plug for 1/16", 1/8" or 3/16" hole, UHMWPE	ea.
A-629	Filter Bottle Cap Plug for 1/16", 1/8" or 3/16" hole, UHMWPE with 20 μm stainless steel frit	ea.

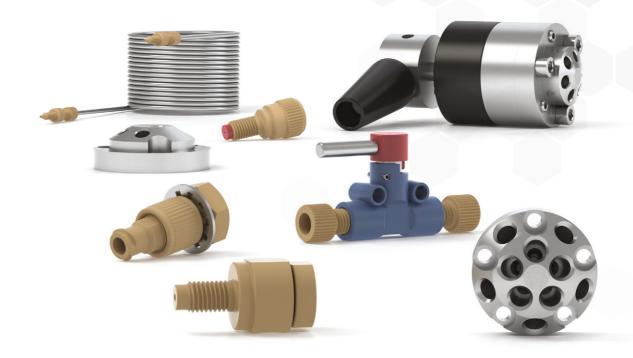

Valves

Our valves are an integral part of advanced fluid-handling solutions for a wide range of analytical instrumentation and clinical diagnostic systems. Our valve options include manual valves for lower frequency use and rotary shear valves that meet the high duty cycle requirements of UHPLC and also come in high and low pressure versions to meet your system requirements. We also offer check valves when there is a need to limit the fluid flow to one direction. Our Back Pressure Regulators products are designed to enhance system performance through outgassing prevention. All of our valve products, components, tools, and accessories are designed keeping our customer's system needs first.

115 VALVE OVERVIEW & FUNCTIONS

115 ROTARY SHEAR VALVES

119 STAND ALONE VALVE PRODUCTS

121 MANUAL VALVES

124 SPARE PARTS

126 VALVE ACCESSORIES

134 FLOW REGULATING VALVES

143 BACK PRESSURE REGULATORS

VALVE OVERVIEW & FUNCTIONS

Valve Module	Flow Configurations	Page
ACTUATED VALVES		
UP TO 15,000 PSI		119
Switching	2-Position, 6-Port 2-Position, 10-Port	
Injection	For Injection, add the appropriately sized Sample Loop to the Switching valves above	
Selection	6-Position, 7-Port	
UP TO 6,000 PSI		119
Switching	2-Position, 6-Port (Analytical and Nano Scale) 2-Position, 10-Port (Analytical and Nano Scale)	
Injection	For Injection, add the appropriately sized Sample Loop to the Switching valves above 2-Position, 6-Port (vertical port)	
Selection	6-Position, 7-Port	
UP TO 125 PSI		119
Switching	2-Position, 6-Port 2-Position, 6-Port (Double 3-Way)	
Selection	6-Position, 7-Port 10-Position, 11-Port	

Valve Module	Flow Configurations	Page
MANUAL VALVES		
UP TO 9,000 PSI		123
Injection	2-Position, 6-Port (Front-Loading, 9,000 psi)	
UP TO 6,000 PSI		123
Switching	2-Position, 6-Port (Analytical and Micro Scale)	
Injection	2-Position, 6-Port	
Selection	6-Position, 7-Port	
UP TO 1,000 PSI		123
Switching	2-Way, Right Angle 4-Position, 4-Port 3-Way, T-Shape 4-Position, 4-Port 4-Way, Diagonal Flow 4-Position, 4 Port	
Injection	2-Position, 6-Port	
Selection	6-Position, 7-Port	

Rotary Shear Valves

Our Rotary Shear Valves were developed in tandem with the evolution of liquid chromatography, where combinations of elevated system pressures, aggressive chemicals, and ever-diminishing fluid volumes continually challenged system manufacturers who required highly precise fluid control and delivery. Today, many other disciplines utilize Rotary Shear Valves for their versatility, reliability, repeatability, long system uptime, and easy preventive maintenance.



Valve Overview and Functions (Cont.)

Choosing a Rotary Shear Valve

Evaluating some simple variables will assist you in choosing the best valve for your needs.

Identify the Operating Pressure of Your Instrument or Application

Valves are designed to repeatedly deliver specific fluids to different locations in a fluidic circuit. Achieving fluidic precision at 15,000 psi requires different valve-design features than those required to achieve fluidic precision at 100 psi. A wide variety of variables such as valve architectures, metals, polymers, coatings, actuation speeds, and manufacturing techniques have been tested to achieve the fluidic accuracy and precision required for the full array of pressure conditions in life science applications. In this catalog, we define four separate pressure groupings:

Up to 15,000 psi (1,035 bar)	UHPLC/Fast Chromatography
Up to 6,000 psi (410 bar)	HPLC
Up to 1,000 psi (69 bar)	Medium Pressure Applications
Up to 125 psi (8.5 bar)	Low Pressure/ Atmospheric Pressure

Identify the Range of Flow Rates in Your System

Because Rotary Shear Valves have been used most often in chromatography systems, certain flow rate ranges have evolved functionally. However, these ranges can apply to any system, not just chromatography:

- Micro/Nano Scale flow rates less than 100 µL per minute
- Analytical Scale flow rates from $100 \, \mu L$ to $10 \, mL$ per minute
- > Prep (or Semi-Prep) Scale flow rates greater than 10 mL per minute

Decide What You Want the Valve to Do

In this chapter Rotary Shear Valves perform three functions:

- > Switching one or more flow paths to a different destination under pressure
- > Injection into a flowing stream under pressure
- > Selection/distribution of a variety of system liquids by means of a common port
- > Read more about valve functions on page 114.

Identify Whether You Want Automated or Manual Control

An automated valve offers more sophisticated functionality. Choose an automated valve if the application requires fast, consistent flow-stream switching. Some other advantages of automated valves include control options (PC- or instrument-triggered), higher torque operation, valve-position feedback, or very small flow paths.

Choose a manual valve if your application involves low frequency of use, demands operator control, or involves injection of smaller sample volumes. (See page 122 for more on Single Mode vs. Dual Mode operation.)

Identify the Chemical Compatibility Requirements Related to Your Fluids

Consulting the chemical compatibility chart in the Technical Resources section at the back of The IDEX Health & Science Laboratory Products catalog helps identify what valve materials to use and avoid — in your application. You can also find Chemical Compatibility information at www.idex-hs.com under Materials and Tools.

Identify Fluidic Connection Requirements in Your System

The rotary shear valves in this catalog accommodate one or more of the following tubing outer diameters: 1/8", 1/16", or 1/32".

Effects of Valves & Tubing on Resolution

The effect of tubing on analytical and microscale analyses can be significant. Since dispersion caused by tubing is proportional to the fourth power of diameter, large bore tubing should be avoided when performing analytical scale or microscale analyses. Tubing ID size ≤ 0.25 mm (0.010") is recommended.

Consider a system with injection and column switching valves and analytical columns with small-bore connecting tubing. The chromatograms below, made using a typical analytical chromatograph, show these effects. Scheme A is the control (injection valve → column → detector) with no valve in the system. In Schemes B and C, two model 7060 Six-Position Switching Valves were placed side by side (injection valve → valve #1 → column → valve #2 → detector).


The injection valve and detector were connected to these valves by the same tubing used in the control. The extra tubing pieces required to connect the valves to the column were a 10 cm length for valve #1-to-column, and a 35 cm length for column-to-valve #2. The diameters of these tubes are indicated in the experimental details, to the right.

Comparison of Observed Column Plates of Analytical and MicroScale Injection Valves

	7725	8125	Δ	
k' = 0.6	2930	5054	72%	
k' = 1.5	4653	6904	48%	
k' = 7.9	7875	8305	5.0%	

UV detector: 1 μ L volume, 4 mm path. Sample volume: 2 μ L, partial-filling method. Column: 2 mm 10×100 mm long, 4 μ m C-18. True plates of column = 11,570.

Effects of Valves and Tubing on Resolution

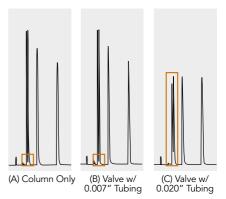

Conclusion: These sequential chromatograms show the effect of adding volume to the flow path through the addition of components.


- (A) Establishes a baseline quality of separation with the minimum volume of liquid in the flow path.


 (B) Adding a valve plus smaller-ID tubing, and thereby increasing the liquid volume only marginally, barely
- affects the separation. However in (C) Adding a valve plus larger-ID tubing, thereby increasing the liquid volume in the flow path to a greater degree, distinctly impairs the quality of the separation and the detectable sample

WHAT IS MAKE-BEFORE-BREAK™, AND WHEN DOES IT MATTER?


Make-Before-Break is a unique design feature of certain dual-mode manual injection valves.

Flow paths of model 7725(i) and 9725(i) with MBB design

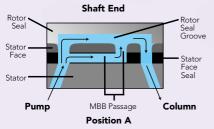

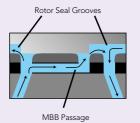


Position B

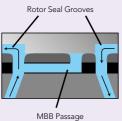


Position C


To maintain a constant, desired highpressure flow, our Make-Before-Break (MBB®) design creates continuous flow between the LOAD and INJECT positions that virtually eliminates pressure transient shock to the system. A passage in the stator face makes a new connection before old connections break. The MBB design an improvement over bypass-style injectors — does not dilute the sample and is easy to maintain and troubleshoot.

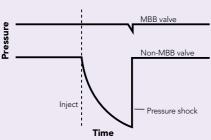

Switching Valves

Switching valves dynamically alternate between two fluid paths without manually disconnecting plumbing. In Chromatography, these valves can be used for column switching, backflushing, sample enrichment, and other techniques. In Diagnostic or Sequencing applications, the switching valve may alternate flow paths to enable back flushing or other fluidic tasks within the instrument.

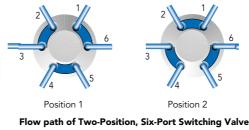

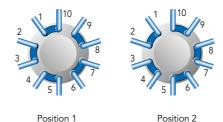


Flow path of Two-Position, Ten-Port Switching Valve

Our switching valves operate between two positions, and may have 6 or 10 ports on the face of the stator (2/6 or 2/10). The flow paths connect ports around the circumference of the stator. The manual switching valves (to 1,000 psi) described on page 123 have different flow path geometry as noted.

Injection Valves

Our injection valves are a form of switching valve. Injection valves can be automated or manual, and they are generally utilized in the two-position, six-port (2/6) configuration and have a sample loop attached.

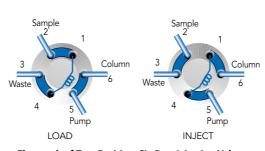

Flow path of Two-Position, Six-Port Injection Valve

The purpose of an Injection valve is to introduce a sample into a flowing stream of liquid. Some Switching valves become Injection valves by the addition of a Sample Loop (a defined length of tubing and fittings configured to match the angle of the valve ports). Sample is loaded and held in the loop until injection is triggered, either manually or automatically.

Injection valves are classified as either Single or Dual Mode based on how the Sample Loop can be filled. A Single Mode Injection valve requires complete filling of the sample loop and is configured for Rear loading, generally in an auto-sample configuration. A Dual Mode Injection valve allows either partial or complete filling of the loop, and introduces sample by syringe through the needle port built into the valve shaft. Complete filling of the sample loop in both the Dual and Single Mode Injection valves provides greater repeatability injection to injection. (See the Application Note, page 131 for greater detail on partial vs. complete loop filling.)

Valve Overview and Functions (Cont.)

Selection Valves

Selection valves enable discrete connections among multiple system liquids (mobile phase, reagents, buffers) by means of a common port (inlet or outlet) connected to a number of different reciprocal ports. In Diagnostic or Sequencing applications, the selection valve alternates between different reagents or sample streams.


Numerous configurations exist among selection valves (e.g., 6-position 7-port, or 10-position 11-port), but these valves typically operate between more than two positions. The ports are usually spaced radially, or outward in some manner around the center port of the stator.



Flow path of Six-Position, Seven-Port Selector Valve



Six column selection using two selection valves.

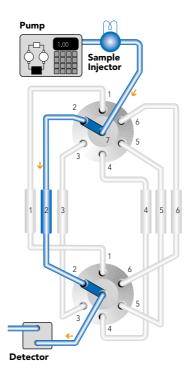



Stand Alone Valve Products

An automated valve offers more sophisticated functionality. Choose an automated valve if the application requires fast, consistent flow-stream switching. Some other advantages of automated valves include control options (PC- or instrument-triggered), higher torque operation, valve-position feedback, or very small flow paths.

MX Series II

MXT to 15,000 psi (1,035 bar)





Add our MX Series II[™] actuated valves to your existing instrument or use in stand-alone lab configurations. MX valves can be controlled remotely or operated manually using the push-button front panel with LED position indicator. MX valves connect to your instrument or PC through contact closure, BCD, serial port, or USB. Commands can be sent to the MX valves using your chromatography software or the included proprietary software for timed-events programmability.

Available flow rates include options for Analytical, Micro/Nano, or Semi-Prep in a range of pressure capabilities. Valve liquid ends are available in materials chosen to be chemically inert and biocompatible. Routine maintenance using authorized RheBuild® kits (page 124) or — for the higher-pressure MXP and MXP valves — the Rapid Replacement Pods™ (page 120) assures optimal performance.

Part No.	Description	Ports, Connections	Wetted Material	Rapid Replacement Pod	Qty.
ACTUATED V	ALVES UP TO 15,000 PSI (1,035 BAR)				
SWITCHING					
MXT715-000	2-Position, 6-Port	10-32 Ports for 1/16" OD Tubing	UltraLife	PD715-000	ea.
MXT715-102	2-Position, 10-Port	10-32 Ports for 1/16" OD Tubing	UltraLife	PD715-102	ea.
INJECTION	For Injection, add the appropriately sized Sar	mple Loop to the Switching valve	s above		
SELECTION					
MXT715-105	6-Position, 7-Port	10-32 Ports for 1/16" OD Tubing	UltraLife	PD715-105	ea.
All of these MXT	valves include a set of 1/16" fittings. Replacement Fitting	gs for MXT valves can be located on pag	e 133.		
ACTUATED V	ALVES UP TO 6,000 PSI (410 BAR)				
SWITCHING					
MXP7900-000	2-Position, 6-Port	10-32 Ports for 1/16" OD Tubing	DuraLife®*	PD7900	ea.
MXP7960-000	2-Position, 10-Port	10-32 Ports for 1/16" OD Tubing	DuraLife	PD7960	ea.
MXP7980-000	2-Position, 6-Port, Nano, 5,000 psi (345 bar)	M4 Ports for 1/32" OD Tubing	DuraLife II	PD7980	ea.
MXP7986-000	2-Position, 10-Port, Nano, 5,000 psi (345 bar)	M4 Ports for 1/32" OD Tubing	DuraLife II	PD7986	ea.
MXP9900-000	2-Position, 6-Port, Biocompatible, 5,000 psi (345 bar)	10-32 Ports for 1/16" OD Tubing	PEEK	PD9900	ea.
MXP9960-000	2-Position, 10-Port, Biocompatible, 5,000 psi (345 bar)	10-32 Ports for 1/16" OD Tubing	PEEK	PD9960	ea.
INJECTION	For Injection, add the appropriately sized Sar	mple Loop to the Switching valve	s above		
MXP7920-000	2-Position, 6-Port, Vertical Port	10-32 Ports for 1/16" OD Tubing	DuraLife	PD7920	ea.
SELECTION					
MXP7970-000	6-Position, 7-Port	10-32 Ports for 1/16" OD Tubing	DuraLife II**	PD7970	ea.
** DuraLife II is	proprietary material combination of SST and an advance a proprietary material combination consisting of Titaniun valves include a set of 1/16" fittings. Replacement Fitting	n and an advanced polymer.	e 133.		
ACTUATED V	ALVES UP TO 125 PSI (8.5 BAR)				
SWITCHING					
MXX777-601	2-Position, 6-Port	Accepts Either 1/16" or 1/8" Tubing	RPC-7*	1/16" and 1/8"	ea.
MXX777-603	2-Position, Double Three Way	Accepts Either 1/16" or 1/8" Tubing	RPC-7	1/16" and 1/8"	ea.
MXX777-612	2-Position, 6-Port, Large Bore	Accepts Either 1/16" or 1/8" Tubing	RPC-7	1/16" and 1/8"	ea.
SELECTION					
MXX777-605	6-Position, 7-Port	Accepts Either 1/16" or 1/8" Tubing	RPC-7	1/16" and 1/8"	ea.
MXX777-616	6-Position, 7-Port, Large Bore	Accepts Either 1/16" or 1/8" Tubing	RPC-7	1/16" and 1/8"	ea.
MXX778-605	10-Position, 11-Port	Accepts Either 1/16" or 1/8" Tubing	RPC-7	1/16" and 1/8"	ea.

^{*} RPC-7 Proprietary Polymer Combination.

All of these MXX valves include a set of 1/16" and 1/8" ferrules. Replacement Fittings for MXX valves can be located on page 133.



Rapid Replacement Pods[®]

For IDEX Health & Science MX Series II Valves

- Zero downtime maintenance
- Improves lab throughput

To help keep your instrument online and performing at maximum precision, select the exact Rapid Replacement Pod for your higher pressure MX Series II valves. Replacement pods are easily exchanged as part of scheduled preventive maintenance, or in an emergency, a pod can be substituted quickly while the original is examined and maintained at your convenience. The pod kit contains complete instructions for removal and replacement.



Part No.	Description	For Valve Part No.	Qty.
TO 15,000 PSI (1,035	5 BAR)		
SWITCHING			
PD715-000	Rapid Replacement Pod	MXT715-000	ea.
PD715-102	Rapid Replacement Pod	MXT715-102	ea.
SELECTION			
PD715-105	Rapid Replacement Pod	MXT715-105	ea.
UP TO 6,000 PSI (410	DBAR)		
SWITCHING			
PD7900	Rapid Replacement Pod	MXP7900-000	ea.
PD7960	Rapid Replacement Pod	MXP7960-000	ea.
PD7980	Rapid Replacement Pod	MXP7980-000	ea.
PD7986	Rapid Replacement Pod	MXP7986-000	ea.
PD9900	Rapid Replacement Pod	MXP9900-000	ea.
PD9960	Rapid Replacement Pod	MXP9960-000	ea.
INJECTION			
PD7920	Rapid Replacement Pod	MXP7920-000	ea.
SELECTION			
PD7970	Rapid Replacement Pod	MXP7970-000	ea.



Manual Valves

Choose a manual valve if your application involves low frequency of use, demands operator control, or involves injection of smaller sample volumes.







7725i

Manual Injection Valve Up to 9,000 psi (600 bar)

7060


Manual Switching Valve Up to 7,000 psi (483 bar)

3725i-038


Manual Switching Valve Up to 7,000 psi (483 bar)




Manual Switching Valve Options

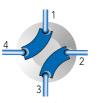

4-Way Diagonal Flow Switching Valves (V-100D, V-101D)

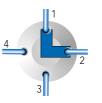


Right Angle Flow Switching Valves (V-100L, V-101L)

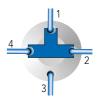

3-Way Flow Switching Valves (V-100T, V-101T)



V-101L Manual Switching Valve Up to 1,000 psi (69 bar)

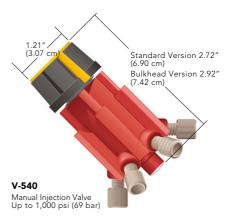



Manual Valves (Cont.)




Switching Valve Applications

Protect sensitive system components (such as a column) during a cleaning cycle with our Diagonal Flow Switching Valve ("D"). This valve eliminates the need to remove, plug and reconnect a low pressure column (see below).

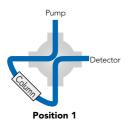


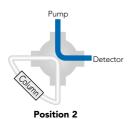

A typical application for a Right Angle Flow Switching Valve ("L") is column switching, allowing two columns to use one detector. Detector switching is another common application for this valve (see below). Plug off the extra port with the included plug.

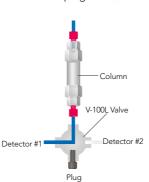


Position 1

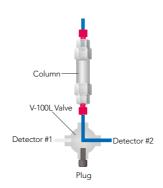

Position 2

Your detector switching application may require the flexibility of routing the column effluent to both detectors simultaneously while retaining the ability to isolate each detector. Use our 3-Way Flow Switching Valve ("T"), plugging off the fourth port with the included plug.

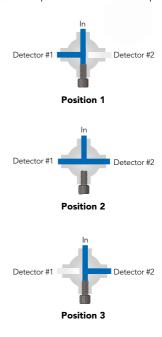

Characteristics of Manual Sample Injection Valves

Scale	Partial Filling Volumes (Range)	Sample Loop Sizes (Range)	Wetted Materials	Max. psi (bar)¹	Max. T (°C)	MBB ²	Model ³
Analytical	1 μL–2.5 mL 1 μL–5.0 mL	2 μL–5.0 mL 2 μL–10 mL	316 SST, Vespel® PEEK, ETFE, ceramic	7,000 (483) 5,000 (340)	80° 50°	Yes Yes	7725, 7725i 9725, 9725i
Micro	0.1 μL–500 μL	5 μL–1.0 mL	316 SST, PEEK, Vespel, ceramic	7,000 (483)	80°	No	8125
Preparative	100 μL–10 mL	2.0 mL-20 mL	316 SST, PEEK PEEK	5,000 (340) 4,000 (276)	50° 50°	Yes Yes	3725(i)-038, 3725i
Analytical	Not Applicable	5 μL–5.0 mL 5 μL–10 mL	316 SST, Vespel PEEK, ETFE, Ceramic	7,000 (483) 5,000 (340)	150° 50°	No No	7010 9010
	Analytical Micro Preparative	Scale Volumes (Range) Analytical 1 μL-2.5 mL 1 μL-5.0 mL Micro 0.1 μL-500 μL Preparative 100 μL-10 mL	Scale Volumes (Range) Sizes (Range) Analytical 1 μL-2.5 mL 1 μL-5.0 mL 2 μL-5.0 mL 2 μL-10 mL Micro 0.1 μL-500 μL 5 μL-1.0 mL Preparative 100 μL-10 mL 2.0 mL-20 mL Analytical Not Applicable 5 μL-5.0 mL	Scale Volumes (Range) Size's (Range) Materials Analytical 1 μL−2.5 mL 1 μL−5.0 mL 2 μL−5.0 mL 2 μL−10 mL 316 SST, Vespel® PEEK, ETFE, ceramic Micro 0.1 μL−500 μL 5 μL−1.0 mL 316 SST, PEEK, Vespel, ceramic Preparative 100 μL−10 mL 2.0 mL−20 mL 316 SST, PEEK PEEK Analytical Not Applicable 5 μL−5.0 mL 316 SST, Vespel	Scale Volumes (Range) Size's (Range) Materials psi (bar)¹ Analytical 1 μL-2.5 mL 1 μL-5.0 mL 2 μL-5.0 mL 2 μL-10 mL 316 SST, Vespel® 7,000 (483) 7,000 (340) Micro 0.1 μL-500 μL 5 μL-1.0 mL 316 SST, PEEK, Vespel, ceramic 7,000 (483) 7,000 (340) Preparative 100 μL-10 mL 2.0 mL-20 mL 316 SST, PEEK PEEK 5,000 (340) 4,000 (276) Application Net Applicable 5 μL-5.0 mL 316 SST, Vespel 7,000 (483)	Scale Volumes (Range) Sizes (Range) Materials psi (bar)¹ T (°C) Analytical 1 μL-2.5 mL 1 μL-5.0 mL 2 μL-5.0 mL 2 μL-10 mL 316 SST, Vespel® PEEK, ETFE, ceramic 7,000 (483) 5,000 (340) 80° 50° Micro 0.1 μL-500 μL 5 μL-1.0 mL 316 SST, PEEK, Vespel, ceramic 7,000 (483) 80° Preparative 100 μL-10 mL 2.0 mL-20 mL 316 SST, PEEK PEEK 5,000 (340) 4,000 (276) 50° 50° Apalytical Not Applicable 5 μL-5.0 mL 316 SST, Vespel 7,000 (483) 150°	Scale Volumes (Range) Sizes (Range) Materials psi (bar)¹ T (°C) MBB² Analytical 1 μL-2.5 mL 1 μL-5.0 mL 2 μL-5.0 mL 2 μL-10 mL 316 SST, Vespel® PEEK, ETFE, ceramic 7,000 (483) 5,000 (340) Vespel, ceramic 80° 7,000 (340) Vespel, ceramic No Preparative 100 μL-10 mL 2.0 mL-20 mL 2.0 mL-20 mL 316 SST, PEEK PEEK 5,000 (340) 4,000 (276) 50° 50° Yes Yes Yes Applytical Not Applicable 5 μL-5.0 mL 5 μL-5.0 mL 316 SST, Vespel 7,000 (483) 7,000 (483) 150° 150° No

This is the maximum pressure to which the valve can be adjusted. Some models are shipped from the factory set for lower pressures.

MBB (Make-Before-Break™) is a design that provides uninterrupted flow when switching between LOAD and INJECT. MBB also greatly reduces transient pressure shocks.

Models with an "i" suffix have a built-in position sensing switch. Models 8125 and 9010 also have a built-in switch.



SPECIFICATIONS & DETAILS

Part No.	Stator Passage Diameter	Factory Set Pressure	Maximum Field Set Pressure	Maximum Temperature (°C)
3000 (PEEK)	1.0 mm (0.040")	3,000 psi (207 bar)	4,000 psi (276 bar)	50°
7000, 7010 (SST)	0.6 mm (0.024")	5,000 psi (340 bar)	7,000 psi (483 bar)	150°
7000L (SST)	1.0 mm (0.040")	3,000 psi (207 bar)	5,000 psi (340 bar)	150°
7030 (SST)	0.6 mm (0.024")	5,000 psi (340 bar)	7,000 psi (483 bar)	150°
7030L (SST)	1.0 mm (0.040")	3,000 psi (207 bar)	5,000 psi (340 bar)	150°
7060 (SST)	0.4 mm (0.016")	5,000 psi (340 bar)	7,000 psi (483 bar)	80°
SST = Stainless Steel				

Manual Valves

Part No.	Description	Tubing/Fitting Size	Wetted Material	Configuration	Qty.
MANUAL VA	ALVES UP TO 9,000 PSI (600 BAR)				
INJECTION					
7725i-188	2-Position, 6-Port, 9,000 psi (600 bar)	10-32 Ports for 1/16" OD Tubing	Stainless Steel, PEEK, Ceramic	Front loading	ea.
MANUAL VA	ALVES UP TO 6,000 PSI (410 BAR)				
SWITCHING	i				
3000	2-Position, 6-Port, Prep Scale	5/16-24 Ports for 1/8" OD Tubing	PEEK	_	ea.
7000	2-Position, 6-Port, Large Bore	10-32 Ports for 1/16" OD Tubing	Stainless Steel & Vespel®	_	ea.
7000L	2-Position, 6-Port, Large Bore	10-32 Ports for 1/16" OD Tubing	Stainless Steel & Vespel	_	ea.
7030	2-Position, 6-Port	10-32 Ports for 1/16" OD Tubing	Stainless Steel & Vespel	Double 3-Way	ea.
7030L	2-Position, 6-Port, Large Bore	10-32 Ports for 1/16" OD Tubing	Stainless Steel & Vespel	Double 3-Way	ea.
INJECTION	*				
7010	2-Position, 6-Port Single Mode	10-32 Ports for 1/16" OD Tubing	Stainless Steel & Vespel	20 μL*	ea.
9010	2-Position, 6-Port Single Mode (Switching, Injection)	10-32 Ports for 1/16" OD Tubing	PEEK, ETFE, Ceramic	20 μL*	ea.
3725-038	2-Position, 6-Port, Prep Scale Dual Mode	5/16-24 Ports for 1/8" Tubing	Stainless Steel & PEEK	10 mL*	ea.
3725i	2-Position, 6-Port, Prep Scale Dual Mode with Switch	5/16-24 Ports for 1/8" Tubing	PEEK	10 mL*	ea.
3725i-038	2-Position, 6-Port, Prep Scale Dual Mode with Switch	5/16-24 Ports for 1/8" Tubing	Stainless Steel & PEEK	10 mL*	ea.
7725	2-Position, 6-Port, Analytical Scale Dual Mode	10-32 Ports for 1/16" OD Tubing	Stainless Steel, Ceramic, Vespel	20 μL*	ea.
7725i	2-Position, 6-Port, Analytical Scale Dual Mode with Switch	10-32 Ports for 1/16" OD Tubing	Stainless Steel, Ceramic, Vespel	20 μL*	ea.
8125**	2-Position, 6-Port, Micro Scale Dual Mode with Switch	10-32 Ports for 0.020" (0.5 mm) or 1/16" Tubing	Stainless Steel, Ceramic, Vespel	5 μL*	ea.
9725	2-Position, 6-Port, Analytical Scale Dual Mode	10-32 Ports for 1/16" OD Tubing	PEEK, ETFE, Ceramic	20 μL*	ea.
9725i	2-Position, 6-Port, Analytical Scale Dual Mode with Switch	10-32 Ports for 1/16" OD Tubing	PEEK, ETFE, Ceramic	20 μL*	ea.
SELECTION					
7060	6-Position, 7-Port	10-32 Ports for 1/16" OD Tubing	Stainless Steel & Vespel	6-Way	ea.
* Cl-::-					

Ships with a sample loop of indicated volume attached to ports 1 and 4.
 The 8125 requires special ferrules for 0.020" (0.5 mm) tubing. 8125-084-0.5 mm ferrule for 8125; 8125-086-0.5 mm ferrule for 8125 — 4-pk.

Part No.	Description	Tubing/Fitting Size	Wetted Material	Configuration	Includes	Qty.
SWITCHING						
V-100D	4-Position, 4-Port, 500 psi (34 bar)	1/4-28 Ports for 1/16" OD Tubing	PEEK, PTFE	Double Diagonal	*	ea.
V-101D	4-Position, 4-Port, Bulkhead, 500 psi (34 bar)	1/4-28 Ports for 1/16" OD Tubing	PEEK, PTFE	Double Diagonal	*	ea.
V-100L	4-Position, 4-Port, 500 psi (34 bar)	1/4-28 Ports for 1/16" OD Tubing	PEEK, PTFE	Right-Angle "L"	**	ea.
V-101L	4-Position, 4-Port, Bulkhead, 500 psi (34 bar)	1/4-28 Ports for 1/16" OD Tubing	PEEK, PTFE	Right-Angle "L"	**	ea.
V-100T	4-Position, 4-Port, 500 psi (34 bar)	1/4-28 Ports for 1/16" OD Tubing	PEEK, PTFE	Single "T"	***	ea.
V-101T	4-Position, 4-Port, Bulkhead, 500 psi (34 bar)	1/4-28 Ports for 1/16" OD Tubing	PEEK, PTFE	Single "T"	***	ea.
INJECTION	For Injection, add the appropriately sized Sample Loop to the Switching valves above					
V-450	2-Position, 6-Port, 1,000 psi (69 bar)	1/4-28 Ports for 1/16" OD Tubing	Polyimide, PTFE	Injection	(6) XP-235	ea.
V-451	2-Position, 6-Port, Bulkhead Version, 1,000 psi (69 bar)	1/4-28 Ports for 1/16" OD Tubing	Polyimide, PTFE	Injection	(6) XP-235	ea.
V-540	2-Position, 6-Port, 1,000 psi (69 bar)	1/4-28 Ports for 1/8" OD Tubing	Polyimide, PTFE	Injection	(6) XP-335	ea.
V-541	2-Position, 6-Port, Bulkhead Version, 1,000 psi (69 bar)	1/4-28 Ports for 1/8" OD Tubing	Polyimide, PTFE	Injection	(6) XP-335	ea.
SELECTION						
V-240	6-Position, 7-Port, 1,000 psi (69 bar)	1/4-28 Ports for 1/16" OD Tubing	Polyimide, PTFE	Multi-port Selection	(6) XP-235	ea.
V-241	6-Position, 7-Port, Bulkhead Version, 1,000 psi (69 bar)	1/4-28 Ports for 1/16" OD Tubing	Polyimide, PTFE	Multi-port Selection	(6) XP-235	ea.
V-340	6-Position, 7-Port, 1,000 psi (69 bar)	1/4-28 Ports for 1/8" OD Tubing	Polyimide, PTFE	Multi-port Selection	(6) XP-335	ea.
V-341	6-Position, 7-Port, Bulkhead Version, 1,000 psi (69 bar)	1/4-28 Ports for 1/8" OD Tubing	Polyimide, PTFE	Multi-port Selection	(6) XP-335	ea.

^{* (4)} P-218BLK, (4) P-240. ** (4) P-218BLK, (4) P-240, (1) P-309. *** (4) P-218BLK, (4) P-240, (1) P-309.

Spare Parts

We offer a full line of genuine spare parts to assist with your use of our valve products. We offer RheBuild® Kits designed for specific valve models. Rotor Seal and Stators are commonly replaceable parts.



Rotor Seals & Stators

The rotor seal is the polymeric disc that makes a high pressure seal against the stator or stator face seal. The seal wears with use and is one of the only parts that may need routine replacement.

Stators are available in 316 stainless steel, PEEK and proprietary materials. Typically, stators need replacement only if the ports or sealing surfaces become damaged. Avoid damage from use of improper injection needles by referring to the "Using Proper Syringe Needles" Application Note on page 130.

Please Note: Rotor seals for MX Series II™ Modules are available in RheBuild® Kits on this page. Stators for MX Series II Modules are available on this page. MX (Series I) Module rotor seals are available in RheBuild Kits on this page.



RheBuild® Kits


RheBuild Kits are available for most valve products. Included in each individualized RheBuild Kit are all parts, tools, and instructions to maintain precision performance of your particular product. RheBuild Kits eliminate individual part ordering.



APPLICATION NOTE

How to Select the Right Rotor Seal

The standard rotor seal in many of our manual valves is made from a Vespel® blend. This polyimide has low wear and high chemical resistance. Vespel tolerates a pH range of 0 to 10. Solutions more basic than pH 10 attack Vespel which damages the rotor seal. If you use any solutions above pH 10, our experts recommend a PEEK blend rotor seal. PEEK offers a high chemical resistance and versatility, and will tolerate the entire pH range from 0 to 14. ETFE blend rotor seals are appropriate for use in applications where PEEK is not generally acceptable, such as when methylene chloride or DMSO in higher concentrations is being used.

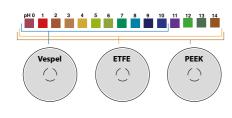



APPLICATION NOTE

How to Avoid Pressure Transients

Air in the sample loop can cause an instantaneous system pressure drop that eventually returns to a normal level. Air causes the pressure to drop when the injector moves from the LOAD to the INJECT position. When large sample loops (\geq 100 μ L) are partially loaded, air present in the needle port tube is pushed into the sample loop (see Figure 1). Air can also enter the sample loop from siphoning which occurs when the vent line is higher than the injection port. In either case, upon injection, the system pressure collapses the air bubble, causing pressure to drop momentarily.

A pressure drop in the system caused by air results in changes in retention time, artifact peaks, and affects column performance.

Avoid pressure drops by removing the air in the needle port tube. Do this by flushing about 1 mL of mobile phase with a luer syringe with needle port cleaner. Keep the needle port tube filled with mobile phase by occasional flushing. Adjust the vent line(s) so the outlet is at the same horizontal level as the needle port (see Figure 2).

For additional injection troubleshooting, refer to our Troubleshooting Guide for HPLC Injection Problems. You may download the Guide from the IDEX Health & Science web site: www.idex-hs.com under Education & Tools.

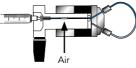


Figure 1 Air present in the needle port tube is pushed by the syringe during loading into the sample loop

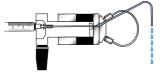


Figure 2 Pathway of the flushing mobile phase using the Needle Port Cleaner, Part # 7125-054 (see page 131) when the injector is in INJECT

VespPEL BIEND ROTOR SEALS Vespel Rotor Seal ea.	Don't No.	Fan Value Madal Na	Description	٠
TODO-016 TODOL, TO40L Vespel Rotor Seal ea. TO10-039 TO10, 7000, 7040 Vespel Rotor Seal ea. TO30-014 7030L Vespel Rotor Seal ea. TO30-014 7030L Vespel Rotor Seal ea. TO60-070 TO60, 7066 Vespel Rotor Seal ea. TO60-070 TO60, 7066 Vespel Rotor Seal ea. TO60-064 TO60L Vespel Rotor Seal ea. TO60-071 TO10, 7010-087, 7000, 7040 ETFE Rotor Seal ea. TO60-071 TO10, 7010-087, 7000, 7040 ETFE Rotor Seal ea. TO60-074 TO30, 9030 ETFE Rotor Seal ea. TO60-074 TO60, 7066, 9020 ETFE Rotor Seal ea. TO60-075 TO60, 7066, 9020 ETFE Rotor Seal ea. TO60-076 TO60, 7066, 9020 PEEK Rotor Seal ea. TO60-076 TO60, 7066, 9020 PEEK Rotor Seal ea. TO60-076 TO60, 7066, 9020 PEEK Rotor Seal ea. TO60-076 TO60, 7072, 7072, 9125, 9725 PEEK Rotor Seal ea. TO60-076 TO60, 7072, 90	Part No.	For Valve Model No.	Description	Qty.
7010-039			Vesnel Rotor Seal	ea
7030-003 7030, 9030 Vespel Rotor Seal ea. 7030-014 7030L Vespel Rotor Seal ea. 7060-070 7060, 7066 Vespel Rotor Seal ea. 7060-064 7060L Vespel Rotor Seal ea. 7125-047 7125, 7725, 7725 Vespel Rotor Seal ea. 7125-047 7125, 7725, 7725 Vespel Rotor Seal ea. 7125-038 7410 Vespel Rotor Seal ea. 7413-013 7413 Vespel Rotor Seal ea. 8125-038 8125 Vespel Rotor Seal ea. 8125-039 7100, 7010, 7010-087, 7000, 7040 ETFE Rotor Seal ea. 8125-079 7125, 7125, 7125, 7125 ETFE Rotor Seal ea. 8125-079 7125, 7125, 7125 ETFE Rotor Seal ea. 910-051 9010 ETFE Rotor Seal ea. 91125-082 9125, 9725 ETFE Rotor Seal ea. 9125-083 3725, 3725-038 PEEK Rotor Seal ea. 9125-081 3725, 3725-038 PEEK Rotor Seal ea. 9125-095 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-096 7000, 7010, 9010 PEEK Rotor Seal ea. 9125-097 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-098 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-098 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-098 MXT15-105 Stator ea. 9123-568 MXT15-100 Stator ea. 9123-568 MXT15-100 Stator ea. 9723-568 MXT715-100 Stator ea. 9724-004 MXP7980-000 Stator ea. 9726-004 MXP7980-000 Stator ea. 9726-004 MXP7980-000 Stator ea. 9726-004 MXP7980-000 Stator ea. 9726-004 MXP7990-000 Stator ea. 9726-004 MXP7990-000 Stator ea. 9726-004 MXP790-000 Stator ea. 9726-004 MXP790-			· ·	
7030-014 7030L Vespel Rotor Seal ea. 7060-070 7060, 7066 Vespel Rotor Seal ea. 7060-064 7060L Vespel Rotor Seal ea. 7125-047 7125, 7725, 9725 Vespel Rotor Seal ea. 7410-038 7410 Vespel Rotor Seal ea. 8125-038 8125 Vespel Rotor Seal ea. 8125-038 8125 Vespel Rotor Seal ea. 8125-038 8125 Vespel Rotor Seal ea. 8125-038 700-071 7010-071 PO10-087, 7000, 7040 ETFE Rotor Seal ea. 8125-097 7125, 7125-081, 7725 ETFE Rotor Seal ea. 8125-097 8125 ETFE Rotor Seal ea. 910-051 9010 ETFE Rotor Seal ea. 9125-092 9125, 9725 ETFE Rotor Seal ea. 9125-093 3725-013 3725-038 PEEK Rotor Seal ea. 9125-095 3725, 3725-038 PEEK Rotor Seal ea. 9125-096 7125, 7725, 9125, 9725				
7060-070 7060, 7060 Vespel Rotor Seal ea. 7060-064 7060L Vespel Rotor Seal ea. 7125-047 7125, 7725, 9725 Vespel Rotor Seal ea. 7410-038 7410 Vespel Rotor Seal ea. 8125-038 8125 Vespel Rotor Seal ea. 8125-038 8125 Vespel Rotor Seal ea. 8125-079 7010, 7010-087, 7000, 7040 ETFE Rotor Seal ea. 7030-015 7030, 9030 ETFE Rotor Seal ea. 8125-079 7125, 7125-081, 7725 ETFE Rotor Seal ea. 8125-079 8125 ETFE Rotor Seal ea. 8125-079 8125 ETFE Rotor Seal ea. 9125-082 9125, 9725 ETFE Rotor Seal ea. 8125-083 3725, 3725-038 PEEK Rotor Seal ea. 8125-19 8125 PEEK Rotor Seal ea. 8125-19 8125 PEEK Rotor Seal ea. 8125-19 8125 PEEK Rotor Seal ea. <td< th=""><th>7030-014</th><th>·</th><th>· ·</th><th></th></td<>	7030-014	·	· ·	
7060-064 7060L Vespel Rotor Seal ea. 7125-047 7125, 7725, 7725 Vespel Rotor Seal ea. 7410-038 7410 Vespel Rotor Seal ea. 7413-013 7413 Vespel Rotor Seal ea. 8125-038 8125 Vespel Rotor Seal ea. 8125-038 8125 Vespel Rotor Seal ea. 7030-011 7010, 7010-087, 7000, 7040 ETFE Rotor Seal ea. 7050-074 7050, 7066, 9060 ETFE Rotor Seal ea. 8125-079 8125 ETFE Rotor Seal ea. 910-051 9010 ETFE Rotor Seal ea. 9125-082 9125, 9725 ETFE Rotor Seal ea. 9125-082 9125, 9725 ETFE Rotor Seal ea. 9125-097 8125 PEEK Rotor Seal ea. 9125-098 9125, 9725 PEEK Rotor Seal ea. 9125-099 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-099 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. </th <th>7060-070</th> <th>7060, 7066</th> <th></th> <th>ea.</th>	7060-070	7060, 7066		ea.
7410-038 7410 Vespel Rotor Seal ea. 7413-013 7413 Vespel Rotor Seal ea. 8125-038 8125 Vespel Rotor Seal ea. 8125-079 7100, 7010-087, 7000, 7040 ETFE Rotor Seal ea. 7060-074 7060, 7066, 9060 ETFE Rotor Seal ea. 7125-079 7125, 7125-081, 7725 ETFE Rotor Seal ea. 8125-097 8125 ETFE Rotor Seal ea. 8125-097 8125 ETFE Rotor Seal ea. 9125-082 9125, 9725 ETFE Rotor Seal ea. 9125-082 9125, 9725 ETFE Rotor Seal ea. 9125-082 9125, 9725 ETFE Rotor Seal ea. 910-065 7000, 7010, 9010 PEEK Rotor Seal ea. 9125-119 8125 PEEK Rotor Seal ea. 9125-119 8125 PEEK Rotor Seal ea. 9125-097 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-119 8125 PEEK Rotor Seal ea. 9125-097 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-098 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-098 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-098 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-098 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-098 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-098 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-098 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-098 9125-090 Stator ea. 9125-090 MX7715-105 Stator ea. 9125-090 MX7715-105 Stator ea. 9125-090 PREVSO-010 Stator ea. 9125-090 PREV	7060-064	7060L		ea.
7413-013 7413 Vespel Rotor Seal ea. 8125-038 8125 Vespel Rotor Seal ea. EFFE BLEND ROTOR SEALS 7010-071 7010, 7010-087, 7000, 7040 ETFE Rotor Seal ea. 7030-015 7030, 9030 ETFE Rotor Seal ea. 7060-074 7060, 7066, 9060 ETFE Rotor Seal ea. 7125-079 7125, 7125-081, 7725 ETFE Rotor Seal ea. 9125-082 9125, 9725 ETFE Rotor Seal ea. 9125-082 9125, 9725 ETFE Rotor Seal ea. 9125-083 3725, 3725-038 PEEK Rotor Seal ea. 9125-085 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-095 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-095 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 7123-584 MXT715-000 Stator ea. 7123-585 MXT715-105 Stator ea. 7123-590 MXP990-000 Stator ea. 7790-229 MXP7920-000 <th>7125-047</th> <th>7125, 7725, 9725</th> <th>Vespel Rotor Seal</th> <th>ea.</th>	7125-047	7125, 7725, 9725	Vespel Rotor Seal	ea.
8125-038 8125 Vespel Rotor Seal ea. ETFE BLEND ROTOR SEALS 7010-071 7010, 7010-087, 7000, 7040 ETFE Rotor Seal ea. 7030-015 7030, 9030 ETFE Rotor Seal ea. 7040-074 7040, 7064, 9060 ETFE Rotor Seal ea. 8125-079 7125, 7125-081, 7725 ETFE Rotor Seal ea. 8125-079 8125 ETFE Rotor Seal ea. 8125-097 8125 ETFE Rotor Seal ea. 91010-051 9010 ETFE Rotor Seal ea. 9125-082 9125, 9725 ETFE Rotor Seal ea. 9125-082 9125, 9725 ETFE Rotor Seal ea. 9125-082 9125, 9725 ETFE Rotor Seal ea. 9125-083 3725, 3725-038 PEEK Rotor Seal ea. 9125-095 7100, 7010, 9010 PEEK Rotor Seal ea. 9125-095 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-095 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9123-548 MXT715-100 Stator ea. 7123-548 MXT715-100 Stator ea. 7123-558 MXT715-102 Stator ea. 7123-558 MXT715-100 Stator ea. 7123-568 MXT715-100 Stator ea. 7770-229 MXP7920-000 Stator ea. 7780-044 MXP7980-000 Stator ea. 7980-004 MXP7980-000 Stator ea. 7980-004 MXP7980-000 Stator ea. 7990-146 MXP990-000 Stator ea. 7990-146 MXP990-000 Stator ea. 7990-148 MXP790-000 Stator ea. 7990-183 MXP7970-000 Stator ea. 7990-183 MXP7970-000 Stator ea. 7990-183 MXP7970-000 Stator ea. 7990-194 MXP7900-000 Stator ea. 7910-194 MXP7900-000 Stator ea. 79123-194 PXFEV50-100 Stator ea. 79123-194 PXFEV50-100 Stator ea. 79123-	7410-038	7410	Vespel Rotor Seal	ea.
### FIFE BLEND ROTOR SEALS 7010-071	7413-013	7413	Vespel Rotor Seal	ea.
7010-071 7010, 7010-087, 7000, 7040 ETFE Rotor Seal ea. 7030-015 7030, 9030 ETFE Rotor Seal ea. ea. 7060-074 7060, 7066, 9060 ETFE Rotor Seal ea. ea. 7125-079 7125, 7125-081, 7725 ETFE Rotor Seal ea. ea. 9125-097 8125 ETFE Rotor Seal ea. ea. 9125-097 8125 ETFE Rotor Seal ea. 9125-082 9125, 9725 ETFE Rotor Seal ea. 9125-082 9125, 9725 ETFE Rotor Seal ea. 9125-082 9125, 9725 ETFE Rotor Seal ea. 9125-083 3725, 3725-038 PEEK Rotor Seal ea. 9125-095 7125, 7725, 9125, 7725 PEEK Rotor Seal ea. 9125-095 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-095 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-095 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-095 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-095 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-095 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-095 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-095 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-095 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-095 PEEK Rotor Seal ea. 9125-0	8125-038	8125	Vespel Rotor Seal	ea.
7030-015 7030, 9030 ETFE Rotor Seal ea. 7060-074 7060, 7066, 9060 ETFE Rotor Seal ea. 7125-079 7125, 7125-081, 7725 ETFE Rotor Seal ea. 8125-097 8125 ETFE Rotor Seal ea. 9010-051 9010 ETFE Rotor Seal ea. 9125-082 9125, 9725 ETFE Rotor Seal ea. 9125-083 3725, 3725-038 PEEK Rotor Seal ea. 9125-095 7000, 7010, 9010 PEEK Rotor Seal ea. 9125-095 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-095 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-095 7123-509 MXT715-000 Stator ea. 7123-548 MXT715-102 Stator ea. 7123-568 MXT715-102 Stator ea. 7123-568 MXT7990-000 Stator ea. 7980-004 MXP7980-000 Stator ea. 7990-146 MXP990-000 Stator ea.	ETFE BLEND RO	TOR SEALS		
7060-074 7060, 7066, 9060 ETFE Rotor Seal ea. 7125-079 7125, 7125-081, 7725 ETFE Rotor Seal ea. 8125-097 8125 ETFE Rotor Seal ea. 9101-051 9010 ETFE Rotor Seal ea. 9125-082 9125, 9725 ETFE Rotor Seal ea. 9125-083 3725, 3725-038 PEEK Rotor Seal ea. 9125-083 3725, 3725-038 PEEK Rotor Seal ea. 9101-065 7000, 7010, 9010 PEEK Rotor Seal ea. 9101-065 7000, 7010, 9010 PEEK Rotor Seal ea. 9125-119 8125 PEEK Rotor Seal ea. 9125-095 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-095 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9123-548 MXT715-100 Stator ea. 7123-548 MXT715-105 Stator ea. 7123-568 MXT715-102 Stator ea. 7770-229 MXP7920-000 Stator ea. 7980-004 MXP7980-000 Stator ea. 7980-004 MXP7980-000 Stator ea. 7990-146 MXP990-000 Stator ea. 7990-146 MXP990-000 Stator ea. 7990-146 MXP990-000 Stator ea. 7990-148 MXP790-000 Stator ea. 7990-183 MXP7970-000 Stator ea. 8- 7760-014 MXP7960-000 Stator ea. 8- 7760-014 MXP7960-000 Stator ea. 8- 8- 8- 8- 8- 8- 8- 8- 8- 8- 8- 8- 8-	7010-071	7010, 7010-087, 7000, 7040	ETFE Rotor Seal	ea.
7125-079 7125,7125-081,7725 ETFE Rotor Seal ea. 8125-097 8125 ETFE Rotor Seal ea. 9010-051 9010 ETFE Rotor Seal ea. 9125-082 9125,9725 ETFE Rotor Seal ea. 9125-081 3725, 3725-038 PEEK Rotor Seal ea. 910-065 7000,7010, 9010 PEEK Rotor Seal ea. 8125-119 8125 PEEK Rotor Seal ea. 9125-095 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-095 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-095 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9123-550 MXT715-000 Stator ea. 7123-548 MXT715-102 Stator ea. 7123-550 MXF7950-000 Stator ea. 788-004 MXP7980-000 Stator ea. 7890-014 MXP990-000 Stator ea. 7900-183 MXP7970-000 Stator ea. 7900-184 </th <th>7030-015</th> <th>7030, 9030</th> <th>ETFE Rotor Seal</th> <th>ea.</th>	7030-015	7030, 9030	ETFE Rotor Seal	ea.
8125-097 8125 ETFE Rotor Seal ea. 9010-051 9010 ETFE Rotor Seal ea. 9125-082 9125, 9725 ETFE Rotor Seal ea. 9125-081 3725, 3725-038 PEEK Rotor Seal ea. 9010-065 7000, 7010, 9010 PEEK Rotor Seal ea. 8125-119 8125 PEEK Rotor Seal ea. 9125-095 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 8125-119 8125 PEEK Rotor Seal ea. 9125-095 7128, 7725, 9125, 9725 PEEK Rotor Seal ea. 8125-119 8125 PEEK Rotor Seal ea. 9125-095 7128, 7725, 9125, 9725 PEEK Rotor Seal ea. 8125-119 8125 MXT715-000 Stator ea. 8125-119 8125 MXT715-000 Stator ea. 8125-119 8125 MXT715-000 Stator ea. 8125-119 MXP7980-000 Stator ea. 77123-129 MXP7980-000 Stator	7060-074	7060, 7066, 9060	ETFE Rotor Seal	ea.
9010-051 9010 ETFE Rotor Seal ea. 9125-082 9125, 9725 ETFE Rotor Seal ea. PEEK BLEND ROTOR SEALS ETFE Rotor Seal ea. 9010-065 7000, 7010, 9010 PEEK Rotor Seal ea. 9010-065 7000, 7010, 9010 PEEK Rotor Seal ea. 9125-095 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-095 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. STATORS FOR MX SERIES II MODULES Stator ea. 7123-548 MXT715-000 Stator ea. 7123-550 MXT715-102 Stator ea. 7123-568 MXT715-102 Stator ea. 7980-004 MXP7980-000 Stator ea. 7980-014 MXP7980-000 Stator ea. 7990-146 MXP7990-000 Stator ea. 7990-183 MXP790-000 Stator ea. 8725-006 3725, 3710-038, 3000-038 Stator ea. 8725-006 3725, 3710-038, 3	7125-079	7125, 7125-081, 7725	ETFE Rotor Seal	ea.
P125-082 9125, 9725 ETFE Rotor Seal ea. 9725-018 3725, 3725-038 PEEK Rotor Seal ea. 9101-065 7000, 7010, 9010 PEEK Rotor Seal ea. 8125-119 8125 PEEK Rotor Seal ea. 9125-095 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-095 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. STATORS FOR MX SERIES II MODULES Stator ea. 7123-550 MXT715-102 Stator ea. 7123-568 MXT715-102 Stator ea. 7980-004 MXP7980-000 Stator ea. 7980-014 MXP7980-000 Stator ea. 7990-146 MXP7980-000 Stator ea. 7990-183 MXP7970-000 Stator ea. 7990-183 MXP7970-000 Stator ea. 7960-014 MXP7960-000 Stator ea. 8725-006 3725, 3710-038, 3000-038 Stator ea. 7010-049 <td< th=""><th>8125-097</th><th>8125</th><th>ETFE Rotor Seal</th><th>ea.</th></td<>	8125-097	8125	ETFE Rotor Seal	ea.
PEEK BLEND ROTOR SEALS 3725-018 3725-018 3725-038 PEEK Rotor Seal ea. 9010-065 7000, 7010, 9010 PEEK Rotor Seal ea. 9125-095 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-095 7126, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-095 7126, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-095 7126, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-095 7126, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-095 7126, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-095 9125, 9725 PEEK Rotor Seal ea. 9125-095 9125, 9710, 9	9010-051	9010	ETFE Rotor Seal	ea.
3725-018 3725, 3725-038 PEEK Rotor Seal ea. 9010-065 7000, 7010, 9010 PEEK Rotor Seal ea. 8125-119 8125 PEEK Rotor Seal ea. 9125-095 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. STATORS FOR MX SERIES II MODULES TURDANIA ea. 7123-548 MXT715-000 Stator ea. 7123-550 MXF715-102 Stator ea. 7770-229 MXP7980-000 Stator ea. 7790-029 MXP7980-000 Stator ea. 7980-004 MXP7980-000 Stator ea. 7990-146 MXP990-000 Stator ea. 7990-141 MXP7900-000 Stator ea. 7990-142 MXP7960-000 Stator ea. 7990-143 MXP7960-000 Stator ea. 3725-006 3725, 3710-038, 3000-038 Stator ea. 3725-006 3725, 3710-038, 3000-038 Stator ea. 7010-049 70000, 1730, 1740L	9125-082	9125, 9725	ETFE Rotor Seal	ea.
9010-065 7000, 7010, 9010 PEEK Rotor Seal ea. 8125-119 8125 PEEK Rotor Seal ea. 9125-095 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. 9125-095 7725, 7725, 9125, 9725 PEEK Rotor Seal ea. STATORS FOR MX SERIES II MODULES T123-584 MXT715-100 Stator ea. 7123-586 MXT715-102 Stator ea. 7770-229 MXP7920-000 Stator ea. 7788-004 MXP7980-000 Stator ea. 7980-004 MXP7980-000 Stator ea. 7990-146 MXP990-000 Stator ea. 7990-146 MXP990-000 Stator ea. 7990-183 MXP7970-000 Stator ea. 7990-183 MXP7970-000 Stator ea. 9760-014 MXP7980-000 Stator ea. 9760-014 MXP7980-000 Stator ea. 9760-015 MXP790-000 Stator ea. 9760-016 MXP9960-000 Stator ea. 9760-017 MXP7900-000 Stator ea. 9760-018 MXP7970-000 Stator ea. 9760-019 MXP7900-000 Stator ea. 9760-010 MXP7980-000 Stator ea. 9760-010 MXP7980-000 Stator ea. 9760-010 MXP7980-000 Stator ea. 9710-040 MXP7980-000 Stator ea. 97110-040 MXP7980-000 Stator ea. 971110-040 MXP7980-000 Stator ea. 9711110-040 MXP7980-000 PR/EVSO0-100 Stator ea. 971110-040 MXP7980-000 STATOR ea. 971110-040 MXP7980-0000 MXP7980-0000 STATOR EA. 971110-040	PEEK BLEND RC	TOR SEALS		
8125-119 8125 PEEK Rotor Seal ea. 9125-095 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. STATORS FOR MX SERIES II MODULES VARIANDES VARIANDES 7123-548 MXT715-000 Stator ea. 7123-568 MXT715-102 Stator ea. 770-229 MXP7920-000 Stator ea. 7980-004 MXP7980-000 Stator ea. 7980-014 MXP990-000 Stator ea. 7900-146 MXP990-000 Stator ea. 7900-179 MXP990-000 Stator ea. 7900-183 MXP790-000 Stator ea. 7960-014 MXP7960-000 Stator ea. 7960-015 MXP9960-000 Stator ea. 3725-006 3725, 3710-038, 3000-038 Stator ea. 3725-006 3725, 3710-038, 3000-038 Stator ea. 7010-040 7010, 7125, 7000, 7030 Stator ea. 7060-039 7060 and 7066 Stator	3725-018	3725, 3725-038	PEEK Rotor Seal	ea.
9125-095 7125, 7725, 9125, 9725 PEEK Rotor Seal ea. STATORS FOR MX SERIES II MODULES 7123-548 MXT715-000 Stator ea. 7123-550 MXT715-102 Stator ea. 7123-568 MXT715-102 Stator ea. 7980-004 MXP7980-000 Stator ea. 7980-004 MXP7980-000 Stator ea. 7990-146 MXP990-000 Stator ea. 7900-146 MXP970-000 Stator ea. 7900-183 MXP790-000 Stator ea. 7960-014 MXP7960-000 Stator ea. 960-018 MXP7960-000 Stator ea. 8725-006 3725, 3710-038, 3000-038 Stator ea. 8725-006 3725, 3710-038, 3000-038 Stator ea. 7010-049 7010, 7125, 7000, 7030 Stator ea. 7060-065 7060, EV501-100 Stator ea. 7123-127 PR/EV500-100 Stator ea.	9010-065	7000, 7010, 9010	PEEK Rotor Seal	ea.
STATORS FOR MX SERIES II MODULES 7123-548 MXT715-000 Stator ea. 7123-550 MXT715-105 Stator ea. 7123-568 MXT715-102 Stator ea. 7790-229 MXP7980-000 Stator ea. 7980-004 MXP7980-000 Stator ea. 7980-014 MXP990-000 Stator ea. 7900-146 MXP9900-000 Stator ea. 7900-183 MXP7970-000 Stator ea. 7900-183 MXP7970-000 Stator ea. 9960-002 MXP9960-000 Stator ea. 9960-002 MXP9960-000 Stator ea. 8725-006 3725, 3710-038, 3000-038 Stator ea. 8725-006 3725, 3710-038, 3000-038 Stator ea. 7010-040 7010, 7125, 7000, 7030 Stator ea. 7010-040 7060 and 7066 Stator ea. 7123-047 PR/EV500-100 Stator ea. 71	8125-119	8125	PEEK Rotor Seal	ea.
7123-548 MXT715-000 Stator ea. 7123-550 MXT715-105 Stator ea. 7123-568 MXT715-102 Stator ea. 7790-229 MXP7980-000 Stator ea. 7980-004 MXP7980-000 Stator ea. 7980-0146 MXP990-000 Stator ea. 7900-146 MXP790-000 Stator ea. 7900-183 MXP790-000 Stator ea. 7900-104 MXP7960-000 Stator ea. 7960-014 MXP7960-000 Stator ea. 7960-014 MXP7960-000 Stator ea. 7960-014 MXP7960-000 Stator ea. 7960-014 MXP7960-000 Stator ea. 7910-002 MXP9960-000 Stator ea. 7910-040 70001, 7301, 7040L Stator ea. 7010-049 70001, 7125, 7000, 7030 Stator ea. 7123-047 PR/EV500-103 Stator ea. <t< th=""><th>9125-095</th><th>7125, 7725, 9125, 9725</th><th>PEEK Rotor Seal</th><th>ea.</th></t<>	9125-095	7125, 7725, 9125, 9725	PEEK Rotor Seal	ea.
7123-550 MXT715-105 Stator ea. 7123-568 MXT715-102 Stator ea. 7770-229 MXP7920-000 Stator ea. 7980-004 MXP7980-000 Stator ea. 7980-014 MXP7980-000 Stator ea. 7900-146 MXP9900-000 Stator ea. 7900-179 MXP7900-000 Stator ea. 7900-183 MXP7970-000 Stator ea. 7960-014 MXP7960-000 Stator ea. 960-002 MXP9960-000 Stator ea. STATORS FOR OTHER IDEX HEALTH & SCIENCE VALVES 3725-006 3725, 3710-038, 3000-038 and 3030-038	STATORS FOR IV	IX SERIES II MODULES		
7123-568 MXT715-102 Stator ea. 7770-229 MXP7920-000 Stator ea. 7980-004 MXP7980-000 Stator ea. 7980-004 MXP7980-000 Stator ea. 7900-146 MXP9900-000 Stator ea. 7900-149 MXP7970-000 Stator ea. 7900-183 MXP7960-000 Stator ea. 7960-014 MXP7960-000 Stator ea. 7960-014 MXP7960-000 Stator ea. 7960-014 MXP7960-000 Stator ea. 7960-014 MXP7960-000 Stator ea. 8725-006 3725, 3710-038, 3000-038 and 3030-038 a	7123-548	MXT715-000	Stator	ea.
7770-229 MXP7920-000 Stator ea. 7980-004 MXP7980-000 Stator ea. 7986-004 MXP7980-000 Stator ea. 7900-146 MXP7900-000 Stator ea. 7900-179 MXP7900-000 Stator ea. 7900-183 MXP7960-000 Stator ea. 7960-014 MXP7960-000 Stator ea. 966-002 MXP960-000 Stator ea. 9760-002 MXP9960-000 Stator ea. STATORS FOR OTHER IDEX HEALTH & SCIENCE VALVES 3725-006 3725, 3710-038, 3000-038 and 3030-038 Stator ea. 7010-049 7000L, 7030L, 7040L Stator ea. 7010-040 7010, 7125, 7000, 7030 and 7040 Stator ea. 7060-039 7060 and 7066 Stator ea. 7123-047 PR/EV500-100 Stator ea. 7123-127 PR/EV500-100 Stator ea. 7123-128 PR/EV500-107 Stator ea.	7123-550	MXT715-105	Stator	ea.
7980-004 MXP7980-000 Stator ea. 7986-004 MXP7986-000 Stator ea. 7900-146 MXP9900-000 Stator ea. 7900-183 MXP7900-000 Stator ea. 7960-014 MXP7960-000 Stator ea. 9960-002 MXP9960-000 Stator ea. STATORS FOR O'THER IDEX HEALTH & SCIENCE VALVES 3725-006 3725, 3710-038, 3000-038 and 3030-038 and 30	7123-568	MXT715-102	Stator	ea.
7986-004 MXP7986-000 Stator ea. 7900-146 MXP9900-000 Stator ea. 7900-179 MXP7900-000 Stator ea. 7900-183 MXP7970-000 Stator ea. 7960-014 MXP7960-000 Stator ea. 9960-002 MXP9960-000 Stator ea. STATORS FOR OTHER IDEX HEALTH & SCIENCE VALVES STATOR OTHER IDEX HEALTH & SCIENCE VALVES <	7770-229	MXP7920-000	Stator	ea.
7900-146 MXP9900-000 Stator ea. 7900-179 MXP7900-000 Stator ea. 7900-183 MXP7970-000 Stator ea. 7960-014 MXP7960-000 Stator ea. 9960-002 MXP9960-000 Stator ea. STATORS FOR OTHER IDEX HEALTH & SCIENCE VALVES 3725-006 3725, 3710-038, 3000-038 and 3030-038 and 3030-038 and 3030-038 Stator ea. 7010-040 7010, 7125, 7000, 7030 and 7040 and 7040 Stator ea. 7060-039 7060 and 7066 Stator ea. 7123-047 PR/EV500-100 Stator ea. 7123-128 PR/EV500-100 Stator ea. 7123-129 PR/EV500-100 Stator ea. 7123-142 PR/EV500-104, EV501-104 Stator ea. 7123-143 PR/EV550-100 Stator ea. 7123-144 PR/EV550-101 Stator ea. 7123-149 PR/EV550-101 Stator ea. 7123-149 PR/EV550-101 <th>7980-004</th> <th>MXP7980-000</th> <th>Stator</th> <th>ea.</th>	7980-004	MXP7980-000	Stator	ea.
7900-179 MXP7900-000 Stator ea. 7900-183 MXP7970-000 Stator ea. 7960-014 MXP7960-000 Stator ea. 9960-002 MXP9960-000 Stator ea. STATORS FOR O'HER IDEX HEALTH & SCIENCE VALVES 3725-3710-038, 300-038 and 3030-038 and 3030-038 and 3030-038 and 3030-038 and 7040 Stator ea. 7010-040 7010, 7125, 7000, 7030 and 7060 Stator ea. 7060-039 7060 and 7066 Stator ea. 7123-047 PR/EV500-100 Stator ea. 7123-127 PR/EV500-100 Stator ea. 7123-128 PR/EV750-107 Stator ea. 7123-142 PR/EV500-104, EV501-104 Stator ea. 7123-143 PR/EV550-100, Stator ea. 7123-144 PR/EV550-100 Stator ea. 7123-148 PR/EV550-101 Stator ea. 7123-149 PR/EV550-101 Stator ea. 7123-180 PR7S3-100 and EV700-105	7986-004	MXP7986-000	Stator	ea.
7900-183 MXP7970-000 Stator ea. 7960-014 MXP7960-000 Stator ea. 9960-002 MXP9960-000 Stator ea. STATORS FOR OTHER IDEX HEALTH & SCIENCE VALVES 3725-006 3725, 3710-038, 3000-038 and 3030-038 Stator ea. 7010-069 7000L, 7030L, 7040L Stator ea. 7060-039 7060 and 7066 Stator ea. 7060-055 7060L, EV501-100 Stator ea. 7123-047 PR/EV500-100 Stator ea. 7123-127 PR/EV500-100 Stator ea. 7123-128 PR/EV750-107 Stator ea. 7123-129 PR/EV500-101 Stator ea. 7123-142 PR/EV500-104, EV501-104 Stator ea. 7123-143 PR/EV550-100 Stator ea. 7123-144 PR/EV550-100 Stator ea. 7123-149 PR/EV550-101 Stator ea. 7123-129 PR/EV550-101 Stator ea. </th <th>7900-146</th> <th>MXP9900-000</th> <th>Stator</th> <th>ea.</th>	7900-146	MXP9900-000	Stator	ea.
7960-014 MXP7960-000 Stator ea. 9960-002 MXP9960-000 Stator ea. STATORS FOR OTHER IDEX HEALTH & SCIENCE VALVES 3725-006 3725, 3710-038, 3000-038 and 3030-038 Stator ea. 7010-069 7000L, 7030L, 7040L Stator ea. 7010-040 37010, 7125, 7000, 7030 and 7040 Stator ea. 7060-039 7060 and 7066 Stator ea. 7060-065 7060L, EV501-100 Stator ea. 7123-047 PR/EV500-100 Stator ea. 7123-127 PR/EV500-107 Stator ea. 7123-128 PR/EV750-107 Stator ea. 7123-142 PR/EV500-104, EV501-104 Stator ea. 7123-143 PR/EV500-104, EV551-104 Stator ea. 7123-1447 PR/EV550-100 Stator ea. 7123-149 PR/EV550-101 Stator ea. 7123-129 PR/EV550-101 Stator ea. 7123-221 PR/S3-100 and EV750-105	7900-179	MXP7900-000	Stator	ea.
9960-002 MXP9960-000 Stator ea. STATORS FOR OTHER IDEX HEALTH & SCIENCE VALVES 3725-006 3725, 3710-038, 3000-038 and 3030-038 Stator ea. 7010-069 7000L, 7030L, 7040L Stator ea. 7010-040 7010, 7125, 7000, 7030 and 7040 Stator ea. 7060-039 7060 and 7066 Stator ea. 7123-047 PR/EV500-100 Stator ea. 7123-127 PR/EV500-107 Stator ea. 7123-128 PR/EV750-107 Stator ea. 7123-129 PR/EV500-104, EV501-104 Stator ea. 7123-141 PR/EV550-100 Stator ea. 7123-142 PR/EV550-104, EV551-104 Stator ea. 7123-143 PR/EV550-100 Stator ea. 7123-149 PR/EV550-101 Stator ea. 7123-149 PR/EV550-101 Stator ea. 7123-221 PR753-100 and EV700-105 Stator ea. 7123-223 PR/EV700-101	7900-183	MXP7970-000	Stator	ea.
STATORS FOR OTHER IDEX HEALTH & SCIENCE VALVES 3725-006 3725, 3710-038, 3000-038 and 3030-038 Stator ea. 7010-069 7000L, 7030L, 7040L Stator ea. 7010-040 7010, 7125, 7000, 7030 and 7040 Stator ea. 7060-039 7060 and 7066 Stator ea. 7060-065 7060L, EV501-100 Stator ea. 7123-047 PR/EV750-107 Stator ea. 7123-128 PR/EV750-107 Stator ea. 7123-128 PR/EV750-107 Stator ea. 7123-142 PR/EV550-104, EV501-104 Stator ea. 7123-143 PR/EV550-100, EV551-104 Stator ea. 7123-144 PR/EV550-101 Stator ea. 7123-148 PR/EV550-101 Stator ea. 7123-149 PR/EV550-101 Stator ea. 7123-129 PR753-100 and EV700-105 Stator ea. 7123-221 PR753-100 and EV750-105 Stator ea. 7520-030 (inlet)			Stator	ea.
3725-006 3725, 3710-038, 3000-038 and 3030-038 and 3030-038 Stator ea. 7010-069 7000L, 7030L, 7040L Stator ea. 7010-040 7010, 7125, 7000, 7030 and 7040 Stator ea. 7060-039 7060 and 7066 Stator ea. 7060-065 7060L, EV501-100 Stator ea. 7123-047 PR/EV750-107 Stator ea. 7123-128 PR/EV750-107 Stator ea. 7123-129 PR/EV500-104, EV501-104 Stator ea. 7123-142 PR/EV500-104, EV501-104 Stator ea. 7123-143 PR/EV550-100, EV551-104 Stator ea. 7123-144 PR/EV550-101 Stator ea. 7123-148 PR/EV550-101 Stator ea. 7123-149 PR/EV550-101 Stator ea. 7123-129 PR/EV50-100 and EV700-105 Stator ea. 7123-221 PR753-100 and EV750-105 Stator ea. 7123-223 PR/EV700-112 Stator ea. <th></th> <th></th> <th></th> <th>ea.</th>				ea.
3725-006 and 3030-038 Stator ea. 7010-069 7000L, 7030L, 7040L Stator ea. 7010-040 7010, 7125, 7000, 7030 and 7040 Stator ea. 7060-039 7060 and 7066 Stator ea. 7060-065 7060L, EV501-100 Stator ea. 7123-047 PR/EV500-100 Stator ea. 7123-128 PR/EV750-107 Stator ea. 7123-128 PR/EV750-104 Stator ea. 7123-128 PR/EV500-104, EV501-104 Stator ea. 7123-142 PR/EV550-104, EV551-104 Stator ea. 7123-143 PR/EV550-100 Stator ea. 7123-144 PR/EV550-101 Stator ea. 7123-148 PR/EV550-101 Stator ea. 7123-149 PR/EV550-101 Stator ea. 7123-129 PR/EV550-101 Stator ea. 7123-221 PR/553-100 and EV750-105 Stator ea. 7520-030 (inlet) 7520	STATORS FOR C		CIENCE VALVES	
7010-069 7000L, 7030L, 7040L Stator ea. 7010-040 7010, 7125, 7000, 7030 and 7040 Stator ea. 7060-039 7060 and 7066 Stator ea. 7060-065 7060L, EV501-100 Stator ea. 7123-047 PR/EV500-100 Stator ea. 7123-127 PR/EV750-107 Stator ea. 7123-128 PR/EV750-107 Stator ea. 7123-142 PR/EV500-104, EV501-104 Stator ea. 7123-143 PR/EV550-100, Ev501-104 Stator ea. 7123-144 PR/EV550-100 Stator ea. 7123-148 PR/EV550-101 Stator ea. 7123-149 PR/EV550-101 Stator ea. 7123-180 PR703-100 and EV700-105 Stator ea. 7123-221 PR753-100 and EV750-105 Stator ea. 7123-223 PR/EV700-112 Stator ea. 7520-030 (inlet) 7520 Stator ea. 7520-035 (outlet)	3725-006		Stator	ea.
7010-040 7010, 7125, 7000, 7030 and 7040 Stator ea. 7060-039 7060 and 7066 Stator ea. 7060-065 7060L, EV501-100 Stator ea. 7123-047 PR/EV500-100 Stator ea. 7123-127 PR/EV500-107 Stator ea. 7123-128 PR/EV750-107 Stator ea. 7123-142 PR/EV500-104, EV501-104 Stator ea. 7123-145 PR/EV550-100, Ev501-104 Stator ea. 7123-147 PR/EV550-100 Stator ea. 7123-148 PR/EV550-101 Stator ea. 7123-149 PR/EV550-101 Stator ea. 7123-180 PR703-100 and EV700-105 Stator ea. 7123-221 PR/EV550-101 Stator ea. 7123-223 PR/EV700-112 Stator ea. 7520-030 (inlet) 7520 Stator ea. 7520-035 (outlet) 7520 Stator ea. 7755-010 7750 <	7010-069		Stator	63
7010-040 and 7040 Stator ea. 7060-039 7060 and 7066 Stator ea. 7060-065 7060L, EV501-100 Stator ea. 7123-047 PR/EV500-100 Stator ea. 7123-127 PR/EV500-107 Stator ea. 7123-128 PR/EV7500-107 Stator ea. 7123-142 PR/EV500-104, EV501-104 Stator ea. 7123-145 PR/EV550-100 Stator ea. 7123-147 PR/EV550-101 Stator ea. 7123-148 PR/EV550-101 Stator ea. 7123-180 PR703-100 and EV700-105 Stator ea. 7123-221 PR753-100 and EV750-105 Stator ea. 7123-223 PR/EV700-112 Stator ea. 7520-030 (inlet) 7520 Stator ea. 7520-035 (outlet) 7520 Stator ea. 7725-010 7750 Stator ea. 7750-070 7750 Stator ea.				
7060-065 7060L, EV501-100 Stator ea. 7123-047 PR/EV500-100 Stator ea. 7123-127 PR/EV750-107 Stator ea. 7123-128 PR/EV750-107 Stator ea. 7123-142 PR/EV500-104, EV501-104 Stator ea. 7123-145 PR/EV550-104, EV551-104 Stator ea. 7123-147 PR/EV550-100 Stator ea. 7123-148 PR/EV550-101 Stator ea. 7123-149 PR/EV550-101 Stator ea. 7123-149 PR/EV550-101 Stator ea. 7123-129 PR/EV550-101 Stator ea. 7123-221 PR753-100 and EV700-105 Stator ea. 7123-223 PR/EV700-112 Stator ea. 7520-030 (inlet) 7520 Stator ea. 7520-035 (outlet) 7520 Stator ea. 755-000 PR/EV700-102 Stator ea. 7750-007 7750 Stator	7010-040		Stator	ea.
7123-047 PR/EV500-100 Stator ea. 7123-127 PR/EV750-107 Stator ea. 7123-128 PR/EV750-107 Stator ea. 7123-142 PR/EV500-104, EV501-104 Stator ea. 7123-145 PR/EV550-104, EV551-104 Stator ea. 7123-147 PR/EV550-100 Stator ea. 7123-148 PR/EV550-101 Stator ea. 7123-149 PR/EV550-101 Stator ea. 7123-149 PR/EV550-101 Stator ea. 7123-129 PR/EV550-101 Stator ea. 7123-221 PR753-100 and EV700-105 Stator ea. 7123-223 PR/EV700-112 Stator ea. 7410-041 7410 and 7413 Stator ea. 7520-030 (inlet) 7520 Stator ea. 7520-035 (outlet) 7520 Stator ea. 7550-002 PR/EV700-102 Stator ea. 7750-038 PR/EV700-100 Stator	7060-039	7060 and 7066	Stator	ea.
7123-127 PR/EV750-107 Stator ea. 7123-128 PR/EV700-107 Stator ea. 7123-142 PR/EV500-104, EV501-104 Stator ea. 7123-145 PR/EV550-100, EV551-104 Stator ea. 7123-147 PR/EV550-100 Stator ea. 7123-148 PR/EV50-101 Stator ea. 7123-149 PR/EV550-101 Stator ea. 7123-180 PR703-100 and EV700-105 Stator ea. 7123-221 PR753-100 and EV750-105 Stator ea. 7123-223 PR/EV700-112 Stator ea. 7410-041 7410 and 7413 Stator ea. 7520-030 (inlet) 7520 Stator ea. 7520-035 (outlet) 7520 Stator ea. 7755-010 7725(j) Stator ea. 7750-038 PR/EV700-100 Stator ea. 7750-038 PR/EV700-100 Stator ea. 8125-098 8125 Stator	7060-065	7060L, EV501-100	Stator	ea.
7123-128 PR/EV700-107 Stator ea. 7123-142 PR/EV500-104, EV501-104 Stator ea. 7123-145 PR/EV550-104, EV551-104 Stator ea. 7123-147 PR/EV550-100 Stator ea. 7123-148 PR/EV550-101 Stator ea. 7123-180 PR703-100 and EV700-105 Stator ea. 7123-221 PR753-100 and EV750-105 Stator ea. 7123-223 PR/EV700-112 Stator ea. 7410-041 7410 and 7413 Stator ea. 7520-030 (inlet) 7520 Stator ea. 7520-035 (outlet) 7520 Stator ea. 7725-010 7725(i) Stator ea. 7750-070 7750 Stator ea. 7750-070 7750 Stator ea. 8125-098 8125 Stator ea. 9650-009 PR/EV750-102 Stator ea.	7123-047	PR/EV500-100	Stator	ea.
7123-142 PR/EV500-104, EV501-104 Stator ea. 7123-145 PR/EV550-104, EV551-104 Stator ea. 7123-147 PR/EV550-100 Stator ea. 7123-148 PR/EV550-101 Stator ea. 7123-149 PR/EV550-101 Stator ea. 7123-180 PR703-100 and EV700-105 Stator ea. 7123-221 PR753-100 and EV750-105 Stator ea. 7123-223 PR/EV700-112 Stator ea. 7410-041 7410 and 7413 Stator ea. 7520-030 (inlet) 7520 Stator ea. 7520-035 (outlet) 7520 Stator ea. 7725-010 7725(i) Stator ea. 7750-030 PR/EV700-102 Stator ea. 7750-038 PR/EV700-100 Stator ea. 8125-098 8125 Stator ea. 9650-009 PR/EV750-102 Stator ea.	7123-127	PR/EV750-107	Stator	ea.
7123-145 PR/EV550-104, EV551-104 Stator ea. 7123-147 PR/EV550-100 Stator ea. 7123-148 PR/EV550-101 Stator ea. 7123-149 PR/EV550-101 Stator ea. 7123-180 PR703-100 and EV700-105 Stator ea. 7123-221 PR753-100 and EV750-105 Stator ea. 7410-041 7410 and 7413 Stator ea. 7520-030 (inlet) 7520 Stator ea. 7520-035 (outlet) 7520 Stator ea. 7752-010 7725(i) Stator ea. 7755-010 7725(i) Stator ea. 7750-038 PR/EV700-100 Stator ea. 8125-098 8125 Stator ea. 9650-009 PR/EV750-102 Stator ea.	7123-128	PR/EV700-107	Stator	ea.
7123-147 PR/EV550-100 Stator ea. 7123-148 PR/EV500-101 Stator ea. 7123-149 PR/EV550-101 Stator ea. 7123-180 PR703-100 and EV700-105 Stator ea. 7123-221 PR753-100 and EV750-105 Stator ea. 7123-223 PR/EV700-112 Stator ea. 7520-030 (inlet) 7520 Stator ea. 7520-035 (outlet) 7520 Stator ea. 7650-002 PR/EV700-102 Stator ea. 7725-010 7725(i) Stator ea. 7750-038 PR/EV700-100 Stator ea. 8125-098 8125 Stator ea. 9125-043 9125, 9010, 9030 and 9725(i) Stator ea. 9650-009 PR/EV750-102 Stator ea.	7123-142	PR/EV500-104, EV501-104	Stator	ea.
7123-148 PR/EV500-101 Stator ea. 7123-149 PR/EV550-101 Stator ea. 7123-180 PR703-100 and EV700-105 Stator ea. 7123-221 PR753-100 and EV750-105 Stator ea. 7123-223 PR/EV700-112 Stator ea. 7410-041 7410 and 7413 Stator ea. 7520-030 (inlet) 7520 Stator ea. 7520-035 (outlet) 7520 Stator ea. 7750-002 PR/EV700-102 Stator ea. 7750-010 7725(i) Stator ea. 7750-038 PR/EV700-100 Stator ea. 8125-098 8125 Stator ea. 9125-043 9125, 9010, 9030 and 9725(i) Stator ea. 9650-009 PR/EV750-102 Stator ea.		PR/EV550-104, EV551-104	Stator	ea.
7123-149 PR/EV550-101 Stator ea. 7123-180 PR703-100 and EV700-105 Stator ea. 7123-221 PR753-100 and EV750-105 Stator ea. 7123-223 PR/EV700-112 Stator ea. 7410-041 7410 and 7413 Stator ea. 7520-030 (inlet) 7520 Stator ea. 7520-035 (outlet) 7520 Stator ea. 7650-002 PR/EV700-102 Stator ea. 7725-010 7725(i) Stator ea. 7750-038 PR/EV700-100 Stator ea. 8125-098 8125 Stator ea. 9125-043 9125, 9010, 9030 and 9725(i) Stator ea. 9650-009 PR/EV750-102 Stator ea.			Stator	ea.
7123-180 PR703-100 and EV700-105 Stator ea. 7123-221 PR753-100 and EV750-105 Stator ea. 7123-223 PR/EV700-112 Stator ea. 7410-041 7410 and 7413 Stator ea. 7520-030 (inlet) 7520 Stator ea. 7520-035 (outlet) 7520 Stator ea. 7650-002 PR/EV700-102 Stator ea. 7725-010 7725(i) Stator ea. 7750-038 PR/EV700-100 Stator ea. 8125-098 8125 Stator ea. 9125-043 9125, 9010, 9030 and 9725(i) Stator ea. 9650-009 PR/EV750-102 Stator ea.				ea.
7123-221 PR753-100 and EV750-105 Stator ea. 7123-223 PR/EV700-112 Stator ea. 7410-041 7410 and 7413 Stator ea. 7520-030 (inlet) 7520 Stator ea. 7550-035 (outlet) 7520 Stator ea. 7650-002 PR/EV700-102 Stator ea. 7755-010 7750() Stator ea. 7750-038 PR/EV700-100 Stator ea. 8125-098 8125 Stator ea. 9125-043 9125, 9010, 9030 and 9725(i) Stator ea. 9650-009 PR/EV750-102 Stator ea.				ea.
7123-223 PR/EV700-112 Stator ea. 7410-041 7410 and 7413 Stator ea. 7520-030 (inlet) 7520 Stator ea. 7520-035 (outlet) 7520 Stator ea. 7650-002 PR/EV700-102 Stator ea. 7725-010 7725(i) Stator ea. 7750-070 7750 Stator ea. 7750-038 PR/EV700-100 Stator ea. 8125-098 8125 Stator ea. 9125-043 9125, 9010, 9030 and 9725(i) Stator ea. 9650-009 PR/EV750-102 Stator ea.				ea.
7410-041 7410 and 7413 Stator ea. 7520-030 (inlet) 7520 Stator ea. 7520-035 (outlet) 7520 Stator ea. 7650-002 PR/EV700-102 Stator ea. 7725-010 7725(i) Stator ea. 7750-070 750 Stator ea. 8125-098 8125 Stator ea. 9125-043 9125, 9010, 9030 and 9725(i) Stator ea. 9650-009 PR/EV750-102 Stator ea.				
7520-030 (inlet) 7520 Stator ea. 7520-035 (outlet) 7520 Stator ea. 7650-002 PR/EV700-102 Stator ea. 7725-010 7725(i) Stator ea. 7750-070 7750 Stator ea. 8125-098 8125 Stator ea. 9125-043 9125, 9010, 9030 and 9725(i) Stator ea. 9650-009 PR/EV750-102 Stator ea.				
7520-035 (outlet) 7520 Stator ea. 7650-002 PR/EV700-102 Stator ea. 7725-010 7725(j) Stator ea. 7750-070 7750 Stator ea. 7750-038 PR/EV700-100 Stator ea. 8125-098 8125 Stator ea. 9125-043 9125, 9010, 9030 and 9725(j) Stator ea. 9650-009 PR/EV750-102 Stator ea.				
7650-002 PR/EV700-102 Stator ea. 7725-010 7725(j) Stator ea. 7750-070 7750 Stator ea. 7750-038 PR/EV700-100 Stator ea. 8125-098 8125 Stator ea. 9125-043 9125, 9010, 9030 and 9725(j) Stator ea. 9650-009 PR/EV750-102 Stator ea.				
7725-010 7725(j) Stator ea. 7750-070 7750 Stator ea. 7750-038 PR/EV700-100 Stator ea. 8125-098 8125 Stator ea. 9125-043 9125, 9010, 9030 and 9725(j) Stator ea. 9650-009 PR/EV750-102 Stator ea.				
7750-070 7750 Stator ea. 7750-038 PR/EV700-100 Stator ea. 8125-098 8125 Stator ea. 9125-043 9125, 9010, 9030 and 9725(i) Stator ea. 9650-009 PR/EV750-102 Stator ea.				
7750-038 PR/EV700-100 Stator ea. 8125-098 8125 Stator ea. 9125-043 9125, 9010, 9030 and 9725(i) Stator ea. 9650-009 PR/EV750-102 Stator ea.		**		
8125-098 8125 Stator ea. 9125-043 9125, 9010, 9030 and 9725(i) Stator ea. 9650-009 PR/EV750-102 Stator ea.				
9125-043 9125, 9010, 9030 and 9725(i) Stator ea. 9650-009 PR/EV750-102 Stator ea.				
9650-009 PR/EV750-102 Stator ea.				
7,00 0E1 117/EV/30-100 Statol ed.			_	
			- Cutton	cu.

RHEBUILD KITS FOR MX SERIES II™ VALVES 7150-999 RheBuild Kit for MXT715-000 (includes 2 rotor seals) 7152-999 RheBuild Kit for MXT715-102 (includes 2 rotor seals) 7155-999 RheBuild Kit for MXT715-105 (includes 2 rotor seals) 7920-999 RheBuild Kit for MXP7920-000 and MXP7900-000 7960-999 RheBuild Kit for MXP9960-000 (includes rotor seal and stator face seal) 7961-999 RheBuild Kit for MXP7960-000 7970-999 RheBuild Kit for MXP7960-000 7970-999 RheBuild Kit for MXP7970-000 79801-999 RheBuild Kit for MXP7980-000 79801-999 RheBuild Kit for MXP7980-000 79801-999 RheBuild Kit for MXP9900-000 (includes rotor seal and stator face seal) RHEBUILD KITS FOR MANUAL VALVES 3725-999 RheBuild Kit for models 3725, 3725i, 3725-038, 3735i-038 7010-997 RheBuild Kit for model S1725, 3725i, 3725-038, 3735i-038 7010-997 RheBuild Kit including Stator for model 7010 7010-998 RheBuild Kit for model 7010 and 7010-type Valves 7125-999 RheBuild Kit for model 7410 7520-999 RheBuild Kit for models 7520 and 7526 (includes inlet stator and seal) 7725-999 RheBuild Kit for model 7725i-188 8125-999 RheBuild Kit for model 7725i-188 8125-999 RheBuild Kit for models 7125 and 7725i 7788-999 RheBuild Kit for models 7125 and 9126 9010-999 RheBuild Kit for models 9125 and 9126 9725-999 RheBuild Kit for models 9125 and 9725i; 7725(i) pH upgrade kit RHEBUILD KITS FOR MX SERIES I™ VALVES 7900-999 RheBuild Kit for models MX7900-000, MX7925-000, MX9900-000, MX9905-000	ea.
7152-999 RheBuild Kit for MXT715-102 (includes 2 rotor seals) 7155-999 RheBuild Kit for MXT715-105 (includes 2 rotor seals) 7920-999 RheBuild Kit for MXP7920-000 and MXP7900-000 7960-999 (includes rotor seal and stator face seal) 7961-999 RheBuild Kit for MXP7960-000 7970-999 RheBuild Kit for MXP7960-000 7970-999 RheBuild Kit for MXP7970-000 79801-999 RheBuild Kit for MXP7980-000 79801-999 RheBuild Kit for MXP7980-000 79801-999 RheBuild Kit for MXP7980-000 7900-999 RheBuild Kit for MXP7980-000 7900-999 RheBuild Kit for MXP7980-000 7910-999 RheBuild Kit for MXP7980-000 7910-999 RheBuild Kit for MXP7980-000 7910-999 RheBuild Kit for models 3725, 3725-338, 3735i-038 7910-997 RheBuild Kit for models 3725, 3725-338, 3735i-038 7910-997 RheBuild Kit for model 7010 and 7010-type Valves 7910-998 RheBuild Kit for model 7010 and 7010-type Valves 79125-999 RheBuild Kit for models 7125 and 7126 7410-999 RheBuild Kit for models 720 and 7526 7410-999 RheBuild Kit for models 7725 and 7725i 7788-999 RheBuild Kit for models 8125 and 8126 7725-999 RheBuild Kit for models 8125 and 8126 7910-999 RheBuild Kit for models 9725 and 9725i 7788-999 RheBuild Kit for models 9725 and 9725i 7725-999 RheBuild Kit for models 9725 and 9725i	ea.
7155-999 RheBuild Kit for MXT715-105 (includes 2 rotor seals) 7920-999 RheBuild Kit for MXP7920-000 and MXP7900-000 7960-999 RheBuild Kit for MXP9960-000 (includes rotor seal and stator face seal) 7961-999 RheBuild Kit for MXP7960-000 7970-999 RheBuild Kit for MXP7980-000 79801-999 RheBuild Kit for MXP7980-000 7980-999 RheBuild Kit for MXP7980-000 7900-999 RheBuild Kit for MXP9900-000 (includes rotor seal and stator face seal) RHEBUILD KITS FOR MANUAL VALVES 3725-999 RheBuild Kit for models 3725, 3725i, 3725-038, 3735i-038 7010-997 RheBuild Kit for models 3725, 3725i, 3725-038, 3735i-038 7010-998 RheBuild Kit for model 7010 and 7010-type Valves 7010-999 RheBuild Kit for model 7010 and 7010-type Valves 7125-999 RheBuild Kit for models 7125 and 7126 7410-999 RheBuild Kit for models 720 and 7526 (includes inlet stator and seal) 7725-999 RheBuild Kit for models 7725 and 7725i 7788-999 RheBuild Kit for models 8125 and 8126 9010-999 RheBuild Kit for models 9125 and 9126 9725-999 RheBuild Kit for models 9725 and 9725i; 7725(i) <	ea.
7920-999 RheBuild Kit for MXP7920-000 and MXP7900-000 7960-999 RheBuild Kit for MXP9960-000 (includes rotor seal and stator face seal) 7961-999 RheBuild Kit for MXP7960-000 7970-999 RheBuild Kit for MXP7980-000 79801-999 RheBuild Kit for MXP7980-000 79861-999 RheBuild Kit for MXP7980-000 7900-999 RheBuild Kit for MXP9900-000 (includes rotor seal and stator face seal) RHEBUILD KITS FOR MANUAL VALVES 3725-999 RheBuild Kit for models 3725, 3725i, 3725-038, 3735i-038 7010-997 RheBuild Kit for models 3725, 3725i, 3725-038, 3735i-038 7010-998 RheBuild Kit, pH Upgrade Kit for model 7000 7010-999 RheBuild Kit for model 7010 and 7010-type Valves 7125-999 RheBuild Kit for models 7125 and 7126 7410-999 RheBuild Kit for models 7220 and 7526 (includes inlet stator and seal) 7725-999 RheBuild Kit for models 7725 and 7725i 7788-999 RheBuild Kit for models 8125 and 8126 8010-999 RheBuild Kit for models 9125 and 9126 9725-999 RheBuild Kit for models 9125 and 9725i; 7725(i) 9725-999 RheBuild Kit for models 9725 and 9725i; 7725(i) PH Up	ea.
7960-999 RheBuild Kit for MXP9960-000 (includes rotor seal and stator face seal) 77961-999 RheBuild Kit for MXP7960-000 7970-999 RheBuild Kit for MXP7980-000 79801-999 RheBuild Kit for MXP7980-000 79801-999 RheBuild Kit for MXP7980-000 7900-999 RheBuild Kit for MXP9900-000 (includes rotor seal and stator face seal) RHEBUILD KITS FOR MANUAL VALVES 3725-999 RheBuild Kit for models 3725, 3725i, 3725-038, 3735i-038 7010-997 RheBuild Kit for models 3725, 3725i, 3725-038, 3735i-038 7010-998 RheBuild Kit for model 7010 and 7010 7010-999 RheBuild Kit, pH Upgrade Kit for model 7000 7010-999 RheBuild Kit for models 7125 and 7126 7410-999 RheBuild Kit for models 7250 and 7526 (includes inlet stator and seal) 7725-999 RheBuild Kit for models 7725 and 7725i 7788-999 RheBuild Kit for models 8125 and 8126 8010-999 RheBuild Kit for models 9125 and 9126 9010-999 RheBuild Kit for models 9725 and 9725i; 7725(i) 9725-999 RheBuild Kit for models 9725 and 9725i; 7725(i) 9725-999 RheBuild Kit for models 9725 and 9725i; 7725(i) 9725-999	ea.
(includes rotor seal and stator face seal) 7861-999 RheBuild Kit for MXP7960-000 78790-999 RheBuild Kit for MXP7980-000 78861-999 RheBuild Kit for MXP7980-000 78901-999 RheBuild Kit for MXP7986-000 78901-999 RheBuild Kit for MXP7986-000 78901-999 RheBuild Kit for MXP7986-000 78901-999 RheBuild Kit for MXP9900-000 78901-999 RheBuild Kit for models 3725, 3725i, 3725-038, 3735i-038 78725-999 RheBuild Kit for models 3725, 3725i, 3725-038, 3735i-038 789101-999 RheBuild Kit for model Stator for model 7010 789101-999 RheBuild Kit for model 7010 and 7010-type Valves 789101-999 RheBuild Kit for models 7125 and 7126 789101-999 RheBuild Kit for models 7520 and 7526 789101-999 RheBuild Kit for models 7725 and 7725i 789101-999 RheBuild Kit for models 7725 and 7725i 789101-999 RheBuild Kit for models 8125 and 8126 789101-999 RheBuild Kit for models 9125 and 9126 789125-999 RheBuild Kit for models 9125 and 9725i; 7725(i) 789101-999 RheBuild Kit for models 9725 and 7725i; 789125-999 RheBuild Kit for models 9125 and 9725i; 789125-999 RheBuild Kit for models 9125 and 9725i; 789125-999 RheBuild Kit for models 9725 and 9725i; 789125-999 RheBuild Kit for models 9725 and 9725i; 789125-999 RheBuild Kit for models 9125 and 9725i; 789125-999 RheB	ea.
7770-999 RheBuild Kit for MXP7970-000 78801-999 RheBuild Kit for MXP7980-000 78861-999 RheBuild Kit for MXP7986-000 7900-999 RheBuild Kit for MXP7986-000 7900-999 RheBuild Kit for MXP7986-000 7900-999 RheBuild Kit for MXP9900-000 (includes rotor seal and stator face seal) 7010-999 RheBuild Kit for models 3725, 3725-338, 3735i-038 7010-997 RheBuild Kit for models 3725, 3725-338, 3735i-038 7010-998 RheBuild Kit, pH Upgrade Kit for model 7000 7010-999 RheBuild Kit for model 7010 and 7010-type Valves 70125-999 RheBuild Kit for models 7125 and 7126 7010-999 RheBuild Kit for models 720 and 7526 (includes inlet stator and seal) 7725-999 RheBuild Kit for models 7725 and 7725i 7788-999 RheBuild Kit for models 8125 and 8126 9010-999 RheBuild Kit for models 8125 and 8126 9010-999 RheBuild Kit for models 9725 and 9725i 9725-999 RheBuild Kit for models 9725 and 9725i 7725-999 RheBuild Kit for models 9725 and 9725i	ea. ea. ea. ea. ea. ea. ea. ea.
79801-999 RheBuild Kit for MXP7980-000 79861-999 RheBuild Kit for MXP7986-000 7900-999 RheBuild Kit for MXP9900-000 RHEBUILD KITS FOR MANUAL VALVES 3725-999 RheBuild Kit for models 3725, 3725i, 3725-038, 3735i-038 7010-997 RheBuild Kit for models 3725, 3725i, 3725-038, 3735i-038 7010-998 RheBuild Kit, pH Upgrade Kit for model 7000 7010-999 RheBuild Kit for model 7010 and 7010-type Valves 7125-999 RheBuild Kit for models 7125 and 7126 7410-999 RheBuild Kit for models 7520 and 7526 7725-999 RheBuild Kit for models 7725 and 7725i 7788-999 RheBuild Kit for models 7725 and 7725i 7788-999 RheBuild Kit for models 8125 and 8126 8010-999 RheBuild Kit for models 9010 9125-999 RheBuild Kit for models 9725 and 9725i; 7725(i) 9725-999 RheBuild Kit for models 9725 and 9725i; 7725(i) PH upgrade kit RheBuild Kit for models MX7900-000, MX7925-000,	ea. ea. ea. ea. ea. ea.
79861-999 RheBuild Kit for MXP7986-000 7900-999 RheBuild Kit for MXP9900-000 (includes rotor seal and stator face seal) RHEBUILD KITS FOR MANUAL VALVES 3725-999 RheBuild Kit for models 3725, 3725i, 3725-038, 3735i-038 7010-997 RheBuild Kit including Stator for model 7010 7010-998 RheBuild Kit, pH Upgrade Kit for model 7000 7010-999 RheBuild Kit for model 7010 and 7010-type Valves 7125-999 RheBuild Kit for models 7125 and 7126 7410-999 RheBuild Kit for models 7520 and 7526 7520-999 RheBuild Kit for models 7725 and 7725i 7788-999 RheBuild Kit for models 7725 and 7725i 7788-999 RheBuild Kit for models 8125 and 8126 8010-999 RheBuild Kit for models 9125 and 9126 9725-999 RheBuild Kit for models 9725 and 9725i; 7725(i) 9725-999 RheBuild Kit for models 9725 and 9725i; 7725(i) PH upgrade kit RheBuild Kit for models MX7900-000, MX7925-000,	ea. ea. ea. ea. ea. ea. ea. ea.
RheBuild Kit for MXP9900-000	ea. ea. ea. ea. ea.
(includes rotor seal and stator face seal) RHEBUILD KITS FOR MANUAL VALVES RheBuild Kit for models 3725, 3725, 3725-038, 3735i-038 RheBuild Kit including Stator for model 7010 RheBuild Kit, pH Upgrade Kit for model 7000 RheBuild Kit for model 7010 and 7010-type Valves RheBuild Kit for model 7125 and 7126 RheBuild Kit for models 7125 and 7126 RheBuild Kit for models 7520 and 7526 (includes inlet stator and seal) RheBuild Kit for model 7725i-188 RheBuild Kit for models 7725i-188 RheBuild Kit for models 8125 and 8126 RheBuild Kit for models 910 RheBuild Kit for models 9125 and 9725i, 7725(i) pH upgrade kit RHEBUILD KITS FOR MX SERIES I** VALVES RheBuild Kit for models MX7900-000, MX7925-000,	ea. ea. ea. ea.
3725-999 RheBuild Kit for models 3725, 3725i, 3725-038, 3735i-038 7010-997 RheBuild Kit including Stator for model 7010 7010-998 RheBuild Kit, pH Upgrade Kit for model 7000 7010-999 RheBuild Kit for model 7010 and 7010-type Valves 7125-999 RheBuild Kit for models 7125 and 7126 7410-999 RheBuild Kit for models 7520 and 7526 (includes inlet stator and seal) 7725-999 RheBuild Kit for models 7725 and 7725i 7788-999 RheBuild Kit for models 7725i-188 8125-999 RheBuild Kit for models 8125 and 8126 9010-999 RheBuild Kit for models 9125 and 9126 9725-999 RheBuild Kit for models 9725 and 9725i; 7725(i) 9725-999 RheBuild Kit for models 9725 and 9725i; 7725(i) PH upgrade kit RHEBUILD KITS FOR MX SERIES I'' VALVES 7000-900 RheBuild Kit for models MX7900-000, MX7925-000,	ea. ea. ea.
RheBuild Kit including Stator for model 7010 RheBuild Kit, pH Upgrade Kit for model 7000 RheBuild Kit, pH Upgrade Kit for model 7000 RheBuild Kit for model 7010 and 7010-type Valves RheBuild Kit for models 7125 and 7126 RheBuild Kit for models 7520 and 7526 (includes inlet stator and seal) RheBuild Kit for models 7725 and 7725i RheBuild Kit for model 7725 and 7725i RheBuild Kit for model 7725i-188 RheBuild Kit for models 8125 and 8126 RheBuild Kit for models 9100 RheBuild Kit for models 9125 and 9126 RheBuild Kit for models 9725 and 9725i; 7725(i) PH upgrade kit RHEBUILD KITS FOR MX SERIES I" VALVES RheBuild Kit for models MX7900-000, MX7925-000,	ea. ea. ea.
7010-998 RheBuild Kit, pH Upgrade Kit for model 7000 7010-999 RheBuild Kit for model 7010 and 7010-type Valves 7125-999 RheBuild Kit for models 7125 and 7126 7410-999 RheBuild Kit for models 7520 and 7526 (includes inlet stator and seal) 7725-999 RheBuild Kit for models 7725 and 7725i 7788-999 RheBuild Kit for model 7725i-188 8125-999 RheBuild Kit for models 8125 and 8126 9010-999 RheBuild Kit for models 9125 and 9126 9725-999 RheBuild Kit for models 9725 and 9725i; 7725(i) pH upgrade kit RHEBUILD KITS FOR MX SERIES I" VALVES 7800-990 RheBuild Kit for models MX7900-000, MX7925-000,	ea. ea. ea.
7010-999 RheBuild Kit for model 7010 and 7010-type Valves 7125-999 RheBuild Kit for models 7125 and 7126 7410-999 RheBuild Kit for models 7520 and 7526 (includes inlet stator and seal) 7725-999 RheBuild Kit for models 7725 and 7725i 7788-999 RheBuild Kit for model 7725i-188 8125-999 RheBuild Kit for models 8125 and 8126 9010-999 RheBuild Kit for model 9010 9125-999 RheBuild Kit for models 9125 and 9126 9725-999 RheBuild Kit for models 9725 and 9725i; 7725(i) pH upgrade kit RHEBUILD KITS FOR MX SERIES I" VALVES 7700-990 RheBuild Kit for models MX7900-000, MX7925-000,	ea. ea.
7125-999 RheBuild Kit for models 7125 and 7126 7410-999 RheBuild Kit for model 7410 7520-999 RheBuild Kit for models 7520 and 7526 (includes inlet stator and seal) 7725-999 RheBuild Kit for models 7725 and 7725i 7788-999 RheBuild Kit for model 7725i-188 8125-999 RheBuild Kit for models 8125 and 8126 9010-999 RheBuild Kit for model 9010 9125-999 RheBuild Kit for models 9125 and 9126 9725-999 RheBuild Kit for models 9725 and 9725i; 7725(i) pH upgrade kit RHEBUILD KITS FOR MX SERIES I" VALVES 7700-990 RheBuild Kit for models MX7900-000, MX7925-000,	ea.
7410-999 RheBuild Kit for model 7410 7520-999 RheBuild Kit for models 7520 and 7526 (includes inlet stator and seal) 7725-999 RheBuild Kit for models 7725 and 7725i 7788-999 RheBuild Kit for model 7725i-188 8125-999 RheBuild Kit for models 8125 and 8126 9010-999 RheBuild Kit for model 9010 9125-999 RheBuild Kit for models 9125 and 9126 9725-999 RheBuild Kit for models 9725 and 9725i; 7725(i) pH upgrade kit RHEBUILD KITS FOR MX SERIES I" VALVES 7700-990 RheBuild Kit for models MX7900-000, MX7925-000,	
RheBuild Kit for models 7520 and 7526 (includes inlet stator and seal) 7725-999 RheBuild Kit for models 7725 and 7725i 7788-999 RheBuild Kit for model 7725i-188 8125-999 RheBuild Kit for models 8125 and 8126 9010-999 RheBuild Kit for model 9010 9125-999 RheBuild Kit for models 9125 and 9126 RheBuild Kit for models 9725 and 9725i; 7725(i) pH upgrade kit RHEBUILD KITS FOR MX SERIES !" VALVES RheBuild Kit for models MX7900-000, MX7925-000,	ea.
(includes inlet stator and seal) 7725-999 RheBuild Kit for models 7725 and 7725i 7788-999 RheBuild Kit for model 7725i-188 8125-999 RheBuild Kit for models 8125 and 8126 9010-999 RheBuild Kit for model 9010 9125-999 RheBuild Kit for models 9125 and 9126 RheBuild Kit for models 9725 and 9725i; 7725(i) pH upgrade kit RHEBUILD KITS FOR MX SERIES !" VALVES RheBuild Kit for models MX7900-000, MX7925-000,	
7788-999 RheBuild Kit for model 7725i-188 8125-999 RheBuild Kit for models 8125 and 8126 9010-999 RheBuild Kit for model 9010 9125-999 RheBuild Kit for models 9125 and 9126 9725-999 RheBuild Kit for models 9725 and 9725i; 7725(i) pH upgrade kit RHEBUILD KITS FOR MX SERIES !" VALVES 7000-990 RheBuild Kit for models MX7900-000, MX7925-000,	ea.
8125-999 RheBuild Kit for models 8125 and 8126 8010-999 RheBuild Kit for model 9010 9125-999 RheBuild Kit for models 9125 and 9126 8725-999 RheBuild Kit for models 9725 and 9725i; 7725(i) pH upgrade kit 8RHEBUILD KITS FOR MX SERIES I** VALVES 8RheBuild Kit for models MX7900-000, MX7925-000,	ea.
9010-999 RheBuild Kit for model 9010 9125-999 RheBuild Kit for models 9125 and 9126 9725-999 RheBuild Kit for models 9725 and 9725i; 7725(i) pH upgrade kit RHEBUILD KITS FOR MX SERIES I" VALVES RheBuild Kit for models MX7900-000, MX7925-000,	ea.
P125-999 RheBuild Kit for models 9125 and 9126 P725-999 RheBuild Kit for models 9725 and 9725i; 7725(i) pH upgrade kit RHEBUILD KITS FOR MX SERIES I" VALVES RheBuild Kit for models MX7900-000, MX7925-000,	ea.
RheBuild Kit for models 9725 and 9725i; 7725(i) pH upgrade kit RHEBUILD KITS FOR MX SERIES I" VALVES RheBuild Kit for models MX7900-000, MX7925-000,	ea.
pH upgrade kit RHEBUILD KITS FOR MX SERIES I'' VALVES RheBuild Kit for models MX7900-000, MX7925-000,	ea.
RheBuild Kit for models MX7900-000, MX7925-000,	ea.
	ea.
7960-999 RheBuild Kit for model MX7960-000	ea.
7980-999 RheBuild Kit for model MX7980-000	ea.
7984-999 RheBuild Kit for model MX7984-000	ea.
7986-999 RheBuild Kit for model MX7986-000	ea.
RHEBUILD KITS FOR LABPRO™ & EV AUTOMATED FLUIDIC INSTRUMENTS	
1006-999 RheBuild Kit for model PR/EV100-106	ea.
5001-999 RheBuild Kit for models PR/EV500-101 and PR/EV550-101	ea.
5100-999 RheBuild Kit for models PR/EV500-100 and PR/EV550-100	ea.
5104-999 RheBuild Kit for models PR/EV500-104 and PR/EV550-104	ea.
7004-999 RheBuild Kit for models PR/EV700-104 and PR/EV750-104	ea.
7112-999 RheBuild Kit for models PR/EV700-112 and PR/EV750-112	ea.
7501-999 RheBuild Kit for models PR/EV700-100 and PR/EV750-100	ea.
7502-999 RheBuild Kit for models PR/EV700-102 and PR/EV750-102	ea.
7507-999 RheBuild Kit for models PR/EV700-107 and PR/EV750-107	ea.
7531-999 RheBuild Kit for models PR703-100 and PR753-100	ea.
TELEVISION NO TO CONTROL TO THE TOTAL TOTAL TO THE TOTAL TOTAL TO THE TOTAL TOTAL TO THE TOTAL TOTAL TOTAL TOTAL TO THE TOTAL TO	ou.







VALVE ACCESSORIES

Our valve accessories include a variety of products that work with and are specific to our valve mechanics. From Sample Loops, driver boards, or mounting brackets we offer a wide array of accessories to meet your system requirements. We also include tools that work specifically with our valves and valve components.

- 127 STAINLESS STEEL SAMPLE LOOPS
- **128** PEEK SAMPLE LOOPS
- 131 SUCTION NEEDLE ADAPTER
- **132** INJECTION PORT ADAPTERS
- 133 WRENCHES, BRACKETS, & REPLACEMENT FITTINGS



Stainless Steel Sample Loops



How to Properly Install Sample Loops: Stainless Steel


Stainless steel sample loops are supplied with fittings that are not swaged onto the tube. It is important that the loop be completely bottomed in the injector port before the ferrule is swaged onto the tube. The depth of the tubing holes may vary slightly from port to port and from valve to valve. A fitting made up in one port may leave a small cavity in another port. The cavity causes high dispersion and peak distortion such as fronting, tailing, or broadening. It is good practice to label loop ends so they will be replaced in the same, respective ports that were used in swaging the ferrules. Hint: swaging ferrules separately on each side, into each respective valve port makes loop installation easier.

These high quality stainless steel sample loops have burr-free, square-cut ends to ensure a flush connection to valve ports. The size designations of loops are nominal. The actual volumes can differ from the theoretical designations because the ID tolerance varies depending on the tubing tolerance of the metal tubing bore. Accuracy of large metal loops (1.0 mm, 0.040" bore) is about $\pm 5\%$, intermediate loops (0.5 mm, 0.020" bore) $\pm 10\%$, and small loops (0.2 mm, 0.007" bore) $\pm 30\%$.

Since both standards and unknowns are usually analyzed using the same sample loop, knowledge of the actual, accurate volume is rarely needed. If the sample loop volume must be known, it is best to calibrate the loop in place on the valve so the flow passages in the valve are also taken into account. An alternative to calibration is to use a dual mode injector and partial-filling method of loading. See the "Sample Loop Loading" Application Note on page 131.

Model 7725 Injector loops are not interchangeable with loops for the model 7125. The port angle for the 7725 is 30° whereas the port angle for the 7125 is 20° requiring the loops to have a different shape. Model 8125 Micro-Scale Sample Injector requires special loops in the 5.0 μ L to 50 μ L range. The 8125 sample loops are made with 0.5 mm (0.020″) OD tubing.

To install the sample loop:

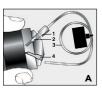

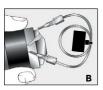




Figure 1 Cut-away view of stainless steel sample loop installation

- Take one end of the loop and place the nut (1) and ferrule (2) onto the tubing (3) with the threaded portion of the nut and tapered portion of the ferrule toward the end. See Figure A. 2 Insert the tubing into port (4). Confirm that the tubing is bottomed in the valve port as shown in Figure A. While firmly pressing down on the tubing, hand-tighten the nut as tight as possible. With the IDEX Wrench (page 51), designed especially for fittings, tighten one quarter turn past finger tight. Remove the loop to confirm the ferrule is swaged onto the tube. Repeat steps a-d with the other end of the loop while the swaged end remains outside
- the valve port. See Figure B.

 6 Reinstall each end of the loop to their respective ports. See Figure C.

Part No.	Volume	Tubing	O+-
		ON VALVES (DO NOT USE FOR 7725)	Qty.
7020	5 μL Sample Loop	0.18 mm (0.007") ID x 1/16" OD	ea.
7021	10 µL Sample Loop	0.30 mm (0.012") ID x 1/16" OD	ea.
7022	20 µL Sample Loop	0.51 mm (0.020") ID x 1/16" OD	ea.
7023	50 µL Sample Loop	0.51 mm (0.020") ID x 1/16" OD	ea.
7024	100 µL Sample Loop	0.51 mm (0.020") ID x 1/16" OD	ea.
7025	200 µL Sample Loop	0.76 mm (0.030") ID x 1/16" OD	ea.
7026	500 µL Sample Loop	0.76 mm (0.030") ID x 1/16" OD	ea.
7027	1.0 mL Sample Loop	0.76 mm (0.030") ID x 1/16" OD	ea.
7028	2.0 mL Sample Loop	1.0 mm (0.040") ID x 1/16" OD	ea.
7029	5.0 mL Sample Loop	1.0 mm (0.040") ID x 1/16" OD	ea.
1876	10 mL Sample Loop	2.0 mm (0.080") ID x 1/8" OD	ea.
STAINLESS STEEL LO	OPS FOR 3725-038, 3725I-038		
3065-018	2.0 mL Sample Loop	2.0 mm (0.080") ID x 1/8" OD	ea.
3065-019	5.0 mL Sample Loop	2.0 mm (0.080") ID x 1/8" OD	ea.
3065-023	10 mL Sample Loop	2.0 mm (0.080") ID x 1/8" OD	ea.
3065-025	20 mL Sample Loop	2.0 mm (0.080") ID x 1/8" OD	ea.
	OPS FOR 7725, 7725I, PR/EV7 DO NOT USE FOR 7125)	00-100, PR/EV703-100, MX MODULE	
7755-020	5 μL Sample Loop	0.18 mm (0.007") ID x 1/16" OD	ea.
7755-021	10 μL Sample Loop	0.30 mm (0.012") ID x 1/16" OD	ea.
7755-022	20 μL Sample Loop	0.30 mm (0.012") ID x 1/16" OD	ea.
7755-023	50 μL Sample Loop	0.51 mm (0.020") ID x 1/16" OD	ea.
7755-024	100 μL Sample Loop	0.51 mm (0.020") ID x 1/16" OD	ea.
7755-025	200 μL Sample Loop	0.76 mm (0.030") ID x 1/16" OD	ea.
7755-026	500 μL Sample Loop	0.76 mm (0.030") ID x 1/16" OD	ea.
7755-027	1.0 mL Sample Loop	0.76 mm (0.030") ID x 1/16" OD	ea.
7755-028	2.0 mL Sample Loop	1.0 mm (0.040") ID x 1/16" OD	ea.
7755-029	5.0 mL Sample Loop	1.0 mm (0.040") ID x 1/16" OD	ea.
1876	10 mL Sample Loop	2.0 mm (0.080") ID x 1/8" OD	ea.
STAINLESS STEEL LOC	OPS FOR 8125 INJECTOR (USE	7755-024 TO 7755-029 FOR VOLUMES >	50 μL)
8020	5 μL Sample Loop	0.20 mm (0.008") ID x 0.020" OD	ea.
8021	10 μL Sample Loop	0.20 mm (0.008") ID x 0.020" OD	ea.
8022	20 μL Sample Loop	0.25 mm (0.010") ID x 0.020" OD	ea.
8023	50 μL Sample Loop	0.30 mm (0.012") ID x 0.020" OD	



PEEK Sample Loops

Flexible PEEK sample loops are alternatives to stainless steel loops. PEEK loop ends are provided with clean, straight cuts for easy valve installation.

PEEK polymer is inert to almost all organic solvents and is biocompatible, giving PEEK loops added versatility. Natural PEEK is used for these sample loops. Like metal loops, the size designations of PEEK loops are nominal. The actual volumes can differ from the theoretical designations because of the tolerance of the tubing bore. Accuracy of large PEEK loops (0.8 mm, 0.030" bore) is about $\pm 14\%$, intermediate loops (0.5 mm, 0.020") $\pm 21\%$, and small loops (0.2 mm, 0.007") $\pm 65\%$.

PEEK loops are also supplied with unswaged RheFlex® fittings but do not require the same swaging precaution. The fittings can reposition along the loop tubing when the fitting is reinserted in the ports for correct loop installation.

Sample Mobile Phase Tube Wall

Figure 1 Schematic of sample flow through mobile phase between tubing walls

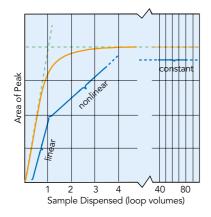


Figure 2 Sample mass (observed peak area) vs. volume of sample dispensed from the syringe, in units of loop volumes injected onto the column from our dual mode injector such as model 7725

APPLICATION NOTE

Fluidic Movement in Tubes

Q: "Why can I load only up to half of the volume of the loop in partial-filling method?"

A: Sample occupies $2 \mu L$ of loop for every $1 \mu L$ loaded from the syringe. For example, $10 \mu L$ of sample spreads out over the entire length of a $20 \mu L$ loop. Any additional sample loaded will overflow the end of the loop and exit out to waste. Reproducibility is poor because the volume of sample in the loop is different from the known volume originally loaded by your syringe.

Fluid spreads in a parabolic shape through a tube instead of moving in one plug because the velocity is different at the center of the tube than at the walls. The velocity at the center of the tube is twice the average velocity, and near the wall the velocity is almost zero, creating a parabolic shape. This fluidic movement is called laminar flow. See Figure 1.

In dual mode injection valves (see "Sample Loop Loading" Application Note on page 131) the sample from the syringe needle loads directly into the sample loop. The sample volume is known since there is no sample waste. The laminar flow phenomenon accounts for the shape of the plot as shown in Figure 2. Note that the plot has three regions:

- 1 Partial-Filling Region. When the volume dispensed is less than half the loop volume, the curve is linear. Sample has not reached the end of the loop. Within this region, performance depends on the syringe and operator.
- 2 Nonlinear Region. When the volume dispensed is between half the loop volume and about two loop volumes, the curve is nonlinear. Sample is lost from the loop, so reproducibility is poor. If you dispense a volume equal to the loop size, you are in this region of poor performance.
- **3** Complete-Filling Region. When the volume of sample dispensed is several loop volumes, the loop contains only pure sample, undiluted by residual mobile phase. Within this region, reproducibility is highest.

In the single mode injection valves the sample must pass through a connecting passage before it reaches the sample loop. Since some of the sample dispensed from the syringe remains in the connecting passageway, an unknown amount enters the sample loop. Therefore, single mode injection valves achieve high reproducibility only by using the complete-filling method.



PEEK Physical Strength Characteristics

Although PEEK material is compatible with virtually all solvents, there are many factors that affect burst pressure of PEEK tubing. Factors such as increases in inner diameter, temperature, exposure time, and concentration of organic solvents affect the degradation of PEEK. Other solvents such a THF, methylene chloride and DMSO cause PEEK tubing to swell while concentrated nitric acid and sulfuric acid weaken the tubing.



How to Find and Fix Common Sample Injector Leaks

Leaks cause valuable sample loss. Nobody wants that. The key to the valve holding pressure is the integrity of the sealing surfaces. If there is a scratch on the sealing surface, or the needle seal in the rotor seal is damaged, a leak may appear. It is also important to realize what appears to be a leak can instead be a result of siphoning. The following are the three most common situations in which fluid leaks occur.

- 1 If fluid leaks out of the needle port only while loading the loop (i.e., while pushing down on the plunger of the syringe), the problem is most likely that the needle seal or the needle port fitting in the loop filler port is not gripping the syringe needle tightly enough. Tighten the needle seal grip by pushing with the eraser end of a pencil on the needle port (See Figure 1). The tightening reduces the hole diameter of the needle seal and port fitting.
- 2 If fluid leaks continuously from the needle port or vent lines and/or from the stator-tostator ring interface, replace the rotor seal and/or stator face assembly. Scratches on the rotor seal or cracks in the stator face assembly allow mobile phase to escape and cause cross port leakage. Genuine IDEX Health & Science RheBuild® Kits are listed on page 124.
- 3 If fluid leaks from the needle port and/or vent lines but eventually stops, the cause is most likely siphoning and not a leak. Siphoning occurs if the vent lines are lower or higher than the needle port. Adjust the vent line(s) so that the outlet is at the same horizontal level as the needle port to prevent siphoning. (See Figure 2).

For other leakage or injection troubleshooting, refer to our Troubleshooting Guide for HPLC Injection Problems. You may download the Guide from our web site: www.idex-hs.com under Education & Tools.

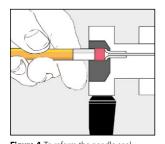


Figure 1 To reform the needle seal, push the eraser end of a pencil against the needle port

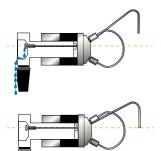
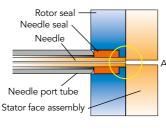


Figure 2 Needle port level compared to the level of vent line outlet:

(A) siphoning occurs when the vent line outlet is above the needle port level

(B) siphoning does not occur if the vent line outlet is the same horizontal level as the needle port

PEEK Sample Loops (Cont.)



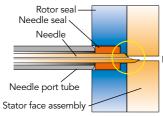


Figure 1 A square cut needle: (A) stops against the stator face assembly; The tip of a pointed needle (B) slips into the stator face and the tip breaks off as the valve rotates



Using Proper Syringe Needles


With front-loading injection valves it is important to use the correct needle when loading the sample loop. An incorrect needle will damage the valve and can cause poor reproducibility. When the needle is too short the tip will not reach the needle seal. When the needle is too small in diameter the seal will not grip tightly enough. Needles with a beveled tip can damage the rotor seal and stator face assembly (see Figure 1). The needle should be #22 gauge (0.028"-0.0285"/0.72 mm), and 90° point style (square cut end). Model 3725i requires a #16 gauge (0.0645"-0.0655"/ 1.65 mm) needle. Never use a beveled, pointed, or tapered needle.

Needle specifications are not critical when using a Loop Filler Port to load the sample loop. However, it is important to tighten the needle port fitting around the needle if using a syringe needle with a slightly smaller diameter than 0.7 mm (0.028").

If the loading method used is complete-filling, a syringe without a needle can be used. A syringe fitted with a Needle Port Cleaner can be used with a front-loading valve (Figure 2A) or with a Loop Filler Port (Figure 2B).

Needle port accessories are listed on page 132.



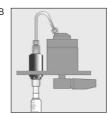

(A) Syringe fitted with Needle Port Cleaner (Part # 7125-054) loading a front-loading valve (model 7725); (B) loading a Loop Filler Port (Part # 7012)

PEEK Sample Loops

Part No.	Volume	Tubing	Valco No.	Qty.
PEEK LOO	PS FOR 3725, 3725I INJECTION VALVES			
3055-018	2.0 mL Sample Loop	1.6 mm (0.062") ID x 1/8" OD	N/A	ea.
3055-019	5.0 mL Sample Loop	1.6 mm (0.062") ID x 1/8" OD	N/A	ea.
3055-023	10 mL Sample Loop	2.0 mm (0.080") ID x 1/8" OD	N/A	ea.
3055-025	20 mL Sample Loop	2.0 mm (0.080") ID x 1/8" OD	N/A	ea.
PEEK LOOP	PS FOR 9725, 9010, PR/EV750-100, PR/EV753-100 INJECTION VALVE			
Part No.	Volume	Bore / Tubing	Valco No.	
9055-020	5.0 µL Sample Loop	0.18 mm (0.007") ID x 1/16" OD	SL5CWPK	ea.
9055-021	10 μL Sample Loop	0.25 mm (0.010") ID x 1/16" OD	SL10WPK	ea.
9055-022	20 μL Sample Loop	0.25 mm (0.010") ID x 1/16" OD	SL20WPK	ea.
9055-023	50 μL Sample Loop	0.51 mm (0.020") ID x 1/16" OD	SL50WPK	ea.
9055-024	100 μL Sample Loop	0.51 mm (0.020") ID x 1/16" OD	SL100WPK	ea.
9055-025	200 μL Sample Loop	0.51 mm (0.020") ID x 1/16" OD	N/A	ea.
9055-026	500 μL Sample Loop	0.76 mm (0.030") ID x 1/16" OD	SL500WPK	ea.
9055-027	1.0 mL Sample Loop	0.76 mm (0.030") ID x 1/16" OD	SL1KCWPK	ea.
9055-028	2.0 mL Sample Loop	0.76 mm (0.030") ID x 1/16" OD	SL2KCWPK	ea.
9055-029	5.0 mL Sample Loop	0.76 mm (0.030") ID x 1/16" OD	N/A	ea.
9055-033	10 mL Sample Loop	0.76 mm (0.030") ID x 1/16" OD	N/A	ea.
PEEK LOO	PS FOR 7725, 7725I, PR/EV700-100			
7123-227	1 μL Sample Loop (models PR/EV700-100 and EV750-100 only)	Internal groove	N/A	ea.
7755-015	2 μL Sample Loop (models 7725, 7725i, and 9725(i) only)	Internal groove	N/A	ea.
REPLACEM	MENT RHEFLEX FITTINGS FOR PEEK LOOPS			
6000-078	Nut/Ferrule Set, Natural PEEK, 5/16-24, for 1/8" OD loops			ea.
6000-079	Ferrules, Natural PEEK, for 1/8" OD loops			5-pk
6000-251	Ferrules, Natural PEEK, for 1/16" OD loops			10-pk
6000-254	Nut/Ferrule Sets, Natural PEEK, 10-32, for 1/16" OD loops			10-pk



Suction Needle Adapter

Our adaptable Loop Filler Ports (Part #7012 and 9012) are used to load sample from syringe needles or luer tips. The Needle Port (Part #9013) conserves sample by minimizing the volume between the needle and the valve.





Dual Mode Sample Loop Loading: Partial-Filling vs. Complete-Filling

Partial-Filling

Use the partial-filling method if you need to conserve sample, or if you want to vary sample volume frequently.

In partial-filling, the syringe sets the volume injected onto the column. There is no sample waste, and the volume injected onto the column is equal to that dispensed from the syringe. Reproducibility is 1.0% relative standard deviation (RSD). The volume of the sample loaded is limited to half the sample loop volume. For example, the most you can load into a 200 μ L sample loop is 100 μ L.

Complete-Filling


Use the complete-filling method if you have plenty of sample, if you do not vary sample volume, or if you need high reproducibility.

In complete-filling, the loop sets the volume loaded onto the column. Use excess sample (two to five loop volumes) to replace all the mobile phase in the loop. See Figure 2. Change the loop to vary the sample volume. Reproducibility is typically 0.1% RSD for loop sizes $\geq 5~\mu L$. Accuracy is limited as loop volumes are nominal.

Q: "Which method should I use and which IDEX Health & Science sample injection valves use this method?"


A: There are two types of injection valves available: dual mode and single mode. Dual mode injection valves allow both partial- and complete-filling whereas single mode injection valves allow only complete-filling. See manual injection valves, page 123.

If you are collecting experimental data, sample is scarce, and/or you want to use different sample volumes, a dual mode injector with a large volume sample loop is appropriate. Only dual mode injection valves allow the partial-filling method for easily varying your volumes (up to half your sample loop volume) by setting the syringe volume. Once you begin routine analysis, and/or you have an abundance of sample, either a dual mode or single mode injector is appropriate. Both types of injection valves allow the complete-filling method in which you overfill the sample loop. Complete-filling maximizes the reproducibility of your results.

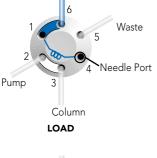

Waste

Flow path for the typical dual mode injector

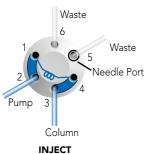

Part No.	Description	Qty.
	DLE ADAPTER & ACCESSORIES	<u> </u>
7012	Stainless Steel Loop Filler Port	ea.
7125-054	Needle Port Cleaner	ea.
9012	PEEK Loop Filler Port	ea.
9013	PEEK Needle Port	ea.
9125-076	Suction Needle Adapter (for Models 7725 and 9725)	ea.




Injection Port Adapters

- > For 360 µm OD tubing
- Mount on bracket or bulkhead


To introduce sample, connect 360 μ m OD capillary tubing to an Injection Port Adapter Assembly. This adapter accepts standard 22 gauge Hamilton-style injection syringe needles. No additional swept volume is added to the fluid pathway by this adapter, as the needle butts directly against the connecting tubing during injections. The adapter can be bulkhead mounted or mounted with the V-447 Kits.

To introduce a sample directly into a 10-32 port, purchase a M-432-03 separately.

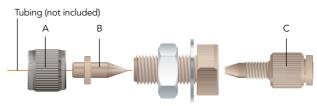

M-432

Micro Injection Port Adapter Assembly

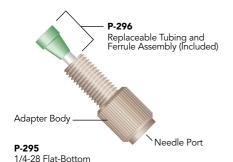

Injection Port Adapter

RELATED PRODUCTS

		_	•
FOR 360 µm OD TUBING			
M-432 and V-447	P-416BLK	F-152	M-432-03

> For use with Injection Valves on page 123.

This simple, biocompatible adapter is designed specifically for the Injection Valves on page 123 and can also convert any 1/4-28 flat-bottom port into a port that can accept a standard 22 gauge HPLC injection needle. This injection port adapter is adjustable, so you can create a snug fit around the needle to prevent any leaking of the analyte. In addition, this product features an internal stop that prevents you from inserting the needle too far, eliminating the possibility of damaging the valve with the needle tip.

Part No.	Description	Qty.
MICRO INJECTIO	ON PORT ADAPTER	
FOR 360 µm OD	TUBING	ea.
F-152	Replacement MicroFerrule for M-432, Natural PEEK	ea.
M-432	Micro Injection Port Adapter Assembly	ea.
M-432-03	Replacement Tubing/Fitting Assembly for M-432 & M-433	ea.
P-416BLK	Replacement Female Nut for M-432, Black PEEK	ea.
V-447	Micro Injection Port Adapter Assembly Actuator Mounting Kit Includes (1) M-432 with mini-actuator bracket and (2) mounting screws	ea.
1/4-28 FLAT-BOT	TTOM INJECTION PORT ADAPTER	
P-295	Adjustable Injection Port Adapter	ea.
P-296	Replacement Tubing/Ferrule Assembly	ea.

Wrenches, Brackets, & Replacement Fittings



Valve Wrenches

- > For convenient wrench-tightening of fittings on high pressure rotary shear valves
- > For removal of knobs on Manual Valves

The smartly designed IDEX Wrench is a double-ended slotted socket wrench that fits over 1/16" and 1/8" OD tubing. It easily loosens and tightens 1/4" and 5/16" hex head stainless steel or PEEK fittings. The "Z" shape of the IDEX Wrench provides ideal leverage for changing sample loops and fittings, and keeps one end from restricting the use of the other.

The V-103 is an Allen (hex-key) wrench designed to remove the knob from our V-101 valves (page 123). The V-104 is an Allen wrench that can be used to remove the knob from our Medium Pressure Selection and Injection Valves (also found on page 123).



Mounting Brackets

Our mounting brackets and panels of different shapes and sizes organize and provide a sturdy support for IDEX Health & Science valves. The Ring Stand Mounting Bracket now allows the valves to mount onto common laboratory equipment.

MXX Replacement Fittings

Use these replacement Ferrules and O-rings for 1/8" and 1/16" tubing with the MXX Series II valves shown on page 119. Please see the part number chart below for a list of individual part numbers.

Part No.	Description	Qty.
VALVE WRENCH	ES	
6810	IDEX Wrench	ea.
MOUNTING BRA	CKET ACCESSORIES	
7160	Mounting Panel	ea.
7160-010	Valve Angle Bracket	ea.
7160-029	Ring Stand Mounting Bracket	ea.
VALVE BRACKET		
M-615-1	Mounting Bracket for IDEX Health & Science Switching Valves	ea.
M-615-2	Mounting Bracket for IDEX Health & Science Injection and Selection Valves	ea.
REPLACEMENT I	FITTINGS	
7770-039	Ferrules for 1/8" OD Tubing	25-pk
7770-040	Ferrules for 1/8" Tubing	50-pk
7770-044	Ferrules for 1/16" OD Tubing	25-pk
7770-124	O-rings for 1/16" OD Tubing	25-pk







FLOW REGULATING VALVES

Our Flow Regulating Valves include specifically designed valves that are used to control or stop the flow of a stream and are ideal for use if your application involves low frequency of use or demands operator control. A variety of types and styles of valves allow you to manage directional flow. In addition, we offer replacement cartridges for all of our flow regulating valves.

- 135 CHECK VALVES
- **140** MICRO-SPLITTER VALVES
- MICRO-METERING VALVES 141
- SHUT-OFF VALVES 142



Inline Cartridge Check Valves

- Low cracking pressures
- Less than 150 μL internal volume
- Materials of construction: PEEK; perfluoroelastomer (CV-3001); gold-plated stainless steel spring (CV-3001); ethylene propylene (CV-3011); and stainless steel spring (CV-3011)

Our cartridge-style Inline Check Valves are designed to limit flow to one direction. These assemblies withstand system pressures of 1,000 psi (69 bar). The cracking pressures for the Inline Check Valve Cartridges are 1.5 psi (0.1 bar) for the CV-3001 and 3 psi (0.2 bar) for the CV-3011. Tolerance on the cracking pressure for CV-3001 is \pm 0.5 psi (0.03 bar) and \pm 1.5 psi (0.1 bar) on CV-3011.



Part No.	Description	Includes	Swept Volume	Qty.
INLINE CART	RIDGE CHECK VALVES			
CV-3000	Inline Check Valve Assembly for 1/16" OD tubing	(1) CV-3001, (2) XP-215	96 μL	ea.
CV-3001	Inline Check Valve Cartridge for CV-3000		91 μL	ea.
CV-3010	Inline Check Valve Assembly for 1/8" OD tubing	(1) CV-3011, (2) XP-315	100 μL	ea.



Add back-flow protection to any 1/4-28 flat-bottom port

- > 15 psi (1 bar) and 3 psi (0.2 bar) cracking pressure versions
- > Excellent chemical resistance
- Materials of construction: PEEK; PCTFE; perfluoroelastomer; PTFE (CV-3301 and CV-3302); stainless steel (CV-3301 and CV-3302); or gold-plated stainless steel (CV-3315 and CV-3316)

Standard 1/4-28 Inline Check Valves

Connect these Inline Check Valves to any 1/4-28 flat-bottom port. Then thread your 1/4-28 flat-bottom fitting into the check valve to connect the tubing. Once installed, the springactuated sealing system eliminates back flow, helping to prevent upstream contamination or damage. In addition, the unique design of this product eliminates the additional tubing cuts and connections required to install conventional inline check valves.









Standard, Outlet 1/4-28 FB Male to 1/4-28 FB Female 3 psi (0.2 bar) cracking pressure



- 1/4-28 Inline Check Valves and Non-Metallic Check Valves with 1/4-28 flat-bottom ports (next page) can be used with any 1/4-28 Flangeless, Super Flangeless™, and VacuTight™ fitting on pages 45, 39, and 42, respectively, of the Fittings Chapter.
- Micro-Volume Inline Check Valves and Non-Metallic Check Valves with 10-32 coned ports (next page) can be used with any 10-32 polymer Fingertight or SealTight™ fitting on page 36. Connect capillary tubing using the optional ferrules listed on page 35 or the NanoTight™ Fittings and Tubing Sleeves on page 37.

Part No.	Description	Cracking Pressure	Qty.
STANDARD 1/4	-28 INLINE CHECK VALVES		
CV-3301	Inlet Check Valve, 1/4-28 FB, M to 1/4-28 FB, F*	15 psi (1 bar)	ea.
CV-3302	Outlet Check Valve, 1/4-28 FB, M to 1/4-28 FB, F*	15 psi (1 bar)	ea.
CV-3315	Inlet Check Valve, 1/4-28 FB, M to 1/4-28 FB, F*	3 psi (0.2 bar)	ea.
CV-3316	Outlet Check Valve, 1/4-28 FB, M to 1/4-28 FB, F*	3 psi (0.2 bar)	ea.

^{*} M = Male (external) threads: F = Female (internal) threads: C = Coned: FB = Flat-Bottom



Non-Metallic Check Valves

- > Cracking pressure of 8 psi (0.6 bar)
- Excellent chemical resistance
- Materials of construction: PEEK and perfluoroelastomer, suitable for biological applications

Our Non-Metallic Check Valves are biocompatible and delivers a low cracking pressure. With a swept volume of only 7.4 μ L, our Check Valve is perfect for applications where low flow path volume is critical, such as delivery to lab-on-a-chip, single-cell analysis and microor nano-LC post-column derivatization. Once installed, this check valve helps prevent back flow and the potential for contamination or damage to sensitive upstream equipment.

10-32 Micro-Volume Inline Check Valves

With a swept volume of only $7.4~\mu$ L, our 10-32 Micro-Volume Inline Check Valves are perfect for applications where low flow path volume is critical, such as delivery to lab-on-achip, single-cell analysis and micro- or nano-LC post-column derivatization. Once installed, this check valve helps prevent back flow and the potential for contamination or damage to sensitive upstream equipment.



Micro-Volume Inline 10-32 C Female to 10-32 C Female



Check valves are specified by:

- Cracking Pressure: the pressure required for the valve to open in the direction of the arrow.
- Maximum Pressure: the maximum pressure the valve can experience in the reverse direction without leaking backwards.
- **Back Pressure Created:** the amount of back pressure generated by the check valve with 50 mL/min room temperature water flowing in the direction of the arrow.



	Swept Volume	Thru-Hole	Max. Pressure Rating	Back Pressure Created	Cracking Pressure Tolerance
TANDARD 1/4	-28 FB				
/-3301, /-3302	20 μL	0.020" (0.50 mm)	2,000 psi (138 bar)	45 psi (3.1 bar)	± 5 psi (0.34 bar)
/-3315, /-3316	16 μL	0.020" (0.50 mm)	2,000 psi (138 bar)	10 psi (0.7 bar)	± 1.5 psi (0.10 bar)
ONMETALLIC	10-32 CONED IV	IICRO-VOLUME			
/-3500	7.4 µL	0.010" (0.25 mm)	3,000 psi (207 bar)	25 psi (1.7 bar)	± 5 psi (0.34 bar)

10-32 Micro-Volume Inline Check Valves

Part No.	Description	Cracking Pressure	Qty.
NONMETALLI	C 10-32 MICRO-VOLUME INLINE CHECK VALVE		
CV-3500	Inlet/Outlet Check Valve, 10-32 C, F to 10-32 C, F*	8 psi (0.6 bar)	ea.
* $M = Male$ (external) threads; $F = Female$ (internal) threads; $C = Coned$; $FB = Flat-Bottom$			

Non-Metallic Check Valves (Cont.)

- > Low cracking pressure of 1 psi (0.07 bar)
- Multiple configurations for different applications
- > Excellent chemical resistance
- Materials of construction: PEEK and perfluoroelastomer

APPLI

APPLICATION NOTE

- The CV-3320 or CV-3321 style can be connected to any 1/4-28 flat-bottom port for trouble-free back flow protection.
- When using a pump after the analytical column, consider placing a CV-3330 Check Valve after the column to prevent fluid from the post-column pump from flowing backwards through the column. This product also serves as an excellent nonmetallic alternative to our CV-3010 (page 135) in sparging applications where the mobile phase may be corrosive to the stainless steel or ethylene propylene components inside the CV-3010 assembly.
- The CV-3335 Inlet and CV-3336 Outlet Check Valves allow tubing larger than 1/16" OD (up to 1/8") to be connected into a 10-32 coned internal port. Use both of these check valves when attaching a larger-volume sample loop to an analytical-scale injection valve. This setup limits the flow of the sample into the loop to one direction, minimizing back flow and sample carry-over.
- The CV-3340 is useful in virtually any high pressure fluid pathway using 1/16" or smaller OD tubing, where limiting the direction of flow is desirable.

1/4-28 & 10-32 Inline Check Valves

Our 1/4-28 & 10-32 Non-Metallic Inline Check Valves provide excellent backflow protection for sensitive equipment along with outstanding chemical resistance guaranteed by the PEEK polymer and perfluoroelastomer construction. Metal-free composition makes these check valves perfect for use with corrosive fluids or biological samples.

These check valves function well up to moderately-high pressure applications. Low internal volume also allows them to be used in areas where flow path volume is important; however, higher flow rates can pass through with minimal pressure drop.





SPECIFICATIONS & DETAILS

	Swept Volume	Max. Pressure Rating	Back Pressure Created	Cracking Pressure Tolerance
CV-3320, CV-3321	37 μL	2,000 psi (138 bar)	30 psi (2.1 bar)	± 0.5 psi (0.03 bar)
CV-3330	34 µL	2,000 psi (138 bar)	30 psi (2.1 bar)	± 0.5 psi (0.03 bar)
CV-3335, CV-3336	49 µL	2,000 psi (138 bar)	30 psi (2.1 bar)	± 0.5 psi (0.03 bar)
CV-3340	34 µL	2,000 psi (138 bar)	30 psi (2.1 bar)	± 0.5 psi (0.03 bar)
CV-3322, CV-3323	49 µL	2,000 psi (138 bar)	30 psi (2.1 bar)	± 0.5 psi (0.03 bar)
CV-3324, CV-3325	182 μL	2,000 psi (138 bar)	30 psi (2.1 bar)	± 0.5 psi (0.03 bar)



NOTE

Upon initial use — or following a period of extended inactivity — the cracking pressure for these check valves may be somewhat higher than the stated cracking pressure.

1/4-28 & 10-32 Inline Check Valves

Part No.	Description	Cracking Pressure	Thru-Hole	Qty.
NONMETA	LLIC 1/4-28 AND 10-32 INLINE CHECK VALVES			
CV-3320	Inlet Check Valve, 1/4-28 FB, M to 1/4-28 FB, F*	1 psi (0.07 bar)	0.020" (0.50 mm)	ea.
CV-3321	Outlet Check Valve, 1/4-28 FB, M to 1/4-28 FB, F*	1 psi (0.07 bar)	0.020" (0.50 mm)	ea.
CV-3322	Inlet Check Valve, 1/4-28 FB, M to 1/4-28 FB, F*	1 psi (0.07 bar)	0.040" (1.0 mm)	ea.
CV-3323	Outlet Check Valve, 1/4-28 FB, M to 1/4-28 FB, F*	1 psi (0.07 bar)	0.040" (1.0 mm)	ea.
CV-3324	Inlet Check Valve, 1/4-28 FB, M to 1/4-28 FB, F*	1 psi (0.07 bar)	0.060" (1.60 mm)	ea.
CV-3325	Outlet Check Valve, 1/4-28 FB, M to 1/4-28 FB, F*	1 psi (0.07 bar)	0.060" (1.60 mm)	ea.
CV-3330	Inlet/Outlet Check Valve, 1/4-28 FB, F to 1/4-28 FB, F*	1 psi (0.07 bar)	0.020" (0.50 mm)	ea.
CV-3335	Inlet Check Valve, 1/4-28 FB, F to 10-32 C, M*	1 psi (0.07 bar)	0.020" (0.50 mm)	ea.
CV-3336	Outlet Check Valve, 1/4-28 FB, F to 10-32 C, M*	1 psi (0.07 bar)	0.020" (0.50 mm)	ea.
CV-3340	Inlet/Outlet Check Valve, 10-32 C, F to 10-32 C, F*	1 psi (0.07 bar)	0.020" (0.50 mm)	ea.

^{*} M = Male (external) threads: F = Female (internal) threads: C = Coned: FB = Flat-Bottom



Quick-Stop Luer Inline Check Valve

- Check valve protection with luer convenience
- > Remains open when engaged
- Materials of construction: PEEK, perfluoroelastomer, and gold-plated stainless steel spring



- > 1/4-28 Inline Check Valves and Non-Metallic Check Valves with 1/4-28 flat-bottom ports (next page) can be used with any 1/4-28 Flangeless, Super Flangeless™, and VacuTight™ fitting on pages 45, 39, and 42, respectively, of the Fittings Chapter.
- Micro-Volume Inline Check Valves and Non-Metallic Check Valves with 10-32 coned ports (page 137) can be used with any 10-32 polymer Fingertight or SealTight™ fitting on page 36. Connect capillary tubing using the optional ferrules listed on page 35 or the NanoTight™ Fittings and Tubing Sleeves on page 37.

The Quick-Stop Luer Check Valve is designed to provide inline luer connect/disconnect convenience without the mess and hazard of spills. Just connect the valve assembly to your inline tubing using standard 1/4-28 flat-bottom fittings (see pages 39 – 46). The check valve is automatically opened once the luer connection is engaged, allowing flow in either direction. Disconnecting the luer union causes the check valve to close. Please see the "Application Note" on this page for specific ideas regarding use of this valve.



APPLICATION NOTE

Inlet Solvent Reservoir:

Quickly change your solvent on the low pressure end of an HPLC system, while preventing potentially hazardous spills! Just install a Quick-Stop Luer Check Valve Assembly between your solvent reservoir and the pump, with the valve towards the bottle. The valve will prevent solvent leakage from the line coming from the reservoir, while the check valves in your pump prevent spills from the line leading to the pump. With both lines still full of solvent, this system also helps reduce the need to reprime your pump.

FIA Sample Injection:

The Quick-Stop Luer Check Valve provides a practical means to introduce a sample into FIA and other low pressure systems, when used in conjunction with a P-612 Pressure Relief Valve Tee (page 148). Simply connect the Tee into the appropriate flow path line with the included fittings and thread the P-697 Quick-Stop Luer Valve onto the 1/4-28 male end of the Tee. Sample can then be introduced conveniently by using a standard luer-tipped syringe. The check valve is automatically opened when the syringe is attached and closed when the syringe is removed.

Post Column Derivitization:

For post-column derivitization, place a CV-3000 Inline Check Valve on the effluent side of your column to prevent derivatizing agents from flowing backwards and poisoning the column. Placement on the post-column reagent line will also prevent mobile phase from contaminating the reagent if the auxiliary pump fails.

Helium Sparging Tank Protection:

Try the CV-3010 Assembly, designed specifically for degassing (sparging) lines to prevent solvent backup if the sparging gas runs out. This check valve will help prevent potential solvent cross-contamination and damage to the gas regulating valve.

Part No.	Description	Includes	Swept Volume	Qty.
QUICK-STOP	LUER CHECK VALVE			
P-696	Quick-Stop Luer Check Valve Assembly	(1) P-697, (1) P-655	127 μL	ea.
P-697	Quick-Stop Luer Check Valve		107 μL	ea.
P-699	Bulkhead Quick-Stop Luer Valve	(1) nut/lock washer set	107 μL	ea.



Micro-Splitter Valves

- > For interfacing LC-MS systems
- Adjustable split stream flow rates
- Versions for up to 800 psi (55 bar) and up to 4,000 psi (276 bar)



APPLICATION NOTE


- With an incoming flow rate of 1 mL/min using room temperature water and equal pressures on both outlet lines, the minimum split flow rate is 2 µL/min for the standard micro-splitter valves and 4.8 µL/min for the high pressure micro-splitter valves.
- All Micro-Splitter Valves have been tested at flow rates to 100 mL/min, with a maximum resulting pressure drop of only 45 psi (3.1 bar) when the valve is fully opened.

Our Micro-Splitter Valves are designed to accurately split and control a low-flow stream off a single incoming supply.

Choose between 1/4-28 flat-bottom and 10-32 coned threaded versions.

The High Pressure Micro-Splitter Valves are designed to operate successfully up to 4,000 psi (276 bar) and the standard Micro-Splitter valves are pressure rated to 800 psi (55 bar).

The Graduated Valve offers many of the benefits and features of Micro-Splitter Valves, plus the ability to adjust and set the split flow to repeatable settings. This allows documentation of settings and the resulting flow rates for easier method development. The graduations also make it easier to employ the valve in a system used to run multiple analyses that require different split flow rates.

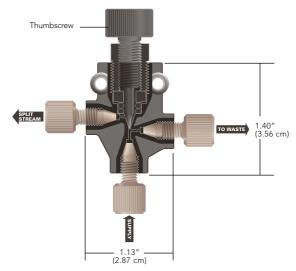

P-451 Standard Micro-Splitter Valve



SPECIFICATIONS & DETAILS

Part No.	Valve Type	Threads	Internal Volume ¹ (closed/fully open)	Max. Operating Pressure
P-450	Standard	1/4-28	2.1 / 4.1 μL	800 psi (55 bar)
P-451	Standard	10-32	1.2 / 2.8 µL	800 psi (55 bar)
P-460S, T	High Pressure	10-32	1.2 / 2.8 µL	4,000 psi (276 bar)
P-470	High Pres. Graduated	10-32	1.2 / 2.8 µL	4,000 psi (276 bar)
4 mm			(a a a a a a a a a a a a a a a a a a a	

¹ The supply and waste port thru-holes have IDs of 0.020" (0.50 mm). The ID for the split-stream port thru-hole is 0.020" (0.50 mm) in standard versions; in capillary versions it is 0.010" (0.25 mm).


Part No.	Description	Includes	Qty.		
MICRO-SPLIT	TTER VALVES				
P-450	Standard, 1/4-28, Biocompatible	(3) XP-235	ea.		
P-451	Standard, 10-32, Biocompatible	(3) F-120	ea.		
P-460S	High Pressure, 10-32, with Stainless Steel Needle	(3) F-120	ea.		
P-460T	High Pressure, 10-32, with Titanium Needle	(3) F-120	ea.		
GRADUATED	GRADUATED MICRO-SPLITTER VALVES				
P-470	High Pressure Graduated, 10-32, with Stainless Steel Needle	(3) F-120	ea.		
* Use with the N	MicroTight Tubing Sleeves, found on page 52.				



Micro-Metering Valves

- > Flow rates as low as 3.5 μL/min*
- ▶ 1/4-28 flat-bottom and 10-32 coned designs available
- **▶ Materials of construction:** PEEK, PTFE
- * At 1.0 mL/min incoming flow rate with room temperature water.

PEEK Micro-Metering Valve

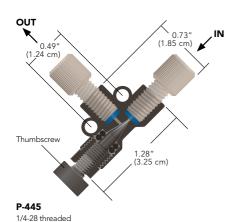

For fine control of fluid flow rates, Micro-Metering Valves can reduce outgoing flow to as low as 3.5 µL/min*. These needle valves are perfect for use with peristaltic pump fluid-transfer applications, mass spectrometry, and fraction collection.

Our Micro-Metering Valves can also be used to regulate gas flow in helium sparging lines and as a flow-dependent variable back pressure regulator. For flow independent regulation of back pressure, please see page 147.

Flow path materials are PEEK polymer and PTFE. All versions of this valve have 0.020'' (0.50 mm) thru-holes.

APPLICATION NOTE

Back Pressure Considerations

The Micro-Splitter Valves are designed to work when both effluent flow path pressures are nearly identical. However, the split flow path will often have higher back pressure than the waste flow path, making it hard to achieve any split flow at all. There are two possible solutions. Place a back pressure regulator (page 143) on the waste flow path that is equal to or slightly greater than the pressure on the split flow path. Or, switch the two effluent pathways such that the split flow pathway is attached to the "waste" port on the valve and the waste flow pathway is attached to the "split" port on the valve. (Please Note: This second method may result in a loss of adjustment sensitivity.)

Multi-Column and Detector Systems

Does your work require analyses with multiple columns and detectors that use the same mobile phase? If so, install one of our High Pressure Micro-Splitter Valves after your injector. A single injection can then be split to two separate columns and detector systems, at two different flow rates. This economical set-up eliminates the need for an additional pump and injector valve, while allowing data to be obtained simultaneously.

Post-Detector Interfacing

Use a Standard Micro-Splitter Valve to route fluid exiting an initial detector to other devices, such as a mass spectrometer and a fraction collector. The valve will split and reduce the flow rate to that required for MS interfacing, while diverting the remainder of the flow to the collector (a back pressure regulator may also be required for this set up, available on page 143).

Other Applications

These valves are also suited for other applications, such as adapting a standard HPLC system to handle microbore analyses. For more information and plumbing diagrams for this application and those listed above, please contact your local distributor or IDEX Health & Science directly.

Part No.	Material	OD Tubing	Thru-hole	Volume*	Includes	Qty.
MICRO-METERII	NG VALVES					
P-445	PEEK, Black	1/16"	0.020" (0.50 mm)	7.7 µL	(2) XP-230	ea.
P-446	PEEK, Black	1/16"	0.020" (0.50 mm)	7.2 µL	(2) F-120	ea.
P-447	PEEK, Black	1/8"	0.020" (0.50 mm)	7.7 µL	(2) XP-330	ea.
* Maximum internal	volume, with valve fully open.					



- > Biocompatible, all-polymer flow path
- Available for 1/16" and 1/8" OD tubing
- > Pressure rated to 500 psi (34 bar)


Shut-Off Valves

Stop a flow stream quickly with IDEX Health & Science biocompatible Shut-Off Valves. The bodies are manufactured from either PEEK or ETFE, and both versions feature a PCTFE rotor, making them highly resistant to chemical attack. The blue colorant used in some valve configurations has proven not to leach out with common HPLC solvents.

Connect semi-rigid or rigid tubing, such as PEEK, stainless steel or fluoropolymer, with the 1/4-28 Flangeless Fittings provided. Soft tubing, such as PharMed® or Tygon® (see page 14), may be connected to these valves using our 1/4-28 barbed adapters, found on page 90.

P-783

ETFE Shut-Off Valve

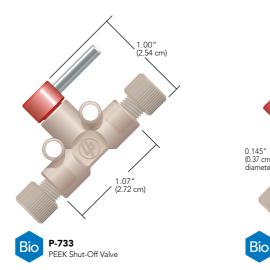

Part No.	Material	OD Tubing	Thru-hole	Internal Volume*	Includes	Qty.
SHUT-OFF VALVES, BIC	COMPATIBLE					
P-721	ETFE, Natural	1/8"	0.040" (1.0 mm)	10.0 μL	(2) P-335, (2) P-300N	ea.
P-732	PEEK, Natural	1/16"	0.020" (0.5 mm)	2.5 µL	(2) XP-235	ea.
P-733	PEEK, Natural	1/8"	0.040" (1.0 mm)	10.0 μL	(2) XP-335	ea.
P-782	ETFE, Blue	1/16"	0.020" (0.5 mm)	2.5 µL	(2) XP-235	ea.
P-783	ETFE, Blue	1/8"	0.040" (1.0 mm)	10.0 μL	(2) XP-335	ea.
* Maximum internal volume,	, with valve fully open.					





BACK PRESSURE REGULATORS

Back Pressure Regulators (BPR) are designed to enhance system performance through outgassing prevention and improved pump check valve efficiency. It includes 5 and 20psi assemblies (replacement cartridges not available), a variety of pressure rated cartridges and assemblies, PEEK and stainless steel BPR holders, high pressure adjustable BPR for pressure between 2000 and 5000psi and ultra-low volume BPRs set to 100 and 500psi.


- 144 ULTRA-LOW VOLUME BACK PRESSURE REGULATOR
- **145** BACK PRESSURE REGULATOR ASSEMBLIES
- 146 BACK PRESSURE REGULATOR HOLDERS
- 147 BACK PRESSURE REGULATOR CARTRIDGES
- **148** PRESSURE RELIEF VALVES



Ultra-Low Volume Back Pressure Regulators (BPRs)

- Wetted flow path materials: PEEK, perfluoroelastomer, and ETFE
- Available pressure settings of 100 or 500 psi (7 or 34 bar)
- Low swept volume of only 6 μL

Our Ultra-Low Volume Back Pressure Regulators (BPRs) were developed to minimize swept volume, which is especially important for multi-detector applications. With a maximum swept volume of only 6 μ L*, it is nearly impossible to detect these BPRs as part of your fluid pathway. To minimize the swept volume added to your flow path, we recommend trimming the length of the attached tubing. And because the flow path is completely polymeric, you are assured of biocompatibility.

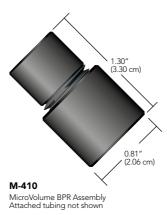

Please Note: Our Ultra-Low Volume Back Pressure Regulators cannot be used as check valves due to their unique internal design. Try our Micro-Volume Inline Check Valve on page 137.

^{*} The maximum internal swept volume listed above is for the back pressure regulator only and does not include the volume of the attached tubing lines



SPECIFICATIONS & DETAILS

	Back Pressure Setting psi (bar)	Flow Rate Recommendations	Recommended Pressure Range psi (bar)	1/16" OD Tubing
M-410	100 2 (7)2	Optimal: 100 µL–1 mL/min Max.: 4 mL/min	40–150 (3–10)	PEEK, 0.010" ID
M-412	500 ² (34) ²	Optimal: 100 µL–1 mL/min Max.: 4 mL/min	250–525 (17–36)	PEEK, 0.010" ID
M-420	100 ³ (7)³	Optimal: 3–8 mL/min Max.: 10 mL/min	40–150 (3–10)	PEEK, 0.020" ID

¹ All data generated using water at room temperature.
² Set at a flow rate of 0.5 mL/min.

Part No.	Description	Pressure Setting	Tubing OD	Includes	Swept Volume	Qty.
ULTRA-LOW VOLUME BPRs						
M-410	Low Flow	100 psi (7 bar)	1/16"	XP-230	6 μL	ea.
M-412	Low Flow	500 psi (34 bar)	1/16"	XP-230	6 μL	ea.
M-420	High Flow	100 psi (7 bar)	1/16"	XP-230	6 μL	ea.

³ Set at a flow rate of 5 mL/min.



Flangeless Nut

Flangeless Ferrule

Cartridge Holder

Color-Coded End Cap


Each BPR Assembly includes a preset BPR Cartridge and IDEX Health & Science fittings for 1/16" OD tubing.


Cartridge


Back Pressure Regulator Assemblies

Choose from our line of Biocompatible and Stainless Steel BPR Assemblies, each complete with a replaceable, factory preset cartridge (except the 5 and 20 psi versions). Our BPR Assemblies create incremental back pressures ranging from 5 to 1,000 psi (0.3

to 69 bar). The Biocompatible BPR Assemblies feature a PEEK holder; polymer-based fittings; biocompatible BPR cartridges and wrenches for tightening. Stainless Steel BPR Assemblies feature the same biocompatible BPR cartridges with a 316 stainless steel holder and polymer fittings.

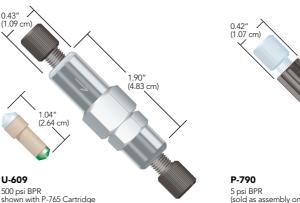

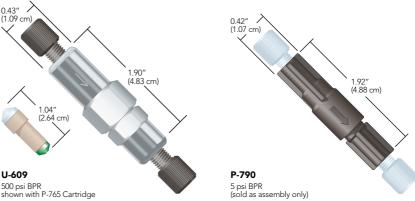

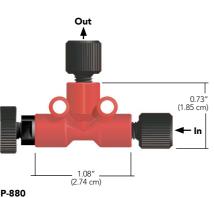

High Pressure Adjustable BPR Includes One-Piece Fingertight Fittings for 1/16" OD tubing

High Pressure Adjustable Back Pressure Regulator

Materials of construction: PEEK, perfluoroelastomer, and PTFE

The biocompatible P-880 High Pressure Adjustable BPR offers the flexibility to adjust your system back pressure between 2,000 and 5,000 psi (138 and 345 bar), independent of the flow. Only 10% fluctuation in pressure generally occurs with flow rates of 0.1–10 mL/min. Lower or higher flow rates will lead to greater fluctuations in pressure. To achieve the desired back pressure setting, simply turn the thumbscrew while monitoring your system pressure. Because this product creates such high back pressure, please check system component specifications prior to using to avoid damaging any sensitive components.

Part No.	Pressure Setting	Holder Material	Includes	Swept Volume	Qty.	
BPR ASSEMBL	IES					
P-790	5 psi (0.3 bar)	PEEK	(2) XP-215	134 µL	ea.	
P-791	20 psi (1.4 bar)	PEEK	(2) XP-215	134 µL	ea.	
P-785	40 psi (2.8 bar)	PEEK	(1) P-761, (2) XP-215	131 μL	ea.	
P-786	75 psi (5.2 bar)	PEEK	(1) P-762, (2) XP-215	131 μL	ea.	
P-787	100 psi (7 bar)	PEEK	(1) P-763, (2) XP-215	131 μL	ea.	
P-788	250 psi (17 bar)	PEEK	(1) P-764, (2) XP-235	102 μL	ea.	
P-789	500 psi (34 bar)	PEEK	(1) P-765, (2) P-250, (2) LT-115	96 μL	ea.	
P-455	1,000 psi (69 bar)	PEEK	(1) P-796, (2) P-250, (2) LT-115	89 µL	ea.	
U-605	40 psi (2.8 bar)	SST	(1) P-761, (2) XP-201	129 µL	ea.	
U-606	75 psi (5.2 bar)	SST	(1) P-762, (2) XP-201	129 µL	ea.	
U-607	100 psi (7 bar)	SST	(1) P-763, (2) XP-201	129 µL	ea.	
U-608	250 psi (17 bar)	SST	(1) P-764, (2) XP-201	99 μL	ea.	
U-609	500 psi (34 bar)	SST	(1) P-765, (2) XP-201	93 μL	ea.	
U-610	750 psi (52 bar)	SST	(1) P-795, (2) P-250, (2) LT-115	91 μL	ea.	
HIGH PRESSU	HIGH PRESSURE ADJUSTABLE BPR ASSEMBLY					
P-880	2,000–5,000 psi (138-345 ba	ar)	(2) F-120BLK	9 μL	ea.	



Back Pressure Regulator Holders

P-465 PEEK and U-469 Stainless Steel BPR Holders work with any of our replacement BPR Cartridges. Each holder comes with fittings for 1/16" OD tubing (see below). The U-469 Holder is surface-treated to prevent galling, a potential problem with large, threaded metal parts.

Please Note: These Back Pressure Regulator Holders are designed to allow each cartridge to operate at its stated pressure setting when tightened to 20 in-lbs. of torque. To approximate this level of torque, first finger tighten the Holder, then tighten an additional 1/8-1/4 turn with the supplied wrenches.





- Using PEEK tubing and supplied fittings. Using stainless steel tubing and supplied fittings

Part No.	Pressure Setting	Holder Material	Includes	Swept Volume	Qty.
BPR HOLDERS					
P-465	Biocompatible BPR	PEEK	(2) P-250, (2) LT-115	7 μL	ea.
U-469	High Pressure BPR	SST	(2) F-300	4 μL	ea.



Back Pressure Regulator Cartridges



APPLICATION NOTE

Small gas bubbles often form as solvent moves from the high pressure of an HPLC column to the low pressure environment leading to the detector. This outgassing can cause erratic baseline readings and loss of sensitivity. Placing a BPR (usually a 40–100 psi) after the detector provides an excellent, low-cost method for reducing this problem by maintaining enough back pressure on the mobile phase to keep gases dissolved in solution.

A back pressure regulator can also be used as a pump preload for low and fluctuating pressure applications. Many of today's pumps require a steady back pressure to function properly. Install an IDEX Health & Science BPR (usually 500–1,000 psi) between the pump and the injector to enhance pump performance.

Caution: Do not exceed the maximum operating pressure of your system — please refer to the operating manuals for your system components before choosing the appropriate BPR.

- Proven outgassing protection
- > Flow-independent pump preload for greater pump efficiency
- ▶ 5 to 1,000 psi cartridges and assemblies available

Back Pressure Regulators are designed to enhance system performance through outgassing prevention and improved pump check valve efficiency.

IDEX Health & Science back pressure regulators include:

- ▶ 5 and 20 psi assemblies (replacement cartridges not available)
- **)** 40, 75, 100, 250, 500, 750, and 1,000 psi cartridges and assemblies
- PEEK and stainless steel BPR holders
- ▶ High pressure adjustable BPR for pressures between 2,000 and 5,000 psi
- Ultra low volume BPRs set to 100 and 500 psi (page 144)
 For flow control options try the Micro-Metering Valves found on page 141.

Back Pressure Regulator Replacement Cartridges

Materials of construction: PEEK, ETFE, perfluoroelastomer, and gold-plated stainless steel

These replacement cartridges will operate in any of the standard BPR holders shown on this page. These cartridges create back pressures from 40 to 1,000 psi (2.8 to 69 bar)—all independent of flow except as noted below.

The recommended operating flow rate range for our BPR Cartridges is 0.1 mL–10 mL/min. Within this range, the amount of back pressure created by the BPR Cartridges and Assemblies will not vary more than $\pm 10\%$. Lower or higher flow rates may result in larger pressure fluctuations.



P-761 40 psi BPR Cartridge



P-796 1,000 psi BPR Cartridge

COLOR CODING

Part No.	Pressure Setting	Body	End-Cap	Swept Volume	Qty.
BPR CARTRIDGES					
P-761	40 psi (2.8 bar)	Tan	Blue	125 μL	ea.
P-762	75 psi (5.2 bar)	Tan	Yellow	125 μL	ea.
P-763	100 psi (7 bar)	Tan	Red	125 μL	ea.
P-764	250 psi (17 bar)	Tan	White	95 μL	ea.
P-765	500 psi (34 bar)	Tan	Green	89 μL	ea.
P-795	750 psi (52 bar)	Black	Blue	87 μL	ea.
P-796	1,000 psi (69 bar)	Black	Green	83 µL	ea.




Pressure Relief Valves

Prevent system over-pressurization

Our Pressure Relief Valves are ideal for preventing system over-pressurization. These products protect system components by diverting fluid flow automatically when inline pressure exceeds the set limit. Choose between preset 100 psi (7 bar) and 5 psi (0.3 bar) assemblies, both shipped with Flangeless Fittings. The 100 psi version is a good, general purpose valve, while the 5 psi version is perfect for protecting syringe and peristaltic pump systems. The void volume of both relief valves is low due to the small 0.020" (0.50 mm) thru-holes in the valve tee body.

If you wish to have the Pressure Relief Valve open at a different pressure than 5 or 100 psi, simply combine one of the other replacement Back Pressure Regulator Assemblies listed on page 145 with the P-612 Pressure Relief Valve Tee. Choose the P-612S for larger bore tubing and higher flow applications.

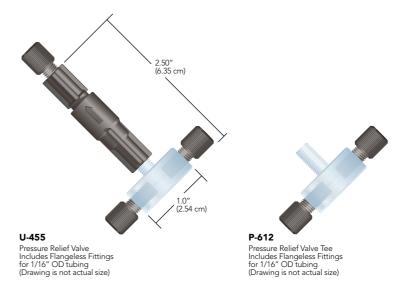

Part No.	Description	Pressure Setting	Tubing OD	Includes	Swept Volume	Qty.
PRESSURE RELIEF VAL	PRESSURE RELIEF VALVES					
U-455	Pressure Relief Assembly	5 psi (0.3 bar)	1/16"	XP-201	148 µL	ea.
U-456	Pressure Relief Assembly	100 psi (7 bar)	1/16"	XP-201, wrenches	139 µL	ea.
P-612	Pressure Relief Tee Only		1/16"	XP-201	14 μL	ea.
P-612S	Pressure Relief Tee Only		3/16"	XP-201	348 µL	ea.

-	
-	
-	
-	



Degassers

Degassers improve fluidic instrument precision and reliability by removing dissolved gases from fluids before they outgas and form problem causing bubbles. Three main types of bubble removing products are available. AF based degassers offer the widest range of chemical compatibility and are used to eliminate retention shifts and baseline fluctuations. Silicone based degassers offer the highest flow rate capabilities for water based systems such as diagnostic and life science instrumentation to improve dispense accuracy and reliability. Poridex based products provide rapid bubble remove for locations where bubble introduction cannot be avoided.

151 DEBUBBLERS

154 DEGASSING SYSTEMS



Debubblers



APPLICATION NOTE

- Liquid handling
- IVD
- > HPLC/UHPLC
- O₂ and CO₂ removal

In medical analyzers, bubbles interfere with critical volumetric reagent dispenses and cause sample failures, wasting time and money. Because bubbles adhere to nearly every part of a dispensing system, high velocity or induced turbulent flow is often used to displace and discharge bubbles from the flow stream and into a waste area. These alternative processes waste reagents and are time consuming, unpredictable, and may additionally require designing the system to recognize bubbles are present. Regardless of how the systems are designed, aqueous systems will always be subject to the laws of physics that cause out-gassing during changes in fluid temperature, pressure, or chemicals mixture. In fluid applications like these, debubblers are the optimal solution to capture and remove formed bubbles to prevent sample dispense inaccuracies, and degassing is ideal to prevent downstream bubble formation from recurring.

Remove Bubbles, Dissolved Gas, or Both!

Dissolved gases and bubbles in system liquids cause dispense volume anomalies in many instruments, negatively affecting both dispense precision and analytical accuracy. Now you have a choice of components for actively removing bubbles with or without also removing dissolved system gases. Online Vacuum Degassing offers operating convenience, high efficiency and low operating costs compared to other common degassing technologies.

Debubbler/Degasser

Combines Vacuum Degassing with Active Bubble Removal

- Improves instrument performance reduces downtime due to bubble formation.
- > Fewer false positives due to reduction of partial reagent dispenses.
- Easily integrates into any pump, degassing tray, or stand-alone degassing application.
- Designed for use with water based solutions with no surfactants. Active degassers are recommended for other solutions.

Active Debubbler

Remove Bubbles in Fluid Stream Before or After the Pump

- Improves instrument performance reduces downtime due to bubble formation.
- > Fewer false positives due to reduction of partial reagent dispenses.
- Easily integrates into any pump, degassing tray, or stand-alone degassing application.

Transfer-Line Degasser

Removes Dissolved Gases During Fluid Transfer


- > Eliminates baseline fluctuations for improved detector sensitivity.
- Coaxial design reduces number of connections, improves reliability.
- > Single lumen design increases degassing reliability.




	ACTIVE DEBUBBLER	DEBUBBLER/ DEGASSER	TRANSFER-LINE DEGASSER
Perfect for applications that require dissolved gas like oxygen for reaction kinetics	~		
Improves dispense precision by capturing and removing bubbles	*	*	
Eliminates false positives and reduces reagent waste by improving instrument performance	~	*	
Easily integrates into fluidic path	*	*	✓
Creates stable instrument performance across system and environmental fluctuations	~	*	✓
Prevents the formation of bubbles downstream of the degasser		*	✓
Eliminates fluctuations for improved detector sensitivity and accuracy by preventing bubble formation		*	✓
Improves fluidic system reliability because coaxial design reduces the number of connections			✓
Flexible design can be implemented as transfer line in new instruments or existing instruments that don't have space available			✓
Minimizes fluidic system internal volumes to reduce reagent cost			~


Debubblers (Cont.)

Overall Dimensions Please note: These drawings are not actual size.

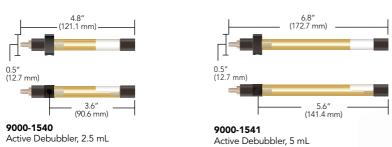

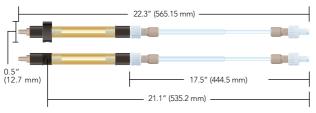

9000-1544 Debubbler / Degasser, 2.5 mL

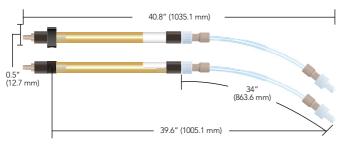


9000-1545 Debubbler/Degasser, 5 mL




9000-1549Transfer-Line Debubbler, 1.1 meter


implementations

Typical Debubbler Implementation

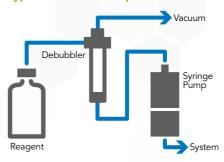

Transfer-Line Degasser Implementation

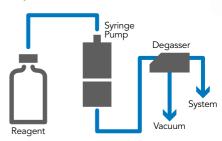

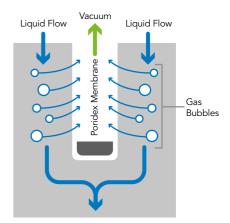


Gas bubbles are actively removed from a flowing liquid stream by vacuum via the PORIDEX membrane.




SPECIFICATIONS (ALL PLATFORMS)

•			
	ACTIVE DEBUBBLERS	DEGASSER/DEBUBBLERS	200 KPA (30 PSI) @ 25 °C
Bubble Removal (volume of air removed/min @ 10 mL/min H₂O)	Up to 30 cc	Up to 30 cc	N/A
Degassing Efficiency [†] @ 1 mL/min H₂O	N/A	$2.5 \text{ mL size: } 36\% \text{ O}_2 \text{ removal, } 5.0 \text{ mL size: } 55\% $ $\text{O}_2 \text{ removal}$	$<$ 4 ppm dissolved O_2 at 5 mL/min
Membrane Material	PORIDEX®	PORIDEX	PORIDEX
Wetted Materials	PORIDEX, Polyolefin, FEP, ETFE, Ultem®	PORIDEX, Polyolefin, FEP, ETFE, Ultem	PORIDEX, Polyolefin, FEP, ETFE
Solvent Compatibility	Solutions > 50% aqueous. Not compatible with detergent concentrations > 0.05%.	Solutions > 50% aqueous. Not compatible with detergent concentrations > 0.05%.	Solutions > 50% aqueous. Not compatible with detergent concentrations > 0.05%.
Standard Bubble Trap Volume	2.5 / 5.0 mL	2.5 / 5.0 mL	N/A
Transfer-Line Volume	N/A	2.5 / 5.0 mL	< 4 mL
Max. Operating Pressure	200 kPa (30 psi) @ 25 °C	200 kPa (30 psi) @ 25 °C	200 kPa (30 psi) @ 25 °C
Max. Operating Temperature	40 °C	40 °C	40 °C
Recommended Vacuum Level	Minimum 16 kPa absolute	Minimum 16 kPa absolute	Minimum 16 kPa absolute
Liquid Connection	1/4-28 fitting system	1/4-28 fitting system	1/4-28 fitting system
Vacuum Connection	Tubing vacuum port(s) for 1/8" (3 mm) ID elastomeric tubing	Tubing vacuum port(s) for 1/8" (3 mm) ID elastomeric tubing	Tubing vacuum port(s) for 1/8" (3 mm) ID elastomeric tubing
Pressure Drop	0.8 mm Hg / mL / min (assumes laminar flow and viscosity of 1 cP)	0.8 mm Hg / mL / min (assumes laminar flow and viscosity of 1 cP)	0.8 mm Hg / mL / min (assumes laminar flow and viscosity of 1 cP)
† Debubbling / degassing efficiency of	can be optimized based on flow rate, fluid to be o	degassed, and gas to be removed.	



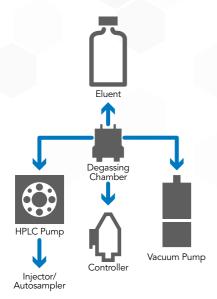

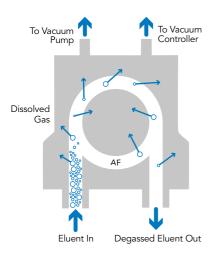

Dissolved gases are actively removed from a flowing liquid stream by vacuum via the IDEX Health & Science AF® membrane.

APPLICATION NOTE

Why Degas Your Mobile Phase?

Dissolved air in HPLC mobile phases can result in flow rate instability and baseline disturbance.

Flow rate instability: Non-degassed mobile phase can outgas in the pump head, causing bubbles to be formed and trapped inside the head or check valves. These bubbles can cause flow disturbances and pressure fluctuations, resulting in flow rate instability.

Baseline disturbance: As the mobile phase passes through the column, it experiences a large pressure drop. Non-degassed mobile phase can outgas due to this pressure differential, causing air bubbles to form. Air bubbles passing through or lodging in the flow cell cause detection disturbances, exhibited as baseline noise.

Why Use a Degassing System?

Helium sparging is a common means of degassing HPLC solvents. This method has its drawbacks, however. Helium tanks are expensive and bulky, and solvent backup and contamination are concerns. In addition, helium sparging can change the composition of a premixed mobile phase over time, due to the difference in the evaporation rates of mobile phase components.

In contrast, the IDEX Health & Science Degassing System has none of these drawbacks, and it is extremely fast and efficient at removing dissolved gases — more efficient than helium sparging or PTFE-based degassing systems.

Tubing Connections

We recommend ETFE tubing (page 27) be used to limit regassing of mobile phase between the degasser and your pump. ETFE is recommended because of its superior impermeability to gases (compared to PTFE, FEP, and PFA tubing). Applicable flangeless fittings for 1/8" OD tubing are found on page 45.

GPC and HFIP Applications

Standard degassing chambers, with PEEK bulkhead unions, are <u>not recommended</u> for GPC applications or for use with HFIP (hexafluoroisopropanol). Special GPC "hardened" versions are available. Please contact us for more information.



Degassing tubing is flexible and therefore can be coiled to shorten the overall length or used to transfer the fluid within an instrument to the next desired location.

Debubblers

Part No.	Description	Standard Bubble Trap Size	Transfer Line Length	Internal Volume	Max Bubble Capacity	Qty.
	ER SERIES – AVAILABLE STA					
9000-1540	2.5 mL Active Debubbler	2.5 mL	_	2.5 mL	2.5 mL	ea.
9000-1541	5 mL Active Debubbler	5 mL	_	5 mL	5 mL	ea.
9000-1544	2.5 mL Debubbler/Degasser	2.5 mL	17.5" (444.5 mm)	2.5 mL in transfer line + 2.5 mL in bubble trap	2.5 mL	ea.
9000-1545	5 mL Debubbler/Degasser	5 mL	34" (863.6 mm)	5 mL in transfer line + 5 mL in bubble trap	5 mL	ea.
9000-1549	1.1 m Transfer-Line Degasser	_	1.1 m (43")	4 mL	N/A	ea.





- Analytical and Prep scale models
- Ultra-high degassing efficiency
- > Low volume, easy to prime
- Patented control eliminates baseline fluctuations
- Inert flow path
- > 5+ year lifetime

Full Stand Alone egassing Systems

No Troubles with Bubbles Anymore!

- prevent bubble formation in your Fluidic System with BIOTECH DEGASi line of Degassers.

Dissolved gasses in a fluidic system can often cause troubles. When the pressure or the temperature changes, the dissolved gasses can form bubbles which affect the accuracy, precision and performance of your equipment. On-line degassing is a very efficient way of removing dissolved gasses from the liquid and preventing bubble formation.



DEGASi® PLUS CLASSIC compared to the old DEGASi® CLASSIC











GENERAL
PURPOSE

GPC ELUENT

LIMITED **VOLUME**


SEMI PREP

Systec AF™ Internal Volume	480 µl	480 µl	100 μΙ	925 µl
No of Channels	F-Î	FΪ	F-Î	F-Î
Biocompatible Flow Path	Υ	N ÃÛ¢ }¢^å	Υ	Υ
Approximate Max Flow Per Channel (ml/min)	3	3	0.5	6
Inner Dimension Flow Path (mm)	1.14	1.14	1.14	1.14



Part Number	Number of Channels	Internal Volume
DEGASi® PLUS CLASSIC		
0003-6351-A	1	480 µl
0003-6352-A	2	480 µl
0003-6353-A	3	480 µl
0003-6354-A	4	480 µl
0003-6355-A	5	480 µl
0003-6356-A	6	480 µl
DEGASi® PLUS GPC		
0003-6621	1	480 µl
0003-6622	2	480 µl
0003-6623	3	480 µl
0003-6624	4	480 µl
0003-6625	5	480 µl
0003-6626	6	480 µl
DEGASi® PLUS MICRO		
0003-6351-S	1	100 µl
0003-6352-S	2	100 µl
0003-6353-S	3	100 µl
0003-6354-S	4	100 µl
0003-6355-S	5	100 µl
0003-6356-S	6	100 µl

Part Number	Number of Channels	Internal Volume
DEGASi® PLUS SEMI-PREF		
0003-6351-L	1	925 µl
0003-6352-L	2	925 µl
0003-6353-L	3	925 µl
0003-6354-L	4	925 µl
0003-6355-L	5	925 µl
0003-6356-L	6	925 µl
DEGASi® PREP		
0001-2053	2	5.3
DEGASi® PREP+		
0001-0120	1	23 ml
0001-0220	2	23 ml
0001-0420	4	23 ml
DEGASi® COMPACT		
0004-2285-W	2	285 µl
0004-2485-W	4	285 µl
0004-2685-W	6	285 µl
DEGASi® HIGH FLOW		60 ml
HF 500-S Stand Alone	1	60 ml
HF 500-A OFM Version	1	60 ml

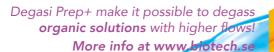




















PREP

ORGANIC SOLVENTS

COMPACT APPLICATION

HIGH FLOW

Systec AF™ Internal Volume	5.3 / 13.8 ml	23 ml	285 µl	60 ml (silicone)
No of Channels	2	1,2 or 4	2, 4, 6	1
Biocompatible Flow Path	N/Y	Υ	Υ	Υ
Approximate Max Flow Per Channel (ml/min)	15 / 50	75-100	2	500
Inner Dimension Flow Path (mm)	1.91	N/A	0.89	N/A

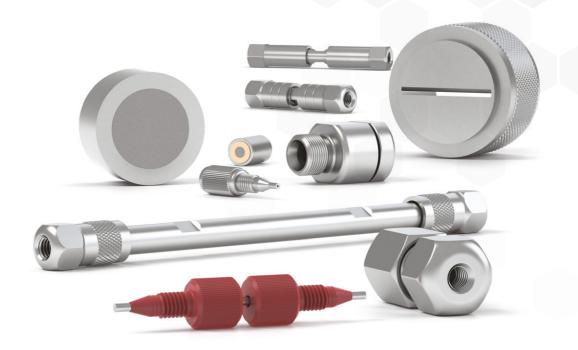

Column Hardware

We offer an extensive line of HPLC and UHPLC Column Hardware that has been optimized to enable selectivity, efficiency, and high-quality separation performance in your flow paths. To browse the full line of our column hardware portfolio please visit: www.idex-hs.com/column-hardware.html

- Biocompatible materials for LC and UHPLC columns
- Accessories for column protection and packing

Our line of column hardware includes protective accessories and connection products that enhance column functionality. Our columns come in a variety of sizes and materials to meet your system requirement.

157

GUARD HARDWARE

158

GUARD COLUMNS



Guard Hardware

Our Guard Hardware portfolio offers a variety of guard cartridges and holders to meet your system requirements. These include guards, guard holders, analytical cartridges, and tools designed to protect your system's valuable columns and help ensure high performance and excellent retention time.



Prep Scale Guard Holders

- 21.2 mm and 30 mm ID column protection
- Improves plate count and symmetry
- New anti-rotation feature aids guard holder assembly

Iso-Prep[™] Guard is a guard cartridge system designed to protect valuable prep columns. It offers superior column protection for adsorptive samples and a proven sample distribution mechanism via a precision machined holder. Iso-Prep Guard is ideal for protecting prep columns with no degradation of peak shape or plate count.

The high-performance guard protects columns in two ways. First, it acts as a filter, trapping particles in the frits. Second, when the guard cartridge is packed with the same material as the prep column, it removes compounds that irreversibly adsorb to the packing material.



Tools


Our tools includes the Iso-Prep Guard Scraper that is specifically designed as a tool used to dress a cartridge upon packing in order to set the frit at the correct depth.



- ▶ 60:40 Acetonitrile:Water
- **▶** 50 mL/min

- > Kromasil 10 µm C18
- ▶ Backpressure: 100 psi (7 bar)

10 μ C18 100 x 21.2 mm, 60:40 Acetonitrile:Water, 20 mL/min

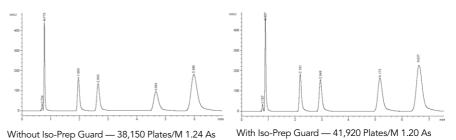

Part No.	Description	Qty.						
PREP SCALE GUARD	PREP SCALE GUARD HOLDERS							
9197-P	lso-Prep Guard Holder	ea.						
9197-P-AR	lso-Prep Guard Holder, 21.2mm, Anti Rotation	ea.						
9197-20	lso-Prep Guard Cartridge, 21.2 mm x 1 cm, 1 Frit	ea.						
9196-P	lso-Prep Guard Holder, 30 mm	ea.						
9196-P-AR	lso-Prep Guard Holder, 30 mm, Anti Rotation	ea.						
9196-20	lso-Prep Guard Cartridge, 30 mm x 1 cm, 2 μm Frit	ea.						
8083-MOD	lso-Prep Guard, 30 mm Frit, 2 µm	ea.						
TOOLS								
9197-S	lso-Prep Guard Finishing Tool	ea.						
9196-S	lso-Prep Guard Finishing Tool, 30 mm	ea.						




Guard Columns

We offer wide selection of pre-packed, cartridge style guard columns in addition to pack-it-yourself hardware. There are options for micro flow applications as well as analytical-scale applications. These guard columns can be immediately implemented into a system for your convenience.



Guard Column Kits

- > 100% biocompatible flow path
- Pressure rated to 4,000 psi (276 bar)
- > Wetted materials are Titanium and PEEK
- Reusable holder complete with fingertight fittings

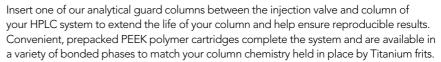

The C-270 Stainless Steel Guard Column Holder is engineered for high-pressure applications to 4,000 psi (276 bar). Each of these holders is surface treated to prevent galling*, a potential problem with threaded metal parts.

The flow path of the C-270 Guard Column Holder is biocompatible. Each comes complete with fittings for 1/16" OD tubing, and can be used with any of the C-28X or C-7XX guard column cartridges listed on the following pages.

* Galling is a form of "cold welding." When two fittings manufactured from the same metal are wrench-tightened too tightly, they can "weld" together, making it virtually impossible to separate the two components.




Microbore Guard Columns

- Ideal for Microbore HPLC
- Easily dry packed (or slurry packed with adapter)
- Made of PEEK polymer and stainless steel

This ultralow volume guard column (1.0 mm ID x 2 cm length) is ideal for narrow-bore chromatography. The unpacked guard column allows you to exactly match the chemistry of your column, resulting in optimum column protection. The total packing volume of 16.2 µL ensures maximum column efficiency and analytical column protection.

Frits often become plugged before a guard column is contaminated. The two $0.5\,\mu m$ frits included with this guard column can be changed in minutes. Optional 2 µm frits may be purchased separately (C-408).

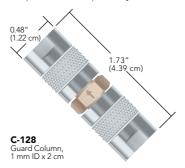



Packing Material Specifications: The cartridges on this page are packed with 5 µm or 10 µm basedeactivated 80 Å spherical silica.




FLUIDICS



Analytical Guard Columns

Easy to pack

The C-130B is our most popular guard column. HPLC users find this column easy to pack and extremely economical. This narrow-bore short column (2.0 mm ID x 2 cm length) creates only a slight pressure increase with virtually no detectable theoretical plate loss when used with a 3 mm ID or larger column. The 2 µm frits are easy to change, prolonging the life of the guard column. With only 62 µL packing volume per guard column, a 3 g bottle of packing material will pack about 30 guard columns.

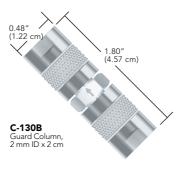




Guard Column Holders

- > 10 mm ID column protection
- Convenient cartridge system
- Easy to pack

The internal volume of our semi-prep guard column is just 780 µL, which only requires approximately 1.50 g of packing material — ideally the same material used in your semi-prep column. The C-1000 Holder will hold to high pressures, and is specially treated to prevent galling*. Use standard 10-32 coned fittings (not included) to connect your 1/16" OD tubing.

* Galling is a form of "cold welding." When two fittings manufactured from the same metal are wrench-tightened too tightly, they can "weld" together, making it virtually impossible to separate the two components.

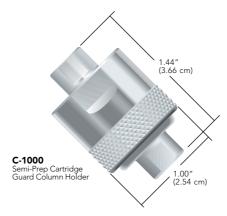



C-1000 Semi-Prep Guard ColumnFittings, tubing, and column shown are not included.

Guard Columns (Cont.)



Accessories

Our accessories include packing adapters and seal, packing funnels, and other useful products to work with your column systems.



Guard Column Cartridges

Our biocompatible Guard Column Cartridges are conveniently prepackaged and are offered in multiple packs. We offer these cartridges in a variety of bonded phases to match your column chemistry. These cartridges are reusable and economical.



APPLICATION NOTE

Why Use A Guard Column?

A guard column can increase the life of your analytical column significantly. Use a guard column with the same packing as your column — it will act as a chemical filter, removing strongly retained materials in your sample that might otherwise contaminate your analytical column. And, it is more economical to replace a guard column cartridge than to buy a new analytical column.



All Guard Columns featured on this page include 10-32 Coned threads. Use any of the 10-32 coned fittings starting on page 32 to connect tubing to these guard columns.



Signs Indicating the Guard Column Needs to be Changed

- System pressure build-up
- > Faster than usual retention times
- Reduced resolution

Guard Columns

Part No.	Description	Includes	Qty.
GUARD COLUMN KITS			
C-281	2.0 mm ID C18 Cartridges (6-pk) with (1) C-270 Assembly		ea.
MICROBORE GUARD COLU	JMNS		
1.0 mm ID x 2 cm Unpacked			
C-128	Guard Column	(2) C-128-31	ea.
C-128-31	0.5 µm Stainless Steel Replacement Frit		ea.
C-408	2 µm Stainless Steel Replacement Frit		ea.
ANALYTICAL GUARD COLU	JMNS		
2.0 mm ID x 2 cm Unpacked			
C-130B	Guard Column	(2) A-100	ea.
C-130-20	Packing Funnel		ea.
A-100	2 µm Stainless Steel Replacement Frit		ea.
A-103	0.5 µm Stainless Steel Replacement Frit		ea.
BIOCOMPATIBLE GUARD C	COLUMN HOLDERS		
10 mm ID x 1 cm			
C-270	High Pressure, Stainless Steel, with (2) F-200 Fittings, Biocc	ompatible	ea.
GUARD COLUMN HOLDERS	S		
10 mm ID x 1 cm			
C-1000	Semi-Prep Cartridge Guard Column Holder		ea.
C-1035	Semi-Prep Cartridge		ea.
ACCESSORIES			
C-128-40	Slurry Packing Adapter		ea.
C-130-40	Slurry Packing Adapter		ea.
C-1030	Threaded Frit Cap with 2 µm Stainless Steel Frit		ea.
GUARD COLUMN CARTRID	GES		
2.0 mm ID x 1 cm, 10 μm Silica			
C-282	Reversed Phase C18		10-pk
4.3 mm ID x 1 cm, 5 µm Silica			
C-751	4.3 mm ID C18 Cartridges (6-pk) with (1) C-270 Assembly		ea.
C-752	Reversed Phase C18		10-pk

Educational Materials

Help is only a click away to access a wealth of product and application support information, including:

Chemical Compatibility

Use our compatibility charts to discover what materials are best suited for your application.

www.idex-hs.com/chemical-compatibility

Chromatography Forum

Browse or join this group that is devoted to the analytical and scientific communities.

www.idex-hs.com/chromatography-forum

Conversion Tools

Find useful Conversion Tools that allow you to easily convert temperatures, measurements, pressure ratings, and other functions.

www.idex-hs.com/conversions

FAOs

Explore our FAQs section to assist in finding the answers to issues you may be facing.

www.idex-hs.com/faqs www.semrock.com/technical-faq.aspx

Fittings Primer

Learn the complexities and functions of nutand-ferrule fittings using our Fittings Primer.

www.idex-hs.com/fittings-primer

Glossary of Terms

Our Glossary of Terms provides a convenient list of definitions and/or explanations to language that is applicable in our industry.

www.idex-hs.com/glossary

HPLC Center

Visit our HPLC Center to learn about high- and ultra-high performance liquid chromatography.

www.idex-hs.com/hplc-center

Injectors Troubleshooting

This guide presents common problems when using manual sample injection valves.

www.idex-hs.com/injectors

IP Rating

Use our IP Rating chart to discover the environmental care of our electrical equipment.

www.idex-hs.com/ip-rating

Materials Guide

Our Materials Guide allows you to view material property and solvent compatibility.

www.idex-hs.com/materials-guide

Optical Filter Plotting Tool

Dur SearchLight™ optical filter plotting tool allows you to easily select elements for your system, visualize the spectral properties, and quickly calculate critical parameters such as signal brightness, autofluorescence levels, and signal-to-noise ratio.

searchlight.semrock.com

Polymer Information

Study the Polymer Information chart to learn about analytical solvent compatibility for HPLC.

www.idex-hs.com/polymer-information

Pump School

Jump over to Pump School, an introduction to rotary positive-displacement (PD) pumps.

www.idex-hs.com/pump-school

Standard Port Drawings

We offer downloadable schematics of our Standard Port Drawings that allow you to easily create your own internal ports using our fittings.

www.idex-hs.com/standard-port-drawings

Tech Tips

Examine various topics that will maximize the precision performance of our products.

www.idex-hs.com/tech-tips

Webinars & Videos

IDEX Health & Science hosts ongoing educational webinars on specific topics, and has also created a product- and process-focused video library.

www.idex-hs.com/webinars-and-videos

Literature & Downloads

Our literature and downloads section provides a wide selection of valuable information about our products and our company. You can access and download a variety of support literature, manuals, designed to assist in the use of our products and services, including:

Operation Manuals & Instruction Cards

Product operation manuals and instruction cards are here for your convenience to download and use as reference.

www.idex-hs.com/manuals

Optical Filter Journals and White Paper Articles

Download and share our extensive collection of white papers and journal reprint articles written by our expert PhD staff. Browse varied topics, from coatings for ultrafast lasers, to polarization, super-resolution microscopy, tunable bandpass filters, and more.

www.semrock.com/tech-library.aspx

Optical Filter Technical Notes

- Visit our General Filter information to learn more about the many optical filters we offer, how to find filter orientation, how to measure OD, and more.
- Learn the basics of Fluorescence Filters, including different types of fluorescence filters, applications, uses, and information.
- Our Raman & Laser Optical System Filters information pages give you an in-depth look into what you need to know and how to avoid issues integrating into your Raman & Laser Optical systems.

www.semrock.com/tech-notes.aspx

Optical Systems, Laser, & Shutter Technical Resources

Have a question or need some technical support? Consult our library of documentation, manuals, data sheets and specifications on our entire family of optical products, including:

- > Technical Note
- White papers
- Journal articles

www.mellesgriot.com/store/technical-library

Support Literature

Download informational IDEX Health & Science product and OEM catalogs, brochures, white papers, data sheets and more.

www.idex-hs.com/support-literature

NDEX

Index by Part Number

0001-6352-A	154		
0001-6353-A	154		
0001-6354-A	154		
0001-6355-A	154	1006-999	125
0001-6622	154	1020-05	107
0001-6623	154	1020-20	107
0001-6624	154	1328	54
0001-6625	154	1356	54
0001-6352-S	154	1531	18
0001-6353-S	154	1531B	18
0001-6354-S	154	1532	18
0001-6355-S	154	1533	18
0001-6352-L	154	1534	18
0001-6353-L	154	1535	18
0001-6354-L	154	1536	18
0001-6355-L	154	1538	18
0001-2053	155	1544	18
0001-6482	155	1560	18
0001-6484	155	1561	18
0001-0120	155	1568	18
0001-0220	155	1569	18
0001-0420	155	1570	18
0004-2285	155	1571	18
0004-4285	155	1572	18
0004-6285	155	1574	18
HF.500-S Stand Alone	155	1576	18
HF.500-A OEM version	155	1577	18

1501	10	7020	1.
1581	18 73	7020 7021	1.
1876	73 127	7022	1.
1958-01	78	7023	1.
		7024	1.
1959-01	76		
3000	123	7025	1.
3055-018	130	7026	10
3055-019	130	7027	•
3055-023	130	7028	1.
3055-025	130	7029	1.
3065-018	127	7030	1.
3065-019	127	7030-003	1.
3065-023	127	7030-014	1.
3065-025	127	7030-015	1.
3255	23	7030L	1.
3505	23	7031-05	10
3725-006	125	7031-20	10
3725-018	125	7060	1.
3725-038	123	7060-039	1.
3725-999	125	7060-064	1.
3725i	123	7060-065	1.
3725i-038	123	7060-070	1.
3755	23	7060-074	1
5001-999	125	7112-999	1.
5100-999	125	7123-047	1.
5104-999	125	7123-127	1.
6000-076	66	7123-128	1:
6000-078	61, 130	7123-142	1.
6000-079	61, 130	7123-145	1:
6000-082	33	7123-147	1:
6000-083	33	7123-148	1:
6000-209		7123-149	1.
6000-210		7123-180	1
6000-211		7123-221	1.
6000-251	61, 130	7123-223	1.
6000-254	61, 130	7123-227	13
6000-255	61	7123-548	1.
6000-262	33	7123-550	1:
6000-282	32	7123-568	1.
6000-360	61 22	7125-047	1.
6255	23	7125-054	1.
6505	23	7125-079	1.
6755	23	7125-999	1.
6810	51, 133	7150-999	1.
7000	123	7152-999	1.
7000-016	125	7155-999	1.
7000L	123	7160	1.
7004-999	125	7160-010	1.
7010	123	7160-029	1.
7010-039	125	7410-038	1.
7010-040	125	7410-041	1.
7010-069	125	7410-999	1.
7010-071	125	7413-013	1.
7010-997	125	7501-999	1:
7010-998	125	7502-999	1:
	125	7507.000	1.
7010-999	125	7507-999	1.

7520-035	125	9000-1549	153	62515	23
7520-999	125	9010	123		23
7531-999	125	9010-051	125	62550	23
7650-002		9010-065			23
7725	123	9010-999			23
7725-010		9012			23
7725-999		9013			23
7725i		9055-020			23
7725i-188		9055-021			23
7750-038		9055-022			23
7750-070		9055-023			23
7755-015		9055-024			23
7755-020		9055-025			125
7755-021		9055-026			125
7755-022		9055-027			23
7755-023		9055-028			23
7755-024		9055-029			23
7755-025		9055-033			23
7755-026		9085-05-10			23
7755-027		9085-20-10			23
7755-028		9125-043			23
7755-029		9125-076			23
7770-039	133	9125-082	125		23
7770-040	133	9125-095	125	610015	23
7770-044	133	9125-999	125	610020	23
7770-124	133	9196-20	157	610050	23
7770-229	125	9196-P	157	615010	23
7788-999	125	9196-P-AR	157	615015	23
7900-146	125	9196-S	157	615020	23
7900-179	125	9197-20	157	615050	23
7900-183	125	9197-P	157	617515	23
7900-999	125	9197-P-AR	157	617520	23
7920-999	125	9197-S	157	617550	23
7960-014	125	9650-009	125	620015	23
7960-999	125	9725	123	620020	23
7961-999	125	9725-999	125	620050	23
7970-999		9725i			23
7980-004		9750-021			23
7980-999		9960-002			23
7984-999		31005			23
7986-004		31505			96, 104, 160
7986-999		32510			96, 106
8020		32515			96, 106
8021		32520			96, 160
8022		32550			96
8023		35010			96
					96
8083-MOD		35015			
8125		35020			96
8125-038		35050			96
8125-097		37510			96, 104
8125-098		37515			102
8125-119		37520			102
8125-999		37550	23		102
9000-1540		61005		A-228	102
9000-1541		61505		A-230A	102
9000-1544	153	62005	23	A-231A	102
9000-1545	153	62510	23	A-242	102

4 040	100	A 700	00	E 440	2.4
	102		98		34 35
	102 102		98 98		35 32
	51		98 98		34
	51		98 98		34
			98 98		34 34
	102 102		98 98		34 34
	102				
			98		34
	102		98		32
	104		110		34
	104		88		35, 37, 80
	106		160		35
	106		96, 160		35
	51		160		35
	28		160		34, 132
	28		160		34
	28		160		41, 93
	107		96	=	34
	96, 107		33, 38, 41		54
	96, 107	C-236x	33		54
A-337	96, 107	C-270	160		54
A-343	96	C-281	160		54
A-350	28		160		54
A-355	106	C-401	96		54
A-356	106	C-402	96		54
A-360	107		96		54
A-410	107	C-408	96, 160		54
A-411	107	C-409	96		36
A-420	96	C-420	96, 110		36
A-427	101	C-425	96, 110	F-195x	36
A-428x	104	C-751	160	F-196x	36
A-429x	104	C-752	160		54
A-430	104	C-1000	160	F-226	54
A-431	104	C-1030	160	F-227	54
A-435	101	C-1035	160	F-228	54
A-437	101	CV-3000	135	F-229	54
A-438	101	CV-3001	135	F-230	54
A-440	101	CV-3010	135	F-231	54
A-441	101	CV-3301	136	F-232	54
A-445	101	CV-3302	136	F-233	54
A-446	101	CV-3315	136	F-234	54
A-451	101	CV-3316	136	F-237	54
A-510	107	CV-3320	138	F-238	54
A-520x	101	CV-3321	138	F-239	54
A-522x	101	CV-3322	138	F-240	54
A-550	101	CV-3323	138	F-241	54
	101	CV-3324	138	F-242	54
A-610	111	CV-3325	138	F-243	54
A-620	111		138	F-244	54
	111		138		54
	112		138		54
	112		138		54
	112		137		54
	112		93		43
	112		93		28
	98, 106		93		36
	98, 106		32		35
,		1 100A		· 000/	

F-330x	35	M-548	109	P-207Sx	43
F-331x	35	M-550	110	P-207x	43
F-333Nx		M-560		P-213	41
F-345		M-572	71	P-217	41
F-346	51	M-615-1	133	P-219	41
F-347		M-615-2		P-229x	
F-350x		M-644-03x		P-232	
F-354x		M-645x		P-240x	
F-356x		M-650x		P-246x	
F-364x		M-652x		P-247x	
F-374		M-653x		P-248x	
F-375		M-655x		P-249x	
F-376		M-660x		P-250x	
F-381		MXP7900-000		P-255x	
F-382		MXP7920-000		P-259	
F-384		MXP7960-000		P-259x	
F-385		MXP7970-000		P-260x	
F-386		MXP7980-000		P-268	
F-387		MXP7986-000		P-270x	
F-388		MXP9900-000		P-272x	
				P-273x	
FS-110		MXP9960-000 MXT715-000		P-274x	
FS-115				P-274x	
		MXT715-102		P-275X	
FS-150		MXT715-105			•
FS-175		MXX777-601		P-278	
FS-315		MXX777-603		P-279	
LT-100x		MXX777-605		P-281	
LT-105x		MXX777-612		P-287	
LT-110x		MXX777-616		P-291	
LT-115x		MXX778-605		P-292	
LT-132x		N-123-02		P-295	
LT-200x		N-290		P-296	
LT-210x	•	N-291		P-297	
LT-300		N-333		P-298	
M-120x		OC-802		P-299	
M-121		OC-803		P-300Nx	
M-125		OC-805	•	P-300x	
M-126		OC-815	•	P-307Sx	
M-128		P-100		P-307x	
M-131		P-116		P-309	
M-133		P-120		P-311	
M-134		P-123		P-312x	
M-140x		P-133Nx		P-314	
M-215		P-133x		P-316	
M-250x		P-134		P-317	
M-410		P-135		P-319	
M-412		P-137x		P-321	
M-420		P-140x	41, 56	P-329x	41
M-432		P-141x		P-331	
M-432-03		P-150		P-336	
M-447		P-154	81	P-337x	
M-500		P-155	81	P-340x	47
M-510	110	P-170	70	P-342x	
M-520	108	P-190	81	P-343x	43
M-525	108	P-191	81	P-344x	
M-540		P-200Nx	47	P-347x	
M-547	109	P-200x	43, 47	P-349x	41

P-350x	41	P-612S	148	P-703	85
P-352x	41	P-618	88		69
P-353x		P-619	88		85
P-355x		P-620	85		86
P-357-2x		P-622			86
P-357x		P-623			86
P-359		P-624			66
P-359x		P-625			88
P-360x		P-626			71
P-363Rx		P-627		= •	142
P-366x		P-628			86
P-372x		P-629			
P-373x		P-630			86
P-374x					70
		P-631			70
P-381		P-632			70
P-387		P-633			142
P-399		P-634			142
P-401x		P-635			69
P-407x		P-642			55
P-416		P-645			55
P-416BLK		P-646	93		55
P-416G	34	P-647	93	P-757	93
P-420	, -	P-648	93	P-760	69
P-430	82	P-650	66	P-761	147
P-432	82	P-652	66	P-762	147
P-440	82	P-654	66	P-763	147
P-441	82	P-655	88	P-764	147
P-441N	82	P-656	88	P-765	147
P-445	141	P-657	88	P-767	93
P-446	141	P-658	88	P-770	78
P-447	141	P-659	88	P-771	71
P-450	140	P-660	88	P-772	71
P-451		P-661	88		68
P-455		P-662			68
P-460S		P-663			69
P-460T		P-668			142
P-465		P-669-01			142
P-470		P-670			145
P-482BLK		P-671			145
P-501x		P-672			145
P-507x		P-673			145
P-512		P-675			145
P-513		P-677			145
P-514		P-678			145
P-515		P-680			93
P-520		P-681			147
P-550		P-683			147
P-551		P-684			93
P-553		P-686			93
P-554		P-689			93
P-555	•	P-691			93
P-556		P-692			93
P-558		P-694			89
P-602	85	P-696			89
P-603		P-697	139	P-831	89
P-604	88	P-699		P-835	89
P-612	148	P-702	85	P-836	89

P-837	89	PK-100x	58	U-160	21
P-838		PK-110x		U-161	
P-840		PK-112		U-162	
P-842x		PK-120BLKx		U-163	
P-844x		PK-126Hx		U-164	
P-845-01		PK-126x		U-165	
P-846x		PK-132x		U-320x	
P-849		PK-152x		U-320x	
P-850				U-321x	
		U-101			• • • • • • • • • • • • • • • • • • • •
P-851	•	U-102		U-400x	
P-852		U-103		U-401x	
P-854		U-104		U-402	
P-857		U-105		U-410x	
P-858		U-106		U-411	
P-859		U-107		U-412	73
P-860	93	U-108	21	U-428	72
P-861	93	U-111	21	U-429	72
P-862	93	U-112	21	U-430	72
P-863	93	U-113	21	U-431	72
P-864	93	U-114	21	U-435	73
P-870	93	U-115	21	U-438	73
P-872	93	U-116		U-450x	
P-874		U-117	21	U-455	
P-875		U-118		U-456	
P-880		U-126		U-466	
P-881		U-127		U-466S	
P-882		U-128		U-467R	
		-			
P-885		U-129		U-469	
P-887		U-130		U-500	
P-888		U-131		U-504	
P-889		U-132		U-510	
P-890	68	U-133	21	U-514	66
P-891	68	U-134	21	U-605	145
P-920-01	66	U-135	21	U-606	145
P-930x	42	U-136	21	U-607	145
P-931x	42	U-137	21	U-608	145
P-938x	42	U-138	21	U-609	145
P-940x	42	U-139	21	U-610	145
P-942x	42	U-140	21	U-650x	56
P-945x	42	U-141	21	U-659	66
P-946x	42	U-142	21	U-660x	56
P-948x		U-143		U-665	
P-982BLKx		U-144		U-800	
P-1000		U-145		U-803	
P-1082BLKx		U-146		U-805	
P-1087x		U-147		U-825	
PD715-000				U-1114	
		U-148			
PD715-102		U-149		U-1115	
PD715-105		U-150		U-1116	
PD7900		U-151		U-1117	
PD7920		U-152		U-1120	
PD7960	120	U-153	21	U-1122	
PD7970	120	U-154	21	U-1125	20
PD7980	120	U-155	21	U-1126	20
PD7986	120	U-156	21	U-1128	20
PD9900	120	U-157	21	U-1130	20
PD9960	120	U-158	21	U-1131	20

U-1132	20
U-1133	
U-1140	
U-1141	
U-1142	
U-1143	
U-1145	
U-1146	
U-1148	
UH-402	
UH-427	
UH-432	
UH-436	
UH-630	78
UH-631-01	78
UH-632	76
UH-634	78
UH-700	74
UH-750	74
UH-752	74
UH-753	78
UH-904x	58
UH-906	78
UP-700	69
UPFP-6025070	
UPFP-6025150	
UPFP-6025250	
UPFP-6025350	
UPFP-6025500	
UPFP-6025600	
UPFP-6050070	
UPFP-6050150	
UPFP-6050250	
UPFP-6050350	
UPFP-6050500	
UPFP-6050600	
UPFP-6075070	
UPFP-6075150	
UPFP-6075250	
UPFP-6075350	
UPFP-6075500	
UPFP-6075600	
UPFP-6100070	
UPFP-6100150	
UPFP-6100250	
UPFP-6100350	
UPFP-6100500	49
UPFP-6100600	49
UPFP-6150070	49
UPFP-6150150	49
UPFP-6150250	49
UPFP-6150350	49
UPFP-6150500	49
UPFP-6150600	49
UPFP-6300070	49
UPFP-6300150	49
UPFP-6300250	

UPFP-630035049	VHP-505	104
UPFP-630050049	VHP-506	104
UPFP-630060049	VHP-550	106
UPFS-610007049	VHP-551	106
UPFS-610015049	VHP-555	106
UPFS-610025049	VHP-556	106
UPFS-610035049	VHP-600	55
UPFS-610050049	VHP-700x	59
UPFS-610060049	VHP-920x	60
UPFS-612507049	VHP-1000	51
UPFS-612515049	VHP-1001	60
UPFS-612525049	VHP-2000	51
UPFS-612535049	VHP-3200	60
UPFS-612550049	VHP-4000	51
UPFS-612560049	VHP-9000	51
UPFS-625407049	XF-358x	47
UPFS-625415049	XF-368x	47
UPFS-625425049	XLT-111x	47
UPFS-625435049	XP-130x	56
UPFS-625450049	XP-131x	56
UPFS-625460049	XP-132x	56
V-100D123	XP-136x	56
V-100L123	XP-137x	56
V-100T123	XP-141x	56
V-101D123	XP-201x	47
V-101L123	XP-202x	47
V-101T123	XP-218x	47
V-240123	XP-230x	47
V-241123	XP-235x	
V-340123	XP-238x	
V-341123	XP-245x	
V-447132	XP-286x	47
V-450123	XP-301x	
V-451123	XP-302x	
V-540123	XP-305x	
V-541123	XP-308x	
VHP-200-01x59	XP-315x	47
VHP-200x59	XP-330x	
VHP-320x60	XP-335x	
VHP-325x60	XU-620x	
VHP-500104	XU-655x	
VHP-501104	XU-662x	56



TRADEMARKS & REGISTERED TRADEMARKS

For a complete list of IDEX Health & Science registered trademarks and their proper attribution, please refer to our trademarks list at: www.idex-hs.com/trademarks







For ordering and technical support, please contact:

EUROPE | BIOTECH AB | +46 (0)300-56 91 80 | info@biotechfluidics.com | www.biotechfluidics.com | USA | BIOTECH USA LLC | 612-703-5718 | sales@biotechfluidics.com | www.biotechfluidics.com | JAPAN | BIONIK INC. | +81 545-38-9125 | info@bionikinc.com | www.bionikinc.com